



Bodystat SDK

Documentation for the Bodystat Software
Developers Kit

Table of Contents

Bodystat API	1
Bodystat API - Main	1
C# .NET Wrapper	56
Bodystat API	57
Files	79
BodystatSDK.cs	80
Visual Basic .NET wrapper	80
Classes	81
BodystatAPI Class	81
BodystatAPI Structures	84
BodystatAPI Fields	114
BodystatAPI Methods	115
Structs, Records, Enums	145
BodystatAPI.BSAnalysisMode Enumeration	145
BodystatAPI.BSDeviceFamily Enumeration	146
BodystatAPI.BSDeviceModel Enumeration	147
BodystatAPI.BSError Enumeration	149
BodystatAPI.BSGender Enumeration	151
FAQs	152
Redistribution	154
Samples	155
Bodystat MFC Sample	155
Classes	155
CBodystatMFCDlg Class	156
CBodystatMFCDlg::CBodystatMFCDlg Constructor	158
CBodystatMFCDlg Data Members	159
CBodystatMFCDlg Methods	165
CBodystatResultsDlg Class	185
CBodystatResultsDlg::~CBodystatResultsDlg Destructor	186
CBodystatResultsDlg::CBodystatResultsDlg Constructor	186
CBodystatResultsDlg Data Members	187
CBodystatResultsDlg Methods	194
Bodystat C# WinForms Sample	204

BodystatCSWinForms Namespace	205
Classes	205
MainForm Class	205
ResultsForm Class	209
Files	215
MainForm.cs	216
ResultsForm.cs	216
Bodystat Visual Basic .NET Sample	216
BodystatVBNetWinFormSample Namespace	217
Classes	217
FrmBodystatMain Class	217
FrmBodystatResults Class	237

Index

a

1 Bodystat API

The Bodystat API contains all the functions, definitions, classes and structures implemented by the Bodystat DLL.

The main API documentation is generic and is relevant to all programming languages. For your convenience, function definitions are provided in the documentation for popular languages (C++, C#, Java, Visual Basic & Matlab). However, you can call the API from any language capable of working with standard Windows DLLs based on C exports.

Language specific wrappers for the API are also provided for some languages. However, we regret we are unable to provide wrappers for all languages - we simply don't have the demand or in-house expertise to deal with them all! If you write a wrapper for a language we don't presently support directly we would be more than happy to include it in our SDK assuming you are willing to share your work with the community.

Modules

Name	Description
Bodystat API - Main (see page 1)	The main Bodystat API containing the functions, definitions, classes and structures implemented by the Bodystat DLL. The main API documentation is generic and is relevant to all programming languages. For your convenience, function definitions are provided in the documentation for popular languages (C++, C#, Java, Visual Basic & Matlab). However, you can call the API from any language capable of working with standard Windows DLLs based on C exports.
C# .NET Wrapper (see page 56)	Bodystat API C# .NET wrapper. Provides a wrapper interface to the Bodystat DLL with C# calls and types you can work with directly in your C# .NET program. The wrapper handles the necessary marshalling and interop issues required to interface with the native unmanaged Bodystat DLL.
Visual Basic .NET wrapper (see page 80)	Bodystat API Visual Basic .NET wrapper. Provides a wrapper interface to the Bodystat DLL with VB calls and types you can work with directly in your Visual Basic .NET program. The wrapper handles the necessary marshalling and interop issues required to interface with the native unmanaged Bodystat DLL.

1.1 Bodystat API - Main

The main Bodystat API containing the functions, definitions, classes and structures implemented by the Bodystat DLL.

The main API documentation is generic and is relevant to all programming languages. For your convenience, function definitions are provided in the documentation for popular languages (C++, C#, Java, Visual Basic & Matlab). However, you can call the API from any language capable of working with standard Windows DLLs based on C exports.

Namespaces

Name	Description
Bodystat (see page 2)	This is the main namespace for the Bodystat API. The entire functions within the Bodystat API are enclosed within the 'Bodystat' namespace.

Files

Name	Description
BodystatSDK.h (see page 55)	<p>Bodystat SDK (C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p> <p>Version 1.36 (09/May/2014) Selected download is now passed through to the results form in the C# sample. Fixed error in VB.Net wrapper (preventing download working) in the VB.Net sample. Fixed some missing normals for MST device. 1.35 (10/Mar/2014) Added Cole methods for ECW/ICW calculations used in the MST device. Fixed a memory pointer issue for the detected bodystat device info structure. 1.34 (12/Nov/2013) Exposed Bluetooth API search timeout values to... more (see page 55)</p>

Macros

Name	Description
BODYSTATSDK_API (see page 53)	The following ifdef block is the standard way of creating macros which make exporting from a DLL simpler. All files within this DLL are compiled with the BODYSTATSDK_EXPORTS symbol defined on the command line. This symbol should not be defined on any project that uses this DLL. This way any other project whose source files include this file see BODYSTATSDK_API functions as being imported from a DLL, whereas this DLL sees symbols defined with this macro as being exported.
M_PI (see page 54)	This is macro M_PI.

1.1.1 Bodystat API

This is the main namespace for the Bodystat API.




















The entire functions within the Bodystat API are enclosed within the 'Bodystat' namespace.









Module

Bodystat API - Main (see page 1)

Functions

	Name	Description
	BSAuthenticateBTDevice (see page 8)	Authenticate (pair) a Bodystat device found during a previous call to BSSearchBTDevices (see page 22). Completely silent implementation.

	BSAutoSetupBT (see page 8)	Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown). This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc)). Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the... more (see page 8)
	BSCalculateNormals (see page 10)	Calculate normal values and normal ranges for a given subject and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from BSCalculateResults (see page 10)).
	BSCalculateResults (see page 10)	Calculate complete result set for a given bio-impedance measurement recorded by the Bodystat device.
	BSCloseComport (see page 11)	Close com port.
	BSConnect (see page 11)	Connect to the Bodystat device and ensure the device is responding.
	BSGetBTBodystatDevice (see page 12)	Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices (see page 22).
	BSGetBTStackInfo (see page 13)	Retrieve information about the Bluetooth stack. Useful for diagnostics.
	BSGetDeviceFamily (see page 13)	Determine the device family of a specific model.
	BSGetDeviceModelName (see page 14)	Determine the full product name of a specific model.
	BSGetDeviceModelNameShort (see page 14)	Determine the short product name of a specific model.
	BSGetSdkLibraryVersion (see page 15)	Determine the Bodystat SDK library (DLL) version.
	BSIsBTAvailable (see page 16)	Check if a supported BT radio is available on this computer. Checks for presence of USB Bluetooth dongles and inbuilt Bluetooth modules in the PC. No communication with Bodystat devices takes place during this query.
	BSOpenComport (see page 16)	Open com port to the device ready for communication.
	BSReadCalibrationTime (see page 17)	Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.
	BSReadCurrentTime (see page 17)	Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).
	BSReadModelVersion (see page 18)	Read device model and firmware version of the connected Bodystat device.
	BSReadPrinterAddress (see page 19)	Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).
	BSReadProtocolInfo (see page 19)	Read protocol information from the connected device. Contains version numbers relating to communication protocols in use, data structures and auxiliary information.
	BSReadSerialNumber (see page 20)	Read the unique serial number of the connected Bodystat device.

	BSReadStoredTestData ( see page 21)	Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged. Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.
	BSReadTestBodystat ( see page 21)	Perform a quick bio-impedance test measurement on the connected bodystat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory. May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion. Useful for diagnostics or performing measurements instigated by software.
	BSSearchBTDevices ( see page 22)	Search BT radio for any Bodystat devices (paired or otherwise). Existing paired devices will be counted even if they are switched off or out of range. New non-paired units will only be counted if switched on, in range and ready for communication (on main startup screen).
	BSUnAuthenticateBTDevices ( see page 23)	Unauthenticate (unpair) ALL Bodystat devices from this PC.
	BSWriteCurrentTime ( see page 24)	Reset the current date/time of the internal clock on the connected device. Date/time is automatically set to match the current date/time of the PC
	BSWriteCurrentTimeAs ( see page 24)	Reset the current date/time of the internal clock on the connected device. Date/time is set to a specific value of your choosing.
	BSWritePrinterAddress ( see page 25)	Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt printer pairing function in the device.
	ESCalculateNormals ( see page 26)	Calculate normal values and normal ranges for a given animal and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from ESCalculateResults ( see page 26)).
	ESCalculateResults ( see page 26)	Calculate complete result set for a given bio-impedance measurement recorded by the Equistat device.
	ESReadStoredTestData ( see page 27)	Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged. Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.
	ESReadTestEquistat ( see page 28)	Perform a quick bio-impedance test measurement on the connected equistat device. Normal test process is skipped (user is not prompted to enter any horse information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory. May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion. Useful for diagnostics or performing measurements instigated by software.

Structs, Records, Enums

	Name	Description
	BSAnalysisMode (see page 30)	BSAnalysisMode: Analysis results are available in several formats: BSAnalysisModeBC= Body composition analysis BSAnalysisModeHN= Hydration analysis BSAnalysisModeBoth= Combined body composition and hydration analysis (default)
	BSDeviceFamily (see page 30)	BSDeviceFamily: Enumerated values that group each device / model variant into a given product family. Devices within a given product family share broad capabilities and features. For details regarding capabilities and features refer to the technical documentation for that product family (http://www.bodystat.com)
	BSDeviceModel (see page 32)	BSDeviceModel: Enumerated values that represent the various different models of the Bodystat devices over the years. Suffix: The suffix after the model indicates specific hardware variations regarding connectivity: OPTO indicates model communicates via opto-isolated serial. DIU indicates model communicates via a data interface unit (induction interface) BT indicates model communicates via Bluetooth serial port.
	BSError (see page 33)	BSError: Bodystat Error codes. A list of Bodystat specific error codes.
	BSGender (see page 35)	BSGender: Enumerated values that represent a person's gender.
	BSMeasurement (see page 36)	BSMeasurement: Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).
	BSNormals (see page 38)	BSNormals: Bodystat Normals structure. Holds the normal values or normal ranges for any given test measurement.
	BSRawData (see page 40)	BSRawData: Bodystat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record starting number.
	BSResults (see page 41)	BSResults: Bodystat Results structure. Holds the calculated results for any given test measurement.
	ESMeasurement (see page 44)	ESMeasurement: Equistat Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the horse, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).
	ESNormals (see page 45)	ESNormals: Equistat Normals structure. Holds the normal values or normal ranges for any given test measurement.
	ESRawData (see page 47)	ESRawData: Equistat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record starting number.
	ESResults (see page 48)	ESResults: Equistat Results structure. Holds the calculated results for any given test measurement.
	ImpData (see page 49)	Represents one frequency measurement

Variables









Name	Description
BS_RAWDATA_ARRAYSIZE (see page 52)	Number of measurement records BSRawData (see page 40) structure should reserved space for (currently fixed)
ES_RAWDATA_ARRAYSIZE (see page 52)	Number of measurement records ESRawData (see page 47) structure should reserved space for (currently fixed)


1.1.1.1 Functions

The following table lists functions in this documentation.

Functions

	Name	Description
≡	BSAuthenticateBTDevice (see page 8)	Authenticate (pair) a Bodystat device found during a previous call to BSSearchBTDevices (see page 22). Completely silent implementation.
≡	BSAutoSetupBT (see page 8)	Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown). This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc)). Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the... more (see page 8)
≡	BSCalculateNormals (see page 10)	Calculate normal values and normal ranges for a given subject and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from BSCalculateResults (see page 10)).
≡	BSCalculateResults (see page 10)	Calculate complete result set for a given bio-impedance measurement recorded by the Bodystat device.
≡	BSCloseComport (see page 11)	Close com port.
≡	BSConnect (see page 11)	Connect to the Bodystat device and ensure the device is responding.
≡	BSGetBTBodystatDevice (see page 12)	Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices (see page 22).
≡	BSGetBTStackInfo (see page 13)	Retrieve information about the Bluetooth stack. Useful for diagnostics.
≡	BSGetDeviceFamily (see page 13)	Determine the device family of a specific model.
≡	BSGetDeviceModelName (see page 14)	Determine the full product name of a specific model.
≡	BSGetDeviceModelNameShort (see page 14)	Determine the short product name of a specific model.
≡	BSGetSdkLibraryVersion (see page 15)	Determine the Bodystat SDK library (DLL) version.
≡	BSIsBTAvailable (see page 16)	Check if a supported BT radio is available on this computer. Checks for presence of USB Bluetooth dongles and inbuilt Bluetooth modules in the PC. No communication with Bodystat devices takes place during this query.
≡	BSOpenComport (see page 16)	Open com port to the device ready for communication.
≡	BSReadCalibrationTime (see page 17)	Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.
≡	BSReadCurrentTime (see page 17)	Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).
≡	BSReadModelVersion (see page 18)	Read device model and firmware version of the connected Bodystat device.
≡	BSReadPrinterAddress (see page 19)	Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).

	BSReadProtocolInfo ( see page 19)	Read protocol information from the connected device. Contains version numbers relating to communication protocols in use, data structures and auxiliary information.
	BSReadSerialNumber ( see page 20)	Read the unique serial number of the connected Bodystat device.
	BSReadStoredTestData ( see page 21)	Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged. Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.
	BSReadTestBodystat ( see page 21)	Perform a quick bio-impedance test measurement on the connected bodystat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory. May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion. Useful for diagnostics or performing measurements instigated by software.
	BSSearchBTDevices ( see page 22)	Search BT radio for any Bodystat devices (paired or otherwise). Existing paired devices will be counted even if they are switched off or out of range. New non-paired units will only be counted if switched on, in range and ready for communication (on main startup screen).
	BSUnAuthenticateBTDevices ( see page 23)	Unauthenticate (unpair) ALL Bodystat devices from this PC.
	BSWriteCurrentTime ( see page 24)	Reset the current date/time of the internal clock on the connected device. Date/time is automatically set to match the current date/time of the PC
	BSWriteCurrentTimeAs ( see page 24)	Reset the current date/time of the internal clock on the connected device. Date/time is set to a specific value of your choosing.
	BSWritePrinterAddress ( see page 25)	Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt printer pairing function in the device.
	ESCalculateNormals ( see page 26)	Calculate normal values and normal ranges for a given animal and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from ESCalculateResults ( see page 26)).
	ESCalculateResults ( see page 26)	Calculate complete result set for a given bio-impedance measurement recorded by the Equistat device.
	ESReadStoredTestData ( see page 27)	Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged. Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.

	ESReadTestEquistat (see page 28)	<p>Perform a quick bio-impedance test measurement on the connected equistat device. Normal test process is skipped (user is not prompted to enter any horse information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.</p> <p>May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.</p> <p>Useful for diagnostics or performing measurements instigated by software.</p>
---	----------------------------------	--

1.1.1.1.1 Bodystat::BSAuthenticateBTDevice Function

Authenticate (pair) a Bodystat device found during a previous call to `BSSearchBTDevices` (see page 22). Completely silent implementation.

C++

```
BODYSTATSDK_API BOOL BSAuthenticateBTDevice(HWND hWndParent = NULL, unsigned short iTimeout = 7);
```

C#

```
BODYSTATSDK_API BOOL BSAuthenticateBTDevice(HWND hWndParent, unsigned short iTimeout);
```

Visual Basic

```
Function BSAuthenticateBTDevice(hWndParent As HWND = NULL, iTimeout As unsigned short = 7) As BODYSTATSDK_API BOOL
```

Java

```
BSAuthenticateBTDevice
```

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSAuthenticateBTDevice(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
HWND hWndParent = NULL	Handle of the parent window.
unsigned short iTimeout = 7	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or timesout.

1.1.1.1.2 Bodystat::BSAutoSetupBT Function

Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device

and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown).

This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc).

Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the device knowing only the com port. It is not necessary to setup the device in this way before each use.

New non-paired units will only be found if switched on, in range and ready for communication (on main startup screen). It is advisable to instruct the user to ensure this is the case before calling.

C++

```
BODYSTATSDK_API BOOL BSAutoSetupBT(LPTSTR lpComPort, int iSize, BOOL bReportErrors, HWND hWndParent = NULL);
```

C#

```
BODYSTATSDK_API BOOL BSAutoSetupBT(LPTSTR lpComPort, int iSize, BOOL bReportErrors, HWND hWndParent);
```

Visual Basic

```
Function BSAutoSetupBT(lpComPort As LPTSTR, iSize As Integer, bReportErrors As BOOL, hWndParent As HWND = NULL) As BODYSTATSDK_API BOOL
```

Java

```
BSAutoSetupBT
```

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSAutoSetupBT( , , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
LPTSTR lpComPort	[out] String buffer to receive the device port the device is assigned to.
int iSize	The size of the buffer.
BOOL bReportErrors	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (still not completely silent).
HWND hWndParent = NULL	Handle of the parent window.

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or timesout.

Errors during the process can be optionally shown/hidden. However, a completely silent approach is not possible using this function (instead use BSSearchBTDevices (see page 22) and BSGetBTBodystatDevice (see page 12) or iterate radio handles manually. Then pass desired handle to BSAuthenticateBTDevice (see page 8) for pairing).

Finally be aware this function is only available for supported locales. Otherwise English messages may be displayed to the user. Adopt silent approach discussed previously if you are targeting a non-supported locale.

1.1.1.1.3 Bodystat::BSCalculateNormals Function

Calculate normal values and normal ranges for a given subject and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from BSCalculateResults (see page 10)).

C++

```
BODYSTATSDK_API BSError BSCalculateNormals(const BSMeasurement * pM, const BSResults * pR,
BSNormals * pN);
```

C#

```
BODYSTATSDK_API BSError BSCalculateNormals(BSMeasurement * pM, BSResults * pR, ref
BSNormals pN);
```

Visual Basic

```
Function BSCalculateNormals(pM As BSMeasurement *, pR As BSResults *, ByRef pN As
BSNormals) As BODYSTATSDK_API BSError
```

Java

```
BSCalculateNormals
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSCalculateNormals(, , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
const BSMeasurement * pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
const BSResults * pR	The results calculated for this measurement previously.
BSNormals * pN	[out] Receives the normal values and normal ranges for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Normals vary by subject attributes, measured bio-impedance values and calculated results.

1.1.1.1.4 Bodystat::BSCalculateResults Function

Calculate complete result set for a given bio-impedance measurement recorded by the Bodystat device.

C++

```
BODYSTATSDK_API BSError BSCalculateResults(const BSMeasurement * pM, BSResults * pR,
BSAnalysisMode iAnalysisMode);
```

C#

```
BODYSTATSDK_API BSError BSCalculateResults(BSMeasurement * pM, ref BSResults pR,
BSAnalysisMode iAnalysisMode);
```

Visual Basic

```
Function BSCalculateResults(pM As BSMeasurement *, ByRef pR As BSResults, iAnalysisMode As
BSAnalysisMode) As BODYSTATSDK_API BSError
```

Java

```
BSCalculateResults
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSCalculateResults( , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
const BSMeasurement * pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
BSResults * pR	[out] Receives the complete calculated results for this measurement.
BSAnalysisMode iAnalysisMode	The desired analysis mode (see BSAnalysisMode (see page 30) for more information).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Results vary by subject attributes and measured bio-impedance values.

1.1.1.1.5 Bodystat::BSCloseComport Function

Close com port.

C++

```
BODYSTATSDK_API void BSCloseComport();
```

C#

```
BODYSTATSDK_API void BSCloseComport();
```

Visual Basic

```
Function BSCloseComport() As BODYSTATSDK_API void
```

Java

```
BSCloseComport
```

MATLAB

```
function [BODYSTATSDK_API void] = BSCloseComport()
```

File

BodystatSDK.h (see page 55)

1.1.1.1.6 Bodystat::BSConnect Function

Connect to the Bodystat device and ensure the device is responding.

C++

```
BODYSTATSDK_API BSError BSConnect();
```

C#

```
BODYSTATSDK_API BSError BSConnect();
```

Visual Basic

```
Function BSConnect() As BODYSTATSDK_API BSError
```

Java

```
BSConnect
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSConnect()
```

File

BodystatSDK.h (see page 55)

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have opened the com port by calling BSOpenComPort previously.

1.1.1.1.7 Bodystat::BSGetBTBodystatDevice Function

Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices (see page 22).

C++

```
BODYSTATSDK_API BOOL BSGetBTBodystatDevice(LPTSTR lpDeviceName, int iDeviceNameBufferSize,
LPTSTR lpBDA, int iBDABufferSize, LPTSTR lpComPort, int iComBufferSize, unsigned short
iTimeout = 7);
```

C#

```
BODYSTATSDK_API BOOL BSGetBTBodystatDevice(LPTSTR lpDeviceName, int iDeviceNameBufferSize,
LPTSTR lpBDA, int iBDABufferSize, LPTSTR lpComPort, int iComBufferSize, unsigned short
iTimeout);
```

Visual Basic

```
Function BSGetBTBodystatDevice(lpDeviceName As LPTSTR, iDeviceNameBufferSize As Integer,
lpBDA As LPTSTR, iBDABufferSize As Integer, lpComPort As LPTSTR, iComBufferSize As Integer,
iTimeout As unsigned short = 7) As BODYSTATSDK_API BOOL
```

Java

BSGetBTBodystatDevice

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSGetBTBodystatDevice(, , , , , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
LPTSTR lpDeviceName	[out] String buffer to receive the product name of the Bodystat device.
int iDeviceNameBufferSize	Size of the device name buffer.
LPTSTR lpBDA	[out] String buffer to receive the Bluetooth hexadecimal address (known as the BDA - unique to each BT module in a device, like a MAC address).
int iBDABufferSize	Size of the BDA buffer.
LPTSTR lpComPort	[out] String buffer to receive device/port the unit is assigned to.
int iComBufferSize	Size of the com port buffer.

unsigned short iTimeout = 7	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).
-----------------------------	--

Returns

TRUE if it succeeds, FALSE if it fails.

1.1.1.1.8 Bodystat::BSGetBTStackInfo Function

Retrieve information about the Bluetooth stack. Useful for diagnostics.

C++

```
BODYSTATSDK_API BOOL BSGetBTStackInfo(LPTSTR lpBuffer, int iBufferSize);
```

C#

```
BODYSTATSDK_API BOOL BSGetBTStackInfo(LPTSTR lpBuffer, int iBufferSize);
```

Visual Basic

```
Function BSGetBTStackInfo(lpBuffer As LPTSTR, iBufferSize As Integer) As BODYSTATSDK_API  
    BOOL
```

Java

```
BSGetBTStackInfo
```

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSGetBTStackInfo(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
LPTSTR lpBuffer	[out] String buffer to receive the information.
int iBufferSize	Size of the buffer.

Returns

true if it succeeds, false if it fails.

Remarks

Note: Information can only be provided for supported BT radios.

1.1.1.1.9 Bodystat::BSGetDeviceFamily Function

Determine the device family of a specific model.

C++

```
BODYSTATSDK_API BSDeviceFamily BSGetDeviceFamily(BSDeviceModel iModel);
```

C#

```
BODYSTATSDK_API BSDeviceFamily BSGetDeviceFamily(BSDeviceModel iModel);
```

Visual Basic

```
Function BSGetDeviceFamily(iModel As BSDeviceModel) As BODYSTATSDK_API BSDeviceFamily
```


Java

BSGetDeviceFamily

MATLAB

```
function [BODYSTATSDK_API BSDeviceFamily] = BSGetDeviceFamily()
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BSDeviceModel iModel	The desired model of interest.

Returns

The device family the model belongs to.

1.1.1.1.10 Bodystat::BSGetDeviceModelName Function

Determine the full product name of a specific model.

C++

```
BODYSTATSDK_API BSError BSGetDeviceModelName(BSDeviceModel iModel, LPTSTR lpBuffer, int iBufferSize);
```

C#

```
BODYSTATSDK_API BSError BSGetDeviceModelName(BSDeviceModel iModel, LPTSTR lpBuffer, int iBufferSize);
```

Visual Basic

```
Function BSGetDeviceModelName(iModel As BSDeviceModel, lpBuffer As LPTSTR, iBufferSize As Integer) As BODYSTATSDK_API BSError
```

Java

BSGetDeviceModelName

MATLAB

```
function [BODYSTATSDK_API BSError] = BSGetDeviceModelName(, , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BSDeviceModel iModel	The desired model of interest.
LPTSTR lpBuffer	[out] String buffer to receive the name.
int iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.1.1.1.11 Bodystat::BSGetDeviceModelNameShort Function

Determine the short product name of a specific model.

C++

```
BODYSTATSDK_API BSError BSGetDeviceModelNameShort(BSDeviceModel iModel, LPTSTR lpBuffer, int iBufferSize);
```

C#

```
BODYSTATSDK_API BSError BSGetDeviceModelNameShort(BSDeviceModel iModel, LPTSTR lpBuffer,
int iBufferSize);
```

Visual Basic

```
Function BSGetDeviceModelNameShort(iModel As BSDeviceModel, lpBuffer As LPTSTR, iBufferSize
As Integer) As BODYSTATSDK_API BSError
```

Java

```
BSGetDeviceModelNameShort
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSGetDeviceModelNameShort(, , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BSDeviceModel iModel	The desired model of interest.
LPTSTR lpBuffer	[out] String buffer to receive the name.
int iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.1.1.1.12 Bodystat::BSGetSdkLibraryVersion Function

Determine the Bodystat SDK library (DLL) version.

C++

```
BODYSTATSDK_API BSError BSGetSdkLibraryVersion(DWORD * pMajorSdkVer, DWORD * pMinorSdkVer);
```

C#

```
BODYSTATSDK_API BSError BSGetSdkLibraryVersion(ref DWORD pMajorSdkVer, ref DWORD
pMinorSdkVer);
```

Visual Basic

```
Function BSGetSdkLibraryVersion(ByRef pMajorSdkVer As DWORD, ByRef pMinorSdkVer As DWORD)
As BODYSTATSDK_API BSError
```

Java

```
BSGetSdkLibraryVersion
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSGetSdkLibraryVersion(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
DWORD * pMajorSdkVer	[out] Receives the major version number of the Bodystat SDK DLL.
DWORD * pMinorSdkVer	[out] Receives the minor version number of the Bodystat SDK DLL.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Caller should check the version number of the DLL if utilising functionality noted only as being available in later revisions of the SDK. Or, alternatively, ensure the SDK DLL is updated appropriately during installation of your application.

1.1.1.1.13 Bodystat::BSIsBTAvailable Function

Check if a supported BT radio is available on this computer. Checks for presence of USB Bluetooth dongles and inbuilt Bluetooth modules in the PC. No communication with Bodystat devices takes place during this query.

C++

```
BODYSTATSDK_API BOOL BSIsBTAvailable();
```

C#

```
BODYSTATSDK_API BOOL BSIsBTAvailable();
```

Visual Basic

```
Function BSIsBTAvailable() As BODYSTATSDK_API BOOL
```

Java

```
BSIsBTAvailable
```

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSIsBTAvailable()
```

File

BodystatSDK.h (see page 55)

Returns

Returns true if supported radio is available and enabled.

Remarks

Note: Only supported BT radios can be queried (presently those compatible with the Microsoft Bluetooth APIs and some Widdcomm models).

Also note, Bluetooth must be enabled on the PC by the user to permit detection. Some PCs, laptops particularly, have physical switches to enable bluetooth/wireless communication.

1.1.1.1.14 Bodystat::BSOpenComport Function

Open com port to the device ready for communication.

C++

```
BODYSTATSDK_API BSError BSOpenComport(LPCTSTR lpComPort, int iSize);
```

C#

```
BODYSTATSDK_API BSError BSOpenComport(LPCTSTR lpComPort, int iSize);
```

Visual Basic

```
Function BSOpenComport(lpComPort As LPCTSTR, iSize As Integer) As BODYSTATSDK_API BSError
```

Java

```
BSOpenComport
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSOpenComport(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
LPCTSTR lpComPort	String buffer containing the device/port to be opened.
int iSize	The size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.1.1.1.15 Bodystat::BSReadCalibrationTime Function

Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.

C++

```
BODYSTATSDK_API BSError BSReadCalibrationTime(__time64_t * pCalibrationTime);
```

C#

```
BODYSTATSDK_API BSError BSReadCalibrationTime(ref __time64_t pCalibrationTime);
```

Visual Basic

```
Function BSReadCalibrationTime(ByRef pCalibrationTime As __time64_t) As BODYSTATSDK_API  
BSError
```

Java

```
BSReadCalibrationTime
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadCalibrationTime()
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
__time64_t * pCalibrationTime	[out] Receives the date the device was last calibrated.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously.

1.1.1.1.16 Bodystat::BSReadCurrentTime Function

Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).

C++

```
BODYSTATSDK_API BSError BSReadCurrentTime(__time64_t * pCurrentTime, int iDST = -1);
```

C#

```
BODYSTATSDK_API BSError BSReadCurrentTime(ref __time64_t pCurrentTime, int iDST);
```

Visual Basic

```
Function BSReadCurrentTime(ByRef pCurrentTime As __time64_t, iDST As Integer = -1) As  
BODYSTATSDK_API BSError
```

Java

```
BSReadCurrentTime
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadCurrentTime(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
__time64_t * pCurrentTime	[out] Receives the current time of the device's internal clock.
int iDST = -1	Daylight Saving Time handling (as per Microsoft time functions) Zero (0) to indicate that standard time is in effect. A value greater than 0 to indicate that daylight saving time is in effect. A value less than zero to have the C run-time library code compute whether standard time or daylight saving time is in effect. (default)

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously.

1.1.1.1.17 Bodystat::BSReadModelVersion Function

Read device model and firmware version of the connected Bodystat device.

C++

```
BODYSTATSDK_API BSError BSReadModelVersion(BSDeviceModel * pModel, BYTE * pMajorVersion,
BYTE * pMinorVersion, BYTE * pPsoc2Version, BYTE * pEepromVersion, BYTE * pBluetoothInfo);
```

C#

```
BODYSTATSDK_API BSError BSReadModelVersion(ref BSDeviceModel pModel, ref BYTE
pMajorVersion, ref BYTE pMinorVersion, ref BYTE pPsoc2Version, ref BYTE pEepromVersion, ref
BYTE pBluetoothInfo);
```

Visual Basic

```
Function BSReadModelVersion(ByRef pModel As BSDeviceModel, ByRef pMajorVersion As BYTE,
ByRef pMinorVersion As BYTE, ByRef pPsoc2Version As BYTE, ByRef pEepromVersion As BYTE,
ByRef pBluetoothInfo As BYTE) As BODYSTATSDK_API BSError
```

Java

```
BSReadModelVersion
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadModelVersion(, , , , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BSDeviceModel * pModel	[out] Receives the device model.
BYTE * pMajorVersion	[out] Receives the major version of the main firmware in the device.
BYTE * pMinorVersion	[out] Receives the minor version of the main firmware in the device.

BYTE * pPsoc2Version	[out] Receives the version of Psoc2 firmware in the device.
BYTE * pEepromVersion	[out] Receives the version of the EEPROM firmware in the device.
BYTE * pBluetoothInfo	[out] Receives the version of the Bluetooth module in the device.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously.

1.1.1.1.18 Bodystat::BSReadPrinterAddress Function

Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).

C++

```
BODYSTATSDK_API BSError BSReadPrinterAddress(LPTSTR lpPrinterAddress, int
iPrinterAddressBufferSize);
```

C#

```
BODYSTATSDK_API BSError BSReadPrinterAddress(LPTSTR lpPrinterAddress, int
iPrinterAddressBufferSize);
```

Visual Basic

```
Function BSReadPrinterAddress(lpPrinterAddress As LPTSTR, iPrinterAddressBufferSize As
Integer) As BODYSTATSDK_API BSError
```

Java

```
BSReadPrinterAddress
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadPrinterAddress(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
LPTSTR lpPrinterAddress	[out] String buffer which receives the printer address.
int iPrinterAddressBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously.

1.1.1.1.19 Bodystat::BSReadProtocolInfo Function

Read protocol information from the connected device. Contains version numbers relating to communication protocols in use, data structures and auxiliary information.

C++

```
BODYSTATSDK_API BSError BSReadProtocolInfo(BYTE * pProtocolVersion, BYTE * pDataVersion,
BYTE * pAuxInfo);
```

C#

```
BODYSTATSDK_API BSError BSReadProtocolInfo(ref BYTE pProtocolVersion, ref BYTE
pDataVersion, ref BYTE pAuxInfo);
```

Visual Basic

```
Function BSReadProtocolInfo(ByRef pProtocolVersion As BYTE, ByRef pDataVersion As BYTE,
ByRef pAuxInfo As BYTE) As BODYSTATSDK_API BSError
```

Java

```
BSReadProtocolInfo
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadProtocolInfo(, , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BYTE * pProtocolVersion	[out] Receives the communication protocol version.
BYTE * pDataVersion	[out] Receives the data structure version.
BYTE * pAuxInfo	[out] Receives auxiliary information.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Generally not required for the majority of SDK users as these details are abstracted by the SDK.

1.1.1.1.20 Bodystat::BSReadSerialNumber Function

Read the unique serial number of the connected Bodystat device.

C++

```
BODYSTATSDK_API BSError BSReadSerialNumber(ULONG * pSerialNumber);
```

C#

```
BODYSTATSDK_API BSError BSReadSerialNumber(ref ULONG pSerialNumber);
```

Visual Basic

```
Function BSReadSerialNumber(ByRef pSerialNumber As ULONG) As BODYSTATSDK_API BSError
```

Java

```
BSReadSerialNumber
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadSerialNumber()
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
ULONG * pSerialNumber	[out] Receives the serial number of the device.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously.

1.1.1.1.21 Bodystat::BSReadStoredTestData Function

Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged.

Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.

C++

```
BODYSTATSDK_API BSError BSReadStoredTestData(BSRawData * pRawData);
```

C#

```
BODYSTATSDK_API BSError BSReadStoredTestData(ref BSRawData pRawData);
```

Visual Basic

```
Function BSReadStoredTestData(ByRef pRawData As BSRawData) As BODYSTATSDK_API BSError
```

Java

```
BSReadStoredTestData
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadStoredTestData()
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BSRawData * pRawData	[out] Pointer to a structure to receive the downloaded data. Caller responsible for allocation and disposing of the structure.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously. Caller must check they are working with a Bodystat (rather than Equistat model, Equistat models will respond without error - but with garbage results).

1.1.1.1.22 Bodystat::BSReadTestBodystat Function

Perform a quick bio-impedance test measurement on the connected bodystat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.

May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.

Useful for diagnostics or performing measurements instigated by software.

C++

```
BODYSTATSDK_API BSError BSReadTestBodystat(UINT * pZ5, UINT * pZ50, UINT * pZ100, UINT * pZ200, UINT * pR50, float * pX50, float * pPA50);
```

C#

```
BODYSTATSDK_API BSError BSReadTestBodystat(ref UINT pZ5, ref UINT pZ50, ref UINT pZ100, ref
```



```
UINT pZ200, ref UINT pR50, ref float pX50, ref float pPA50);
```

Visual Basic

```
Function BSReadTestBodystat(ByRef pZ5 As UINT, ByRef pZ50 As UINT, ByRef pZ100 As UINT,
ByRef pZ200 As UINT, ByRef pR50 As UINT, ByRef pX50 As float, ByRef pPA50 As float) As
BODYSTATSDK_API BSError
```

Java

```
BSReadTestBodystat
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSReadTestBodystat(, , , , , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
UINT * pZ5	[out] Receives the impedance measured at 5 kHz (result in ohms).
UINT * pZ50	[out] Receives the impedance measured at 50 kHz (result in ohms).
UINT * pZ100	[out] Receives the impedance measured at 100 kHz (result in ohms).
UINT * pZ200	[out] Receives the impedance measured at 200 kHz (result in ohms).
UINT * pR50	[out] Receives the resistance measured at 50 kHz (result in ohms).
float * pX50	[out] Receives the reactance measured at 50 kHz (result in ohms).
float * pPA50	[out] Receives the phase angle measured at 50 kHz (result in degrees).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously. Caller must check they are working with a Bodystat (rather than Equistat model, Equistat models will respond without error - but with garbage results).

1.1.1.1.23 Bodystat::BSSearchBTDevices Function

Search BT radio for any Bodystat devices (paired or otherwise). Existing paired devices will be counted even if they are switched off or out of range. New non-paired units will only be counted if switched on, in range and ready for communication (on main startup screen).

C++

```
BODYSTATSDK_API BOOL BSSearchBTDevices(int & iNewDevices, int & iAuthenticatedDevices, BOOL
bReportErrors = TRUE, HWND hWndParent = NULL, unsigned short iTimeout = 7);
```

C#

```
BODYSTATSDK_API BOOL BSSearchBTDevices(ref int iNewDevices, ref int iAuthenticatedDevices,
BOOL bReportErrors, HWND hWndParent, unsigned short iTimeout);
```

Visual Basic

```
Function BSSearchBTDevices(ByRef iNewDevices As Integer, ByRef iAuthenticatedDevices As
Integer, bReportErrors As BOOL = TRUE, hWndParent As HWND = NULL, iTimeout As unsigned
short = 7) As BODYSTATSDK_API BOOL
```

Java

```
BSSearchBTDevices
```

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSSearchBTDevices( , , , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
int & iNewDevices	[out] The number of new Bodystat devices found (not paired).
int & iAuthenticatedDevices	[out] The number of authenticated Bodystat devices (paired).
BOOL bReportErrors = TRUE	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (completely silent).
HWND hWndParent = NULL	Handle of the parent window.
unsigned short iTimeout = 7	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or timesout.

1.1.1.1.24 Bodystat::BSUnAuthenticateBTDevices Function

Unauthenticate (unpair) ALL Bodystat devices from this PC.

C++

```
BODYSTATSDK_API BOOL BSUnAuthenticateBTDevices(BOOL bReportErrors = TRUE, HWND hWndParent = NULL, unsigned short iTimeout = 7);
```

C#

```
BODYSTATSDK_API BOOL BSUnAuthenticateBTDevices(BOOL bReportErrors, HWND hWndParent, unsigned short iTimeout);
```

Visual Basic

```
Function BSUnAuthenticateBTDevices(bReportErrors As BOOL = TRUE, hWndParent As HWND = NULL, iTimeout As unsigned short = 7) As BODYSTATSDK_API BOOL
```

Java

```
BSUnAuthenticateBTDevices
```

MATLAB

```
function [BODYSTATSDK_API BOOL] = BSUnAuthenticateBTDevices( , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
BOOL bReportErrors = TRUE	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user.
HWND hWndParent = NULL	Handle of the parent window.
unsigned short iTimeout = 7	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or timesout.

1.1.1.1.25 Bodystat::BSWriteCurrentTime Function

Reset the current date/time of the internal clock on the connected device. Date/time is automatically set to match the current date/time of the PC

C++

```
BODYSTATSDK_API BSError BWriteCurrentTime();
```

C#

```
BODYSTATSDK_API BSError BWriteCurrentTime();
```

Visual Basic

```
Function BWriteCurrentTime() As BODYSTATSDK_API BSError
```

Java

```
BWriteCurrentTime
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BWriteCurrentTime()
```

File

BodystatSDK.h (see page 55)

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously. Locale Daylight Saving Time (DST) rules are observed. Use sister function BSWriteCurrentTime(__time64_t CurrentTime) to set a specific date/time which differs from the PC clock or override DST behaviour.

1.1.1.1.26 Bodystat::BSWriteCurrentTimeAs Function

Reset the current date/time of the internal clock on the connected device. Date/time is set to a specific value of your choosing.

C++

```
BODYSTATSDK_API BSError BSWriteCurrentTimeAs(__time64_t CurrentTime);
```

C#

```
BODYSTATSDK_API BSError BSWriteCurrentTimeAs(__time64_t CurrentTime);
```

Visual Basic

```
Function BSWriteCurrentTimeAs(CurrentTime As __time64_t) As BODYSTATSDK_API BSError
```

Java

```
BSWriteCurrentTimeAs
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSWriteCurrentTimeAs()
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
__time64_t CurrentTime	Date/time to apply to the internal clock of the device (should be current!).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously. Use sister function (omit date/time parameter entirely) BSWriteCurrentTime (see page 24)() to automatically set to the current date/time of the PC.

1.1.1.1.27 Bodystat::BSWritePrinterAddress Function

Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt printer pairing function in the device.

C++

```
BODYSTATSDK_API BSError BSWritePrinterAddress(LPCTSTR lpPrinterAddress, int iBufferSize);
```

C#

```
BODYSTATSDK_API BSError BSWritePrinterAddress(LPCTSTR lpPrinterAddress, int iBufferSize);
```

Visual Basic

```
Function BSWritePrinterAddress(lpPrinterAddress As LPCTSTR, iBufferSize As Integer) As  
BODYSTATSDK_API BSError
```

Java

```
BSWritePrinterAddress
```

MATLAB

```
function [BODYSTATSDK_API BSError] = BSWritePrinterAddress(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
LPCTSTR lpPrinterAddress	String buffer containing the Bluetooth address (BDA) of the printer.
int iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Note not all device models support this function (not all support direct printing).

1.1.1.1.28 Bodystat::ESCalculateNormals Function

Calculate normal values and normal ranges for a given animal and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from EScalculateResults (see page 26)).

C++

```
BODYSTATSDK_API BSError EScalculateNormals(const ESMeasurement * pM, const ESResults * pR,
ESNormals * pN);
```

C#

```
BODYSTATSDK_API BSError EScalculateNormals(ESMeasurement * pM, ESResults * pR, ref
ESNormals pN);
```

Visual Basic

```
Function EScalculateNormals(pM As ESMeasurement *, pR As ESResults *, ByRef pN As
ESNormals) As BODYSTATSDK_API BSError
```

Java

```
ESCalculateNormals
```

MATLAB

```
function [BODYSTATSDK_API BSError] = EScalculateNormals( , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
const ESMeasurement * pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
const ESResults * pR	The results calculated for this measurement previously.
ESNormals * pN	[out] Receives the normal values and normal ranges for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Normals vary by input attributes, measured bio-impedance values and calculated results.

1.1.1.1.29 Bodystat::ESCalculateResults Function

Calculate complete result set for a given bio-impedance measurement recorded by the Equistat device.

C++

```
BODYSTATSDK_API BSError EScalculateResults(const ESMeasurement * pM, ESResults * pR);
```

C#

```
BODYSTATSDK_API BSError EScalculateResults(ESMeasurement * pM, ref ESResults pR);
```

Visual Basic

```
Function ESCalculateResults(pM As ESMeasurement *, ByRef pR As ESResults) As BODYSTATSDK_API BSError
```

Java

```
ESCalculateResults
```

MATLAB

```
function [BODYSTATSDK_API BSError] = ESCalculateResults(, )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
const ESMeasurement * pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
ESResults * pR	[out] Receives the complete calculated results for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Results vary by input attributes and measured bio-impedance values.

1.1.1.1.30 Bodystat::ESReadStoredTestData Function

Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged.

Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.

C++

```
BODYSTATSDK_API BSError ESReadStoredTestData(ESRawData * pRawData);
```

C#

```
BODYSTATSDK_API BSError ESReadStoredTestData(ref ESRawData pRawData);
```

Visual Basic

```
Function ESReadStoredTestData(ByRef pRawData As ESRawData) As BODYSTATSDK_API BSError
```

Java

```
ESReadStoredTestData
```

MATLAB

```
function [BODYSTATSDK_API BSError] = ESReadStoredTestData( )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
ESRawData * pRawData	[out] Pointer to a structure to receive the downloaded data. Caller responsible for allocation and disposing of the structure.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 11) previously. Caller must check they are working with an Equistat (rather than Bodystat model, Bodystat models will respond without error - but with garbage results).

1.1.1.1.31 Bodystat::ESReadTestEquistat Function

Perform a quick bio-impedance test measurement on the connected equistat device. Normal test process is skipped (user is not prompted to enter any horse information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.

May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.

Useful for diagnostics or performing measurements instigated by software.

C++

```
BODYSTATSDK_API BSError ESReadTestEquistat(UINT * pZ5, UINT * pZ16, UINT * pZ24, UINT *
pZ50, UINT * pZ140, UINT * pZ200, UINT * pZ280, float * pPA50);
```

C#

```
BODYSTATSDK_API BSError ESReadTestEquistat(ref UINT pZ5, ref UINT pZ16, ref UINT pZ24, ref
UINT pZ50, ref UINT pZ140, ref UINT pZ200, ref UINT pZ280, ref float pPA50);
```

Visual Basic

```
Function ESReadTestEquistat(ByRef pZ5 As UINT, ByRef pZ16 As UINT, ByRef pZ24 As UINT,
ByRef pZ50 As UINT, ByRef pZ140 As UINT, ByRef pZ200 As UINT, ByRef pZ280 As UINT, ByRef
pPA50 As float) As BODYSTATSDK_API BSError
```

Java

```
ESReadTestEquistat
```

MATLAB

```
function [BODYSTATSDK_API BSError] = ESReadTestEquistat( , , , , , , , )
```

File

BodystatSDK.h (see page 55)

Parameters

Parameters	Description
UINT * pZ5	[out] Receives the impedance measured at 5 kHz (result in ohms).
UINT * pZ16	[out] Receives the impedance measured at 16 kHz (result in ohms).
UINT * pZ24	[out] Receives the impedance measured at 24 kHz (result in ohms).
UINT * pZ50	[out] Receives the impedance measured at 50 kHz (result in ohms).
UINT * pZ140	[out] Receives the impedance measured at 140 kHz (result in ohms).
UINT * pZ200	[out] Receives the impedance measured at 200 kHz (result in ohms).
UINT * pZ280	[out] Receives the impedance measured at 280 kHz (result in ohms).
float * pPA50	[out] Receives the phase angle measured at 50 kHz (result in degrees) (ES Pro Only).

Returns

Returns `BSError::NoError` if successful. Specific error code otherwise.






Remarks

You must have connected to the device by calling `BSConnect` (see page 11) previously. Caller must check they are working with an Equistat (rather than Bodystat model, Bodystat models will respond without error - but with garbage results).







1.1.1.2 Structs, Records, Enums




The following table lists structs, records, enums in this documentation.

Enumerations

	Name	Description
	<code>BSAnalysisMode</code> (see page 30)	<code>BSAnalysisMode</code> : Analysis results are available in several formats: <code>BSAnalysisModeBC</code> = Body composition analysis <code>BSAnalysisModeHN</code> = Hydration analysis <code>BSAnalysisModeBoth</code> = Combined body composition and hydration analysis (default)
	<code>BSDeviceFamily</code> (see page 30)	<code>BSDeviceFamily</code> : Enumerated values that group each device / model variant into a given product family. Devices within a given product family share broad capabilities and features. For details regarding capabilities and features refer to the technical documentation for that product family (http://www.bodystat.com)
	<code>BSDeviceModel</code> (see page 32)	<code>BSDeviceModel</code> : Enumerated values that represent the various different models of the Bodystat devices over the years. Suffix: The suffix after the model indicates specific hardware variations regarding connectivity: <code>OPTO</code> indicates model communicates via opto-isolated serial. <code>DIU</code> indicates model communicates via a data interface unit (induction interface) <code>BT</code> indicates model communicates via Bluetooth serial port.
	<code>BSError</code> (see page 33)	<code>BSError</code> : Bodystat Error codes. A list of Bodystat specific error codes.
	<code>BSGender</code> (see page 35)	<code>BSGender</code> : Enumerated values that represent a person's gender.

Structures

	Name	Description
	<code>BSMeasurement</code> (see page 36)	<code>BSMeasurement</code> : Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).
	<code>BSNormals</code> (see page 38)	<code>BSNormals</code> : Bodystat Normals structure. Holds the normal values or normal ranges for any given test measurement.
	<code>BSRawData</code> (see page 40)	<code>BSRawData</code> : Bodystat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record starting number.
	<code>BSResults</code> (see page 41)	<code>BSResults</code> : Bodystat Results structure. Holds the calculated results for any given test measurement.
	<code>ESMeasurement</code> (see page 44)	<code>ESMeasurement</code> : Equistat Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the horse, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).
	<code>ESNormals</code> (see page 45)	<code>ESNormals</code> : Equistat Normals structure. Holds the normal values or normal ranges for any given test measurement.

	ESRawData (see page 47)	ESRawData: Equistat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record starting number.
	ESResults (see page 48)	ESResults: Equistat Results structure. Holds the calculated results for any given test measurement.
	ImpData (see page 49)	Represents one frequency measurement

1.1.1.2.1 Bodystat::BSAnalysisMode Enumeration

BSAnalysisMode: Analysis results are available in several formats:

BSAnalysisModeBC= Body composition analysis

BSAnalysisModeHN= Hydration analysis

BSAnalysisModeBoth= Combined body composition and hydration analysis (default)

C++

```
enum BSAnalysisMode {
    BSAnalysisModeNone = 0,
    BSAnalysisModeBC = 1,
    BSAnalysisModeHN = 2,
    BSAnalysisModeBoth = 3
};
```

C#

```
public enum BSAnalysisMode {
    BSAnalysisModeNone = 0,
    BSAnalysisModeBC = 1,
    BSAnalysisModeHN = 2,
    BSAnalysisModeBoth = 3
}
```

Visual Basic

```
Public Enum BSAnalysisMode
    BSAnalysisModeNone = 0
    BSAnalysisModeBC = 1
    BSAnalysisModeHN = 2
    BSAnalysisModeBoth = 3
End Enum
```

Java

```
public class BSAnalysisMode;
```

MATLAB

```
BSAnalysisMode
```

File

BodystatSDK.h ([see page 55](#))

Members

Members	Description
BSAnalysisModeNone = 0	No analysis mode specified (raw data only)
BSAnalysisModeBC = 1	Body composition enabled
BSAnalysisModeHN = 2	Hydration enabled
BSAnalysisModeBoth = 3	Both body composition and hydration (default)

1.1.1.2.2 Bodystat::BSDeviceFamily Enumeration

BSDeviceFamily: Enumerated values that group each device / model variant into a given product family. Devices within a given product family share broad capabilities and features. For details regarding capabilities and features refer to the

technical documentation for that product family (<http://www.bodystat.com>)

C++

```
enum BSDeviceFamily {
    BSUnknown_Family = 0,
    BS1500_Family = -1,
    BSMDD_Family = -2,
    BSDualScan_Family = -3,
    BSQuadScan_Family = -4,
    BSMultiScan_Family = -5,
    BSEquistatLite_Family = -6,
    BSEquistatPro_Family = -7
};
```

C#

```
public enum BSDeviceFamily {
    BSUnknown_Family = 0,
    BS1500_Family = -1,
    BSMDD_Family = -2,
    BSDualScan_Family = -3,
    BSQuadScan_Family = -4,
    BSMultiScan_Family = -5,
    BSEquistatLite_Family = -6,
    BSEquistatPro_Family = -7
}
```

Visual Basic

```
Public Enum BSDeviceFamily
    BSUnknown_Family = 0
    BS1500_Family = -1
    BSMDD_Family = -2
    BSDualScan_Family = -3
    BSQuadScan_Family = -4
    BSMultiScan_Family = -5
    BSEquistatLite_Family = -6
    BSEquistatPro_Family = -7
End Enum
```

Java

```
public class BSDeviceFamily;
```

MATLAB

```
BSDeviceFamily
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
BSUnknown_Family = 0	Unknown product family
BS1500_Family = -1	Bodystat 1500 family
BSMDD_Family = -2	Bodystat MDD family
BSDualScan_Family = -3	Bodystat DualScan family (legacy)
BSQuadScan_Family = -4	Bodystat QuadScan family
BSMultiScan_Family = -5	Bodystat Multiscan family
BSEquistatLite_Family = -6	Equistat Lite family (not supported by this SDK)
BSEquistatPro_Family = -7	Equistat Pro family (not supported by this SDK)

Remarks

Note that not all families are supported by this SDK. Some are listed for historic reference only.

1.1.1.2.3 Bodystat::BSDeviceModel Enumeration

BSDeviceModel: Enumerated values that represent the various different models of the Bodystat devices over the years.

Suffix: The suffix after the model indicates specific hardware variations regarding connectivity:

OPTO indicates model communicates via opto-isolated serial.

DIU indicates model communicates via a data interface unit (induction interface)

BT indicates model communicates via Bluetooth serial port.

C++

```
enum BSDeviceModel {
    BSUnknown = 0,
    BS1500_OPTO = 1,
    BSMDD_OPTO = 2,
    BS1500_DIU = 3,
    BSMDD_DIU = 4,
    BSDualScan_OPTO = 5,
    BSDualScan_DIU = 6,
    BSQuadScan_DIU = 7,
    BSMultiScan_DIU = 8,
    BS1500_BT = 9,
    BSMDD_BT = 10,
    BSQuadScan_BT = 11,
    BSEquistatPro_BT = 12,
    BSEquistatLite_BT = 13,
    BSQuadScanTouch = 14,
    BSMultiScanTouch = 15
};
```

C#

```
public enum BSDeviceModel {
    BSUnknown = 0,
    BS1500_OPTO = 1,
    BSMDD_OPTO = 2,
    BS1500_DIU = 3,
    BSMDD_DIU = 4,
    BSDualScan_OPTO = 5,
    BSDualScan_DIU = 6,
    BSQuadScan_DIU = 7,
    BSMultiScan_DIU = 8,
    BS1500_BT = 9,
    BSMDD_BT = 10,
    BSQuadScan_BT = 11,
    BSEquistatPro_BT = 12,
    BSEquistatLite_BT = 13,
    BSQuadScanTouch = 14,
    BSMultiScanTouch = 15
}
```

Visual Basic

```
Public Enum BSDeviceModel
    BSUnknown = 0
    BS1500_OPTO = 1
    BSMDD_OPTO = 2
    BS1500_DIU = 3
    BSMDD_DIU = 4
    BSDualScan_OPTO = 5
    BSDualScan_DIU = 6
    BSQuadScan_DIU = 7
    BSMultiScan_DIU = 8
    BS1500_BT = 9
    BSMDD_BT = 10
    BSQuadScan_BT = 11
    BSEquistatPro_BT = 12
    BSEquistatLite_BT = 13
    BSQuadScanTouch = 14

```

```
BSMultiScanTouch = 15
End Enum
```

Java

```
public class BSDeviceModel;
```

MATLAB

```
BSDeviceModel
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
BSUnknown = 0	Unknown model
BS1500_OPTO = 1	Bodystat 1500 with opto-isolated interface (legacy)
BSMDD_OPTO = 2	Bodystat MDD with opto-isolated interface (legacy)
BS1500_DIU = 3	Bodystat 1500 with data interface unit (legacy)
BSMDD_DIU = 4	Bodystat MDD with data interface unit (legacy)
BSDualScan_OPTO = 5	Bodystat DualScan with opto-isolated interface (legacy)
BSDualScan_DIU = 6	Bodystat DualScan with data interface unit (legacy)
BSQuadScan_DIU = 7	Bodystat QuadScan with data interface unit (legacy)
BSMultiScan_DIU = 8	Bodystat MultiScan with data interface unit (legacy)
BS1500_BT = 9	Bodystat 1500 with Bluetooth interface
BSMDD_BT = 10	Bodystat 1500MDD with Bluetooth interface
BSQuadScan_BT = 11	Bodystat QuadScan 4000 with Bluetooth interface
BSEquistatPro_BT = 12	Equistat Pro with Bluetooth interface (not supported by this SDK)
BSEquistatLite_BT = 13	Equistat Lite with Bluetooth interface (not supported by this SDK)
BSQuadScanTouch = 14	Bodystat QuadScan 4000 Touch with Wifi/Bluetooth interface
BSMultiScanTouch = 15	Bodystat Multiscan with Wifi/Bluetooth interface

Remarks

Note that not all models are supported by this SDK. Some are listed for historic reference only.

1.1.1.2.4 Bodystat::BSError Enumeration

BSError: Bodystat Error codes. A list of Bodystat specific error codes.

C++

```
enum BSError {
    NoError = 0,
    InvalidDeviceModel = -1,
    InvalidDeviceFamily = -2,
    InvalidDeviceMode = -3,
    InvalidGender = -4,
    InvalidAge = -5,
    InvalidHeight = -6,
    InvalidWeight = -7,
    InvalidActivity = -8,
    InvalidWaist = -9,
    InvalidHip = -10,
    InvalidImpedanceZ5 = -11,
    InvalidImpedanceZ50 = -12,
    InvalidImpedanceZ100 = -13,
    InvalidImpedanceZ200 = -14,
    InvalidResistanceR50 = -15,
    InvalidReactanceX50 = -16,
```

```

InvalidPhaseAnglePA50 = -17,
InvalidParamSupplied = -18,
ComErrorOpeningPort = -50,
ComErrorConfiguringPort = -51,
ComErrorConnecting = -52,
ComErrorUnsupportedCommand = -53,
ComErrorProcessingCommand = -54,
ComErrorReading = -55,
ComErrorWriting = -56,
ComErrorUserAborted = -57
};

```

C#

```

public enum BSError {
    NoError = 0,
    InvalidDeviceModel = -1,
    InvalidDeviceFamily = -2,
    InvalidDeviceMode = -3,
    InvalidGender = -4,
    InvalidAge = -5,
    InvalidHeight = -6,
    InvalidWeight = -7,
    InvalidActivity = -8,
    InvalidWaist = -9,
    InvalidHip = -10,
    InvalidImpedanceZ5 = -11,
    InvalidImpedanceZ50 = -12,
    InvalidImpedanceZ100 = -13,
    InvalidImpedanceZ200 = -14,
    InvalidResistanceR50 = -15,
    InvalidReactanceX50 = -16,
    InvalidPhaseAnglePA50 = -17,
    InvalidParamSupplied = -18,
    ComErrorOpeningPort = -50,
    ComErrorConfiguringPort = -51,
    ComErrorConnecting = -52,
    ComErrorUnsupportedCommand = -53,
    ComErrorProcessingCommand = -54,
    ComErrorReading = -55,
    ComErrorWriting = -56,
    ComErrorUserAborted = -57
}

```

Visual Basic

```

Public Enum BSError
    NoError = 0
    InvalidDeviceModel = -1
    InvalidDeviceFamily = -2
    InvalidDeviceMode = -3
    InvalidGender = -4
    InvalidAge = -5
    InvalidHeight = -6
    InvalidWeight = -7
    InvalidActivity = -8
    InvalidWaist = -9
    InvalidHip = -10
    InvalidImpedanceZ5 = -11
    InvalidImpedanceZ50 = -12
    InvalidImpedanceZ100 = -13
    InvalidImpedanceZ200 = -14
    InvalidResistanceR50 = -15
    InvalidReactanceX50 = -16
    InvalidPhaseAnglePA50 = -17
    InvalidParamSupplied = -18
    ComErrorOpeningPort = -50
    ComErrorConfiguringPort = -51
    ComErrorConnecting = -52
    ComErrorUnsupportedCommand = -53
    ComErrorProcessingCommand = -54
    ComErrorReading = -55
    ComErrorWriting = -56

```

```
ComErrorUserAborted = -57
End Enum
```

Java

```
public class BSError;
```

MATLAB

```
BSError
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
NoError = 0	No error occurred
InvalidDeviceModel = -1	Invalid model specified
InvalidDeviceFamily = -2	Invalid family specified
InvalidDeviceMode = -3	Invalid analysis mode
InvalidGender = -4	Invalid BSGender (see page 35) specified
InvalidAge = -5	Invalid age specified
InvalidHeight = -6	Invalid height measurement specified
InvalidWeight = -7	Invalid weight measurement specified
InvalidActivity = -8	Invalid activity level specified
InvalidWaist = -9	Invalid waist measurement specified
InvalidHip = -10	Invalid hip measurement specified
InvalidImpedanceZ5 = -11	Invalid impedance measurement specified
InvalidImpedanceZ50 = -12	Invalid impedance measurement specified
InvalidImpedanceZ100 = -13	Invalid impedance measurement specified
InvalidImpedanceZ200 = -14	Invalid impedance measurement specified
InvalidResistanceR50 = -15	Invalid resistance measurement specified
InvalidReactanceX50 = -16	Invalid reactance measurement specified
InvalidPhaseAnglePA50 = -17	Invalid phase angle measurement specified
InvalidParamSupplied = -18	Invalid parameter supplied
ComErrorOpeningPort = -50	Error opening the com port
ComErrorConfiguringPort = -51	Error configuring the com port parameters (required baud rate, buffers, stop bits, etc)
ComErrorConnecting = -52	Error connecting to the device
ComErrorUnsupportedCommand = -53	The requested command is not supported by the model being asked to perform it
ComErrorProcessingCommand = -54	Error occurred whilst processing the requested command
ComErrorReading = -55	Error occurred whilst reading data from the com port
ComErrorWriting = -56	Error occurred whilst writing data to the com port
ComErrorUserAborted = -57	The user aborted the process

1.1.1.2.5 Bodystat::BSGender Enumeration

BSGender: Enumerated values that represent a person's gender.

C++

```
enum BSGender {
    BSFemale = 0,
    BSMale = 1
};
```

C#

```
public enum BSGender {
    BSFemale = 0,
    BSMale = 1
}
```

Visual Basic

```
Public Enum BSGender
    BSFemale = 0
    BSMale = 1
End Enum
```

Java

```
public class BSGender;
```

MATLAB

```
BSGender
```

File

```
BodystatSDK.h (see page 55)
```

Members

Members	Description
BSFemale = 0	Female
BSMale = 1	Male

1.1.1.2.6 Bodystat::BSMeasurement Structure

BSMeasurement: Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).

C++

```
struct BSMeasurement {
    BSDeviceModel iDeviceModel;
    ULONG ulDeviceSerialNumber;
    __time64_t tTestDate;
    BSGender iGender;
    int iAge;
    int iHeight;
    float fWeight;
    int iActivity;
    int iWaist;
    int iHip;
    int iZ_5kHz;
    int iZ_50kHz;
    int iZ_100kHz;
    int iZ_200kHz;
    int iR_50kHz;
    float fX_50kHz;
    float fPA_50kHz;
    int iFrequencies;
    ImpData* pMultifreqData;
};
```

C#

```
public struct BSMeasurement {
    public BSDeviceModel iDeviceModel;
    public ULONG ulDeviceSerialNumber;
    public __time64_t tTestDate;
    public BSGender iGender;
    public int iAge;
    public int iHeight;
    public float fWeight;
```

```

public int iActivity;
public int iWaist;
public int iHip;
public int iZ_5kHz;
public int iZ_50kHz;
public int iZ_100kHz;
public int iZ_200kHz;
public int iR_50kHz;
public float fX_50kHz;
public float fPA_50kHz;
public int iFrequencies;
public ImpData* pMultifreqData;
}

```

Visual Basic

```

Public Structure BSMeasurement
    Public iDeviceModel As BSDeviceModel
    Public ulDeviceSerialNumber As ULONG
    Public tTestDate As __time64_t
    Public iGender As BSGender
    Public iAge As Integer
    Public iHeight As Integer
    Public fWeight As float
    Public iActivity As Integer
    Public iWaist As Integer
    Public iHip As Integer
    Public iZ_5kHz As Integer
    Public iZ_50kHz As Integer
    Public iZ_100kHz As Integer
    Public iZ_200kHz As Integer
    Public iR_50kHz As Integer
    Public fX_50kHz As float
    Public fPA_50kHz As float
    Public iFrequencies As Integer
    Public pMultifreqData As ImpData*
End Structure

```

Java

```
public class BSMeasurement;
```

MATLAB

```
BSMeasurement
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
BSDeviceModel iDeviceModel;	Device model which performed the measurement
ULONG ulDeviceSerialNumber;	Unique serial number of the device which performed the measurement
__time64_t tTestDate;	Date & time the test measurement was performed (Windows 64-bit time structure) (MDD/Quad only)
BSGender iGender;	Inputted gender of the subject
int iAge;	Inputted age of the subject (years)
int iHeight;	Inputted height of the subject (centimeters)
float fWeight;	Inputted weight of the subject (kilograms)
int iActivity;	Inputted physical activity level of the subject (range 1 to 5: 1=very low 2=low/medium 3=medium 4=medium/high 5=very high)
int iWaist;	Inputted waist size of the subject (centimeters)
int iHip;	Inputted hip size of the subject (centimeters)
int iZ_5kHz;	Impedance measured at 5 kHz (result in ohms) (MDD/Quad only)

int iZ_50kHz;	Impedance measured at 50 kHz (result in ohms)
int iZ_100kHz;	Impedance measured at 100 kHz (result in ohms) (Quad only)
int iZ_200kHz;	Impedance measured at 200 kHz (result in ohms) (Quad only)
int iR_50kHz;	Resistance measured at 50 kHz (result in ohms) (MDD/Quad only)
float fX_50kHz;	Reactance measured at 50 kHz (result in ohms) (MDD/Quad only)
float fPA_50kHz;	Phase angle measured at 50 kHz (result in degrees) (MDD/Quad only)
int iFrequencies;	number of frequencies used 11/50
ImpData* pMultifreqData;	pointer to ImpData (see page 49) array

1.1.1.2.7 Bodystat::BSNormals Structure

BSNormals: Bodystat Normals structure. Holds the normal values or normal ranges for any given test measurement.

C++

```
struct BSNormals {
    int iFatPerc_L;
    int iFatPerc_H;
    int iFatKg_L;
    int iFatKg_H;
    int iLeanPerc_L;
    int iLeanPerc_H;
    int iLeanKg_L;
    int iLeanKg_H;
    int iTTotalWeight_L;
    int iTTotalWeight_H;
    int iTTotalWeightMethod;
    int iTBW_Perc_L;
    int iTBW_Perc_H;
    int iTBW_L;
    int iTBW_H;
    int iECW_Perc_Norm;
    int iICW_Perc_Norm;
    int iBFMI_L;
    int iBFMI_H;
    int iFFMI_L;
    int iFFMI_H;
    float fNutrition_Norm;
    float fIllness_L;
    float fIllness_H;
    int iBMI_L;
    int iBMI_H;
    float fWaistHip_Norm;
    float fWellness_L;
    float fWellness_H;
};
```

C#

```
public struct BSNormals {
    public int iFatPerc_L;
    public int iFatPerc_H;
    public int iFatKg_L;
    public int iFatKg_H;
    public int iLeanPerc_L;
    public int iLeanPerc_H;
    public int iLeanKg_L;
    public int iLeanKg_H;
    public int iTTotalWeight_L;
    public int iTTotalWeight_H;
    public int iTTotalWeightMethod;
    public int iTBW_Perc_L;
```

```

public int iTBWPerC_H;
public int iTBW_L;
public int iTBW_H;
public int iECWPerC_Norm;
public int iICWPerC_Norm;
public int iBFMI_L;
public int iBFMI_H;
public int iFFMI_L;
public int iFFMI_H;
public float fNutrition_Norm;
public float fIllness_L;
public float fIllness_H;
public int iBMI_L;
public int iBMI_H;
public float fWaistHip_Norm;
public float fWellness_L;
public float fWellness_H;
}

```

Visual Basic

```

Public Structure BSNormals
    Public iFatPerC_L As Integer
    Public iFatPerC_H As Integer
    Public iFatKg_L As Integer
    Public iFatKg_H As Integer
    Public iLeanPerC_L As Integer
    Public iLeanPerC_H As Integer
    Public iLeanKg_L As Integer
    Public iLeanKg_H As Integer
    Public iTotalWeight_L As Integer
    Public iTotalWeight_H As Integer
    Public iTotalWeightMethod As Integer
    Public iTBWPerC_L As Integer
    Public iTBWPerC_H As Integer
    Public iTBW_L As Integer
    Public iTBW_H As Integer
    Public iECWPerC_Norm As Integer
    Public iICWPerC_Norm As Integer
    Public iBFMI_L As Integer
    Public iBFMI_H As Integer
    Public iFFMI_L As Integer
    Public iFFMI_H As Integer
    Public fNutrition_Norm As float
    Public fIllness_L As float
    Public fIllness_H As float
    Public iBMI_L As Integer
    Public iBMI_H As Integer
    Public fWaistHip_Norm As float
    Public fWellness_L As float
    Public fWellness_H As float
End Structure

```

Java

```
public class BSNormals;
```

MATLAB

```
BSNormals
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
int iFatPerC_L;	Normal range for fat percentage - (low end)
int iFatPerC_H;	Normal range for fat percentage - (high end)
int iFatKg_L;	Normal range for fat (in kg) - (low end)
int iFatKg_H;	Normal range for fat (in kg) - (high end)

int iLeanPerc_L;	Normal range for lean percentage - (low end)
int iLeanPerc_H;	Normal range for lean percentage - (high end)
int iLeanKg_L;	Normal range for lean (in kg) - (low end)
int iLeanKg_H;	Normal range for lean (in kg) - (high end)
int iTotalWeight_L;	Normal range for Total Weight (in kg) - (low end)
int iTotalWeight_H;	Normal range for Total Weight (in kg) - (high end)
int iTotalWeightMethod;	Method used when calculating Total Weight normal range (0 = Composition Method, 1 = BMI Method)
int iTBW_Perc_L;	Normal range for Total Body Water percentage - (low end)
int iTBW_Perc_H;	Normal range for Total Body Water percentage - (high end)
int iTBW_L;	Normal range for Total Body Water (in litres) - (low end)
int iTBW_H;	Normal range for Total Body Water (in litres) - (high end)
int iECW_Perc_Norm;	Normal value for ECW Percentage (Quad only)
int iICW_Perc_Norm;	Normal value for ICW Percentage (Quad only)
int iBFMI_L;	Normal range for Body Fat Mass Index (MDD/Quad only) - (low end)
int iBFMI_H;	Normal range for Body Fat Mass Index (MDD/Quad only) - (high end)
int iFFMI_L;	Normal range for Fat Free Mass Index (MDD/Quad only) - (low end)
int iFFMI_H;	Normal range for Fat Free Mass Index (MDD/Quad only) - (high end)
float fNutrition_Norm;	Normal value for Nutrition Index (Quad only)
float fIllness_L;	Normal range for Prediction Marker (TM) - formerly known as Illness Marker (TM) - (Quad only) - (low end)
float fIllness_H;	Normal range for Prediction Marker (TM) - formerly known as Illness Marker (TM) - (Quad only) - (high end)
int iBMI_L;	Normal range for Body Mass Index - (low end)
int iBMI_H;	Normal range for Body Mass Index - (high end)
float fWaistHip_Norm;	Normal value for Waist/Hip ratio
float fWellness_L;	Normal range for Wellness Marker (TM) (MDD only) - (low end)
float fWellness_H;	Normal range for Wellness Marker (TM) (MDD only) - (high end)

1.1.1.2.8 Bodystat::BSRawData Structure

BSRawData: Bodystat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record starting number.

C++

```
struct BSRawData {
    int iRecordArraySize;
    BSMeasurement record[BS_RAWDATA_ARRAYSIZE];
    int iTotalNumRecs;
    int ulFirstTestNum;
};
```

C#

```
public struct BSRawData {
    public int iRecordArraySize;
    public BSMeasurement record;
    public int iTotalNumRecs;
    public int ulFirstTestNum;
}
```

Visual Basic

```
Public Structure BSRawData
    Public iRecordArraySize As Integer
    Public record As BSMeasurement
    Public iTotalNumRecs As Integer
    Public ulFirstTestNum As Integer
End Structure
```

Java

```
public class BSRawData;
```

MATLAB

```
BSRawData
```

File

```
BodystatSDK.h (see page 55)
```

Members

Members	Description
int iRecordArraySize;	Maximum number of records this structure has reserved space for (currently fixed)
BSMeasurement record[BS_RAWDATA_ARRAYSIZE];	Array of test measurement data
int iTotalNumRecs;	Number of populated records in the array
int ulFirstTestNum;	Test number of the first record stored in the array.

Remarks

Most models index 1000 records and store the last 100 at any given time.

1.1.1.2.9 Bodystat::BSResults Structure

BSResults: Bodystat Results structure. Holds the calculated results for any given test measurement.

C++

```
struct BSResults {
    float fFatPerc;
    float fFatKg;
    float fLeanPerc;
    float fLeanKg;
    float fTotalWeight;
    float fDryLW;
    float fTBWPerc;
    float fTBW;
    float fECWPerc;
    float fECW;
    float fICWPerc;
    float fICW;
    float fBCM;
    float fThirdSpace;
    float fNutrition;
    float fIllness;
    float fBMR;
    float fBMRkg;
    float fEstAvg;
    float fBMI;
    float fBFMI;
    float fFFMI;
    float fWaistHip;
    float fWellness;
    float fECW_Legacy;
    float fTBW_Legacy;
    float fOHY;
    float fSkMuscle;
    float fCm;
```

```

float fRext;
float fRint;
float fFc;
float fAlpha;
};

```

C#

```

public struct BSResults {
    public float fFatPerc;
    public float fFatKg;
    public float fLeanPerc;
    public float fLeanKg;
    public float fTotalWeight;
    public float fDryLW;
    public float fTBWPerc;
    public float fTBW;
    public float fECWPerc;
    public float fECW;
    public float fICWPerc;
    public float fICW;
    public float fBCM;
    public float fThirdSpace;
    public float fNutrition;
    public float fIllness;
    public float fBMR;
    public float fBMRkg;
    public float fEstAvg;
    public float fBMI;
    public float fBFMI;
    public float fFFMI;
    public float fWaistHip;
    public float fWellness;
    public float fECW_Legacy;
    public float fTBW_Legacy;
    public float fOHY;
    public float fSkMuscle;
    public float fCm;
    public float fRext;
    public float fRint;
    public float fFc;
    public float fAlpha;
}

```

Visual Basic

```

Public Structure BSResults
    Public fFatPerc As float
    Public fFatKg As float
    Public fLeanPerc As float
    Public fLeanKg As float
    Public fTotalWeight As float
    Public fDryLW As float
    Public fTBWPerc As float
    Public fTBW As float
    Public fECWPerc As float
    Public fECW As float
    Public fICWPerc As float
    Public fICW As float
    Public fBCM As float
    Public fThirdSpace As float
    Public fNutrition As float
    Public fIllness As float
    Public fBMR As float
    Public fBMRkg As float
    Public fEstAvg As float
    Public fBMI As float
    Public fBFMI As float
    Public fFFMI As float
    Public fWaistHip As float
    Public fWellness As float
    Public fECW_Legacy As float
    Public fTBW_Legacy As float

```

```

Public fOHY As float
Public fSkMuscle As float
Public fCm As float
Public fRext As float
Public fRint As float
Public fFc As float
Public fAlpha As float
End Structure

```

Java

```
public class BSResults;
```

MATLAB

```
BSResults
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
float fFatPerc;	Fat percentage
float fFatKg;	Fat (in kg)
float fLeanPerc;	Lean percentage
float fLeanKg;	Lean (in kg)
float fTotalWeight;	Total Weight (in kg)
float fDryLW;	Dry Lean Weight (in kg)
float fTBWPerc;	Water percentage or Total Body Water percentage
float fTBW;	Water (in litres) or Total Body Water (in litres)
float fECWPerc;	Extra Cellular Water percentage (MDD/Quad only)
float fECW;	Extra Cellular Water (in litres) (MDD/Quad only)
float fICWPerc;	Intra Cellular Water percentage (MDD/Quad only)
float fICW;	Intra Cellular Water (in litres) (MDD/Quad only)
float fBCM;	Body cell mass (Quad only)
float fThirdSpace;	3rd space water (in litres) (Quad only)
float fNutrition;	Nutrition index (Quad only)
float fIllness;	Prediction Marker (TM) - formerly known as Illness Marker (TM) (Quad only)
float fBMR;	Basal Metabolic Rate (in kcal)
float fBMRkg;	Basal Metabolic Rate per kilogram in (kcal/kg)
float fEstAvg;	Estimated Average Requirement (in kcal)
float fBMI;	Body Mass Index
float fBFMI;	Body Fat Mass Index (MDD/Quad only)
float fFFMI;	Fat Free Mass Index (MDD/Quad only)
float fWaistHip;	Waist/Hip ratio
float fWellness;	Wellness Marker (TM) (MDD only)
float fECW_Legacy;	Legacy ECW calculation (QuanScan mode)
float fOHY;	Over hydration
float fSkMuscle;	Skeletal muscle mass
float fCm;	Cell membrane capacitance
float fRext;	R extracellular
float fRint;	R intracellular
float fFc;	Characteristic frequency
float fAlpha;	Alpha angle

1.1.1.2.10 Bodystat::ESMeasurement Structure

ESMeasurement: Equistat Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the horse, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).

C++

```
struct ESMeasurement {
    BSDeviceModel iDeviceModel;
    ULONG ulDeviceSerialNumber;
    __time64_t tTestDate;
    int iHeight;
    int iBCScore;
    int iZ_5kHz;
    int iZ_16kHz;
    int iZ_24kHz;
    int iZ_50kHz;
    int iZ_140kHz;
    int iZ_200kHz;
    int iZ_280kHz;
    float fPA_50kHz;
};
```

C#

```
public struct ESMeasurement {
    public BSDeviceModel iDeviceModel;
    public ULONG ulDeviceSerialNumber;
    public __time64_t tTestDate;
    public int iHeight;
    public int iBCScore;
    public int iZ_5kHz;
    public int iZ_16kHz;
    public int iZ_24kHz;
    public int iZ_50kHz;
    public int iZ_140kHz;
    public int iZ_200kHz;
    public int iZ_280kHz;
    public float fPA_50kHz;
}
```

Visual Basic

```
Public Structure ESMeasurement
    Public iDeviceModel As BSDeviceModel
    Public ulDeviceSerialNumber As ULONG
    Public tTestDate As __time64_t
    Public iHeight As Integer
    Public iBCScore As Integer
    Public iZ_5kHz As Integer
    Public iZ_16kHz As Integer
    Public iZ_24kHz As Integer
    Public iZ_50kHz As Integer
    Public iZ_140kHz As Integer
    Public iZ_200kHz As Integer
    Public iZ_280kHz As Integer
    Public fPA_50kHz As float
End Structure
```

Java

```
public class ESMeasurement;
```

MATLAB

```
ESMeasurement
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
BSDDeviceModel iDeviceModel;	Device model which performed the measurement
ULONG ulDeviceSerialNumber;	Unique serial number of the device which performed the measurement
__time64_t tTestDate;	Date & time the test measurement was performed (Windows 64-bit time structure) (ES Pro only)
int iHeight;	Inputted height of the subject (centimeters)
int iBCScore;	Inputted BC score
int iZ_5kHz;	Impedance measured at 5 kHz (result in ohms)
int iZ_16kHz;	Impedance measured at 16 kHz (result in ohms)
int iZ_24kHz;	Impedance measured at 24 kHz (result in ohms)
int iZ_50kHz;	Impedance measured at 50 kHz (result in ohms)
int iZ_140kHz;	Impedance measured at 140 kHz (result in ohms)
int iZ_200kHz;	Impedance measured at 200 kHz (result in ohms)
int iZ_280kHz;	Impedance measured at 280 kHz (result in ohms)
float fPA_50kHz;	Phase angle measured at 50 kHz (result in degrees) (ES Pro only)

1.1.1.2.11 Bodystat::ESNormals Structure

ESNormals: Equistat Normals structure. Holds the normal values or normal ranges for any given test measurement.

C++

```
struct ESNormals {
    int iFatKg_L;
    int iFatKg_H;
    int iFatKg_PBMPerc_L;
    int iFatKg_PBMPerc_H;
    int iTFV_L;
    int iTFV_H;
    int iTFVPerc_L;
    int iTFVPerc_H;
    int iECFV_L;
    int iECFV_H;
    int iECFV_PBM_L;
    int iECFV_PBM_H;
    int iECFV_TFVPerc_L;
    int iECFV_TFVPerc_H;
    int iICFV_L;
    int iICFV_H;
    int iICFV_PBM_L;
    int iICFV_PBM_H;
    int iICFV_TFVPerc_L;
    int iICFV_TFVPerc_H;
    int iPV_L;
    int iPV_H;
    int iPV_ECFVPerc_L;
    int iPV_ECFVPerc_H;
};
```

C#

```
public struct ESNormals {
    public int iFatKg_L;
    public int iFatKg_H;
    public int iFatKg_PBMPerc_L;
    public int iFatKg_PBMPerc_H;
    public int iTFV_L;
    public int iTFV_H;
    public int iTFVPerc_L;
    public int iTFVPerc_H;
```



```

public int iECFV_L;
public int iECFV_H;
public int iECFV_PBM_L;
public int iECFV_PBM_H;
public int iECFV_TFVPerc_L;
public int iECFV_TFVPerc_H;
public int iICFV_L;
public int iICFV_H;
public int iICFV_PBM_L;
public int iICFV_PBM_H;
public int iICFV_TFVPerc_L;
public int iICFV_TFVPerc_H;
public int iPV_L;
public int iPV_H;
public int iPV_ECFVPerc_L;
public int iPV_ECFVPerc_H;
}

```

Visual Basic

```

Public Structure ESNormals
    Public iFatKg_L As Integer
    Public iFatKg_H As Integer
    Public iFatKg_PBMPerc_L As Integer
    Public iFatKg_PBMPerc_H As Integer
    Public iTFV_L As Integer
    Public iTFV_H As Integer
    Public iTFVPerc_L As Integer
    Public iTFVPerc_H As Integer
    Public iECFV_L As Integer
    Public iECFV_H As Integer
    Public iECFV_PBM_L As Integer
    Public iECFV_PBM_H As Integer
    Public iECFV_TFVPerc_L As Integer
    Public iECFV_TFVPerc_H As Integer
    Public iICFV_L As Integer
    Public iICFV_H As Integer
    Public iICFV_PBM_L As Integer
    Public iICFV_PBM_H As Integer
    Public iICFV_TFVPerc_L As Integer
    Public iICFV_TFVPerc_H As Integer
    Public iPV_L As Integer
    Public iPV_H As Integer
    Public iPV_ECFVPerc_L As Integer
    Public iPV_ECFVPerc_H As Integer
End Structure

```

Java

```
public class ESNormals;
```

MATLAB

```
ESNormals
```

File

BodystatSDK.h (📄 see page 55)

Members

Members	Description
int iFatKg_L;	Normal range for fat (in kg) - (low end)
int iFatKg_H;	Normal range for fat (in kg) - (high end)
int iFatKg_PBMPerc_L;	Normal range for fat kg pbm percentage - (low end)
int iFatKg_PBMPerc_H;	Normal range for fat kg pbm percentage - (high end)
int iTFV_L;	Normal range for TFV - (low end)
int iTFV_H;	Normal range for TFV - (high end)
int iTFVPerc_L;	Normal range for TFV percentage - (low end)
int iTFVPerc_H;	Normal range for TFV percentage - (high end)

int iECFV_L;	Normal range for ECFV - (low end)
int iECFV_H;	Normal range for ECFV - (high end)
int iECFV_PBM_L;	Normal range for ECFV PBM - (low end)
int iECFV_PBM_H;	Normal range for ECFV PBM - (high end)
int iECFV_TFVPerc_L;	Normal range for ECFV TFV percentage - (low end)
int iECFV_TFVPerc_H;	Normal range for ECFV TFV percentage - (high end)
int iICFV_L;	Normal range for ICFV - (low end)
int iICFV_H;	Normal range for ICFV - (high end)
int iICFV_PBM_L;	Normal range for ICFV PBM - (low end)
int iICFV_PBM_H;	Normal range for ICFV PBM - (high end)
int iICFV_TFVPerc_L;	Normal range for ICFV TFV percentage - (low end)
int iICFV_TFVPerc_H;	Normal range for ICFV TFV percentage - (high end)
int iPV_L;	Normal range for PV - (low end)
int iPV_H;	Normal range for PV - (high end)
int iPV_ECFVPerc_L;	Normal range for PV ECFV percentage - (low end)
int iPV_ECFVPerc_H;	Normal range for PV ECFV percentage - (high end)

1.1.1.2.12 Bodystat::ESRawData Structure

ESRawData: Equistat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record starting number.

C++

```
struct ESRawData {
    int iRecordArraySize;
    ESMeasurement record[ES_RAWDATA_ARRAYSIZE];
    int iTotalNumRecs;
    int ulFirstTestNum;
};
```

C#

```
public struct ESRawData {
    public int iRecordArraySize;
    public ESMeasurement record;
    public int iTotalNumRecs;
    public int ulFirstTestNum;
}
```

Visual Basic

```
Public Structure ESRawData
    Public iRecordArraySize As Integer
    Public record As ESMeasurement
    Public iTotalNumRecs As Integer
    Public ulFirstTestNum As Integer
End Structure
```

Java

```
public class ESRawData;
```

MATLAB

```
ESRawData
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
int iRecordArraySize;	Maximum number of records this structure has reserved space for (currently fixed)
ESMeasurement record[ES_RAWDATA_ARRAYSIZE];	Array of test measurement data
int iTotalNumRecs;	Number of populated records in the array
int ulFirstTestNum;	Test number of the first record stored in the array.

Remarks

Most models index 1000 records and store the last 100 at any given time.

1.1.1.2.13 Bodystat::ESResults Structure

ESResults: Equistat Results structure. Holds the calculated results for any given test measurement.

C++

```
struct ESResults {
    float fFatKg;
    float fPBM;
    float fIBM;
    float fFatKg_PBMperc;
    float fTFV;
    float fTFVperc;
    float fECFV;
    float fECFV_PBM;
    float fECFV_TFVperc;
    float fICFV;
    float fICFV_PBM;
    float fICFV_TFVperc;
    float fPV;
    float fPV_ECFVperc;
    float fDehydration;
};
```

C#

```
public struct ESResults {
    public float fFatKg;
    public float fPBM;
    public float fIBM;
    public float fFatKg_PBMperc;
    public float fTFV;
    public float fTFVperc;
    public float fECFV;
    public float fECFV_PBM;
    public float fECFV_TFVperc;
    public float fICFV;
    public float fICFV_PBM;
    public float fICFV_TFVperc;
    public float fPV;
    public float fPV_ECFVperc;
    public float fDehydration;
}
```

Visual Basic

```
Public Structure ESResults
    Public fFatKg As float
    Public fPBM As float
    Public fIBM As float
    Public fFatKg_PBMperc As float
    Public fTFV As float
    Public fTFVperc As float
    Public fECFV As float
    Public fECFV_PBM As float
    Public fECFV_TFVperc As float
```

```

Public fICFV As float
Public fICFV_PBM As float
Public fICFV_TFVPerc As float
Public fPV As float
Public fPV_ECFVPerc As float
Public fDehydration As float
End Structure

```

Java

```
public class ESResults;
```

MATLAB

```
ESResults
```

File

BodystatSDK.h (see page 55)

Members

Members	Description
float fFatKg;	Fat (in kg)
float fPBM;	PBM
float fIBM;	IBM
float fFatKg_PBMperc;	Fat Kg PBM percentage
float fTFV;	TFV
float fTFVperc;	TFV Percentage
float fECFV;	ECFV
float fECFV_PBM;	ECFV PBM
float fECFV_TFVperc;	ECFV TFV percentage
float fICFV;	ICFV
float fICFV_PBM;	ICFV PBM
float fICFV_TFVperc;	ICFV TFV percentage
float fPV;	PV
float fPV_ECFVperc;	PV ECFV percentage
float fDehydration;	Dehydration

1.1.1.2.14 ImpData Structure

Represents one frequency measurement

C++

```

struct ImpData {
    float fFrequency;
    float fImpedance;
    float fPhaseAngle;
};

```

C#

```

public struct ImpData {
    public float fFrequency;
    public float fImpedance;
    public float fPhaseAngle;
}

```

Visual Basic

```

Public Structure ImpData
    Public fFrequency As float
    Public fImpedance As float
    Public fPhaseAngle As float
End Structure

```

Java

```
public class ImpData;
```




MATLAB

```
ImpData
```



File

```
BodystatSDK.h (see page 55)
```

ImpData Data Members

	Name	Description
	fFrequency (see page 50)	Frequency in kHz
	fImpedance (see page 50)	Absolute value of impedance measured at given frequency (result in ohms)
	fPhaseAngle (see page 51)	Phase angle (result in degrees)

ImpData Methods

	Name	Description
	ImZ (see page 51)	Imaginary part of impedance, also reactance $\text{Im } Z = Z \cdot \sin(\text{Phase})$
	ReZ (see page 51)	Real part of impedance, also resistance $\text{Re } Z = Z \cdot \cos(\text{Phase})$

1.1.1.2.14.1 ImpData Data Members

1.1.1.2.14.1.1 ImpData::fFrequency Data Member

C++

```
float fFrequency;
```

C#

```
public float fFrequency;
```

Visual Basic

```
Public fFrequency As float
```

Java

```
public float fFrequency;
```

MATLAB

```
fFrequency
```

Description

```
Frequency in kHz
```

1.1.1.2.14.1.2 ImpData::fImpedance Data Member

C++

```
float fImpedance;
```

C#

```
public float fImpedance;
```

Visual Basic

```
Public fImpedance As float
```

Java

```
public float fImpedance;
```

MATLAB

```
fImpedance
```

Description

Absolute value of impedance measured at given frequency (result in ohms)

1.1.1.2.14.1.3 ImpData::fPhaseAngle Data Member**C++**

```
float fPhaseAngle;
```

C#

```
public float fPhaseAngle;
```

Visual Basic

```
Public fPhaseAngle As float
```

Java

```
public float fPhaseAngle;
```

MATLAB

```
fPhaseAngle
```

Description

Phase angle (result in degrees)

1.1.1.2.14.2 ImpData Methods**1.1.1.2.14.2.1 ImpData::ImZ Method**

Imaginary part of impedance, also reactance $\text{Im } Z = |Z| \cdot \sin(\text{Phase})$

C++

```
float ImZ() const;
```

C#

```
public float ImZ();
```

Visual Basic

```
Public Function ImZ() As float
```

Java

```
public float ImZ();
```

MATLAB

```
ImZ
```

Body Source

```
float ImZ() const
{
    return fImpedance * sin(fPhaseAngle * M_PI / 180);
}
```

1.1.1.2.14.2.2 ImpData::ReZ Method

Real part of impedance, also resistance $\text{Re } Z = |Z| \cdot \cos(\text{Phase})$

C++

```
float ReZ() const;
```

C#

```
public float ReZ();
```

Visual Basic

```
Public Function ReZ() As float
```

Java

```
public float ReZ();
```

MATLAB

```
ReZ
```

Body Source

```
float ReZ() const
{
    return fImpedance * cos(fPhaseAngle * M_PI / 180);
}
```

1.1.1.3 Variables

The following table lists variables in this documentation.

Variables

Name	Description
BS_RAWDATA_ARRAYSIZE (see page 52)	Number of measurement records BSRawData (see page 40) structure should reserved space for (currently fixed)
ES_RAWDATA_ARRAYSIZE (see page 52)	Number of measurement records ESRawData (see page 47) structure should reserved space for (currently fixed)

1.1.1.3.1 Bodystat::BS_RAWDATA_ARRAYSIZE Variable

C++

```
const int BS_RAWDATA_ARRAYSIZE = 100;
```

C#

```
BS_RAWDATA_ARRAYSIZE
```

Visual Basic

```
BS_RAWDATA_ARRAYSIZE
```

Java

```
BS_RAWDATA_ARRAYSIZE
```

MATLAB

```
BS_RAWDATA_ARRAYSIZE
```

File

```
BodystatSDK.h (see page 55)
```

Description

Number of measurement records BSRawData (see page 40) structure should reserved space for (currently fixed)

1.1.1.3.2 Bodystat::ES_RAWDATA_ARRAYSIZE Variable

C++

```
const int ES_RAWDATA_ARRAYSIZE = 100;
```

C#

```
ES_RAWDATA_ARRAYSIZE
```

Visual Basic

```
ES_RAWDATA_ARRAYSIZE
```

Java

ES_RAWDATA_ARRAYSIZE

MATLAB

ES_RAWDATA_ARRAYSIZE

File

BodystatSDK.h (see page 55)

Description

Number of measurement records ESRawData (see page 47) structure should reserved space for (currently fixed)

1.1.2 Macros

The following table lists macros in this documentation.

Macros

Name	Description
BODYSTATSDK_API (see page 53)	The following ifdef block is the standard way of creating macros which make exporting from a DLL simpler. All files within this DLL are compiled with the BODYSTATSDK_EXPORTS symbol defined on the command line. This symbol should not be defined on any project that uses this DLL. This way any other project whose source files include this file see BODYSTATSDK_API functions as being imported from a DLL, whereas this DLL sees symbols defined with this macro as being exported.
M_PI (see page 54)	This is macro M_PI.

Module

Bodystat API - Main (see page 1)

1.1.2.1 BODYSTATSDK_API Macro

C++

```
#define BODYSTATSDK_API extern "C" __declspec(dllimport)
```

C#

```
#define BODYSTATSDK_API extern "C" __declspec(dllimport)
```

Visual Basic

```
#Const BODYSTATSDK_API = extern "C" __declspec(dllimport)
```

Java

BODYSTATSDK_API

MATLAB

BODYSTATSDK_API

File

BodystatSDK.h (see page 55)

Description

The following ifdef block is the standard way of creating macros which make exporting from a DLL simpler. All files within this DLL are compiled with the BODYSTATSDK_EXPORTS symbol defined on the command line. This symbol should not be defined on any project that uses this DLL. This way any other project whose source files include this file see BODYSTATSDK_API functions as being imported from a DLL, whereas this DLL sees symbols defined with this macro as

being exported.

1.1.2.2 M_PI Macro

C++

```
#define M_PI 3.14159265358979323846
```

C#

```
#define M_PI 3.14159265358979323846
```

Visual Basic

```
#Const M_PI = 3.14159265358979323846
```

Java

```
M_PI
```

MATLAB

```
M_PI
```

File

BodystatSDK.h (see page 55)

Description

This is macro M_PI.

1.1.3 Files

The following table lists files in this documentation.

Files

Name	Description
BodystatSDK.h (see page 55)	<div><div>Bodystat SDK</div><div>(C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved.</div><div>Use subject to license.</div><div>Web: http://www.bodystat.com</div><div>Email: info@bodystat.com</div><div>Tel: +44 (0)1624 629 571</div><div>Fax: +44 (0)1624 611 544</div><div>Version</div><div>1.36 (09/May/2014)</div><div>Selected download is now passed through to the results form in the C# sample.</div><div>Fixed error in VB.Net wrapper (preventing download working) in the VB.Net sample.</div><div>Fixed some missing normals for MST device.</div><div>1.35 (10/Mar/2014)</div><div>Added Cole methods for ECW/ICW calculations used in the MST device.</div><div>Fixed a memory pointer issue for the detected bodystat device info structure.</div><div>1.34 (12/Nov/2013)</div><div>Exposed Bluetooth API search timeout values to... more (see page 55)</div></div>

Module

Bodystat API - Main (see page 1)

1.1.3.1 BodystatSDK.h

Bodystat SDK

(C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved.

Use subject to license.

Web: <http://www.bodystat.com>

Email: info@bodystat.com

Tel: +44 (0)1624 629 571

Fax: +44 (0)1624 611 544

Version

1.36 (09/May/2014)

Selected download is now passed through to the results form in the C# sample.

Fixed error in VB.Net wrapper (preventing download working) in the VB.Net sample.

Fixed some missing normals for MST device.

1.35 (10/Mar/2014)

Added Cole methods for ECW/ICW calculations used in the MST device.

Fixed a memory pointer issue for the detected bodystat device info structure.

1.34 (12/Nov/2013)

Exposed Bluetooth API search timeout values to the SDK

1.33 (28/Sep/2012)

Removed overly aggressive constraints on measurement parameters in CalculateResults and CalculateNormals()

Added ECW and ICW to MDD results.

Updated MFC/VB/C# samples to support ECW and ICW for MDD results.

Renamed Illness Maker(TM) to Prediction Marker(TM) - naming/labeling change only, code left unchanged to avoid breaking change.

1.32 (08/Aug/2011)

Corrected a bug in which GetBodystatDevice could erroneously return success with an invalid device even though no device was found

Corrected message box title when pairing new devices

Updated C# wrapper with get/set accessors for databinding [Breaking Change: Minor refactoring will be required to accommodate new accessor names]

Updated C# wrapper with additional wrappers for calls with StringBuilder parameters (now callable directly with string parameters)

Updated C# wrapper to force imports as Cdecl types

Fixed bug in C# wrapper causing invalid parameter error during download

Updated VB wrapper to force imports as Cdecl types

Updated C# WinForms sample to use new string wrappers

1.31 (26/Oct/2010)

Added C# wrapper and C# WinForms sample

Added initial support for Equistat devices

1.30 (29/Aug/2010)

Added VB.Net wrapper and VB WinForms sample

1.29 (30/Jun/2010)

Initial Release

Macros

Name	Description
BODYSTATSDK_API (see page 53)	The following ifdef block is the standard way of creating macros which make exporting from a DLL simpler. All files within this DLL are compiled with the BODYSTATSDK_EXPORTS symbol defined on the command line. This symbol should not be defined on any project that uses this DLL. This way any other project whose source files include this file see BODYSTATSDK_API functions as being imported from a DLL, whereas this DLL sees symbols defined with this macro as being exported.
M_PI (see page 54)	This is macro M_PI.

Namespaces

Name	Description
Bodystat (see page 2)	This is the main namespace for the Bodystat API. The entire functions within the Bodystat API are enclosed within the 'Bodystat' namespace.

1.2 C# .NET Wrapper

Bodystat API C# .NET wrapper.

Provides a wrapper interface to the Bodystat DLL with C# calls and types you can work with directly in your C# .NET program. The wrapper handles the necessary marshalling and interop issues required to interface with the native unmanaged Bodystat DLL.

Namespaces

Name	Description
BodystatAPI (see page 57)	This is the Bodystat API namespace for our Bodystat DLL wrapper. See main BodystatAPI class members for more information.

Files

Name	Description
BodystatSDK.cs (see page 80)	<p>Bodystat SDK - C# Wrapper (C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p> <p>Version For version information please refer to main BodystatSDK API</p>

1.2.1 Bodystat API

This is the Bodystat API namespace for our Bodystat DLL wrapper.

See main BodystatAPI class members for more information.

Description

Main documentation for the functions, definitions, classes and structures provided by the Bodystat API.

Module

C# .NET Wrapper (see page 56)

1.2.1.1 BodystatAPI.BodystatAPI.BSAutoSetupBT@out string@bool@IntPtr

Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown).

This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc).

Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the device knowing only the com port. It is not necessary to setup the device in this way before each use.

New non-paired units will only be found if switched on, in range and ready for communication (on main startup screen). It is advisable to instruct the user to ensure this is the case before calling.

Parameters

Parameters	Description
lpComPort	[out] String buffer to receive the device port the device is assigned to.
bReportErrors	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (still not completely silent).
hWndParent	Handle of the parent window.

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or times out.

Errors during the process can be optionally shown/hidden. However, a completely silent approach is not possible using this function (instead use `BSSearchBTDevices` and `BSGetBTBodystatDevice` or iterate radio handles manually. Then pass desired handle to `BSAuthenticateBTDevice` for pairing).

Finally be aware this function is only available for supported locales. Otherwise English messages may be displayed to the user. Adopt silent approach discussed previously if you are targeting a non-supported locale.

1.2.1.2

BodystatAPI.BodystatAPI.BSAutoSetupBT@StringBuilder@int@bool@IntPtr

Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown).

This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc).

Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the device knowing only the com port. It is not necessary to setup the device in this way before each use.

New non-paired units will only be found if switched on, in range and ready for communication (on main startup screen). It is advisable to instruct the user to ensure this is the case before calling.

Parameters

Parameters	Description
lpComPort	[out] String buffer to receive the device port the device is assigned to.
iSize	The size of the buffer.
bReportErrors	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (still not completely silent).
hWndParent	Handle of the parent window.

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or times out.

Errors during the process can be optionally shown/hidden. However, a completely silent approach is not possible using this function (instead use `BSSearchBTDevices` and `BSGetBTBodystatDevice` or iterate radio handles manually. Then pass desired handle to `BSAuthenticateBTDevice` for pairing).

Finally be aware this function is only available for supported locales. Otherwise English messages may be displayed to the user. Adopt silent approach discussed previously if you are targeting a non-supported locale.

1.2.1.3

BodystatAPI.BodystatAPI.BSCalculateNormals @BSMeasurement @BSResults @out BSNormals

Calculate normal values and normal ranges for a given subject and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from BSCalculateResults).

Parameters

Parameters	Description
pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
pR	The results calculated for this measurement previously.
pN	[out] Receives the normal values and normal ranges for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Normals vary by subject attributes, measured bio-impedance values and calculated results.

1.2.1.4

BodystatAPI.BodystatAPI.BSCalculateResults @BSMeasurement @out BSResults @BSAnalysisMode

Calculate complete result set for a given bio-impedance measurement recorded by the Bodystat device.

Parameters

Parameters	Description
pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
pR	[out] Receives the complete calculated results for this measurement.
iAnalysisMode	The desired analysis mode (see BSAnalysisMode for more information).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Results vary by subject attributes and measured bio-impedance values.

1.2.1.5 BodystatAPI.BodystatAPI.BSCloseComport

Close com port.

1.2.1.6 BodystatAPI.BodystatAPI.BSConnect

Connect to the Bodystat device and ensure the device is responding.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have opened the com port by calling BSOpenComPort previously. This connect command must be issued before sending other any other commands.

1.2.1.7 BodystatAPI.BodystatAPI.BSGetBTStackInfo@out string

Retrieve information about the Bluetooth stack. Useful for diagnostics.

Parameters

Parameters	Description
lpBuffer	[out] String buffer to receive the information.

Returns

true if it succeeds, false if it fails.

1.2.1.8

BodystatAPI.BodystatAPI.BSGetBTStackInfo@StringBuilder@int

Retrieve information about the Bluetooth stack. Useful for diagnostics.

Parameters

Parameters	Description
lpBuffer	[out] String buffer to receive the information.
iBufferSize	Size of the buffer.

Returns

true if it succeeds, false if it fails.

Remarks

Note: Information can only be provided for supported BT radios.

1.2.1.9

BodystatAPI.BodystatAPI.BSGetDeviceModelName@BSDeviceModel@out string

Determine the full product name of a specific model.

Parameters

Parameters	Description
iModel	The desired model of interest.
lpBuffer	[out] String buffer to receive the name.
iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.2.1.10**BodystatAPI.BodystatAPI.BSGetDeviceModelName@BSDeviceModel@StringBuilder@int**

Determine the full product name of a specific model.

Parameters

Parameters	Description
iModel	The desired model of interest.
lpBuffer	[out] String buffer to receive the name.
iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.2.1.11**BodystatAPI.BodystatAPI.BSGetDeviceModelNameShort@BSDeviceModel@out string**

Determine the short product name of a specific model.

Parameters

Parameters	Description
iModel	The desired model of interest.
lpBuffer	[out] String buffer to receive the name.
iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.2.1.12**BodystatAPI.BodystatAPI.BSGetDeviceModelNameShort@BSDeviceModel@StringBuilder@int**

Determine the short product name of a specific model.

Parameters

Parameters	Description
iModel	The desired model of interest.
lpBuffer	[out] String buffer to receive the name.
iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.2.1.13 BodystatAPI.BodystatAPI.BSGetSdkLibraryVersion@out int@out int

Determine the Bodystat SDK library (DLL) version.

Parameters

Parameters	Description
pMajorSdkVer	[out] Receives the major version number of the Bodystat SDK DLL.
pMinorSdkVer	[out] Receives the minor version number of the Bodystat SDK DLL.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Caller should check the version number of the DLL if utilising functionality noted only as being available in later revisions of the SDK. Or, alternatively, ensure the SDK DLL is updated appropriately during installation of your application.

1.2.1.14 BodystatAPI.BodystatAPI.BSIsBTAvailable

Check if a supported BT radio is available on this computer. Checks for presence of USB Bluetooth dongles and inbuilt Bluetooth modules in the PC. No communication with Bodystat devices takes place during this query.

Returns

Returns true if supported radio is available and enabled.

Remarks

Note: Only supported BT radios can be queried (presently those compatible with the Microsoft Bluetooth APIs and some Widdcomm models).

Also note, Bluetooth must be enabled on the PC by the user to permit detection. Some PCs, laptops particularly, have physical switches to enable bluetooth/wireless communication.

1.2.1.15 BodystatAPI.BodystatAPI.BSNormals.AerobicCapacity_H

High for Average grade in Aerobic Capacity normals (mlsO2/kg/min)

1.2.1.16 BodystatAPI.BodystatAPI.BSNormals.AerobicCapacity_L

Low for Average grade in Aerobic Capacity normals (mlsO2/kg/min)

1.2.1.17 BodystatAPI.BodystatAPI.BSNormals.BFMI_H

Normal range for Body Fat Mass Index (MDD/Quad only) - (high end)

1.2.1.18 BodystatAPI.BodystatAPI.BSNormals.BFMI_L

Normal range for Body Fat Mass Index (MDD/Quad only) - (low end)

1.2.1.19 BodystatAPI.BodystatAPI.BSNormals.BMI_H

Normal range for Body Mass Index - (high end)

1.2.1.20 BodystatAPI.BodystatAPI.BSNormals.BMI_L

Normal range for Body Mass Index - (low end)

1.2.1.21 BodystatAPI.BodystatAPI.BSNormals.Cholesterol_HDL_H

Cholesterol HDL high range in mg/dl

1.2.1.22 BodystatAPI.BodystatAPI.BSNormals.Cholesterol_HDL_L

Cholesterol HDL low range in mg/dl

1.2.1.23 BodystatAPI.BodystatAPI.BSNormals.Cholesterol_LDL_H

Cholesterol VDL high range in mg/dl

1.2.1.24 BodystatAPI.BodystatAPI.BSNormals.Cholesterol_LDL_L

Cholesterol VDL low range in mg/dl

1.2.1.25

BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Total_H

Cholesterol total high range in mg/dl

1.2.1.26

BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Total_L

Cholesterol total low range in mg/dl

1.2.1.27

BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Triglycerides_H

Triglycerides high range in mg/dl

1.2.1.28

BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Triglycerides_L

Triglycerides low range in mg/dl

1.2.1.29

BodystatAPI.BodystatAPI.BSNormals.Cholesterol_VLDL_H

Cholesterol VLDL high range in mg/dl

1.2.1.30

BodystatAPI.BodystatAPI.BSNormals.Cholesterol_VLDL_L

Cholesterol VLDL low range in mg/dl

1.2.1.31 BodystatAPI.BodystatAPI.BSNormals.ECWPerC_Norm

Normal value for ECW Percentage (MDD/Quad only)

1.2.1.32 BodystatAPI.BodystatAPI.BSNormals.FatKg_H

Normal range for fat (in kg) - (high end)

1.2.1.33 BodystatAPI.BodystatAPI.BSNormals.FatKg_L

Normal range for fat (in kg) - (low end)

1.2.1.34 BodystatAPI.BodystatAPI.BSNormals.FatPerc_H

Normal range for fat percentage - (high end)

1.2.1.35 BodystatAPI.BodystatAPI.BSNormals.FatPerc_L

Normal range for fat percentage - (low end)

1.2.1.36 BodystatAPI.BodystatAPI.BSNormals.FFMI_H

Normal range for Fat Free Mass Index (MDD/Quad only) - (high end)

1.2.1.37 BodystatAPI.BodystatAPI.BSNormals.FFMI_L

Normal range for Fat Free Mass Index (MDD/Quad only) - (low end)

1.2.1.38 BodystatAPI.BodystatAPI.BSNormals.Flexibility_H

Flexibility high range of normal in cm

1.2.1.39 BodystatAPI.BodystatAPI.BSNormals.Flexibility_L

Flexibility low range of normal in cm

1.2.1.40 BodystatAPI.BodystatAPI.BSNormals.GlucoseFasting_H

Glucose (mg/dL)

1.2.1.41 BodystatAPI.BodystatAPI.BSNormals.GlucoseFasting_L

Glucose (mg/dL)

1.2.1.42

BodystatAPI.BodystatAPI.BSNormals.GlucoseNonFasting_H

Glucose (mg/dL)

1.2.1.43

BodystatAPI.BodystatAPI.BSNormals.GlucoseNonFasting_L

Glucose (mg/dL)

1.2.1.44 BodystatAPI.BodystatAPI.BSNormals.GripLeft

Grip left normal value in kg

1.2.1.45 BodystatAPI.BodystatAPI.BSNormals.GripRight

Grip right normal value in kg

1.2.1.46 BodystatAPI.BodystatAPI.BSNormals.ICWPerc_Norm

Normal value for ICW Percentage (MDD/Quad only)

1.2.1.47 BodystatAPI.BodystatAPI.BSNormals.Illness_H

Normal range for Illness Marker (TM) (Quad only) - (high end)

1.2.1.48 BodystatAPI.BodystatAPI.BSNormals.Illness_L

Normal range for Illness Marker (TM) (Quad only) - (low end)

1.2.1.49 BodystatAPI.BodystatAPI.BSNormals.LeanKg_H

Normal range for lean (in kg) - (high end)

1.2.1.50 BodystatAPI.BodystatAPI.BSNormals.LeanKg_L

Normal range for lean (in kg) - (low end)

1.2.1.51 BodystatAPI.BodystatAPI.BSNormals.LeanPerc_H

Normal range for lean percentage - (high end)

1.2.1.52 BodystatAPI.BodystatAPI.BSNormals.LeanPerc_L

Normal range for lean percentage - (low end)

1.2.1.53 BodystatAPI.BodystatAPI.BSNormals.Nutrition_Norm

Normal value for Nutrition Index (Quad only)

1.2.1.54 BodystatAPI.BodystatAPI.BSNormals.TBW_H

Normal range for Total Body Water (in litres) - (high end)

1.2.1.55 BodystatAPI.BodystatAPI.BSNormals.TBW_L

Normal range for Total Body Water (in litres) - (low end)

1.2.1.56 BodystatAPI.BodystatAPI.BSNormals.TBWPerc_H

Normal range for Total Body Water percentage - (high end)

1.2.1.57 BodystatAPI.BodystatAPI.BSNormals.TBWPerc_L

Normal range for Total Body Water percentage - (low end)

1.2.1.58 BodystatAPI.BodystatAPI.BSNormals.TotalWeight_H

Normal range for Total Weight (in kg) - (high end)

1.2.1.59 BodystatAPI.BodystatAPI.BSNormals.TotalWeight_L

Normal range for Total Weight (in kg) - (low end)

1.2.1.60 BodystatAPI.BodystatAPI.BSNormals.TotalWeightMethod

Method used when calculating Total Weight normal range (0 = Composition Method, 1 = BMI Method)

1.2.1.61 BodystatAPI.BodystatAPI.BSNormals.WaistHip_Norm

Normal value for Waist/Hip ratio

1.2.1.62 BodystatAPI.BodystatAPI.BSNormals.Wellness_H

Normal range for Wellness Marker (TM) (MDD only) - (high end)

1.2.1.63 BodystatAPI.BodystatAPI.BSNormals.Wellness_L

Normal range for Wellness Marker (TM) (MDD only) - (low end)

1.2.1.64 BodystatAPI.BodystatAPI.BSOpenComport@string

Open com port to the device ready for communication.

Parameters

Parameters	Description
lpComPort	String buffer containing the device/port to be opened.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

The port must be closed once you are finished sending commands. To save battery power on the device, it is recommended you chain command sequences together and open/close the port immediately before and after the chain sequence.

1.2.1.65

BodystatAPI.BodystatAPI.BSOpenComport@StringBuilder@int

Open com port to the device ready for communication.

Parameters

Parameters	Description
lpComPort	String buffer containing the device/port to be opened.
iSize	The size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

The port must be closed once you are finished sending commands. To save battery power on the device, it is recommended you chain command sequences together and open/close the port immediately before and after the chain sequence.

1.2.1.66 BodystatAPI.BodystatAPI.BSRawData.BSRawData@int

Override the constructor to initialise required variables in the structure. You must ensure you allocate the structure in a manner which calls this function.

Parameters

Parameters	Description
	The parameter is not required by the wrapper and is ignored.
iDontCare	=0

1.2.1.67 BodystatAPI.BodystatAPI.BSReadCalibrationTime@out DateTime

Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.

Parameters

Parameters	Description
dtCalibrationTime	[out] Receives the date the device was last calibrated.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.68 BodystatAPI.BodystatAPI.BSReadCalibrationTime@out Int64

Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.

Parameters

Parameters	Description
pCalibrationTime	[out] Receives the date the device was last calibrated.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.69 BodystatAPI.BodystatAPI.BSReadCurrentTime@out DateTime@int

Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).

Parameters

Parameters	Description
dtCurrentTime	[out] Receives the current time of the device's internal clock.
iDST	Daylight Saving Time handling (as per Microsoft time functions) Zero (0) to indicate that standard time is in effect. A value greater than 0 to indicate that daylight saving time is in effect. A value less than zero to have the C run-time library code compute whether standard time or daylight saving time is in effect. (default) = -1

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.70 BodystatAPI.BodystatAPI.BSReadCurrentTime@out Int64@int

Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).

Parameters

Parameters	Description
pCurrentTime	[out] Receives the current time of the device's internal clock.
iDST	Daylight Saving Time handling (as per Microsoft time functions) Zero (0) to indicate that standard time is in effect. A value greater than 0 to indicate that daylight saving time is in effect. A value less than zero to have the C run-time library code compute whether standard time or daylight saving time is in effect. (default) = -1

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.71 BodystatAPI.BodystatAPI.BSReadModelVersion@out BSDeviceModel@out byte@out byte@out byte@out byte@out byte

Read device model and firmware version of the connected Bodystat device.

Parameters

Parameters	Description
pModel	[out] Receives the device model.
pMajorVersion	[out] Receives the major version of the main firmware in the device.
pMinorVersion	[out] Receives the minor version of the main firmware in the device.
pPsoc2Version	[out] Receives the version of Psoc2 firmware in the device.
pEepromVersion	[out] Receives the version of the EEPROM firmware in the device.
pBluetoothInfo	[out] Receives the version of the Bluetooth module in the device.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.72 BodystatAPI.BodystatAPI.BSReadPrinterAddress@out string

Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).

Parameters

Parameters	Description
lpPrinterAddress	[out] String buffer containing the Bluetooth address (BDA) of the printer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.2.1.73

BodystatAPI.BodystatAPI.BSReadPrinterAddress@StringBuilder@int

Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).

Parameters

Parameters	Description
lpPrinterAddress	[out] String buffer which receives the printer address.

iPrinterAddressBufferSize	Size of the buffer.
---------------------------	---------------------

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.74 BodystatAPI.BodystatAPI.BSReadProtocolInfo@out byte@out byte@out byte

Read protocol information from the connected device. Contains version numbers relating to communication protocols in use, data structures and auxiliary information.

Parameters

Parameters	Description
pProtocolVersion	[out] Receives the communication protocol version.
pDataVersion	[out] Receives the data structure version.
pAuxInfo	[out] Receives auxiliary information.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Generally not required for the majority of SDK users as these details are abstracted by the SDK. You must have connected to the device by calling BSConnect previously.

1.2.1.75 BodystatAPI.BodystatAPI.BSReadSerialNumber@out uint

Read the unique serial number of the connected Bodystat device.

Parameters

Parameters	Description
pSerialNumber	[out] Receives the serial number of the device.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously.

1.2.1.76 BodystatAPI.BodystatAPI.BSReadTestBodystat@out int@out int@out int@out int@out int@out float@out float

Perform a quick bio-impedance test measurement on the connected bodystat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.

May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.

Useful for diagnostics or performing measurements instigated by software.

Parameters

Parameters	Description
pZ5	[out] Receives the impedance measured at 5 kHz (result in ohms).
pZ50	[out] Receives the impedance measured at 50 kHz (result in ohms).
pZ100	[out] Receives the impedance measured at 100 kHz (result in ohms).
pZ200	[out] Receives the impedance measured at 200 kHz (result in ohms).
pR50	[out] Receives the resistance measured at 50 kHz (result in ohms).
pX50	[out] Receives the reactance measured at 50 kHz (result in ohms).
pPA50	[out] Receives the phase angle measured at 50 kHz (result in degrees).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously. Caller must check they are working with a Bodystat (rather than Equistat model, Equistat models will respond without error - but with garbage results).

1.2.1.77 BodystatAPI.BodystatAPI.BSResults.BCM

Body cell mass (Quad only)

1.2.1.78 BodystatAPI.BodystatAPI.BSResults.BFMI

Body Fat Mass Index (MDD/Quad only)

1.2.1.79 BodystatAPI.BodystatAPI.BSResults.BMI

Body Mass Index

1.2.1.80 BodystatAPI.BodystatAPI.BSResults.BMR

Basal Metabolic Rate (in kcal)

1.2.1.81 BodystatAPI.BodystatAPI.BSResults.BMRkg

Basal Metabolic Rate per kilogram in (kcal/kg)

1.2.1.82 BodystatAPI.BodystatAPI.BSResults.DryLW

Dry Lean Weight (in kg)

1.2.1.83 BodystatAPI.BodystatAPI.BSResults.ECW

Extra Cellular Water (in litres) (MDD/Quad only)

1.2.1.84 BodystatAPI.BodystatAPI.BSResults.ECWPerc

Extra Cellular Water percentage (MDD/Quad only)

1.2.1.85 BodystatAPI.BodystatAPI.BSResults.EstAvg

Estimated Average Requirement (in kcal)

1.2.1.86 BodystatAPI.BodystatAPI.BSResults.FatKg

Fat (in kg)

1.2.1.87 BodystatAPI.BodystatAPI.BSResults.FatPerc

Fat percentage

1.2.1.88 BodystatAPI.BodystatAPI.BSResults.FFMI

Fat Free Mass Index (MDD/Quad only)

1.2.1.89 BodystatAPI.BodystatAPI.BSResults.ICW

Intra Cellular Water (in litres) (MDD/Quad only)

1.2.1.90 BodystatAPI.BodystatAPI.BSResults.ICWPerc

Intra Cellular Water percentage (MDD/Quad only)

1.2.1.91 BodystatAPI.BodystatAPI.BSResults.Illness

Prediction Marker (TM) - formerly known as Illness Marker (TM) - (Quad only)

1.2.1.92 BodystatAPI.BodystatAPI.BSResults.LeanKg

Lean (in kg)

1.2.1.93 BodystatAPI.BodystatAPI.BSResults.LeanPerc

Lean percentage

1.2.1.94 BodystatAPI.BodystatAPI.BSResults.Nutrition

Nutrition index (Quad only)

1.2.1.95 BodystatAPI.BodystatAPI.BSResults.TBW

Water (in litres) or Total Body Water (in litres)

1.2.1.96 BodystatAPI.BodystatAPI.BSResults.TBWPerc

Water percentage or Total Body Water percentage

1.2.1.97 BodystatAPI.BodystatAPI.BSResults.ThirdSpace

3rd space water (in litres) (Quad only)

1.2.1.98 BodystatAPI.BodystatAPI.BSResults.TotalWeight

Total Weight (in kg)

1.2.1.99 BodystatAPI.BodystatAPI.BSResults.WaistHip

Waist/Hip ratio

1.2.1.100 BodystatAPI.BodystatAPI.BSResults.Wellness

Wellness Marker (TM) (MDD only)

1.2.1.101 BodystatAPI.BodystatAPI.BSWriteCurrentTime

Reset the current date/time of the internal clock on the connected device. Date/time is automatically set to match the current date/time of the PC

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously. Locale Daylight Saving Time (DST) rules are observed. Use sister function BSWriteCurrentTime(__time64_t CurrentTime) to set a specific date/time which differs from the PC clock or override DST behaviour.

1.2.1.102

BodystatAPI.BodystatAPI.BSWritePrinterAddress@string@int

Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt

printer pairing function in the device.

Parameters

Parameters	Description
lpPrinterAddress	String buffer containing the Bluetooth address (BDA) of the printer.
iPrinterAddressBufferSize	Size of the printer address buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

1.2.1.103

BodystatAPI.BodystatAPI.BSWritePrinterAddress@StringBuilder@int

Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt printer pairing function in the device.

Parameters

Parameters	Description
lpPrinterAddress	String buffer containing the Bluetooth address (BDA) of the printer.
iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Note not all device models support this function (not all support direct printing).

1.2.1.104

BodystatAPI.BodystatAPI.ESCalculateNormals@ESMeasurement@ESResults@outESNormals

Calculate normal values and normal ranges for a given animal and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from ESCalculateResults).

Parameters

Parameters	Description
pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
pR	The results calculated for this measurement previously.
pN	[out] Receives the normal values and normal ranges for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Normals vary by input attributes, measured bio-impedance values and calculated results.

1.2.1.105**BodystatAPI.BodystatAPI.ESCalculateResults@ESMeasurement@out ESResults**

Calculate complete result set for a given bio-impedance measurement recorded by the Equistat device.

Parameters

Parameters	Description
pM	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
pR	[out] Receives the complete calculated results for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Results vary by input attributes and measured bio-impedance values.

1.2.1.106 BodystatAPI.BodystatAPI.ESRawData.ESRawData@int

Override the constructor to initialise required variables in the structure. You must ensure you allocate the structure in a manner which calls this function.

Parameters

Parameters	Description
iDontCare	=0

Description

The parameter is not required by the wrapper and is ignored.

1.2.1.107 BodystatAPI.BodystatAPI.ESReadTestEquistat@out int@out int@out int@out int@out int@out int@out float

Perform a quick bio-impedance test measurement on the connected Equistat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.

May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.

Useful for diagnostics or performing measurements instigated by software.

Parameters

Parameters	Description
pZ5	[out] Receives the impedance measured at 5 kHz (result in ohms).
pZ16	[out] Receives the impedance measured at 16 kHz (result in ohms).
pZ24	[out] Receives the impedance measured at 24 kHz (result in ohms).
pZ50	[out] Receives the impedance measured at 50 kHz (result in ohms).
pZ140	[out] Receives the impedance measured at 140 kHz (result in ohms).
pZ200	[out] Receives the impedance measured at 200 kHz (result in ohms).
pZ280	[out] Receives the impedance measured at 280 kHz (result in ohms).
pPA50	[out] Receives the phase angle measured at 50 kHz (result in degrees) (ES Pro Only).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect previously. Caller must check they are working with an Equistat (rather than Bodystat model, Bodystat models will respond without error - but with garbage results).

1.2.1.108 BodystatAPI.BodystatAPI.BSGetBTBodystatDevice@out string@out string@out string@UInt16

Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices.

Parameters

Parameters	Description
lpDeviceName	[out] String buffer to receive the product name of the Bodystat device.
lpBDA	[out] String buffer to receive the Bluetooth hexadecimal address (known as the BDA - unique to each BT module in a device, like a MAC address).
lpComPort	[out] String buffer to receive the device port the device is assigned to.
iTimeout	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds). = 7

Returns

TRUE if it succeeds, FALSE if it fails.

1.2.1.109

BodystatAPI.BodystatAPI.BSGetBTBodystatDevice@StringBuilder@int@StringBuilder@int@StringBuilder@int@UInt16

Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices.

Parameters

Parameters	Description
lpDeviceName	[out] String buffer to receive the product name of the Bodystat device.
iDeviceNameBufferSize	Size of the device name buffer.
lpBDA	[out] String buffer to receive the Bluetooth hexadecimal address (known as the BDA - unique to each BT module in a device, like a MAC address).
iBDABufferSize	Size of the BDA buffer.
lpComPort	[out] String buffer to receive device/port the unit is assigned to.
iComBufferSize	Size of the com port buffer.
iTimeout	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds). = 7

Returns

TRUE if it succeeds, FALSE if it fails.

1.2.1.110 BodystatAPI.BodystatAPI.BSSearchBTDevices@out int@out int@bool@IntPtr@UInt16

Search BT radio for any Bodystat devices (paired or otherwise). Existing paired devices will be counted even if they are switched off or out of range. New non-paired units will only be counted if switched on, in range and ready for communication (on main startup screen).

Parameters

Parameters	Description
iNewDevices	[out] The number of new Bodystat devices found (not paired).
iAuthenticatedDevices	[out] The number of authenticated Bodystat devices (paired).
bReportErrors	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (completely silent).
hWndParent	Handle of the parent window.

iTimeout	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds). = 7
----------	--

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or times out.

1.2.1.111

BodystatAPI.BodystatAPI.BSUnAuthenticateBTDevices@bool@IntPtr@UInt16

Unauthenticate (unpair) ALL Bodystat devices from this PC.

Parameters

Parameters	Description
bReportErrors	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user.
hWndParent	Handle of the parent window.
iTimeout	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds). = 7

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or times out.

1.2.2 Files

The following table lists files in this documentation.

Files

Name	Description
BodystatSDK.cs (📄 see page 80)	<p>Bodystat SDK - C# Wrapper (C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p> <p>Version For version information please refer to main BodystatSDK API</p>

Module

C# .NET Wrapper (📄 see page 56)

1.2.2.1 BodystatSDK.cs

Bodystat SDK - C# Wrapper

(C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved.

Use subject to license.

Web: <http://www.bodystat.com>

Email: info@bodystat.com

Tel: +44 (0)1624 629 571

Fax: +44 (0)1624 611 544

Version

For version information please refer to main BodystatSDK API

Namespaces


Name	Description
BodystatAPI (📄 see page 57)	<p>This is the Bodystat API namespace for our Bodystat DLL wrapper.</p> <p>See main BodystatAPI class members for more information.</p>

1.3 Visual Basic .NET wrapper






Bodystat API Visual Basic .NET wrapper.

Provides a wrapper interface to the Bodystat DLL with VB calls and types you can work with directly in your Visual Basic .NET program. The wrapper handles the necessary marshalling and interop issues required to interface with the native unmanaged Bodystat DLL.

Classes

	Name	Description
	BodystatAPI (see page 81)	This is the main class for our Bodystat DLL wrapper.


Structs, Records, Enums

	Name	Description
	BSAnalysisMode (see page 145)	BSAnalysisMode: Analysis results are available in several formats: BSAnalysisModeBC= Body composition analysis BSAnalysisModeHN= Hydration analysis BSAnalysisModeBoth= Combined body composition and hydration analysis (default)
	BSDeviceFamily (see page 146)	BSDeviceFamily: Enumerated values that group each device / model variant into a given product family. Devices within a given product family share broad capabilities and features. For details regarding capabilities and features refer to the technical documentation for that product family (http://www.bodystat.com)
	BSDeviceModel (see page 147)	BSDeviceModel: Enumerated values that represent the various different models of the Bodystat devices over the years. Suffix: The suffix after the model indicates specific hardware variations regarding connectivity: OPTO indicates model communicates via opto-isolated serial. DIU indicates model communicates via a data interface unit (induction interface) BT indicates model communicates via Bluetooth serial port.
	BSError (see page 149)	BSError: Bodystat Error codes. A list of Bodystat specific error codes.
	BSGender (see page 151)	BSGender: Enumerated values that represent a person's gender.

1.3.1 Classes

The following table lists classes in this documentation.

Classes

	Name	Description
	BodystatAPI (see page 81)	This is the main class for our Bodystat DLL wrapper.

Module

Visual Basic .NET wrapper ([see page 80](#))

1.3.1.1 BodystatAPI Class

This is the main class for our Bodystat DLL wrapper.

Class Hierarchy

```
BodystatAPI
└─ BodystatAPI
```

C++

```
class BodystatAPI;
```

C#

```
public class BodystatAPI;
```

Visual Basic

```
Public Class BodystatAPI
```

Java

```
public class BodystatAPI;
```



MATLAB

```
BodystatAPI
```



























File

















```
BodystatSDK.vb
```

BodystatAPI Fields

	Name	Description
	BS_RAWDATA_ARRAYSIZE (see page 114)	Number of measurement records BSRawData (see page 96) structure should reserved space for (currently fixed)
	ES_RAWDATA_ARRAYSIZE (see page 115)	Number of measurement records ESRawData (see page 109) structure should reserved space for (currently fixed)

BodystatAPI Methods

	Name	Description
 	BSAuthenticateBTDevice (see page 115)	Authenticate (pair) a Bodystat device found during a previous call to BSSearchBTDevices (see page 136). Completely silent implementation.
 	BSAutoSetupBT (see page 116)	Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown). This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc). Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the... more (see page 116)
 	BSCalculateNormals (see page 117)	Calculate normal values and normal ranges for a given subject and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from BSCalculateResults (see page 118)).
 	BSCalculateResults (see page 118)	Calculate complete result set for a given bio-impedance measurement recorded by the Bodystat device.
 	BSCloseComport (see page 119)	Close com port.
 	BSConnect (see page 120)	Connect to the Bodystat device and ensure the device is responding.
 	BSGetBTBodystatDevice (see page 120)	Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices (see page 136).
 	BSGetBTStackInfo (see page 121)	Retrieve information about the Bluetooth stack. Useful for diagnostics.
 	BSGetDeviceFamily (see page 122)	Determine the device family of a specific model.
 	BSGetDeviceModelName (see page 123)	Determine the full product name of a specific model.
 	BSGetDeviceModelNameShort (see page 124)	Determine the short product name of a specific model.
 	BSGetSdkLibraryVersion (see page 124)	Determine the Bodystat SDK library (DLL) version.
 	BSIsBTAvailable (see page 125)	Check if a supported BT radio is available on this computer. Checks for presence of USB Bluetooth dongles and inbuilt Bluetooth modules in the PC. No communication with Bodystat devices takes place during this query.

	BSOpenComport (see page 126)	Open com port to the device ready for communication.
	BSReadCalibrationTime (see page 127)	Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.
	BSReadCurrentTime (see page 128)	Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).
	BSReadModelVersion (see page 130)	Read device model and firmware version of the connected Bodystat device.
	BSReadPrinterAddress (see page 131)	Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).
	BSReadProtocolInfo (see page 132)	Read protocol information from the connected device. Contains version numbers relating to communication protocols in use, data structures and auxiliary information.
	BSReadSerialNumber (see page 133)	Read the unique serial number of the connected Bodystat device.
	BSReadStoredTestData (see page 133)	Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged. Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.
	BSReadTestBodystat (see page 134)	Perform a quick bio-impedance test measurement on the connected bodystat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory. May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion. Useful for diagnostics or performing measurements instigated by software.
	BSSearchBTDevices (see page 136)	Search BT radio for any Bodystat devices (paired or otherwise). Existing paired devices will be counted even if they are switched off or out of range. New non-paired units will only be counted if switched on, in range and ready for communication (on main startup screen).
	BSUnAuthenticateBTDevices (see page 137)	Unauthenticate (unpair) ALL Bodystat devices from this PC.
	BSWriteCurrentTime (see page 138)	Reset the current date/time of the internal clock on the connected device. Date/time is automatically set to match the current date/time of the PC
	BSWriteCurrentTimeAs (see page 138)	Reset the current date/time of the internal clock on the connected device. Date/time is set to a specific value of your choosing.
	BSWritePrinterAddress (see page 140)	Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt printer pairing function in the device.
	ESCalculateNormals (see page 141)	Calculate normal values and normal ranges for a given animal and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from ESCalculateResults (see page 142)).
	ESCalculateResults (see page 142)	Calculate complete result set for a given bio-impedance measurement recorded by the Equistat device.

	ESReadStoredTestData (see page 142)	Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged. Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.
	ESReadTestEquistat (see page 143)	Perform a quick bio-impedance test measurement on the connected Equistat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory. May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion. Useful for diagnostics or performing measurements instigated by software.

BodystatAPI Structures

	Name	Description
	BSMeasurement (see page 84)	BSMeasurement: Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).
	BSNormals (see page 93)	BSNormals: Bodystat Normals structure. Holds the normal values or normal ranges for any given test measurement.
	BSRawData (see page 96)	BSRawData: Bodystat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record (see page 98) starting number.
	BSResults (see page 99)	BSResults: Bodystat Results structure. Holds the calculated results for any given test measurement.
	ESMeasurement (see page 101)	ESMeasurement: Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).
	ESNormals (see page 108)	ESNormals: Equistat Normals structure. Holds the normal values or normal ranges for any given test measurement.
	ESRawData (see page 109)	ESRawData: Equistat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record (see page 112) starting number.
	ESResults (see page 112)	ESResults: Equistat Results structure. Holds the calculated results for any given test measurement.
	ImpData (see page 114)	Represents one frequency measurement

1.3.1.1.1 BodystatAPI Structures

1.3.1.1.1.1 BodystatAPI.BSMeasurement Structure

BSMeasurement: Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).

C++

```
[System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)]
struct BSMeasurement {
    BSDeviceModel iDeviceModel;
    UInteger ulDeviceSerialNumber;
    Int64 llTestDate;
    BSGender iGender;
    Integer iAge;
    Integer iHeight;
    Single fWeight;
    Integer iActivity;
    Integer iWaist;
    Integer iHip;
    Integer iZ_5kHz;
    Integer iZ_50kHz;
    Integer iZ_100kHz;
    Integer iZ_200kHz;
    Integer iR_50kHz;
    Single fX_50kHz;
    Single fPA_50kHz;
    Integer iFrequencies;
    ImpData pMultifreqData;
};
```

C#

```
[System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)]
public struct BSMeasurement {
    public BSDeviceModel iDeviceModel;
    public UInteger ulDeviceSerialNumber;
    public Int64 llTestDate;
    public BSGender iGender;
    public Integer iAge;
    public Integer iHeight;
    public Single fWeight;
    public Integer iActivity;
    public Integer iWaist;
    public Integer iHip;
    public Integer iZ_5kHz;
    public Integer iZ_50kHz;
    public Integer iZ_100kHz;
    public Integer iZ_200kHz;
    public Integer iR_50kHz;
    public Single fX_50kHz;
    public Single fPA_50kHz;
    public Integer iFrequencies;
    public ImpData() pMultifreqData;
}
```

Visual Basic

```
<System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)>
Public Structure BSMeasurement
    Public iDeviceModel As BSDeviceModel
    Public ulDeviceSerialNumber As UInteger
    Public llTestDate As Int64
    Public iGender As BSGender
    Public iAge As Integer
    Public iHeight As Integer
    Public fWeight As Single
    Public iActivity As Integer
    Public iWaist As Integer
    Public iHip As Integer
    Public iZ_5kHz As Integer
    Public iZ_50kHz As Integer
    Public iZ_100kHz As Integer
    Public iZ_200kHz As Integer
    Public iR_50kHz As Integer
    Public fX_50kHz As Single
```



```
Public fPA_50kHz As Single
Public iFrequencies As Integer
Public pMultifreqData As ImpData()
End Structure
```

Java

```
[System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)]
public class BSMeasurement;
```




















MATLAB

```
BSMeasurement
```


File

```
BodystatSDK.vb
```

BSMeasurement Fields

	Name	Description
	fPA_50kHz (see page 86)	Phase angle measured at 50 kHz (result in degrees) (MDD/Quad only)
	fWeight (see page 87)	Inputted weight of the subject (kilograms)
	fX_50kHz (see page 87)	Reactance measured at 50 kHz (result in ohms) (MDD/Quad only)
	iActivity (see page 87)	Inputted physical activity level of the subject (range 1 to 5: 1=very low 2=low/medium 3=medium 4=medium/high 5=very high)
	iAge (see page 88)	Inputted age of the subject (years)
	iDeviceModel (see page 88)	Device model which performed the measurement
	iFrequencies (see page 88)	Number of frequencies used 11/50
	iGender (see page 89)	Inputted gender of the subject
	iHeight (see page 89)	Inputted height of the subject (centimeters)
	iHip (see page 89)	Inputted hip size of the subject (centimeters)
	iR_50kHz (see page 90)	Resistance measured at 50 kHz (result in ohms) (MDD/Quad only)
	iWaist (see page 90)	Inputted waist size of the subject (centimeters)
	iZ_100kHz (see page 90)	Impedance measured at 100 kHz (result in ohms) (Quad only)
	iZ_200kHz (see page 91)	Impedance measured at 200 kHz (result in ohms) (Quad only)
	iZ_50kHz (see page 91)	Impedance measured at 50 kHz (result in ohms)
	iZ_5kHz (see page 91)	Impedance measured at 5 kHz (result in ohms) (MDD/Quad only)
	IITestDate (see page 92)	Date & time the test measurement was performed (Windows 64-bit time structure) (MDD/Quad only)
	pMultifreqData (see page 92)	Pointer to ImpData (see page 114) array
	ulDeviceSerialNumber (see page 92)	Unique serial number of the device which performed the measurement

BSMeasurement Properties

	Name	Description
	tTestDate (see page 93)	Convert underlying 64-bit time structure to/from Visual Basic Date/Time

1.3.1.1.1.1.1 BSMeasurement Fields**1.3.1.1.1.1.1.1 BodystatAPI.BSMeasurement.fPA_50kHz Field****C++**

```
Single fPA_50kHz;
```

C#

```
public Single fPA_50kHz;
```

Visual Basic

```
Public fPA_50kHz As Single
```

Java

```
public Single fPA_50kHz;
```

MATLAB

```
fPA_50kHz
```

Description

Phase angle measured at 50 kHz (result in degrees) (MDD/Quad only)

1.3.1.1.1.1.2 BodystatAPI.BSMMeasurement.fWeight Field**C++**

```
Single fWeight;
```

C#

```
public Single fWeight;
```

Visual Basic

```
Public fWeight As Single
```

Java

```
public Single fWeight;
```

MATLAB

```
fWeight
```

Description

Inputted weight of the subject (kilograms)

1.3.1.1.1.1.3 BodystatAPI.BSMMeasurement.fX_50kHz Field**C++**

```
Single fX_50kHz;
```

C#

```
public Single fX_50kHz;
```

Visual Basic

```
Public fX_50kHz As Single
```

Java

```
public Single fX_50kHz;
```

MATLAB

```
fX_50kHz
```

Description

Reactance measured at 50 kHz (result in ohms) (MDD/Quad only)

1.3.1.1.1.1.4 BodystatAPI.BSMMeasurement.iActivity Field**C++**

```
Integer iActivity;
```

C#

```
public Integer iActivity;
```

Visual Basic

```
Public iActivity As Integer
```

Java

```
public Integer iActivity;
```

MATLAB

```
iActivity
```

Description

Inputted physical activity level of the subject (range 1 to 5: 1=very low 2=low/medium 3=medium 4=medium/high 5=very high)

1.3.1.1.1.1.1.5 BodystatAPI.BSMMeasurement.iAge Field**C++**

```
Integer iAge;
```

C#

```
public Integer iAge;
```

Visual Basic

```
Public iAge As Integer
```

Java

```
public Integer iAge;
```

MATLAB

```
iAge
```

Description

Inputted age of the subject (years)

1.3.1.1.1.1.1.6 BodystatAPI.BSMMeasurement.iDeviceModel Field**C++**

```
BSDeviceModel iDeviceModel;
```

C#

```
public BSDeviceModel iDeviceModel;
```

Visual Basic

```
Public iDeviceModel As BSDeviceModel
```

Java

```
public BSDeviceModel iDeviceModel;
```

MATLAB

```
iDeviceModel
```

Description

Device model which performed the measurement

1.3.1.1.1.1.1.7 BodystatAPI.BSMMeasurement.iFrequencies Field**C++**

```
Integer iFrequencies;
```

C#

```
public Integer iFrequencies;
```

Visual Basic

```
Public iFrequencies As Integer
```

Java

```
public Integer iFrequencies;
```

MATLAB

```
iFrequencies
```

Description

Number of frequencies used 11/50

1.3.1.1.1.1.1.8 BodystatAPI.BSMMeasurement.iGender Field**C++**

```
BSGender iGender;
```

C#

```
public BSGender iGender;
```

Visual Basic

```
Public iGender As BSGender
```

Java

```
public BSGender iGender;
```

MATLAB

```
iGender
```

Description

Inputted gender of the subject

1.3.1.1.1.1.1.9 BodystatAPI.BSMMeasurement.iHeight Field**C++**

```
Integer iHeight;
```

C#

```
public Integer iHeight;
```

Visual Basic

```
Public iHeight As Integer
```

Java

```
public Integer iHeight;
```

MATLAB

```
iHeight
```

Description

Inputted height of the subject (centimeters)

1.3.1.1.1.1.1.10 BodystatAPI.BSMMeasurement.iHip Field**C++**

```
Integer iHip;
```

C#

```
public Integer iHip;
```

Visual Basic

```
Public iHip As Integer
```

Java

```
public Integer iHip;
```

MATLAB

```
iHip
```

Description

Inputted hip size of the subject (centimeters)

1.3.1.1.1.1.11 BodystatAPI.BSMMeasurement.iR_50kHz Field**C++**

```
Integer iR_50kHz;
```

C#

```
public Integer iR_50kHz;
```

Visual Basic

```
Public iR_50kHz As Integer
```

Java

```
public Integer iR_50kHz;
```

MATLAB

```
iR_50kHz
```

Description

Resistance measured at 50 kHz (result in ohms) (MDD/Quad only)

1.3.1.1.1.1.1.12 BodystatAPI.BSMMeasurement.iWaist Field**C++**

```
Integer iWaist;
```

C#

```
public Integer iWaist;
```

Visual Basic

```
Public iWaist As Integer
```

Java

```
public Integer iWaist;
```

MATLAB

```
iWaist
```

Description

Inputted waist size of the subject (centimeters)

1.3.1.1.1.1.1.13 BodystatAPI.BSMMeasurement.iZ_100kHz Field**C++**

```
Integer iZ_100kHz;
```

C#

```
public Integer iZ_100kHz;
```

Visual Basic

```
Public iZ_100kHz As Integer
```

Java

```
public Integer iZ_100kHz;
```

MATLAB

```
iZ_100kHz
```

Description

Impedance measured at 100 kHz (result in ohms) (Quad only)

1.3.1.1.1.1.14 BodystatAPI.BSMMeasurement.iZ_200kHz Field**C++**

```
Integer iZ_200kHz;
```

C#

```
public Integer iZ_200kHz;
```

Visual Basic

```
Public iZ_200kHz As Integer
```

Java

```
public Integer iZ_200kHz;
```

MATLAB

```
iZ_200kHz
```

Description

Impedance measured at 200 kHz (result in ohms) (Quad only)

1.3.1.1.1.1.15 BodystatAPI.BSMMeasurement.iZ_50kHz Field**C++**

```
Integer iZ_50kHz;
```

C#

```
public Integer iZ_50kHz;
```

Visual Basic

```
Public iZ_50kHz As Integer
```

Java

```
public Integer iZ_50kHz;
```

MATLAB

```
iZ_50kHz
```

Description

Impedance measured at 50 kHz (result in ohms)

1.3.1.1.1.1.16 BodystatAPI.BSMMeasurement.iZ_5kHz Field**C++**

```
Integer iZ_5kHz;
```

C#

```
public Integer iZ_5kHz;
```

Visual Basic

```
Public iZ_5kHz As Integer
```

Java

```
public Integer iZ_5kHz;
```

MATLAB

```
iZ_5kHz
```

Description

Impedance measured at 5 kHz (result in ohms) (MDD/Quad only)

1.3.1.1.1.1.17 BodystatAPI.BSMMeasurement.llTestDate Field**C++**

```
Int64 llTestDate;
```

C#

```
public Int64 llTestDate;
```

Visual Basic

```
Public llTestDate As Int64
```

Java

```
public Int64 llTestDate;
```

MATLAB

```
llTestDate
```

Description

Date & time the test measurement was performed (Windows 64-bit time structure) (MDD/Quad only)

1.3.1.1.1.1.18 BodystatAPI.BSMMeasurement.pMultifreqData Field**C++**

```
ImpData pMultifreqData;
```

C#

```
public ImpData() pMultifreqData;
```

Visual Basic

```
Public pMultifreqData As ImpData()
```

Java

```
public ImpData() pMultifreqData;
```

MATLAB

```
pMultifreqData
```

Description

Pointer to ImpData (see page 114) array

1.3.1.1.1.1.19 BodystatAPI.BSMMeasurement.ulDeviceSerialNumber Field**C++**

```
UInteger ulDeviceSerialNumber;
```

C#

```
public UInteger ulDeviceSerialNumber;
```

Visual Basic

```
Public ulDeviceSerialNumber As UInteger
```

Java

```
public UInteger ulDeviceSerialNumber;
```

MATLAB

```
ulDeviceSerialNumber
```

Description

Unique serial number of the device which performed the measurement

1.3.1.1.1.1.2 BSMeasurement Properties**1.3.1.1.1.1.2.1 BodystatAPI.BSMMeasurement.tTestDate Property****C++**

```
__property DateTime tTestDate;
```

C#

```
public DateTime tTestDate;
```

Visual Basic

```
Public Property tTestDate() As DateTime
```

Java

```
tTestDate
```

MATLAB

```
tTestDate
```

Description

Convert underlying 64-bit time structure to/from Visual Basic Date/Time

1.3.1.1.1.2 BodystatAPI.BodystatAPI.BSNormals Structure

BSNormals: Bodystat Normals structure. Holds the normal values or normal ranges for any given test measurement.

C++

```
[StructLayout(LayoutKind.Sequential)]
struct BSNormals {
    Integer iFatPerc_L;
    Integer iFatPerc_H;
    Integer iFatKg_L;
    Integer iFatKg_H;
    Integer iLeanPerc_L;
    Integer iLeanPerc_H;
    Integer iLeanKg_L;
    Integer iLeanKg_H;
    Integer iTotalWeight_L;
    Integer iTotalWeight_H;
    Integer iTotalWeightMethod;
    Integer iTBW_Perc_L;
    Integer iTBW_Perc_H;
    Integer iTBW_L;
    Integer iTBW_H;
    Integer iECW_Perc_Norm;
    Integer iICW_Perc_Norm;
    Integer iBFMI_L;
    Integer iBFMI_H;
    Integer iFFMI_L;
    Integer iFFMI_H;
    Single fNutrition_Norm;
    Single fIllness_L;
    Single fIllness_H;
    Integer iBMI_L;
    Integer iBMI_H;
```



```

    Single fWaistHip_Norm;
    Single fWellness_L;
    Single fWellness_H;
};

```

C#

```

[StructLayout(LayoutKind.Sequential)]
public struct BSNormals {
    public Integer iFatPerc_L;
    public Integer iFatPerc_H;
    public Integer iFatKg_L;
    public Integer iFatKg_H;
    public Integer iLeanPerc_L;
    public Integer iLeanPerc_H;
    public Integer iLeanKg_L;
    public Integer iLeanKg_H;
    public Integer iTotalWeight_L;
    public Integer iTotalWeight_H;
    public Integer iTotalWeightMethod;
    public Integer iTBW_Perc_L;
    public Integer iTBW_Perc_H;
    public Integer iTBW_L;
    public Integer iTBW_H;
    public Integer iECW_Perc_Norm;
    public Integer iICW_Perc_Norm;
    public Integer iBFMI_L;
    public Integer iBFMI_H;
    public Integer iFFMI_L;
    public Integer iFFMI_H;
    public Single fNutrition_Norm;
    public Single fIllness_L;
    public Single fIllness_H;
    public Integer iBMI_L;
    public Integer iBMI_H;
    public Single fWaistHip_Norm;
    public Single fWellness_L;
    public Single fWellness_H;
}

```

Visual Basic

```

<StructLayout(LayoutKind.Sequential)>
Public Structure BSNormals
    Public iFatPerc_L As Integer
    Public iFatPerc_H As Integer
    Public iFatKg_L As Integer
    Public iFatKg_H As Integer
    Public iLeanPerc_L As Integer
    Public iLeanPerc_H As Integer
    Public iLeanKg_L As Integer
    Public iLeanKg_H As Integer
    Public iTotalWeight_L As Integer
    Public iTotalWeight_H As Integer
    Public iTotalWeightMethod As Integer
    Public iTBW_Perc_L As Integer
    Public iTBW_Perc_H As Integer
    Public iTBW_L As Integer
    Public iTBW_H As Integer
    Public iECW_Perc_Norm As Integer
    Public iICW_Perc_Norm As Integer
    Public iBFMI_L As Integer
    Public iBFMI_H As Integer
    Public iFFMI_L As Integer
    Public iFFMI_H As Integer
    Public fNutrition_Norm As Single
    Public fIllness_L As Single
    Public fIllness_H As Single
    Public iBMI_L As Integer
    Public iBMI_H As Integer
    Public fWaistHip_Norm As Single
    Public fWellness_L As Single
    Public fWellness_H As Single

```

End Structure

Java

```
[StructLayout(LayoutKind.Sequential)]
public class BSNormals;
```

MATLAB

BSNormals

File

BodystatSDK.vb

Members

Members	Description
public Integer iFatPerc_L;	Normal range for fat percentage - (low end)
public Integer iFatPerc_H;	Normal range for fat percentage - (high end)
public Integer iFatKg_L;	Normal range for fat (in kg) - (low end)
public Integer iFatKg_H;	Normal range for fat (in kg) - (high end)
public Integer iLeanPerc_L;	Normal range for lean percentage - (low end)
public Integer iLeanPerc_H;	Normal range for lean percentage - (high end)
public Integer iLeanKg_L;	Normal range for lean (in kg) - (low end)
public Integer iLeanKg_H;	Normal range for lean (in kg) - (high end)
public Integer iTotalWeight_L;	Normal range for Total Weight (in kg) - (low end)
public Integer iTotalWeight_H;	Normal range for Total Weight (in kg) - (high end)
public Integer iTotalWeightMethod;	Method used when calculating Total Weight normal range (0 = Composition Method, 1 = BMI Method)
public Integer iTBW_Perc_L;	Normal range for Total Body Water percentage - (low end)
public Integer iTBW_Perc_H;	Normal range for Total Body Water percentage - (high end)
public Integer iTBW_L;	Normal range for Total Body Water (in litres) - (low end)
public Integer iTBW_H;	Normal range for Total Body Water (in litres) - (high end)
public Integer iECW_Perc_Norm;	Normal value for ECW Percentage (MDD/Quad only)
public Integer iICW_Perc_Norm;	Normal value for ICW Percentage (MDD/Quad only)
public Integer iBFMI_L;	Normal range for Body Fat Mass Index (MDD/Quad only) - (low end)
public Integer iBFMI_H;	Normal range for Body Fat Mass Index (MDD/Quad only) - (high end)
public Integer iFFMI_L;	Normal range for Fat Free Mass Index (MDD/Quad only) - (low end)
public Integer iFFMI_H;	Normal range for Fat Free Mass Index (MDD/Quad only) - (high end)
public Single fNutrition_Norm;	Normal value for Nutrition Index (Quad only)
public Single fIllness_L;	Normal range for Prediction Marker (TM) - formerly known as Illness Marker (TM) - (Quad only) - (low end)
public Single fIllness_H;	Normal range for Prediction Marker (TM) - formerly known as Illness Marker (TM) - (Quad only) - (high end)
public Integer iBMI_L;	Normal range for Body Mass Index - (low end)
public Integer iBMI_H;	Normal range for Body Mass Index - (high end)
public Single fWaistHip_Norm;	Normal value for Waist/Hip ratio
public Single fWellness_L;	Normal range for Wellness Marker (TM) (MDD only) - (low end)
public Single fWellness_H;	Normal range for Wellness Marker (TM) (MDD only) - (high end)

1.3.1.1.1.3 BodystatAPI.BSRawData Structure

BSRawData: Bodystat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record (see page 98) starting number.

C++

```
[StructLayout(LayoutKind.Sequential)]
struct BSRawData {
    Integer iRecordArraySize;
    [System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
BSMeasurement record;
    Integer iTotalNumRecs;
    Integer ulFirstTestNum;
};
```

C#

```
[StructLayout(LayoutKind.Sequential)]
public struct BSRawData {
    public Integer iRecordArraySize;
    [System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
    public BSMeasurement() record;
    public Integer iTotalNumRecs;
    public Integer ulFirstTestNum;
}
```

Visual Basic

```
<StructLayout(LayoutKind.Sequential)>
Public Structure BSRawData
    Public iRecordArraySize As Integer
    <System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)>
    Public record As BSMeasurement()
    Public iTotalNumRecs As Integer
    Public ulFirstTestNum As Integer
End Structure
```

Java

```
[StructLayout(LayoutKind.Sequential)]
public class BSRawData;
```

MATLAB

BSRawData


File

BodystatSDK.vb





Remarks

Most models index 1000 records and store the last 100 at any given time.

Methods

	Name	Description
	New (see page 97)	Override the constructor to initialise required variables in the structure. You must ensure you allocate the structure in a manner which calls this function.

BSRawData Fields

	Name	Description
	iRecordArraySize (see page 97)	Maximum number of records this structure has reserved space for (currently fixed)
	iTotalNumRecs (see page 98)	Number of populated records in the array
	record (see page 98)	Array of test measurement data
	ulFirstTestNum (see page 99)	Test number of the first record (see page 98) stored in the array.

1.3.1.1.3.1 BodystatAPI.BSRawData.New Constructor

Override the constructor to initialise required variables in the structure. You must ensure you allocate the structure in a manner which calls this function.

C++

```
BSRawData (ByVal Integer dontcare);
```

C#

```
public New (ByVal Integer dontcare);
```

Visual Basic

```
Public Sub New (ByVal dontcare As Integer)
```

Java

```
public BSRawData (ByVal Integer dontcare);
```

MATLAB

```
New
```

Parameters

Parameters	Description
ByVal dontcare As Integer	Ignored: VB requires New() constructors to specify a non-optional parameter for structures. The parameter is not required by the wrapper and is ignored.

Body Source

```
Public Sub New (ByVal dontcare As Integer)
    Me.iRecordArraySize = BS_RAWDATA_ARRAYSIZE ' Value must be populated with size of
array (currently fixed)
    Me.ulFirstTestNum = 0
    Me.iTotalNumRecs = 0
End Sub
```

1.3.1.1.3.2 BSRawData Fields

1.3.1.1.3.2.1 BodystatAPI.BSRawData.iRecordArraySize Field

C++

```
Integer iRecordArraySize;
```

C#

```
public Integer iRecordArraySize;
```

Visual Basic

```
Public iRecordArraySize As Integer
```

Java

```
public Integer iRecordArraySize;
```

MATLAB

```
iRecordArraySize
```

Description

Maximum number of records this structure has reserved space for (currently fixed)

1.3.1.1.3.2.2 BodystatAPI.BSRawData.iTotalNumRecs Field**C++**

```
Integer iTotalNumRecs;
```

C#

```
public Integer iTotalNumRecs;
```

Visual Basic

```
Public iTotalNumRecs As Integer
```

Java

```
public Integer iTotalNumRecs;
```

MATLAB

```
iTotalNumRecs
```

Description

Number of populated records in the array

1.3.1.1.3.2.3 BodystatAPI.BSRawData.record Field**C++**

```
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
BSMeasurement record;
```

C#

```
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
public BSMeasurement() record;
```

Visual Basic

```
<System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)>
Public record As BSMeasurement()
```

Java

```
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=BS_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
public BSMeasurement() record;
```

MATLAB

```
record
```

Description

Array of test measurement data

1.3.1.1.1.3.2.4 BodystatAPI.BSRawData.ulFirstTestNum Field

C++

```
Integer ulFirstTestNum;
```

C#

```
public Integer ulFirstTestNum;
```

Visual Basic

```
Public ulFirstTestNum As Integer
```

Java

```
public Integer ulFirstTestNum;
```

MATLAB

```
ulFirstTestNum
```

Description

Test number of the first record (see page 98) stored in the array.

1.3.1.1.1.4 BodystatAPI.BodystatAPI.BSResults Structure

BSResults: Bodystat Results structure. Holds the calculated results for any given test measurement.

C++

```
[StructLayout(LayoutKind.Sequential)]
struct BSResults {
    Single fFatPerc;
    Single fFatKg;
    Single fLeanPerc;
    Single fLeanKg;
    Single fTotalWeight;
    Single fDryLW;
    Single fTBWPerc;
    Single fTBW;
    Single fECWPerc;
    Single fECW;
    Single fICWPerc;
    Single fICW;
    Single fBCM;
    Single fThirdSpace;
    Single fNutrition;
    Single fIllness;
    Single fBMR;
    Single fBMRkg;
    Single fEstAvg;
    Single fBMI;
    Single fBFMI;
    Single fFFMI;
    Single fWaistHip;
    Single fWellness;
    Single fECW_Legacy;
    Single fTBW_Legacy;
    Single fOHY;
    Single fSkMuscle;
    Single fCm;
    Single fRext;
    Single fRint;
    Single fFc;
    Single fAlpha;
};
```

C#

```
[StructLayout(LayoutKind.Sequential)]
public struct BSResults {
    public Single fFatPerc;
```

```

public Single fFatKg;
public Single fLeanPerc;
public Single fLeanKg;
public Single fTotalWeight;
public Single fDryLW;
public Single fTBWPerc;
public Single fTBW;
public Single fECWPerc;
public Single fECW;
public Single fICWPerc;
public Single fICW;
public Single fBCM;
public Single fThirdSpace;
public Single fNutrition;
public Single fIllness;
public Single fBMR;
public Single fBMRkg;
public Single fEstAvg;
public Single fBMI;
public Single fBFMI;
public Single fFFMI;
public Single fWaistHip;
public Single fWellness;
public Single fECW_Legacy;
public Single fTBW_Legacy;
public Single fOHY;
public Single fSkMuscle;
public Single fCm;
public Single fRext;
public Single fRint;
public Single fFc;
public Single fAlpha;
}

```

Visual Basic

```

<StructLayout(LayoutKind.Sequential)>
Public Structure BSResults
    Public fFatPerc As Single
    Public fFatKg As Single
    Public fLeanPerc As Single
    Public fLeanKg As Single
    Public fTotalWeight As Single
    Public fDryLW As Single
    Public fTBWPerc As Single
    Public fTBW As Single
    Public fECWPerc As Single
    Public fECW As Single
    Public fICWPerc As Single
    Public fICW As Single
    Public fBCM As Single
    Public fThirdSpace As Single
    Public fNutrition As Single
    Public fIllness As Single
    Public fBMR As Single
    Public fBMRkg As Single
    Public fEstAvg As Single
    Public fBMI As Single
    Public fBFMI As Single
    Public fFFMI As Single
    Public fWaistHip As Single
    Public fWellness As Single
    Public fECW_Legacy As Single
    Public fTBW_Legacy As Single
    Public fOHY As Single
    Public fSkMuscle As Single
    Public fCm As Single
    Public fRext As Single
    Public fRint As Single
    Public fFc As Single
    Public fAlpha As Single
End Structure

```

Java

```
[StructLayout(LayoutKind.Sequential)]
public class BSResults;
```

MATLAB

```
BSResults
```

File

```
BodystatSDK.vb
```

Members

Members	Description
public Single fFatPerc;	Fat percentage
public Single fFatKg;	Fat (in kg)
public Single fLeanPerc;	Lean percentage
public Single fLeanKg;	Lean (in kg)
public Single fTotalWeight;	Total Weight (in kg)
public Single fDryLW;	Dry Lean Weight (in kg)
public Single fTBWPerc;	Water percentage or Total Body Water percentage
public Single fTBW;	Water (in litres) or Total Body Water (in litres)
public Single fECWPerc;	Extra Cellular Water percentage (MDD/Quad only)
public Single fECW;	Extra Cellular Water (in litres) (MDD/Quad only)
public Single fICWPerc;	Intra Cellular Water percentage (MDD/Quad only)
public Single fICW;	Intra Cellular Water (in litres) (MDD/Quad only)
public Single fBCM;	Body cell mass (Quad only)
public Single fThirdSpace;	3rd space water (in litres) (Quad only)
public Single fNutrition;	Nutrition index (Quad only)
public Single fIllness;	Prediction Marker (TM) - formerly known as Illness Marker (TM) (Quad only)
public Single fBMR;	Basal Metabolic Rate (in kcal)
public Single fBMRkg;	Basal Metabolic Rate per kilogram in (kcal/kg)
public Single fEstAvg;	Estimated Average Requirement (in kcal)
public Single fBMI;	Body Mass Index
public Single fBFMI;	Body Fat Mass Index (MDD/Quad only)
public Single fFFMI;	Fat Free Mass Index (MDD/Quad only)
public Single fWaistHip;	Waist/Hip ratio
public Single fWellness;	Wellness Marker (TM) (MDD only)
public Single fECW_Legacy;	Legacy ECW calculation (QuanScan mode)
public Single fOHY;	Over hydration
public Single fSkMuscle;	Skeletal muscle mass
public Single fCm;	Cell membrane capacitance
public Single fRext;	R extracellular
public Single fRint;	R intracellular
public Single fFc;	Characteristic frequency
public Single fAlpha;	Alpha angle

1.3.1.1.1.5 BodystatAPI.ESMeasurement Structure

ESMeasurement: Measurement structure. Each test measurement recorded in the device is stored one of these structures. Holds the user inputted information about the subject, together with bio-impedance readings measured by the device (raw electrical readings only, not full results).

C++

```
[System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)]
struct ESMeasurement {
    BSDeviceModel iDeviceModel;
    UInteger ulDeviceSerialNumber;
    Int64 llTestDate;
    Integer iHeight;
    Integer iBCScore;
    Integer iZ_5kHz;
    Integer iZ_16kHz;
    Integer iZ_24kHz;
    Integer iZ_50kHz;
    Integer iZ_140kHz;
    Integer iZ_200kHz;
    Integer iZ_280kHz;
    Single fPA_50kHz;
};
```

C#

```
[System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)]
public struct ESMeasurement {
    public BSDeviceModel iDeviceModel;
    public UInteger ulDeviceSerialNumber;
    public Int64 llTestDate;
    public Integer iHeight;
    public Integer iBCScore;
    public Integer iZ_5kHz;
    public Integer iZ_16kHz;
    public Integer iZ_24kHz;
    public Integer iZ_50kHz;
    public Integer iZ_140kHz;
    public Integer iZ_200kHz;
    public Integer iZ_280kHz;
    public Single fPA_50kHz;
}
```

Visual Basic

```
<System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)>
Public Structure ESMeasurement
    Public iDeviceModel As BSDeviceModel
    Public ulDeviceSerialNumber As UInteger
    Public llTestDate As Int64
    Public iHeight As Integer
    Public iBCScore As Integer
    Public iZ_5kHz As Integer
    Public iZ_16kHz As Integer
    Public iZ_24kHz As Integer
    Public iZ_50kHz As Integer
    Public iZ_140kHz As Integer
    Public iZ_200kHz As Integer
    Public iZ_280kHz As Integer
    Public fPA_50kHz As Single
End Structure
```

Java

```
[System.Runtime.InteropServices.StructLayoutAttribute(System.Runtime.InteropServices.LayoutKind.Sequential)]
public class ESMeasurement;
```














MATLAB

```
ESMeasurement
```


File

```
BodystatSDK.vb
```

ESMeasurement Fields

	Name	Description
	fPA_50kHz (see page 103)	Phase angle measured at 50 kHz (result in degrees) (ES Pro only)
	iBCScore (see page 103)	Inputted BC score
	iDeviceModel (see page 104)	Device model which performed the measurement
	iHeight (see page 104)	Inputted height of the subject (centimeters)
	iZ_140kHz (see page 104)	Impedance measured at 140 kHz (result in ohms)
	iZ_16kHz (see page 105)	Impedance measured at 16 kHz (result in ohms)
	iZ_200kHz (see page 105)	Impedance measured at 200 kHz (result in ohms)
	iZ_24kHz (see page 105)	Impedance measured at 24 kHz (result in ohms)
	iZ_280kHz (see page 106)	Impedance measured at 280 kHz (result in ohms)
	iZ_50kHz (see page 106)	Impedance measured at 50 kHz (result in ohms)
	iZ_5kHz (see page 106)	Impedance measured at 5 kHz (result in ohms)
	lITestDate (see page 107)	Date & time the test measurement was performed (Windows 64-bit time structure) (MDD/Quad only)
	ulDeviceSerialNumber (see page 107)	Unique serial number of the device which performed the measurement

ESMeasurement Properties

	Name	Description
	tTestDate (see page 107)	Convert underlying 64-bit time structure to/from Visual Basic Date/Time

1.3.1.1.1.5.1 ESMeasurement Fields**1.3.1.1.1.5.1.1 BodystatAPI.ESMeasurement.fPA_50kHz Field****C++**

```
Single fPA_50kHz;
```

C#

```
public Single fPA_50kHz;
```

Visual Basic

```
Public fPA_50kHz As Single
```

Java

```
public Single fPA_50kHz;
```

MATLAB

```
fPA_50kHz
```

Description

Phase angle measured at 50 kHz (result in degrees) (ES Pro only)

1.3.1.1.1.5.1.2 BodystatAPI.ESMeasurement.iBCScore Field**C++**

```
Integer iBCScore;
```

C#

```
public Integer iBCScore;
```

Visual Basic

```
Public iBCScore As Integer
```

Java

```
public Integer iBCScore;
```

MATLAB

```
iBCScore
```

Description

Inputted BC score

1.3.1.1.1.5.1.3 BodystatAPI.ESMeasurement.iDeviceModel Field**C++**

```
BSDeviceModel iDeviceModel;
```

C#

```
public BSDeviceModel iDeviceModel;
```

Visual Basic

```
Public iDeviceModel As BSDeviceModel
```

Java

```
public BSDeviceModel iDeviceModel;
```

MATLAB

```
iDeviceModel
```

Description

Device model which performed the measurement

1.3.1.1.1.5.1.4 BodystatAPI.ESMeasurement.iHeight Field**C++**

```
Integer iHeight;
```

C#

```
public Integer iHeight;
```

Visual Basic

```
Public iHeight As Integer
```

Java

```
public Integer iHeight;
```

MATLAB

```
iHeight
```

Description

Inputted height of the subject (centimeters)

1.3.1.1.1.5.1.5 BodystatAPI.ESMeasurement.iZ_140kHz Field**C++**

```
Integer iZ_140kHz;
```

C#

```
public Integer iZ_140kHz;
```

Visual Basic

```
Public iZ_140kHz As Integer
```

Java

```
public Integer iZ_140kHz;
```

MATLAB

```
iZ_140kHz
```

Description

Impedance measured at 140 kHz (result in ohms)

1.3.1.1.5.1.6 BodystatAPI.ESMeasurement.iZ_16kHz Field**C++**

```
Integer iZ_16kHz;
```

C#

```
public Integer iZ_16kHz;
```

Visual Basic

```
Public iZ_16kHz As Integer
```

Java

```
public Integer iZ_16kHz;
```

MATLAB

```
iZ_16kHz
```

Description

Impedance measured at 16 kHz (result in ohms)

1.3.1.1.5.1.7 BodystatAPI.ESMeasurement.iZ_200kHz Field**C++**

```
Integer iZ_200kHz;
```

C#

```
public Integer iZ_200kHz;
```

Visual Basic

```
Public iZ_200kHz As Integer
```

Java

```
public Integer iZ_200kHz;
```

MATLAB

```
iZ_200kHz
```

Description

Impedance measured at 200 kHz (result in ohms)

1.3.1.1.5.1.8 BodystatAPI.ESMeasurement.iZ_24kHz Field**C++**

```
Integer iZ_24kHz;
```

C#

```
public Integer iZ_24kHz;
```

Visual Basic

```
Public iZ_24kHz As Integer
```

Java

```
public Integer iZ_24kHz;
```

MATLAB

```
iZ_24kHz
```

Description

Impedance measured at 24 kHz (result in ohms)

1.3.1.1.1.5.1.9 BodystatAPI.ESMeasurement.iZ_280kHz Field**C++**

```
Integer iZ_280kHz;
```

C#

```
public Integer iZ_280kHz;
```

Visual Basic

```
Public iZ_280kHz As Integer
```

Java

```
public Integer iZ_280kHz;
```

MATLAB

```
iZ_280kHz
```

Description

Impedance measured at 280 kHz (result in ohms)

1.3.1.1.1.5.1.10 BodystatAPI.ESMeasurement.iZ_50kHz Field**C++**

```
Integer iZ_50kHz;
```

C#

```
public Integer iZ_50kHz;
```

Visual Basic

```
Public iZ_50kHz As Integer
```

Java

```
public Integer iZ_50kHz;
```

MATLAB

```
iZ_50kHz
```

Description

Impedance measured at 50 kHz (result in ohms)

1.3.1.1.1.5.1.11 BodystatAPI.ESMeasurement.iZ_5kHz Field**C++**

```
Integer iZ_5kHz;
```

C#

```
public Integer iZ_5kHz;
```

Visual Basic

```
Public iZ_5kHz As Integer
```

Java

```
public Integer iZ_5kHz;
```

MATLAB

```
iZ_5kHz
```

Description

Impedance measured at 5 kHz (result in ohms)

1.3.1.1.1.5.1.12 BodystatAPI.ESMeasurement.llTestDate Field**C++**

```
Int64 llTestDate;
```

C#

```
public Int64 llTestDate;
```

Visual Basic

```
Public llTestDate As Int64
```

Java

```
public Int64 llTestDate;
```

MATLAB

```
llTestDate
```

Description

Date & time the test measurement was performed (Windows 64-bit time structure) (MDD/Quad only)

1.3.1.1.1.5.1.13 BodystatAPI.ESMeasurement.ulDeviceSerialNumber Field**C++**

```
UInteger ulDeviceSerialNumber;
```

C#

```
public UInteger ulDeviceSerialNumber;
```

Visual Basic

```
Public ulDeviceSerialNumber As UInteger
```

Java

```
public UInteger ulDeviceSerialNumber;
```

MATLAB

```
ulDeviceSerialNumber
```

Description

Unique serial number of the device which performed the measurement

1.3.1.1.1.5.2 ESMeasurement Properties**1.3.1.1.1.5.2.1 BodystatAPI.ESMeasurement.tTestDate Property****C++**

```
__property DateTime tTestDate;
```

C#

```
public DateTime tTestDate;
```

Visual Basic

```
Public Property tTestDate() As DateTime
```

Java

```
tTestDate
```

MATLAB

```
tTestDate
```

Description

Convert underlying 64-bit time structure to/from Visual Basic Date/Time

1.3.1.1.1.6 BodystatAPI.BodystatAPI.ESNormals Structure

ESNormals: Equistat Normals structure. Holds the normal values or normal ranges for any given test measurement.

C++

```
[StructLayout(LayoutKind.Sequential)]
struct ESNormals {
    Integer iFatKg_L;
    Integer iFatKg_H;
    Integer iFatKg_PBMPerC_L;
    Integer iFatKg_PBMPerC_H;
    Integer iTFV_L;
    Integer iTFV_H;
    Integer iTFVPerC_L;
    Integer iTFVPerC_H;
    Integer iECFV_L;
    Integer iECFV_H;
    Integer iECFV_PBM_L;
    Integer iECFV_PBM_H;
    Integer iECFV_TFVPerC_L;
    Integer iECFV_TFVPerC_H;
    Integer iICFV_L;
    Integer iICFV_H;
    Integer iICFV_PBM_L;
    Integer iICFV_PBM_H;
    Integer iICFV_TFVPerC_L;
    Integer iICFV_TFVPerC_H;
    Integer iPV_L;
    Integer iPV_H;
    Integer iPV_ECFVPerC_L;
    Integer iPV_ECFVPerC_H;
};
```

C#

```
[StructLayout(LayoutKind.Sequential)]
public struct ESNormals {
    public Integer iFatKg_L;
    public Integer iFatKg_H;
    public Integer iFatKg_PBMPerC_L;
    public Integer iFatKg_PBMPerC_H;
    public Integer iTFV_L;
    public Integer iTFV_H;
    public Integer iTFVPerC_L;
    public Integer iTFVPerC_H;
    public Integer iECFV_L;
    public Integer iECFV_H;
    public Integer iECFV_PBM_L;
    public Integer iECFV_PBM_H;
    public Integer iECFV_TFVPerC_L;
    public Integer iECFV_TFVPerC_H;
    public Integer iICFV_L;
    public Integer iICFV_H;
    public Integer iICFV_PBM_L;
    public Integer iICFV_PBM_H;
    public Integer iICFV_TFVPerC_L;
};
```

```

    public Integer iICFV_TFVPerc_H;
    public Integer iPV_L;
    public Integer iPV_H;
    public Integer iPV_ECFVPerc_L;
    public Integer iPV_ECFVPerc_H;
}

```

Visual Basic

```

<StructLayout(LayoutKind.Sequential)>
Public Structure ESNormals
    Public iFatKg_L As Integer
    Public iFatKg_H As Integer
    Public iFatKg_PBMPerc_L As Integer
    Public iFatKg_PBMPerc_H As Integer
    Public iTFV_L As Integer
    Public iTFV_H As Integer
    Public iTFVPerc_L As Integer
    Public iTFVPerc_H As Integer
    Public iECFV_L As Integer
    Public iECFV_H As Integer
    Public iECFV_PBM_L As Integer
    Public iECFV_PBM_H As Integer
    Public iECFV_TFVPerc_L As Integer
    Public iECFV_TFVPerc_H As Integer
    Public iICFV_L As Integer
    Public iICFV_H As Integer
    Public iICFV_PBM_L As Integer
    Public iICFV_PBM_H As Integer
    Public iICFV_TFVPerc_L As Integer
    Public iICFV_TFVPerc_H As Integer
    Public iPV_L As Integer
    Public iPV_H As Integer
    Public iPV_ECFVPerc_L As Integer
    Public iPV_ECFVPerc_H As Integer
End Structure

```

Java

```

[StructLayout(LayoutKind.Sequential)]
public class ESNormals;

```

MATLAB

```
ESNormals
```

File

```
BodystatSDK.vb
```

1.3.1.1.1.7 BodystatAPI.ESRawData Structure

ESRawData: Equistat Raw Data structure. Used when downloading test measurements from the device. Contains a fixed array of measurement records together with information about the number of records occupied within the array and test record (see page 112) starting number.

C++

```

[StructLayout(LayoutKind.Sequential)]
struct ESRawData {
    Integer iRecordArraySize;
    [System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
    SizeConst:=ES_RAWDATA_ARRAYSIZE,
    ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
    ESMeasurement record;
    Integer iTotalNumRecs;
    Integer ulFirstTestNum;
};

```

C#

```

[StructLayout(LayoutKind.Sequential)]
public struct ESRawData {

```



```
public Integer iRecordArraySize;
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.Unmanage
dType.ByValArray,
SizeConst:=ES_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
public ESMeasurement() record;
public Integer iTotalNumRecs;
public Integer ulFirstTestNum;
}
```

Visual Basic

```
<StructLayout(LayoutKind.Sequential)>
Public Structure ESRawData
    Public iRecordArraySize As Integer
    <System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.Unmanage
dType.ByValArray,
SizeConst:=ES_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)>
    Public record As ESMeasurement()
    Public iTotalNumRecs As Integer
    Public ulFirstTestNum As Integer
End Structure
```

Java

```
[StructLayout(LayoutKind.Sequential)]
public class ESRawData;
```

MATLAB

```
ESRawData
```


File

```
BodystatSDK.vb
```





Remarks

Most models index 1000 records and store the last 100 at any given time.

Methods

	Name	Description
	New (see page 110)	Override the constructor to initialise required variables in the structure. You must ensure you allocate the structure in a manner which calls this function.

ESRawData Fields

	Name	Description
	iRecordArraySize (see page 111)	Maximum number of records this structure has reserved space for (currently fixed)
	iTotalNumRecs (see page 111)	Number of populated records in the array
	record (see page 112)	Array of test measurement data
	ulFirstTestNum (see page 112)	Test number of the first record (see page 112) stored in the array.

1.3.1.1.1.7.1 BodystatAPI.ESRawData.New Constructor

Override the constructor to initialise required variables in the structure. You must ensure you allocate the structure in a manner which calls this function.

C++

```
ESRawData(ByVal Integer dontcare);
```

C#

```
public New(ByVal Integer dontcare);
```

Visual Basic

```
Public Sub New(ByVal dontcare As Integer)
```

Java

```
public ESRawData(ByVal Integer dontcare);
```

MATLAB

```
New
```

Parameters

Parameters	Description
ByVal dontcare As Integer	Ignored: VB requires New() constructors to specify a non-optional parameter for structures. The parameter is not required by the wrapper and is ignored.

Body Source

```
Public Sub New(ByVal dontcare As Integer)
    Me.iRecordArraySize = ES_RAWDATA_ARRAYSIZE ' Value must be populated with size of
array (currently fixed)
End Sub
```

1.3.1.1.1.7.2 ESRawData Fields

1.3.1.1.1.7.2.1 BodystatAPI.ESRawData.iRecordArraySize Field

C++

```
Integer iRecordArraySize;
```

C#

```
public Integer iRecordArraySize;
```

Visual Basic

```
Public iRecordArraySize As Integer
```

Java

```
public Integer iRecordArraySize;
```

MATLAB

```
iRecordArraySize
```

Description

Maximum number of records this structure has reserved space for (currently fixed)

1.3.1.1.1.7.2.2 BodystatAPI.ESRawData.iTotalNumRecs Field

C++

```
Integer iTotalNumRecs;
```

C#

```
public Integer iTotalNumRecs;
```

Visual Basic

```
Public iTotalNumRecs As Integer
```

Java

```
public Integer iTotalNumRecs;
```

MATLAB

```
iTotalNumRecs
```

Description

Number of populated records in the array

1.3.1.1.1.7.2.3 BodystatAPI.ESRawData.record Field

C++

```
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=ES_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
ESMeasurement record;
```

C#

```
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=ES_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
public ESMeasurement() record;
```

Visual Basic

```
<System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=ES_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)>
Public record As ESMeasurement()
```

Java

```
[System.Runtime.InteropServices.MarshalAsAttribute(System.Runtime.InteropServices.UnmanagedType.ByValArray,
SizeConst:=ES_RAWDATA_ARRAYSIZE,
ArraySubType:=System.Runtime.InteropServices.UnmanagedType.Struct)]
public ESMeasurement() record;
```

MATLAB

```
record
```

Description

Array of test measurement data

1.3.1.1.1.7.2.4 BodystatAPI.ESRawData.ulFirstTestNum Field

C++

```
Integer ulFirstTestNum;
```

C#

```
public Integer ulFirstTestNum;
```

Visual Basic

```
Public ulFirstTestNum As Integer
```

Java

```
public Integer ulFirstTestNum;
```

MATLAB

```
ulFirstTestNum
```

Description

Test number of the first record (see page 112) stored in the array.

1.3.1.1.1.8 BodystatAPI.BodystatAPI.ESResults Structure

ESResults: Equistst Results structure. Holds the calculated results for any given test measurement.

C++

```
[StructLayout(LayoutKind.Sequential)]
struct ESResults {
```

```

Single fFatKg;
Single fPBM;
Single fIBM;
Single fFatKg_PBMperc;
Single fTFV;
Single fTFVperc;
Single fECFV;
Single fECFV_PBM;
Single fECFV_TFVperc;
Single fICFV;
Single fICFV_PBM;
Single fICFV_TFVperc;
Single fPV;
Single fPV_ECFVperc;
Single fDehydration;
};

```

C#

```

[StructLayout(LayoutKind.Sequential)]
public struct ESResults {
    public Single fFatKg;
    public Single fPBM;
    public Single fIBM;
    public Single fFatKg_PBMperc;
    public Single fTFV;
    public Single fTFVperc;
    public Single fECFV;
    public Single fECFV_PBM;
    public Single fECFV_TFVperc;
    public Single fICFV;
    public Single fICFV_PBM;
    public Single fICFV_TFVperc;
    public Single fPV;
    public Single fPV_ECFVperc;
    public Single fDehydration;
}

```

Visual Basic

```

<StructLayout(LayoutKind.Sequential)>
Public Structure ESResults
    Public fFatKg As Single
    Public fPBM As Single
    Public fIBM As Single
    Public fFatKg_PBMperc As Single
    Public fTFV As Single
    Public fTFVperc As Single
    Public fECFV As Single
    Public fECFV_PBM As Single
    Public fECFV_TFVperc As Single
    Public fICFV As Single
    Public fICFV_PBM As Single
    Public fICFV_TFVperc As Single
    Public fPV As Single
    Public fPV_ECFVperc As Single
    Public fDehydration As Single
End Structure

```

Java

```

[StructLayout(LayoutKind.Sequential)]
public class ESResults;

```

MATLAB

```

ESResults

```

File

```

BodystatSDK.vb

```

1.3.1.1.1.9 BodystatAPI.BodystatAPI.ImpData Structure

Represents one frequency measurement

C++

```
[StructLayout(LayoutKind.Sequential)]
struct ImpData {
    Single fFrequency;
    Single fImpedance;
    Single fPhaseAngle;
};
```

C#

```
[StructLayout(LayoutKind.Sequential)]
public struct ImpData {
    public Single fFrequency;
    public Single fImpedance;
    public Single fPhaseAngle;
}
```

Visual Basic

```
<StructLayout(LayoutKind.Sequential)>
Public Structure ImpData
    Public fFrequency As Single
    Public fImpedance As Single
    Public fPhaseAngle As Single
End Structure
```

Java

```
[StructLayout(LayoutKind.Sequential)]
public class ImpData;
```

MATLAB

```
ImpData
```

File

```
BodystatSDK.vb
```

Members

Members	Description
public Single fFrequency;	Frequency in kHz
public Single fImpedance;	Absolute value of impedance measured at given frequency (result in ohms)
public Single fPhaseAngle;	Phase angle (result in degrees)

1.3.1.1.2 BodystatAPI Fields

1.3.1.1.2.1 BodystatAPI.BS_RAWDATA_ARRAYSIZE Field

C++

```
Integer BS_RAWDATA_ARRAYSIZE = 100;
```

C#

```
public const Integer BS_RAWDATA_ARRAYSIZE = 100;
```

Visual Basic

```
Public Const BS_RAWDATA_ARRAYSIZE As Integer = 100
```

Java

```
public Integer BS_RAWDATA_ARRAYSIZE = 100;
```

MATLAB

BS_RAWDATA_ARRAYSIZE

Description

Number of measurement records BSRawData (see page 96) structure should reserved space for (currently fixed)

1.3.1.1.2.2 BodystatAPI.ES_RAWDATA_ARRAYSIZE Field**C++**

Integer **ES_RAWDATA_ARRAYSIZE** = 100;

C#

public const Integer **ES_RAWDATA_ARRAYSIZE** = 100;

Visual Basic

Public Const **ES_RAWDATA_ARRAYSIZE** **As Integer** = 100

Java

public Integer **ES_RAWDATA_ARRAYSIZE** = 100;

MATLAB

ES_RAWDATA_ARRAYSIZE

Description

Number of measurement records ESRawData (see page 109) structure should reserved space for (currently fixed)

1.3.1.1.3 BodystatAPI Methods**1.3.1.1.3.1 BodystatAPI.BSAuthenticateBTDevice Method**

Authenticate (pair) a Bodystat device found during a previous call to BSSearchBTDevices (see page 136). Completely silent implementation.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
Integer BSAuthenticateBTDevice(ByVal IntPtr hWndParent, Optional ByVal UInt16 iTimeout = 7
_);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSAuthenticateBTDevice(ByVal IntPtr hWndParent, Optional ByVal UInt16
iTimeout);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSAuthenticateBTDevice(ByVal hWndParent As IntPtr, Optional ByVal
iTimeout As UInt16 = 7 _) As Integer
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSAuthenticateBTDevice(ByVal IntPtr hWndParent, Optional ByVal UInt16
iTimeout);
```

MATLAB

BSAuthenticateBTDevice

Parameters

Parameters	Description
ByVal hWndParent As IntPtr	Handle of the parent window.
Optional ByVal iTimeout As UInt16 = 7 _	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or times out.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSAuthenticateBTDevice( _
ByVal hWndParent As IntPtr _
, Optional ByVal iTimeout As UInt16 = 7 _
) As Integer
End Function
```

1.3.1.1.3.2 BodystatAPI.BSAutoSetupBT Method

Automatic Bluetooth setup. Automatically searches available Bluetooth devices, exchanges passkeys with Bodystat device and enables serial port service required for communication with the device. If multiple Bodystat devices are detected the user is prompted to chose the desired device (model and BT BDA address shown).

This function is typically called as a one time setup for the device (e.g. after first installation, or from the hardware configuration page of your application's settings, or when the user obtains a replacement device (upgrades/replacement/etc).

Typically, once setup you need only store the device com port returned. You can then subsequently communicate with the device knowing only the com port. It is not necessary to setup the device in this way before each use.

New non-paired units will only be found if switched on, in range and ready for communication (on main startup screen). It is advisable to instruct the user to ensure this is the case before calling.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
Integer BSAutoSetupBT(ByVal StringBuilder lpComPort, ByVal Integer iComPortBufferSize,
ByVal Integer bReportErrors, ByVal IntPtr hWndParent);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSAutoSetupBT(ByVal StringBuilder lpComPort, ByVal Integer
iComPortBufferSize, ByVal Integer bReportErrors, ByVal IntPtr hWndParent);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSAutoSetupBT(ByVal lpComPort As StringBuilder, ByVal
iComPortBufferSize As Integer, ByVal bReportErrors As Integer, ByVal hWndParent As IntPtr)
As Integer
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSAutoSetupBT(ByVal StringBuilder lpComPort, ByVal Integer
iComPortBufferSize, ByVal Integer bReportErrors, ByVal IntPtr hWndParent);
```

MATLAB

```
BSAutoSetupBT
```

Parameters

Parameters	Description
ByVal lpComPort As StringBuilder	[out] String buffer to receive the device port the device is assigned to.
ByVal bReportErrors As Integer	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (still not completely silent).
ByVal hWndParent As IntPtr	Handle of the parent window.
iSize	The size of the buffer.

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or timesout.

Errors during the process can be optionally shown/hidden. However, a completely silent approach is not possible using this function (instead use BSSearchBTDevices (see page 136) and BSGetBTBodystatDevice (see page 120) or iterate radio handles manually. Then pass desired handle to BSAuthenticateBTDevice (see page 115) for pairing).

Finally be aware this function is only available for supported locales. Otherwise English messages may be displayed to the user. Adopt silent approach discussed previously if you are targeting a non-supported locale.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSAutoSetupBT( _
ByVal lpComPort As StringBuilder _
, ByVal iComPortBufferSize As Integer _
, ByVal bReportErrors As Integer _
, ByVal hWndParent As IntPtr _
) As Integer
End Function
```

1.3.1.1.3.3 BodystatAPI.BSCalculateNormals Method

Calculate normal values and normal ranges for a given subject and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from BSCalculateResults (see page 118)).

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSCalculateNormals(ByRef BSMeasurement pM, ByRef BSResults pR, ByRef BSNormals pN);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSCalculateNormals(ByRef BSMeasurement pM, ByRef BSResults pR, ByRef
```



```
BSNormals pN);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSCalculateNormals(ByRef pM As BSMeasurement, ByRef pR As BSResults,
ByRef pN As BSNormals) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSCalculateNormals(ByRef BSMeasurement pM, ByRef BSResults pR, ByRef
BSNormals pN);
```

MATLAB

```
BSCalculateNormals
```

Parameters

Parameters	Description
ByRef pM As BSMeasurement	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
ByRef pR As BSResults	The results calculated for this measurement previously.
ByRef pN As BSNormals	[out] Receives the normal values and normal ranges for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Normals vary by subject attributes, measured bio-impedance values and calculated results.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSCalculateNormals( _
ByRef pM As BSMeasurement _
, ByRef pR As BSResults _
, ByRef pN As BSNormals _
) As BSError
End Function
```

1.3.1.1.3.4 BodystatAPI.BSCalculateResults Method

Calculate complete result set for a given bio-impedance measurement recorded by the Bodystat device.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSCalculateResults(ByRef BSMeasurement pM, ByRef BSResults pR, ByVal BSAnalysisMode
iAnalysisMode);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSCalculateResults(ByRef BSMeasurement pM, ByRef BSResults pR, ByVal
BSAnalysisMode iAnalysisMode);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSCalculateResults(ByRef pM As BSMeasurement, ByRef pR As BSResults,
ByVal iAnalysisMode As BSAnalysisMode) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSCalculateResults(ByRef BSMeasurement pM, ByRef BSResults pR, ByVal
BSAnalysisMode iAnalysisMode);
```

MATLAB

```
BSCalculateResults
```

Parameters

Parameters	Description
ByRef pM As BSMeasurement	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
ByRef pR As BSResults	[out] Receives the complete calculated results for this measurement.
ByVal iAnalysisMode As BSAnalysisMode	The desired analysis mode (see BSAnalysisMode (see page 145) for more information).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Results vary by subject attributes and measured bio-impedance values.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSCalculateResults( _
ByRef pM As BSMeasurement _
, ByRef pR As BSResults _
, ByVal iAnalysisMode As BSAnalysisMode _
) As BSError
End Function
```

1.3.1.1.3.5 BodystatAPI.BSCloseComport Method

Close com port.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSCloseComport(_);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSCloseComport(_);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Sub BSCloseComport( As _)
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSCloseComport(_ );
```

MATLAB

```
BSCloseComport
```

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Sub BSCloseComPort( _
)
End Sub
```

1.3.1.1.3.6 BodystatAPI.BSConnect Method

Connect to the Bodystat device and ensure the device is responding.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSConnect(_);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSConnect(_);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSConnect( As _ ) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSConnect(_ );
```

MATLAB

```
BSConnect
```

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have opened the com port by calling BSOpenComPort previously.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSConnect( _
) As BSError
End Function
```

1.3.1.1.3.7 BodystatAPI.BSGetBTBodystatDevice Method

Used to get information about any Bodystat device found after a previous call to BSSearchBTDevices (see page 136).

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
Integer BSGetBTBodystatDevice(ByVal StringBuilder lpDeviceName, ByVal Integer
iDeviceNameBufferSize, ByVal StringBuilder lpBDA, ByVal Integer iBDABufferSize, ByVal
StringBuilder lpComPort, ByVal Integer iComPortBufferSize, Optional ByVal UInt16 iTimeout =
7 _);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSGetBTBodystatDevice(ByVal StringBuilder lpDeviceName, ByVal Integer
iDeviceNameBufferSize, ByVal StringBuilder lpBDA, ByVal Integer iBDABufferSize, ByVal
StringBuilder lpComPort, ByVal Integer iComPortBufferSize, Optional ByVal UInt16 iTimeout);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSGetBTBodystatDevice(ByVal lpDeviceName As StringBuilder, ByVal
iDeviceNameBufferSize As Integer, ByVal lpBDA As StringBuilder, ByVal iBDABufferSize As
Integer, ByVal lpComPort As StringBuilder, ByVal iComPortBufferSize As Integer, Optional
ByVal iTimeout As UInt16 = 7 _) As Integer
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSGetBTBodystatDevice(ByVal StringBuilder lpDeviceName, ByVal Integer
iDeviceNameBufferSize, ByVal StringBuilder lpBDA, ByVal Integer iBDABufferSize, ByVal
StringBuilder lpComPort, ByVal Integer iComPortBufferSize, Optional ByVal UInt16 iTimeout);
```

MATLAB

```
BSGetBTBodystatDevice
```

Parameters

Parameters	Description
ByVal lpDeviceName As StringBuilder	[out] String buffer to receive the product name of the Bodystat device.
ByVal iDeviceNameBufferSize As Integer	Size of the device name buffer.
ByVal lpBDA As StringBuilder	[out] String buffer to receive the Bluetooth hexadecimal address (known as the BDA - unique to each BT module in a device, like a MAC address).
ByVal iBDABufferSize As Integer	Size of the BDA buffer.
ByVal lpComPort As StringBuilder	[out] String buffer to receive device/port the unit is assigned to.
Optional ByVal iTimeout As UInt16 = 7 _	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).
iComBufferSize	Size of the com port buffer.

Returns

TRUE if it succeeds, FALSE if it fails.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSGetBTBodystatDevice( _
ByVal lpDeviceName As StringBuilder _
, ByVal iDeviceNameBufferSize As Integer _
, ByVal lpBDA As StringBuilder _
, ByVal iBDABufferSize As Integer _
, ByVal lpComPort As StringBuilder _
, ByVal iComPortBufferSize As Integer _
, Optional ByVal iTimeout As UInt16 = 7 _
) As Integer
End Function
```

1.3.1.1.3.8 BodystatAPI.BSGetBTStackInfo Method

Retrieve information about the Bluetooth stack. Useful for diagnostics.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
```

```
Integer BSGetBTStackInfo(ByVal StringBuilder lpBuffer, ByVal Integer iBufferSize);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSGetBTStackInfo(ByVal StringBuilder lpBuffer, ByVal Integer
iBufferSize);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSGetBTStackInfo(ByVal lpBuffer As StringBuilder, ByVal iBufferSize
As Integer) As Integer
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSGetBTStackInfo(ByVal StringBuilder lpBuffer, ByVal Integer
iBufferSize);
```

MATLAB

```
BSGetBTStackInfo
```

Parameters

Parameters	Description
ByVal lpBuffer As StringBuilder	[out] String buffer to receive the information.
ByVal iBufferSize As Integer	Size of the buffer.

Returns

true if it succeeds, false if it fails.

Remarks

Note: Information can only be provided for supported BT radios.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSGetBTStackInfo( _
ByVal lpBuffer As StringBuilder _
, ByVal iBufferSize As Integer _
) As Integer
End Function
```

1.3.1.1.3.9 BodystatAPI.BSGetDeviceFamily Method

Determine the device family of a specific model.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSDeviceFamily BSGetDeviceFamily(ByVal BSDeviceModel iModel);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSDeviceFamily BSGetDeviceFamily(ByVal BSDeviceModel iModel);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSGetDeviceFamily(ByVal iModel As BSDeviceModel) As BSDeviceFamily
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
```

```
public Shared BSDeviceFamily BSGetDeviceFamily(ByVal BSDeviceModel iModel);
```

MATLAB

```
BSGetDeviceFamily
```

Parameters

Parameters	Description
ByVal iModel As BSDeviceModel	The desired model of interest.

Returns

The device family the model belongs to.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSGetDeviceFamily( _
ByVal iModel As BSDeviceModel _
) As BSDeviceFamily
End Function
```

1.3.1.1.3.10 BodystatAPI.BSGetDeviceModelName Method

Determine the full product name of a specific model.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSGetDeviceModelName(ByVal BSDeviceModel iModel, ByVal StringBuilder lpBuffer,
ByVal Integer iBufferSize);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSGetDeviceModelName(ByVal BSDeviceModel iModel, ByVal StringBuilder
lpBuffer, ByVal Integer iBufferSize);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSGetDeviceModelName(ByVal iModel As BSDeviceModel, ByVal lpBuffer
As StringBuilder, ByVal iBufferSize As Integer) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSGetDeviceModelName(ByVal BSDeviceModel iModel, ByVal StringBuilder
lpBuffer, ByVal Integer iBufferSize);
```

MATLAB

```
BSGetDeviceModelName
```

Parameters

Parameters	Description
ByVal iModel As BSDeviceModel	The desired model of interest.
ByVal lpBuffer As StringBuilder	[out] String buffer to receive the name.
ByVal iBufferSize As Integer	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSGetDeviceModelName( _
```

```

ByVal iModel As BSDeviceModel _
, ByVal lpBuffer As StringBuilder _
, ByVal iBufferSize As Integer _
) As BSError
End Function

```

1.3.1.1.3.11 BodystatAPI.BSGetDeviceModelNameShort Method

Determine the short product name of a specific model.

C++

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSGetDeviceModelNameShort(ByVal BSDeviceModel iModel, ByVal StringBuilder lpBuffer,
ByVal Integer iBufferSize);

```

C#

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSGetDeviceModelNameShort(ByVal BSDeviceModel iModel, ByVal
StringBuilder lpBuffer, ByVal Integer iBufferSize);

```

Visual Basic

```

<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSGetDeviceModelNameShort(ByVal iModel As BSDeviceModel, ByVal
lpBuffer As StringBuilder, ByVal iBufferSize As Integer) As BSError

```

Java

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSGetDeviceModelNameShort(ByVal BSDeviceModel iModel, ByVal
StringBuilder lpBuffer, ByVal Integer iBufferSize);

```

MATLAB

```
BSGetDeviceModelNameShort
```

Parameters

Parameters	Description
ByVal iModel As BSDeviceModel	The desired model of interest.
ByVal lpBuffer As StringBuilder	[out] String buffer to receive the name.
ByVal iBufferSize As Integer	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Body Source

```

<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSGetDeviceModelNameShort( _
ByVal iModel As BSDeviceModel _
, ByVal lpBuffer As StringBuilder _
, ByVal iBufferSize As Integer _
) As BSError
End Function

```

1.3.1.1.3.12 BodystatAPI.BSGetSdkLibraryVersion Method

Determine the Bodystat SDK library (DLL) version.

C++

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSGetSdkLibraryVersion(ByRef Integer pMajorSdkVer, ByRef Integer pMinorSdkVer);

```

```
C#
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSGetSdkLibraryVersion(ByRef Integer pMajorSdkVer, ByRef Integer
pMinorSdkVer);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSGetSdkLibraryVersion(ByRef pMajorSdkVer As Integer, ByRef
pMinorSdkVer As Integer) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSGetSdkLibraryVersion(ByRef Integer pMajorSdkVer, ByRef Integer
pMinorSdkVer);
```

MATLAB

```
BSGetSdkLibraryVersion
```

Parameters

Parameters	Description
ByRef pMajorSdkVer As Integer	[out] Receives the major version number of the Bodystat SDK DLL.
ByRef pMinorSdkVer As Integer	[out] Receives the minor version number of the Bodystat SDK DLL.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Caller should check the version number of the DLL if utilising functionality noted only as being available in later revisions of the SDK. Or, alternatively, ensure the SDK DLL is updated appropriately during installation of your application.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSGetSdkLibraryVersion( _
ByRef pMajorSdkVer As Integer _
, ByRef pMinorSdkVer As Integer _
) As BSError
End Function
```

1.3.1.1.3.13 BodystatAPI.BSIsBTAvailable Method

Check if a supported BT radio is available on this computer. Checks for presence of USB Bluetooth dongles and inbuilt Bluetooth modules in the PC. No communication with Bodystat devices takes place during this query.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
Integer BSIsBTAvailable(_);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSIsBTAvailable(_);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSIsBTAvailable( As _) As Integer
```


Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSIsBTAvailable(_ );
```

MATLAB

```
BSIsBTAvailable
```

Returns

Returns true if supported radio is available and enabled.

Remarks

Note: Only supported BT radios can be queried (presently those compatible with the Microsoft Bluetooth APIs and some Widdcomm models).

Also note, Bluetooth must be enabled on the PC by the user to permit detection. Some PCs, laptops particularly, have physical switches to enable bluetooth/wireless communication.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSIsBTAvailable( _
) As Integer
End Function
```

1.3.1.1.3.14 BodystatAPI.BSOpenComport Method

Open com port to the device ready for communication.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSOpenComport(ByVal StringBuilder lpComPort, ByVal Integer iComPortBufferSize);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSOpenComport(ByVal StringBuilder lpComPort, ByVal Integer
iComPortBufferSize);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSOpenComport(ByVal lpComPort As StringBuilder, ByVal
iComPortBufferSize As Integer) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSOpenComport(ByVal StringBuilder lpComPort, ByVal Integer
iComPortBufferSize);
```

MATLAB

```
BSOpenComport
```

Parameters

Parameters	Description
ByVal lpComPort As StringBuilder	String buffer containing the device/port to be opened.
iSize	The size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSOpenComport( _
ByVal lpComPort As StringBuilder _
, ByVal iComPortBufferSize As Integer _
) As BSError
End Function
```

1.3.1.1.3.15 BSReadCalibrationTime Method**1.3.1.1.3.15.1 BodystatAPI.BSReadCalibrationTime Method (DateTime)**

Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.

C++

```
BSError BSReadCalibrationTime(ByRef DateTime dtCalibrationTime);
```

C#

```
public Shared BSError BSReadCalibrationTime(ByRef DateTime dtCalibrationTime);
```

Visual Basic

```
Public Shared Function BSReadCalibrationTime(ByRef dtCalibrationTime As DateTime) As BSError
```

Java

```
public Shared BSError BSReadCalibrationTime(ByRef DateTime dtCalibrationTime);
```

MATLAB

```
BSReadCalibrationTime
```

Parameters

Parameters	Description
ByRef dtCalibrationTime As DateTime	[out] Receives the date the device was last calibrated.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
Public Shared Function BSReadCalibrationTime( _
ByRef dtCalibrationTime As DateTime _
) As BSError

Dim llCalibrationDateTime As Int64
Dim bsErr As BSError = BSReadCalibrationTime(llCalibrationDateTime)

dtCalibrationTime = "#1/1/1970#"
dtCalibrationTime = dtCalibrationTime.AddSeconds(llCalibrationDateTime)

Return bsErr
End Function
```

1.3.1.1.3.15.2 BodystatAPI.BSReadCalibrationTime Method (Int64)

Read the calibration date from the connected device (when the unit was last calibrated). Only the date element of this value should be shown. Exact time value is not relevant.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
```

```
BSError BSReadCalibrationTime(ByRef Int64 pCalibrationTime);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadCalibrationTime(ByRef Int64 pCalibrationTime);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadCalibrationTime(ByRef pCalibrationTime As Int64) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadCalibrationTime(ByRef Int64 pCalibrationTime);
```

MATLAB

```
BSReadCalibrationTime
```

Parameters

Parameters	Description
ByRef pCalibrationTime As Int64	[out] Receives the date the device was last calibrated.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadCalibrationTime( _
ByRef pCalibrationTime As Int64 _
) As BSError
End Function
```

1.3.1.1.3.16 BSReadCurrentTime Method

1.3.1.1.3.16.1 BodystatAPI.BSReadCurrentTime Method (DateTime, Integer)

Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).

C++

```
BSError BSReadCurrentTime(ByRef DateTime dtCurrentTime, Optional ByVal Integer iDST = -1 _);
```

C#

```
public Shared BSError BSReadCurrentTime(ByRef DateTime dtCurrentTime, Optional ByVal
Integer iDST);
```

Visual Basic

```
Public Shared Function BSReadCurrentTime(ByRef dtCurrentTime As DateTime, Optional ByVal
iDST As Integer = -1 _) As BSError
```

Java

```
public Shared BSError BSReadCurrentTime(ByRef DateTime dtCurrentTime, Optional ByVal
Integer iDST);
```

MATLAB

```
BSReadCurrentTime
```

Parameters

Parameters	Description
ByRef dtCurrentTime As DateTime	[out] Receives the current time of the device's internal clock.
Optional ByVal iDST As Integer = -1 _	Daylight Saving Time handling (as per Microsoft time functions) Zero (0) to indicate that standard time is in effect. A value greater than 0 to indicate that daylight saving time is in effect. A value less than zero to have the C run-time library code compute whether standard time or daylight saving time is in effect. (default)

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
Public Shared Function BSReadCurrentTime( _
    ByRef dtCurrentTime As DateTime _
    , Optional ByVal iDST As Integer = -1 _
) As BSError

    Dim llCurrentDateTime As Int64
    Dim bsErr As BSError = BSReadCurrentTime(llCurrentDateTime, iDST)

    dtCurrentTime = "#1/1/1970#"
    dtCurrentTime = dtCurrentTime.AddSeconds(llCurrentDateTime)
    If (iDST <> 0) Then
        dtCurrentTime = dtCurrentTime.ToLocalTime()
    End If

    Return bsErr
End Function
```

1.3.1.1.3.16.2 BodystatAPI.BSReadCurrentTime Method (Int64, Integer)

Read the current time from the internal clock of the connected device. Note not all device models support this function (not all have an internal clock).

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadCurrentTime(ByRef Int64 pCurrentTime, Optional ByVal Integer iDST = -1 _);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadCurrentTime(ByRef Int64 pCurrentTime, Optional ByVal Integer iDST);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadCurrentTime(ByRef pCurrentTime As Int64, Optional ByVal iDST
As Integer = -1 _) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadCurrentTime(ByRef Int64 pCurrentTime, Optional ByVal Integer iDST);
```

MATLAB

```
BSReadCurrentTime
```

Parameters

Parameters	Description
ByRef pCurrentTime As Int64	[out] Receives the current time of the device's internal clock.
Optional ByVal iDST As Integer = -1 _	Daylight Saving Time handling (as per Microsoft time functions) Zero (0) to indicate that standard time is in effect. A value greater than 0 to indicate that daylight saving time is in effect. A value less than zero to have the C run-time library code compute whether standard time or daylight saving time is in effect. (default)

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadCurrentTime( _
ByRef pCurrentTime As Int64 _
, Optional ByVal iDST As Integer = -1 _
) As BSError
End Function
```

1.3.1.1.3.17 BodystatAPI.BSReadModelVersion Method

Read device model and firmware version of the connected Bodystat device.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadModelVersion(ByRef BSDeviceModel pModel, ByRef Byte pMajorVersion, ByRef Byte
pMinorVersion, ByRef Byte pPsoc2Version, ByRef Byte pEepromVersion, ByRef Byte
pBluetoothInfo);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadModelVersion(ByRef BSDeviceModel pModel, ByRef Byte
pMajorVersion, ByRef Byte pMinorVersion, ByRef Byte pPsoc2Version, ByRef Byte
pEepromVersion, ByRef Byte pBluetoothInfo);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadModelVersion(ByRef pModel As BSDeviceModel, ByRef
pMajorVersion As Byte, ByRef pMinorVersion As Byte, ByRef pPsoc2Version As Byte, ByRef
pEepromVersion As Byte, ByRef pBluetoothInfo As Byte) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadModelVersion(ByRef BSDeviceModel pModel, ByRef Byte
pMajorVersion, ByRef Byte pMinorVersion, ByRef Byte pPsoc2Version, ByRef Byte
pEepromVersion, ByRef Byte pBluetoothInfo);
```

MATLAB

BSReadModelVersion

Parameters

Parameters	Description
ByRef pModel As BSDeviceModel	[out] Receives the device model.

ByRef pMajorVersion As Byte	[out] Receives the major version of the main firmware in the device.
ByRef pMinorVersion As Byte	[out] Receives the minor version of the main firmware in the device.
ByRef pPsoc2Version As Byte	[out] Receives the version of Psoc2 firmware in the device.
ByRef pEepromVersion As Byte	[out] Receives the version of the EEPROM firmware in the device.
ByRef pBluetoothInfo As Byte	[out] Receives the version of the Bluetooth module in the device.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadModelVersion( _
ByRef pModel As BSDeviceModel _
, ByRef pMajorVersion As Byte _
, ByRef pMinorVersion As Byte _
, ByRef pPsoc2Version As Byte _
, ByRef pEepromVersion As Byte _
, ByRef pBluetoothInfo As Byte _
) As BSError
End Function
```

1.3.1.1.3.18 BodystatAPI.BSReadPrinterAddress Method

Read the printer address the connected device is paired with (for direct printing). Not all device models support this function (not all support direct printing).

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadPrinterAddress(ByVal StringBuilder lpPrinterAddress, ByVal Integer
iPrinterAddressBufferSize);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadPrinterAddress(ByVal StringBuilder lpPrinterAddress, ByVal
Integer iPrinterAddressBufferSize);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadPrinterAddress(ByVal lpPrinterAddress As StringBuilder, ByVal
iPrinterAddressBufferSize As Integer) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadPrinterAddress(ByVal StringBuilder lpPrinterAddress, ByVal
Integer iPrinterAddressBufferSize);
```

MATLAB

```
BSReadPrinterAddress
```

Parameters

Parameters	Description
ByVal lpPrinterAddress As StringBuilder	[out] String buffer which receives the printer address.
ByVal iPrinterAddressBufferSize As Integer	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadPrinterAddress( _
ByVal lpPrinterAddress As StringBuilder _
, ByVal iPrinterAddressBufferSize As Integer _
) As BSError
End Function
```

1.3.1.1.3.19 BodystatAPI.BSReadProtocolInfo Method

Read protocol information from the connected device. Contains version numbers relating to communication protocols in use, data structures and auxiliary information.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadProtocolInfo(ByRef Byte pProtocolVersion, ByRef Byte pDataVersion, ByRef Byte
pAuxInfo);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadProtocolInfo(ByRef Byte pProtocolVersion, ByRef Byte
pDataVersion, ByRef Byte pAuxInfo);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadProtocolInfo(ByRef pProtocolVersion As Byte, ByRef
pDataVersion As Byte, ByRef pAuxInfo As Byte) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadProtocolInfo(ByRef Byte pProtocolVersion, ByRef Byte
pDataVersion, ByRef Byte pAuxInfo);
```

MATLAB

```
BSReadProtocolInfo
```

Parameters

Parameters	Description
ByRef pProtocolVersion As Byte	[out] Receives the communication protocol version.
ByRef pDataVersion As Byte	[out] Receives the data structure version.
ByRef pAuxInfo As Byte	[out] Receives auxiliary information.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Generally not required for the majority of SDK users as these details are abstracted by the SDK.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadProtocolInfo( _
ByRef pProtocolVersion As Byte _
, ByRef pDataVersion As Byte _
, ByRef pAuxInfo As Byte _
) As BSError
End Function
```

1.3.1.1.3.20 BodystatAPI.BSReadSerialNumber Method

Read the unique serial number of the connected Bodystat device.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadSerialNumber(ByRef Integer pSerialNumber);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadSerialNumber(ByRef Integer pSerialNumber);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadSerialNumber(ByRef pSerialNumber As Integer) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadSerialNumber(ByRef Integer pSerialNumber);
```

MATLAB

```
BSReadSerialNumber
```

Parameters

Parameters	Description
ByRef pSerialNumber As Integer	[out] Receives the serial number of the device.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadSerialNumber( _
ByRef pSerialNumber As Integer _
) As BSError
End Function
```

1.3.1.1.3.21 BodystatAPI.BSReadStoredTestData Method

Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged.

Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test

records before overwriting the oldest records automatically.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadStoredTestData(ByRef BSRawData pRawData);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadStoredTestData(ByRef BSRawData pRawData);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadStoredTestData(ByRef pRawData As BSRawData) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadStoredTestData(ByRef BSRawData pRawData);
```

MATLAB

```
BSReadStoredTestData
```

Parameters

Parameters	Description
ByRef pRawData As BSRawData	[out] Pointer to a structure to receive the downloaded data. Caller responsible for allocation and disposing of the structure.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Caller must check they are working with a Bodystat (rather than Equistat model, Equistat models will respond without error - but with garbage results).

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSReadStoredTestData( _
ByRef pRawData As BSRawData _
) As BSError
End Function
```

1.3.1.1.3.22 BodystatAPI.BSReadTestBodystat Method

Perform a quick bio-impedance test measurement on the connected bodystat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.

May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.

Useful for diagnostics or performing measurements instigated by software.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSReadTestBodystat(ByRef Integer pZ5, ByRef Integer pZ50, ByRef Integer pZ100,
ByRef Integer pZ200, ByRef Integer pR50, ByRef Single pX50, ByRef Single pPA50);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadTestBodystat(ByRef Integer pZ5, ByRef Integer pZ50, ByRef
Integer pZ100, ByRef Integer pZ200, ByRef Integer pR50, ByRef Single pX50, ByRef Single
pPA50);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadTestBodystat(ByRef pZ5 As Integer, ByRef pZ50 As Integer,
ByRef pZ100 As Integer, ByRef pZ200 As Integer, ByRef pR50 As Integer, ByRef pX50 As
Single, ByRef pPA50 As Single) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSReadTestBodystat(ByRef Integer pZ5, ByRef Integer pZ50, ByRef
Integer pZ100, ByRef Integer pZ200, ByRef Integer pR50, ByRef Single pX50, ByRef Single
pPA50);
```

MATLAB

```
BSReadTestBodystat
```

Parameters

Parameters	Description
ByRef pZ5 As Integer	[out] Receives the impedance measured at 5 kHz (result in ohms).
ByRef pZ50 As Integer	[out] Receives the impedance measured at 50 kHz (result in ohms).
ByRef pZ100 As Integer	[out] Receives the impedance measured at 100 kHz (result in ohms).
ByRef pZ200 As Integer	[out] Receives the impedance measured at 200 kHz (result in ohms).
ByRef pR50 As Integer	[out] Receives the resistance measured at 50 kHz (result in ohms).
ByRef pX50 As Single	[out] Receives the reactance measured at 50 kHz (result in ohms).
ByRef pPA50 As Single	[out] Receives the phase angle measured at 50 kHz (result in degrees).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Caller must check they are working with a Bodystat (rather than Equistat model, Equistat models will respond without error - but with garbage results).

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSReadTestBodystat( _
ByRef pZ5 As Integer _
, ByRef pZ50 As Integer _
, ByRef pZ100 As Integer _
, ByRef pZ200 As Integer _
, ByRef pR50 As Integer _
, ByRef pX50 As Single _
, ByRef pPA50 As Single _
) As BSError
End Function
```

1.3.1.1.3.23 BodystatAPI.BSSearchBTDevices Method

Search BT radio for any Bodystat devices (paired or otherwise). Existing paired devices will be counted even if they are switched off or out of range. New non-paired units will only be counted if switched on, in range and ready for communication (on main startup screen).

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
Integer BSSearchBTDevices(ByRef Integer iNewDevices, ByRef Integer iAuthenticatedDevices,
ByVal Boolean bReportErrors, ByVal IntPtr hWndParent, Optional ByVal UInt16 iTimeout = 7 _);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSSearchBTDevices(ByRef Integer iNewDevices, ByRef Integer
iAuthenticatedDevices, ByVal Boolean bReportErrors, ByVal IntPtr hWndParent, Optional ByVal
UInt16 iTimeout);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSSearchBTDevices(ByRef iNewDevices As Integer, ByRef
iAuthenticatedDevices As Integer, ByVal bReportErrors As Boolean, ByVal hWndParent As
IntPtr, Optional ByVal iTimeout As UInt16 = 7 _) As Integer
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSSearchBTDevices(ByRef Integer iNewDevices, ByRef Integer
iAuthenticatedDevices, ByVal Boolean bReportErrors, ByVal IntPtr hWndParent, Optional ByVal
UInt16 iTimeout);
```

MATLAB

BSSearchBTDevices

Parameters

Parameters	Description
ByRef iNewDevices As Integer	[out] The number of new Bodystat devices found (not paired).
ByRef iAuthenticatedDevices As Integer	[out] The number of authenticated Bodystat devices (paired).
ByVal bReportErrors As Boolean	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user (completely silent).
ByVal hWndParent As IntPtr	Handle of the parent window.
Optional ByVal iTimeout As UInt16 = 7 _	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task completes or timesout.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSSearchBTDevices( _
ByRef iNewDevices As Integer _
, ByRef iAuthenticatedDevices As Integer _
, ByVal bReportErrors As Boolean _
, ByVal hWndParent As IntPtr _
, Optional ByVal iTimeout As UInt16 = 7 _
) As Integer
End Function
```

1.3.1.1.3.24 BodystatAPI.BSUnAuthenticateBTDevices Method

Unauthenticate (unpair) ALL Bodystat devices from this PC.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
Integer BSUnAuthenticateBTDevices(ByVal Boolean bReportErrors, ByVal IntPtr hWndParent,
Optional ByVal UInt16 iTimeout = 7 _);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSUnAuthenticateBTDevices(ByVal Boolean bReportErrors, ByVal IntPtr
hWndParent, Optional ByVal UInt16 iTimeout);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSUnAuthenticateBTDevices(ByVal bReportErrors As Boolean, ByVal
hWndParent As IntPtr, Optional ByVal iTimeout As UInt16 = 7 _) As Integer
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared Integer BSUnAuthenticateBTDevices(ByVal Boolean bReportErrors, ByVal IntPtr
hWndParent, Optional ByVal UInt16 iTimeout);
```

MATLAB

BSUnAuthenticateBTDevices

Parameters

Parameters	Description
ByVal bReportErrors As Boolean	TRUE to report errors to the user (shown as message boxes, locale constraints apply). FALSE to hide errors from the user.
ByVal hWndParent As IntPtr	Handle of the parent window.
Optional ByVal iTimeout As UInt16 = 7 _	Timeout multiplier supplied to the Bluetooth radio search API. Controls the amount of time to spend searching for the device. This function will block for this duration so is better called from a thread. Each multiplier increment equates to about 1.28 seconds. Thus a timeout value of 10 = 12.8s. Maximum time permitted by the API is a multiplier of 45 (just under 1 minute), values higher than this will be capped at 45. Default value is 7 (about 10 seconds).

Returns

TRUE if it succeeds, FALSE if it fails.

Remarks

Note this process can be time consuming to complete and executes synchronously. Either call from a threaded function or illustrate operation in progress to the user e.g. by means of an indefinite progress bar or busy mouse cursor, until the task

completes or timesout.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSUnAuthenticateBTDevices( _
ByVal bReportErrors As Boolean _
, ByVal hWndParent As IntPtr _
, Optional ByVal iTimeout As UInt16 = 7 _
) As Integer
End Function
```

1.3.1.1.3.25 BodystatAPI.BSWriteCurrentTime Method

Reset the current date/time of the internal clock on the connected device. Date/time is automatically set to match the current date/time of the PC

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSWriteCurrentTime(_);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSWriteCurrentTime(_);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSWriteCurrentTime( As _) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSWriteCurrentTime(_ );
```

MATLAB

```
BSWriteCurrentTime
```

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Locale Daylight Saving Time (DST) rules are observed. Use sister function BSWriteCurrentTime(__time64_t CurrentTime) to set a specific date/time which differs from the PC clock or override DST behaviour.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSWriteCurrentTime( _
) As BSError
End Function
```

1.3.1.1.3.26 BSWriteCurrentTimeAs Method

1.3.1.1.3.26.1 BodystatAPI.BSWriteCurrentTimeAs Method (DateTime)

Reset the current date/time of the internal clock on the connected device. Date/time is set to a specific value of your choosing.

C++

```
BSError BSWriteCurrentTimeAs(ByVal DateTime dtCurrentTime);
```

C#

```
public Shared BSError BSWriteCurrentTimeAs(ByVal DateTime dtCurrentTime);
```

Visual Basic

```
Public Shared Function BSWriteCurrentTimeAs(ByVal dtCurrentTime As DateTime) As BSError
```

Java

```
public Shared BSError BSWriteCurrentTimeAs(ByVal DateTime dtCurrentTime);
```

MATLAB

```
BSWriteCurrentTimeAs
```

Parameters

Parameters	Description
ByVal dtCurrentTime As DateTime	Date/time to apply to the internal clock of the device (should be current!).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Use sister function (omit date/time parameter entirely) BSWriteCurrentTime (see page 138)() to automatically set to the current date/time of the PC.

Body Source

```
Public Shared Function BSWriteCurrentTimeAs( _
    ByVal dtCurrentTime As DateTime _
) As BSError

    Dim dtBaseTime As DateTime = "#1/1/1970#"
    Dim elapsed = dtCurrentTime - dtBaseTime
    Dim llCurrentTime As Int64
    llCurrentTime = elapsed.TotalSeconds

    Return BSWriteCurrentTimeAs(llCurrentTime)
End Function
```

1.3.1.1.3.26.2 BodystatAPI.BSWriteCurrentTimeAs Method (Int64)

Reset the current date/time of the internal clock on the connected device. Date/time is set to a specific value of your choosing.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
    CharSet:=CharSet.Unicode)]
BSError BSWriteCurrentTimeAs(ByVal Int64 CurrentTime);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
    CharSet:=CharSet.Unicode)]
public Shared BSError BSWriteCurrentTimeAs(ByVal Int64 CurrentTime);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
    CharSet:=CharSet.Unicode)>
Public Shared Function BSWriteCurrentTimeAs(ByVal CurrentTime As Int64) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
    CharSet:=CharSet.Unicode)]
public Shared BSError BSWriteCurrentTimeAs(ByVal Int64 CurrentTime);
```

MATLAB

```
BSWriteCurrentTimeAs
```

Parameters

Parameters	Description
ByVal CurrentTime As Int64	Date/time to apply to the internal clock of the device (should be current!).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Use sister function (omit date/time parameter entirely) BSWriteCurrentTime (see page 138)() to automatically set to the current date/time of the PC.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSWriteCurrentTimeAs( _
ByVal CurrentTime As Int64 _
) As BSError
End Function
```

1.3.1.1.3.27 BodystatAPI.BSWritePrinterAddress Method

Set the address of the printer the connected device should use for direct printing. Programmatic implementation of the inbuilt printer pairing function in the device.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError BSWritePrinterAddress(ByVal StringBuilder lpPrinterAddress, ByVal Integer
iPrinterAddressBufferSize);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSWritePrinterAddress(ByVal StringBuilder lpPrinterAddress, ByVal
Integer iPrinterAddressBufferSize);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function BSWritePrinterAddress(ByVal lpPrinterAddress As StringBuilder, ByVal
iPrinterAddressBufferSize As Integer) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError BSWritePrinterAddress(ByVal StringBuilder lpPrinterAddress, ByVal
Integer iPrinterAddressBufferSize);
```

MATLAB

BSWritePrinterAddress

Parameters

Parameters	Description
ByVal lpPrinterAddress As StringBuilder	String buffer containing the Bluetooth address (BDA) of the printer.
iBufferSize	Size of the buffer.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Note not all device models support this function (not all support direct printing).

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function BSWritePrinterAddress( _
ByVal lpPrinterAddress As StringBuilder _
, ByVal iPrinterAddressBufferSize As Integer _
) As BSError
End Function
```

1.3.1.1.3.28 BodystatAPI.ESCalculateNormals Method

Calculate normal values and normal ranges for a given animal and their bio-impedance measurement together with the matching results structure for that measurement. You must supply a populated results structure which pertains to the measurement being supplied (obtained from `ESCalculateResults` (see page 142)).

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError ESCalculateNormals(ByRef ESMeasurement pM, ByRef ESResults pR, ByRef ESNormals pN);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESCalculateNormals(ByRef ESMeasurement pM, ByRef ESResults pR, ByRef
ESNormals pN);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function ESCalculateNormals(ByRef pM As ESMeasurement, ByRef pR As ESResults,
ByRef pN As ESNormals) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESCalculateNormals(ByRef ESMeasurement pM, ByRef ESResults pR, ByRef
ESNormals pN);
```

MATLAB

```
ESCalculateNormals
```

Parameters

Parameters	Description
ByRef pM As ESMeasurement	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
ByRef pR As ESResults	The results calculated for this measurement previously.
ByRef pN As ESNormals	[out] Receives the normal values and normal ranges for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Normals vary by input attributes, measured bio-impedance values and calculated results.

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function ESCalculateNormals( _
```



```

ByRef pM As ESMeasurement _
, ByRef pR As ESResults _
, ByRef pN As ESNormals _
) As BSError
End Function

```

1.3.1.1.3.29 BodystatAPI.ESCalculateResults Method

Calculate complete result set for a given bio-impedance measurement recorded by the Equistat device.

C++

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError ESCalculateResults(ByRef ESMeasurement pM, ByRef ESResults pR);

```

C#

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESCalculateResults(ByRef ESMeasurement pM, ByRef ESResults pR);

```

Visual Basic

```

<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function ESCalculateResults(ByRef pM As ESMeasurement, ByRef pR As ESResults)
As BSError

```

Java

```

[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESCalculateResults(ByRef ESMeasurement pM, ByRef ESResults pR);

```

MATLAB

```
ESCalculateResults
```

Parameters

Parameters	Description
ByRef pM As ESMeasurement	The measurement details downloaded from the device (includes model, user inputted parameters and bio-impedance measurement readings).
ByRef pR As ESResults	[out] Receives the complete calculated results for this measurement.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

Results vary by input attributes and measured bio-impedance values.

Body Source

```

<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function ESCalculateResults( _
ByRef pM As ESMeasurement _
, ByRef pR As ESResults _
) As BSError
End Function

```

1.3.1.1.3.30 BodystatAPI.ESReadStoredTestData Method

Download test data from the connected device. Retrieves the last 100 test measurements (records) stored in the memory of the device. Downloading data does not affect the results stored in the memory - they remain unchanged.

Note: you should store downloaded data locally if it is required permanently. Most devices store only the last 100 test records before overwriting the oldest records automatically.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError ESReadStoredTestData(ByRef ESRawData pRawData);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESReadStoredTestData(ByRef ESRawData pRawData);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function ESReadStoredTestData(ByRef pRawData As ESRawData) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESReadStoredTestData(ByRef ESRawData pRawData);
```

MATLAB

```
ESReadStoredTestData
```

Parameters

Parameters	Description
ByRef pRawData As ESRawData	[out] Pointer to a structure to receive the downloaded data. Caller responsible for allocation and disposing of the structure.

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Caller must check they are working with an Equistat (rather than Bodystat model, Bodystat models will respond without error - but with garbage results).

Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function ESReadStoredTestData( _
ByRef pRawData As ESRawData _
) As BSError
End Function
```

1.3.1.1.3.31 BodystatAPI.ESReadTestEquistat Method

Perform a quick bio-impedance test measurement on the connected Equistat device. Normal test process is skipped (user is not prompted to enter any subject information). Normal test protocol is disregarded (measurement proceeds regardless of lead status). Test measurement is not stored in device memory.

May give wildly erratic readings if leads are not connected properly. Caller should sanity check results to detect for such scenarios. Results are intentionally not sanity checked to provide a true raw reading and provide caller with total discretion.

Useful for diagnostics or performing measurements instigated by software.

C++

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
BSError ESReadTestEquistat(ByRef Integer pZ5, ByRef Integer pZ16, ByRef Integer pZ24, ByRef
Integer pZ50, ByRef Integer pZ140, ByRef Integer pZ200, ByRef Integer pZ280, ByRef Single
pPA50);
```

C#

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
```

```
CharSet:=CharSet.Unicode)]
public Shared BSError ESReadTestEquistat(ByRef Integer pZ5, ByRef Integer pZ16, ByRef
Integer pZ24, ByRef Integer pZ50, ByRef Integer pZ140, ByRef Integer pZ200, ByRef Integer
pZ280, ByRef Single pPA50);
```

Visual Basic

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)>
Public Shared Function ESReadTestEquistat(ByRef pZ5 As Integer, ByRef pZ16 As Integer,
ByRef pZ24 As Integer, ByRef pZ50 As Integer, ByRef pZ140 As Integer, ByRef pZ200 As
Integer, ByRef pZ280 As Integer, ByRef pPA50 As Single) As BSError
```

Java

```
[DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)]
public Shared BSError ESReadTestEquistat(ByRef Integer pZ5, ByRef Integer pZ16, ByRef
Integer pZ24, ByRef Integer pZ50, ByRef Integer pZ140, ByRef Integer pZ200, ByRef Integer
pZ280, ByRef Single pPA50);
```

MATLAB

```
ESReadTestEquistat
```

Parameters

Parameters	Description
ByRef pZ5 As Integer	[out] Receives the impedance measured at 5 kHz (result in ohms).
ByRef pZ16 As Integer	[out] Receives the impedance measured at 16 kHz (result in ohms).
ByRef pZ24 As Integer	[out] Receives the impedance measured at 24 kHz (result in ohms).
ByRef pZ50 As Integer	[out] Receives the impedance measured at 50 kHz (result in ohms).
ByRef pZ140 As Integer	[out] Receives the impedance measured at 140 kHz (result in ohms).
ByRef pZ200 As Integer	[out] Receives the impedance measured at 200 kHz (result in ohms).
ByRef pZ280 As Integer	[out] Receives the impedance measured at 280 kHz (result in ohms).
ByRef pPA50 As Single	[out] Receives the phase angle measured at 50 kHz (result in degrees) (ES Pro Only).

Returns

Returns BSError::NoError if successful. Specific error code otherwise.

Remarks

You must have connected to the device by calling BSConnect (see page 120) previously. Caller must check they are working with an Equistat (rather than Bodystat model, Bodystat models will respond without error - but with garbage results).






Body Source

```
<DllImport("BodystatSDK.dll", CallingConvention:=CallingConvention.Cdecl,
CharSet:=CharSet.Unicode)> _
Public Shared Function ESReadTestEquistat( _
ByRef pZ5 As Integer _
, ByRef pZ16 As Integer _
, ByRef pZ24 As Integer _
, ByRef pZ50 As Integer _
, ByRef pZ140 As Integer _
, ByRef pZ200 As Integer _
, ByRef pZ280 As Integer _
, ByRef pPA50 As Single _
) As BSError
End Function
```

1.3.2 Structs, Records, Enums

The following table lists structs, records, enums in this documentation.

Enumerations

	Name	Description
	BSAnalysisMode (see page 145)	BSAnalysisMode: Analysis results are available in several formats: BSAnalysisModeBC= Body composition analysis BSAnalysisModeHN= Hydration analysis BSAnalysisModeBoth= Combined body composition and hydration analysis (default)
	BSDeviceFamily (see page 146)	BSDeviceFamily: Enumerated values that group each device / model variant into a given product family. Devices within a given product family share broad capabilities and features. For details regarding capabilities and features refer to the technical documentation for that product family (http://www.bodystat.com)
	BSDeviceModel (see page 147)	BSDeviceModel: Enumerated values that represent the various different models of the Bodystat devices over the years. Suffix: The suffix after the model indicates specific hardware variations regarding connectivity: OPTO indicates model communicates via opto-isolated serial. DIU indicates model communicates via a data interface unit (induction interface) BT indicates model communicates via Bluetooth serial port.
	BSError (see page 149)	BSError: Bodystat Error codes. A list of Bodystat specific error codes.
	BSGender (see page 151)	BSGender: Enumerated values that represent a person's gender.

Module

Visual Basic .NET wrapper ([see page 80](#))

1.3.2.1 BodystatAPI.BSAnalysisMode Enumeration

BSAnalysisMode: Analysis results are available in several formats:

BSAnalysisModeBC= Body composition analysis

BSAnalysisModeHN= Hydration analysis

BSAnalysisModeBoth= Combined body composition and hydration analysis (default)

C++

```
enum BSAnalysisMode {
    BSAnalysisModeNone = 0,
    BSAnalysisModeBC = 1,
    BSAnalysisModeHN = 2,
    BSAnalysisModeBoth = 3
};
```

C#

```
public enum BSAnalysisMode {
    BSAnalysisModeNone = 0,
    BSAnalysisModeBC = 1,
    BSAnalysisModeHN = 2,
    BSAnalysisModeBoth = 3
}
```

Visual Basic

```
Public Enum BSAnalysisMode
    BSAnalysisModeNone = 0
    BSAnalysisModeBC = 1
    BSAnalysisModeHN = 2
    BSAnalysisModeBoth = 3
End Enum
```

Java

```
public class BSAnalysisMode;
```

MATLAB

```
BSAnalysisMode
```

File

```
BodystatSDK.vb
```

Members

Members	Description
BSAnalysisModeNone = 0	No analysis mode specified (raw data only)
BSAnalysisModeBC = 1	Body composition enabled
BSAnalysisModeHN = 2	Hydration enabled
BSAnalysisModeBoth = 3	Both body composition and hydration (default)

1.3.2.2 BodystatAPI.BSDeviceFamily Enumeration

BSDeviceFamily: Enumerated values that group each device / model variant into a given product family. Devices within a given product family share broad capabilities and features. For details regarding capabilities and features refer to the technical documentation for that product family (<http://www.bodystat.com>)

C++

```
enum BSDeviceFamily {
    BSUnknown_Family = 0,
    BS1500_Family = -1,
    BSMDD_Family = -2,
    BSDualScan_Family = -3,
    BSQuadScan_Family = -4,
    BSMultiScan_Family = -5,
    BSEquistatLite_Family = -6,
    BSEquistatPro_Family = -7
};
```

C#

```
public enum BSDeviceFamily {
    BSUnknown_Family = 0,
    BS1500_Family = -1,
    BSMDD_Family = -2,
    BSDualScan_Family = -3,
    BSQuadScan_Family = -4,
    BSMultiScan_Family = -5,
    BSEquistatLite_Family = -6,
    BSEquistatPro_Family = -7
}
```

Visual Basic

```
Public Enum BSDeviceFamily
    BSUnknown_Family = 0
    BS1500_Family = -1
    BSMDD_Family = -2
    BSDualScan_Family = -3
    BSQuadScan_Family = -4
    BSMultiScan_Family = -5
```

```
BSEquistatLite_Family = -6
BSEquistatPro_Family = -7
End Enum
```

Java

```
public class BSDeviceFamily;
```

MATLAB

```
BSDeviceFamily
```

File

```
BodystatSDK.vb
```

Members

Members	Description
BSUnknown_Family = 0	Unknown product family
BS1500_Family = -1	Bodystat 1500 family
BSMDD_Family = -2	Bodystat MDD family
BSDualScan_Family = -3	Bodystat DualScan family (legacy)
BSQuadScan_Family = -4	Bodystat QuadScan family
BSMultiScan_Family = -5	Bodystat Multiscan family
BSEquistatLite_Family = -6	Equistat Lite family (not supported by this SDK)
BSEquistatPro_Family = -7	Equistat Pro family (not supported by this SDK)

Remarks

Note that not all families are supported by this SDK. Some are listed for historic reference only.

1.3.2.3 BodystatAPI.BSDeviceModel Enumeration

BSDeviceModel: Enumerated values that represent the various different models of the Bodystat devices over the years.

Suffix: The suffix after the model indicates specific hardware variations regarding connectivity:

OPTO indicates model communicates via opto-isolated serial.

DIU indicates model communicates via a data interface unit (induction interface)

BT indicates model communicates via Bluetooth serial port.

C++

```
enum BSDeviceModel {
    BSUnknown = 0,
    BS1500_OPTO = 1,
    BSMDD_OPTO = 2,
    BS1500_DIU = 3,
    BSMDD_DIU = 4,
    BSDualScan_OPTO = 5,
    BSDualScan_DIU = 6,
    BSQuadScan_DIU = 7,
    BSMultiScan_DIU = 8,
    BS1500_BT = 9,
    BSMDD_BT = 10,
    BSQuadScan_BT = 11,
    BSEquistatPro_BT = 12,
    BSEquistatLite_BT = 13,
    BSQuadScanTouch = 14,
    BSMultiScanTouch = 15
};
```

C#

```
public enum BSDeviceModel {
    BSUnknown = 0,
```

```

    BS1500_OPTO = 1,
    BSMDD_OPTO = 2,
    BS1500_DIU = 3,
    BSMDD_DIU = 4,
    BSDualScan_OPTO = 5,
    BSDualScan_DIU = 6,
    BSQuadScan_DIU = 7,
    BSMultiScan_DIU = 8,
    BS1500_BT = 9,
    BSMDD_BT = 10,
    BSQuadScan_BT = 11,
    BSEquistatPro_BT = 12,
    BSEquistatLite_BT = 13,
    BSQuadScanTouch = 14,
    BSMultiScanTouch = 15
}

```

Visual Basic

```

Public Enum BSDeviceModel
    BSUnknown = 0
    BS1500_OPTO = 1
    BSMDD_OPTO = 2
    BS1500_DIU = 3
    BSMDD_DIU = 4
    BSDualScan_OPTO = 5
    BSDualScan_DIU = 6
    BSQuadScan_DIU = 7
    BSMultiScan_DIU = 8
    BS1500_BT = 9
    BSMDD_BT = 10
    BSQuadScan_BT = 11
    BSEquistatPro_BT = 12
    BSEquistatLite_BT = 13
    BSQuadScanTouch = 14
    BSMultiScanTouch = 15
End Enum

```

Java

```
public class BSDeviceModel;
```

MATLAB

```
BSDeviceModel
```

File

```
BodystatSDK.vb
```

Members

Members	Description
BSUnknown = 0	Unknown model
BS1500_OPTO = 1	Bodystat 1500 with opto-isolated interface (legacy)
BSMDD_OPTO = 2	Bodystat MDD with opto-isolated interface (legacy)
BS1500_DIU = 3	Bodystat 1500 with data interface unit (legacy)
BSMDD_DIU = 4	Bodystat MDD with data interface unit (legacy)
BSDualScan_OPTO = 5	Bodystat DualScan with opto-isolated interface (legacy)
BSDualScan_DIU = 6	Bodystat DualScan with data interface unit (legacy)
BSQuadScan_DIU = 7	Bodystat QuadScan with data interface unit (legacy)
BSMultiScan_DIU = 8	Bodystat MultiScan with data interface unit (legacy)
BS1500_BT = 9	Bodystat 1500 with Bluetooth interface
BSMDD_BT = 10	Bodystat 1500MDD with Bluetooth interface
BSQuadScan_BT = 11	Bodystat QuadScan 4000 with Bluetooth interface
BSEquistatPro_BT = 12	Equistat Pro with Bluetooth interface (not supported by this SDK)

BSEquistatLite_BT = 13	Equistat Lite with Bluetooth interface (not supported by this SDK)
BSQuadScanTouch = 14	Bodystat QuadScan 4000 Touch with WiFi/Bluetooth interface
BSMultiScanTouch = 15	Bodystat Multiscan with WiFi/Bluetooth interface

Remarks

Note that not all models are supported by this SDK. Some are listed for historic reference only.

1.3.2.4 BodystatAPI.BSError Enumeration

BSError: Bodystat Error codes. A list of Bodystat specific error codes.

C++

```
enum BSError {
    NoError = 0,
    InvalidDeviceModel = -1,
    InvalidDeviceFamily = -2,
    InvalidDeviceMode = -3,
    InvalidGender = -4,
    InvalidAge = -5,
    InvalidHeight = -6,
    InvalidWeight = -7,
    InvalidActivity = -8,
    InvalidWaist = -9,
    InvalidHip = -10,
    InvalidImpedanceZ5 = -11,
    InvalidImpedanceZ50 = -12,
    InvalidImpedanceZ100 = -13,
    InvalidImpedanceZ200 = -14,
    InvalidResistanceR50 = -15,
    InvalidReactanceX50 = -16,
    InvalidPhaseAnglePA50 = -17,
    InvalidParamSupplied = -18,
    ComErrorOpeningPort = -50,
    ComErrorConfiguringPort = -51,
    ComErrorConnecting = -52,
    ComErrorUnsupportedCommand = -53,
    ComErrorProcessingCommand = -54,
    ComErrorReading = -55,
    ComErrorWriting = -56,
    ComErrorUserAborted = -57
};
```

C#

```
public enum BSError {
    NoError = 0,
    InvalidDeviceModel = -1,
    InvalidDeviceFamily = -2,
    InvalidDeviceMode = -3,
    InvalidGender = -4,
    InvalidAge = -5,
    InvalidHeight = -6,
    InvalidWeight = -7,
    InvalidActivity = -8,
    InvalidWaist = -9,
    InvalidHip = -10,
    InvalidImpedanceZ5 = -11,
    InvalidImpedanceZ50 = -12,
    InvalidImpedanceZ100 = -13,
    InvalidImpedanceZ200 = -14,
    InvalidResistanceR50 = -15,
    InvalidReactanceX50 = -16,
    InvalidPhaseAnglePA50 = -17,
    InvalidParamSupplied = -18,
    ComErrorOpeningPort = -50,
```



```

ComErrorConfiguringPort = -51,
ComErrorConnecting = -52,
ComErrorUnsupportedCommand = -53,
ComErrorProcessingCommand = -54,
ComErrorReading = -55,
ComErrorWriting = -56,
ComErrorUserAborted = -57
}

```

Visual Basic

```

Public Enum BSError
    NoError = 0
    InvalidDeviceModel = -1
    InvalidDeviceFamily = -2
    InvalidDeviceMode = -3
    InvalidGender = -4
    InvalidAge = -5
    InvalidHeight = -6
    InvalidWeight = -7
    InvalidActivity = -8
    InvalidWaist = -9
    InvalidHip = -10
    InvalidImpedanceZ5 = -11
    InvalidImpedanceZ50 = -12
    InvalidImpedanceZ100 = -13
    InvalidImpedanceZ200 = -14
    InvalidResistanceR50 = -15
    InvalidReactanceX50 = -16
    InvalidPhaseAnglePA50 = -17
    InvalidParamSupplied = -18
    ComErrorOpeningPort = -50
    ComErrorConfiguringPort = -51
    ComErrorConnecting = -52
    ComErrorUnsupportedCommand = -53
    ComErrorProcessingCommand = -54
    ComErrorReading = -55
    ComErrorWriting = -56
    ComErrorUserAborted = -57
End Enum

```

Java

```
public class BSError;
```

MATLAB

```
BSError
```

File

```
BodystatSDK.vb
```

Members

Members	Description
NoError = 0	No error occurred
InvalidDeviceModel = -1	Invalid model specified
InvalidDeviceFamily = -2	Invalid family specified
InvalidDeviceMode = -3	Invalid analysis mode
InvalidGender = -4	Invalid BSGender (see page 151) specified
InvalidAge = -5	Invalid age specified
InvalidHeight = -6	Invalid height measurement specified
InvalidWeight = -7	Invalid weight measurement specified
InvalidActivity = -8	Invalid activity level specified
InvalidWaist = -9	Invalid waist measurement specified
InvalidHip = -10	Invalid hip measurement specified
InvalidImpedanceZ5 = -11	Invalid impedance measurement specified

InvalidImpedanceZ50 = -12	Invalid impedance measurement specified
InvalidImpedanceZ100 = -13	Invalid impedance measurement specified
InvalidImpedanceZ200 = -14	Invalid impedance measurement specified
InvalidResistanceR50 = -15	Invalid resistance measurement specified
InvalidReactanceX50 = -16	Invalid reactance measurement specified
InvalidPhaseAnglePA50 = -17	Invalid phase angle measurement specified
InvalidParamSupplied = -18	Invalid parameter supplied
ComErrorOpeningPort = -50	Error opening the com port
ComErrorConfiguringPort = -51	Error configuring the com port parameters (required baud rate, buffers, stop bits, etc)
ComErrorConnecting = -52	Error connecting to the device
ComErrorUnsupportedCommand = -53	The requested command is not supported by the model being asked to perform it
ComErrorProcessingCommand = -54	Error occurred whilst processing the requested command
ComErrorReading = -55	Error occurred whilst reading data from the com port
ComErrorWriting = -56	Error occurred whilst writing data to the com port
ComErrorUserAborted = -57	The user aborted the process

1.3.2.5 BodystatAPI.BSGender Enumeration

BSGender: Enumerated values that represent a person's gender.

C++

```
enum BSGender {
    BSFemale = 0,
    BSMale = 1
};
```

C#

```
public enum BSGender {
    BSFemale = 0,
    BSMale = 1
}
```

Visual Basic

```
Public Enum BSGender
    BSFemale = 0
    BSMale = 1
End Enum
```

Java

```
public class BSGender;
```

MATLAB

```
BSGender
```

File

```
BodystatSDK.vb
```

2 FAQs

Frequently Asked Questions

- C++ Linking Errors:

If you receive errors regarding unresolved symbols when building/linking your project you have most likely forgotten to add the BodystatSDK.lib file to your project.

You can do this by adding a #pragma lib directive to one of your source files like those shown below. An example using this method is also shown in the MFC sample (BodystatMFCDlg.cpp).

```
#if defined _M_X64
#pragma comment ( lib, "..\\..\\BodystatSDK\\Lib\\x64\\BodystatSDK.lib" )
#pragma message ( "BodystatMFCDlg.cpp::Linking to 64-bit BodystatSDK.lib" )
#else
#pragma comment ( lib, "..\\..\\BodystatSDK\\Lib\\x86\\BodystatSDK.lib" )
#pragma message ( "BodystatMFCDlg.cpp::Linking to 32-bit BodystatSDK.lib" )
#endif
```

Alternatively, you can reference the Bodystat.lib dependency in the additional dependencies section of your project settings (Project Settings->Configuration Properties->Linker->Input->Additional Dependencies). You will need to specify the file using a relative path to the project or add the path to your VC Library directories.

A simpler alternative still, available to projects in VS2010 (or later), is to simply add the Bodystat.lib file to your solution/project. VS2010 understands .lib files and will do the rest for you.

If you continue to encounter unresolved symbols relating to the Bodystat API, and you have referenced the Bodystat.lib file, another possibility is that your project is using an unsupported character set. The Bodystat SDK is designed to use the Unicode character set. Check your project settings are configured to use Unicode (Project Settings->Configuration Properties->General->Character Set). You may use another character set in your project if you prefer, but ensure your calls to the Bodystat SDK work with Unicode strings.

Once your project builds properly, whichever method you chose, remember you must also copy the BodystatSDK.dll dependency alongside your program.exe or you will see a runtime error. The Bodystat MFC sample automatically copies the dll alongside the app as a pre-build event command line. Ultimately you will need to redistribute the BodystatSDK.DLL in your setup routine for installation along with your application. Remember to reference the appropriate DLL binary for your target platform (e.g. x86). See the Redistribution (see page 154) section for more information.

- .NET Runtime Error: "BadImageFormatException"

If during runtime, your program complains it is unable to load the BodystatSDK.dll because it is an incorrect format you should check the CPU target settings in your configuration. Be careful if you target 'AnyCPU' as this is not presently supported by our SDK (though you can might be able to work around this with careful mappings in your assembly). Instead you simply specify your specific target (x86/x64) then ensure you are referencing the corresponding x86/x64 version of the

BodystatSDK.dll.

3 Redistribution

Redistribution

Once your application is complete and ready for shipping you will need to redistribute the appropriate BodystatSDK.dll together with your application.

The BodystatSDK.DLL can be found in the Lib\x86 or Lib\x64 folder of the Bodystat SDK installation, for 32-bit and 64-bit installations respectively. Simply configure your installation script to install the appropriate DLL alongside your application binaries. Note that 64-bit support is planned and not available at time of writing (please contact us if you have a need for it).

Note: You need only redistribute the BodystatSDK.DLL file. The other .exp and .lib files are only required by your development environment.

The BodystatSDK.dll has no dependencies other than the Windows kernel.

Redistribution of our BodystatSDK.dll remains subject to the terms of our license agreement.

Finally, we recommend that you should always develop and test your software using the latest version of the Bodystat SDK.

4 Samples

Sample projects included within the SDK. The sample projects provide detailed examples which utilise the various functions within the Bodystat API. Samples are provided in a variety of popular programming languages (though we regret we are unable to provide examples for all languages). Please reference the sample in the language you are most familiar with.



Modules

Name	Description
Bodystat MFC Sample (see page 155)	A sample application written in Visual C++ based upon the Microsoft Foundation Classes (MFC). Solutions are provided for Visual Studio 2008 and 2010. You must have Visual C++ 2008 or 2010 installed to work with this sample.
Bodystat C# WinForms Sample (see page 204)	A sample application written in C# .NET based upon WinForms within the Microsoft .NET Framework. Solutions are provided for Visual Studio 2008 and 2010. You must have Visual Studio 2008 or 2010 installed to work with this sample.
Bodystat Visual Basic .NET Sample (see page 216)	A sample application written in Visual Basic .NET based upon the Microsoft .NET Framework. Solutions are provided for Visual Basic .NET 2008 and 2010. You must have Visual Basic .NET 2008 or 2010 installed to work with this sample.

4.1 Bodystat MFC Sample

A sample application written in Visual C++ based upon the Microsoft Foundation Classes (MFC). Solutions are provided for Visual Studio 2008 and 2010. You must have Visual C++ 2008 or 2010 installed to work with this sample.



Classes

	Name	Description
	CBodystatMFCDlg (see page 156)	CBodystatMFCDlg. Dialog for viewing downloaded Bodystat test measurements records and testing key functions within the Bodystat API.
	CBodystatResultsDlg (see page 185)	CBodystatResultsDlg. Dialog for viewing the Bodystat test measurement results and normals.

4.1.1 Classes

The following table lists classes in this documentation.

Classes

	Name	Description
	CBodystatMFCDlg (see page 156)	CBodystatMFCDlg. Dialog for viewing downloaded Bodystat test measurements records and testing key functions within the Bodystat API.
	CBodystatResultsDlg (see page 185)	CBodystatResultsDlg. Dialog for viewing the Bodystat test measurement results and normals.

Module

Bodystat MFC Sample (see page 155)

4.1.1.1 CBodystatMFCDlg Class

CBodystatMFCDlg. Dialog for viewing downloaded Bodystat test measurements records and testing key functions within the Bodystat API.

Class Hierarchy



C++

```
class CBodystatMFCDlg : public CDialogEx;
```

C#

```
public class CBodystatMFCDlg : CDialogEx;
```

Visual Basic

```
Public Class CBodystatMFCDlg
Inherits CDialogEx
```

Java

```
public class CBodystatMFCDlg CDialogEx;
```

MATLAB

```
CBodystatMFCDlg
```

File







BodystatMFCDlg.h

Methods
















	Name	Description
	CBodystatMFCDlg (see page 158)	The main CBodystatMFCDlg constructor











CBodystatMFCDlg Data Members

	Name	Description
	m_bsRawData (see page 159)	Used to hold the downloaded test measurements from the device
	m_btnAuthenticate (see page 159)	The button authenticate
	m_btnAutoSetup (see page 159)	The button automatic setup
	m_btnConnect (see page 160)	The button connect
	m_btnDisconnect (see page 160)	The button disconnect
	m_btnReadCalibrationTime (see page 160)	The button read calibration time
	m_btnReadCurrentTime (see page 161)	The button read current time
	m_btnReadModelVersion (see page 161)	The button read model version
	m_btnReadPrinterAddr (see page 161)	The button read printer address
	m_btnReadProtocolInfo (see page 162)	The button read protocol info
	m_btnReadSerialNumber (see page 162)	The button read serial number
	m_btnReadStoredTestData (see page 162)	The button read stored test data
	m_btnReadTestMeasurement (see page 162)	The button read test measurement

	m_btnUnauthenticate (see page 163)	The button unauthenticate
	m_btnWriteCurrentTime (see page 163)	The button write current time
	m_hlcon (see page 163)	Application icon
	m_iTimeout (see page 164)	Bluetooth timeout multiplier. Initial value of 7 (equates to about 10 seconds). This is incremented upon each Bluetooth search. A maximum value of 45 is permitted (about 1 minute).
	m_strComPort (see page 164)	The string com port
	m_TestMeasurementListCtrl (see page 164)	The test measurement list control

CBodystatMFCDlg Methods

	Name	Description
	DoDataExchange (see page 165)	DDX/DDV support Bind the form controls to member variables and set basic input range limits
	OnClickedBtAuthenticateButton (see page 165)	Executes the Bluetooth Authenticate button action. Illustrates use of BSAuthenticateBTDevice (see page 8) from the Bodystat API
	OnClickedBtAutoSetupButton (see page 166)	Executes the Bluetooth Auto Setup button action. Illustrates use of BSAutoSetupBT (see page 8) from the Bodystat API
	OnClickedBtDeviceInfoButton (see page 167)	Executes the Bluetooth Device Info button action. Illustrates use of BSGetBTBodystatDevice (see page 12) from the Bodystat API
	OnClickedBtSearchButton (see page 167)	Executes the Bluetooth Search button action. Illustrates use of BSSearchBTDevices (see page 22) from the Bodystat API
	OnClickedBtStackInfoButton (see page 168)	Executes the Bluetooth Get Stack Info button action. Illustrates use of BSGetBTStackInfo (see page 13) from the Bodystat API
	OnClickedBtStatusButton (see page 169)	Executes the Bluetooth Get Status button action. Illustrates use of BSIsBTAvailable (see page 16) from the Bodystat API
	OnClickedBtUnauthenticateButton (see page 169)	Executes the Bluetooth Unauthenticate button action. Illustrates use of BSUnAuthenticateBTDevices (see page 23) from the Bodystat API
	OnClickedComConnectButton (see page 170)	Executes the communication Connect button action. Illustrates use of BSOpenComport (see page 16) from the Bodystat API
	OnClickedComReadCalibrationTimeButton (see page 170)	Executes the communication Read Calibration Time button action. Illustrates use of BSReadCalibrationTime (see page 17) from the Bodystat API
	OnClickedComReadCurrentTimeButton (see page 171)	Executes the communication Read Current Time button action. Illustrates use of BSReadCurrentTime (see page 17) from the Bodystat API
	OnClickedComReadModelVersionButton (see page 172)	Executes the communication Read Model Version button action. Illustrates use of BSReadModelVersion (see page 18) from the Bodystat API
	OnClickedComReadPrinterAddrButton (see page 173)	Executes the communication Read Printer Address button action. Illustrates use of BSReadPrinterAddress (see page 19) from the Bodystat API
	OnClickedComReadProtocolInfoButton (see page 174)	Executes the communication ReadProtocolInfo button action. Illustrates use of BSReadProtocolInfo (see page 19) from the Bodystat API
	OnClickedComReadSerialNumberButton (see page 175)	Executes the communication Read Serial Number button action. Illustrates use of BSReadSerialNumber (see page 20) from the Bodystat API

	OnClickedComReadStoredTestDataButton (see page 175)	Executes the communication Read Stored Test Data button action. Illustrates use of BSReadStoredTestData (see page 21) from the Bodystat API to download data
	OnClickedComReadTestMeasurementButton (see page 176)	Executes the communication Read Test Measurement button action. Performs an immediate test measurement on the device using the BSReadTestBodystat (see page 21) function from the Bodystat API
	OnClickedComWriteCurrentTimeButton (see page 177)	Executes the communication Write Current Time button action. Illustrates use of BSWriteCurrentTime (see page 24) from the Bodystat API
	OnClickedResultsButton (see page 178)	Executes the Results button action. Displays the results dialog with the selected test measurement or a sample result if none is selected.
	OnDbclckTestMeasurementList (see page 179)	Executes when an item in the measurement list is double clicked. Displays the results dialog with the selected test measurement (or a sample result if none is selected.)
 	OnInitDialog (see page 179)	Generated message map functions CBodystatMFCDlg message handlers
	OnPaint (see page 180)	If you add a minimize button to your dialog, you will need the code below to draw the icon. For MFC applications using the document/view model, this is automatically done for you by the framework.
	OnQueryDragIcon (see page 181)	The system calls this function to obtain the cursor to display while the user drags the minimized window.
	OnSysCommand (see page 181)	This is OnSysCommand, a member of class CBodystatMFCDlg.
	PopulateMeasurementList (see page 182)	Populate the measurement list control with test measurement data from the API. Illustrates iteration through a BSRawData (see page 40) structure, working with data in a BSMeasurement (see page 36) structure and helper functions from the Bodystat API like BSGetDeviceModelName (see page 14)
	ShowResultsDialog (see page 184)	Displays the results dialog with the selected test measurement (or a sample result if none is selected.)

4.1.1.1.1 CBodystatMFCDlg::CBodystatMFCDlg Constructor

The main CBodystatMFCDlg constructor

C++

```
CBodystatMFCDlg(CWnd* pParent = NULL);
```

C#

```
public CBodystatMFCDlg(ref CWnd pParent);
```

Visual Basic

```
Public Sub New(ByRef pParent As CWnd = NULL)
```

Java

```
public CBodystatMFCDlg(CWnd pParent);
```

MATLAB

```
CBodystatMFCDlg
```

Parameters

Parameters	Description
CWnd* pParent = NULL	=NULL

Description

standard constructor

Body Source

```

CBodystatMFCDlg::CBodystatMFCDlg(CWnd* pParent /*=NULL*/)
: CDialogEx(CBodystatMFCDlg::IDD, pParent)
, m_strComPort(_T(""))
{
    m_hIcon = AfxGetApp()->LoadIcon(IDR_MAINFRAME);
    m_iTimeout = 7; // Initial Bluetooth timeout multiplier of 7 (equates to about 10
seconds). This is incremented upon each Bluetooth search. A maximum value of 45 is
permitted (about 1 minute).

}

```

4.1.1.1.2 CBodystatMFCDlg Data Members**4.1.1.1.2.1 CBodystatMFCDlg::m_bsRawData Data Member****C++**

```
BSRawData m_bsRawData;
```

C#

```
protected BSRawData m_bsRawData;
```

Visual Basic

```
Protected m_bsRawData As BSRawData
```

Java

```
protected BSRawData m_bsRawData;
```

MATLAB

```
m_bsRawData
```

Description

Used to hold the downloaded test measurements from the device

4.1.1.1.2.2 CBodystatMFCDlg::m_btnAuthenticate Data Member

The button authenticate

C++

```
CButton m_btnAuthenticate;
```

C#

```
protected CButton m_btnAuthenticate;
```

Visual Basic

```
Protected m_btnAuthenticate As CButton
```

Java

```
protected CButton m_btnAuthenticate;
```

MATLAB

```
m_btnAuthenticate
```

4.1.1.1.2.3 CBodystatMFCDlg::m_btnAutoSetup Data Member

The button automatic setup

C++

```
CButton m_btnAutoSetup;
```

C#

```
protected CButton m_btnAutoSetup;
```

Visual Basic

```
Protected m_btnAutoSetup As CButton
```

Java

```
protected CButton m_btnAutoSetup;
```

MATLAB

```
m_btnAutoSetup
```

4.1.1.1.2.4 CBodystatMFCDlg::m_btnConnect Data Member

The button connect

C++

```
CButton m_btnConnect;
```

C#

```
protected CButton m_btnConnect;
```

Visual Basic

```
Protected m_btnConnect As CButton
```

Java

```
protected CButton m_btnConnect;
```

MATLAB

```
m_btnConnect
```

4.1.1.1.2.5 CBodystatMFCDlg::m_btnDisconnect Data Member

The button disconnect

C++

```
CButton m_btnDisconnect;
```

C#

```
protected CButton m_btnDisconnect;
```

Visual Basic

```
Protected m_btnDisconnect As CButton
```

Java

```
protected CButton m_btnDisconnect;
```

MATLAB

```
m_btnDisconnect
```

4.1.1.1.2.6 CBodystatMFCDlg::m_btnReadCalibrationTime Data Member

The button read calibration time

C++

```
CButton m_btnReadCalibrationTime;
```

C#

```
protected CButton m_btnReadCalibrationTime;
```

Visual Basic

```
Protected m_btnReadCalibrationTime As CButton
```

Java

```
protected CButton m_btnReadCalibrationTime;
```

MATLAB

```
m_btnReadCalibrationTime
```

4.1.1.1.2.7 CBodystatMFCDlg::m_btnReadCurrentTime Data Member

The button read current time

C++

```
CButton m_btnReadCurrentTime;
```

C#

```
protected CButton m_btnReadCurrentTime;
```

Visual Basic

```
Protected m_btnReadCurrentTime As CButton
```

Java

```
protected CButton m_btnReadCurrentTime;
```

MATLAB

```
m_btnReadCurrentTime
```

4.1.1.1.2.8 CBodystatMFCDlg::m_btnReadModelVersion Data Member

The button read model version

C++

```
CButton m_btnReadModelVersion;
```

C#

```
protected CButton m_btnReadModelVersion;
```

Visual Basic

```
Protected m_btnReadModelVersion As CButton
```

Java

```
protected CButton m_btnReadModelVersion;
```

MATLAB

```
m_btnReadModelVersion
```

4.1.1.1.2.9 CBodystatMFCDlg::m_btnReadPrinterAddr Data Member

The button read printer address

C++

```
CButton m_btnReadPrinterAddr;
```

C#

```
protected CButton m_btnReadPrinterAddr;
```

Visual Basic

```
Protected m_btnReadPrinterAddr As CButton
```

Java

```
protected CButton m_btnReadPrinterAddr;
```

MATLAB

```
m_btnReadPrinterAddr
```

4.1.1.1.2.10 CBodystatMFCDlg::m_btnReadProtocolInfo Data Member

The button read protocol info

C++

```
CButton m_btnReadProtocolInfo;
```

C#

```
protected CButton m_btnReadProtocolInfo;
```

Visual Basic

```
Protected m_btnReadProtocolInfo As CButton
```

Java

```
protected CButton m_btnReadProtocolInfo;
```

MATLAB

```
m_btnReadProtocolInfo
```

4.1.1.1.2.11 CBodystatMFCDlg::m_btnReadSerialNumber Data Member

The button read serial number

C++

```
CButton m_btnReadSerialNumber;
```

C#

```
protected CButton m_btnReadSerialNumber;
```

Visual Basic

```
Protected m_btnReadSerialNumber As CButton
```

Java

```
protected CButton m_btnReadSerialNumber;
```

MATLAB

```
m_btnReadSerialNumber
```

4.1.1.1.2.12 CBodystatMFCDlg::m_btnReadStoredTestData Data Member

The button read stored test data

C++

```
CButton m_btnReadStoredTestData;
```

C#

```
protected CButton m_btnReadStoredTestData;
```

Visual Basic

```
Protected m_btnReadStoredTestData As CButton
```

Java

```
protected CButton m_btnReadStoredTestData;
```

MATLAB

```
m_btnReadStoredTestData
```

4.1.1.1.2.13 CBodystatMFCDlg::m_btnReadTestMeasurement Data Member

The button read test measurement

C++

```
CButton m_btnReadTestMeasurement;
```

C#

```
protected CButton m_btnReadTestMeasurement;
```

Visual Basic

```
Protected m_btnReadTestMeasurement As CButton
```

Java

```
protected CButton m_btnReadTestMeasurement;
```

MATLAB

```
m_btnReadTestMeasurement
```

4.1.1.1.2.14 CBodystatMFCDlg::m_btnUnauthenticate Data Member

The button unauthenticate

C++

```
CButton m_btnUnauthenticate;
```

C#

```
protected CButton m_btnUnauthenticate;
```

Visual Basic

```
Protected m_btnUnauthenticate As CButton
```

Java

```
protected CButton m_btnUnauthenticate;
```

MATLAB

```
m_btnUnauthenticate
```

4.1.1.1.2.15 CBodystatMFCDlg::m_btnWriteCurrentTime Data Member

The button write current time

C++

```
CButton m_btnWriteCurrentTime;
```

C#

```
protected CButton m_btnWriteCurrentTime;
```

Visual Basic

```
Protected m_btnWriteCurrentTime As CButton
```

Java

```
protected CButton m_btnWriteCurrentTime;
```

MATLAB

```
m_btnWriteCurrentTime
```

4.1.1.1.2.16 CBodystatMFCDlg::m_hIcon Data Member**C++**

```
HICON m_hIcon;
```

C#

```
protected HICON m_hIcon;
```

Visual Basic

```
Protected m_hIcon As HICON
```

Java

```
protected HICON m_hIcon;
```

MATLAB

m_hIcon

Description

Application icon

4.1.1.1.2.17 CBodystatMFCDlg::m_iTimeout Data Member**C++**

```
unsigned short m_iTimeout;
```

C#

```
protected unsigned short m_iTimeout;
```

Visual Basic

```
Protected m_iTimeout As unsigned short
```

Java

```
protected unsigned short m_iTimeout;
```

MATLAB

m_iTimeout

Description

Bluetooth timeout multiplier. Initial value of 7 (equates to about 10 seconds). This is incremented upon each Bluetooth search. A maximum value of 45 is permitted (about 1 minute).

4.1.1.1.2.18 CBodystatMFCDlg::m_strComPort Data Member

The string com port

C++

```
CString m_strComPort;
```

C#

```
protected CString m_strComPort;
```

Visual Basic

```
Protected m_strComPort As CString
```

Java

```
protected CString m_strComPort;
```

MATLAB

m_strComPort

4.1.1.1.2.19 CBodystatMFCDlg::m_TestMeasurementListCtrl Data Member

The test measurement list control

C++

```
CListCtrl m_TestMeasurementListCtrl;
```

C#

```
protected CListCtrl m_TestMeasurementListCtrl;
```

Visual Basic

```
Protected m_TestMeasurementListCtrl As CListCtrl
```

Java

```
protected CListCtrl m_TestMeasurementListCtrl;
```

MATLAB

```
m_TestMeasurementListCtrl
```

4.1.1.1.3 CBodystatMFCDlg Methods**4.1.1.1.3.1 CBodystatMFCDlg::DoDataExchange Method****C++**

```
virtual void DoDataExchange(CDataExchange* pDX);
```

C#

```
protected DoDataExchange(ref CDataExchange pDX);
```

Visual Basic

```
Protected Sub DoDataExchange(ByRef pDX As CDataExchange)
```

Java

```
protected DoDataExchange(CDataExchange pDX);
```

MATLAB

```
DoDataExchange
```

Description

DDX/DDV support

Bind the form controls to member variables and set basic input range limits

Body Source

```
void CBodystatMFCDlg::DoDataExchange(CDataExchange* pDX)
{
    CDialogEx::DoDataExchange(pDX);

    DDX_Control(pDX, IDC_BT_AUTOSSETUP_BUTTON, m_btnAutoSetup);
    DDX_Control(pDX, IDC_BT_AUTHENTICATE_BUTTON, m_btnAuthenticate);
    DDX_Control(pDX, IDC_BT_UNAUTHENTICATE_BUTTON, m_btnUnauthenticate);

    DDX_Text(pDX, IDC_COMPORT_EDIT, m_strComPort);
    DDX_Control(pDX, IDC_COM_CONNECT_BUTTON, m_btnConnect);
    DDX_Control(pDX, IDC_COM_READCALIBRATIONTIME_BUTTON, m_btnReadCalibrationTime);
    DDX_Control(pDX, IDC_COM_READCURRENTTIME_BUTTON, m_btnReadCurrentTime);
    DDX_Control(pDX, IDC_COM_READMODELVERSION_BUTTON, m_btnReadModelVersion);
    DDX_Control(pDX, IDC_COM_READPRINTERADDR_BUTTON, m_btnReadPrinterAddr);
    DDX_Control(pDX, IDC_COM_READPROTOCOLINFO_BUTTON, m_btnReadProtocolInfo);
    DDX_Control(pDX, IDC_COM_READSERIALNUMBER_BUTTON, m_btnReadSerialNumber);
    DDX_Control(pDX, IDC_COM_READSTOREDTESTDATA_BUTTON, m_btnReadStoredTestData);
    DDX_Control(pDX, IDC_COM_READTESTMEASUREMENT_BUTTON, m_btnReadTestMeasurement);
    DDX_Control(pDX, IDC_COM_WRITECURRENTTIME_BUTTON, m_btnWriteCurrentTime);
    DDX_Control(pDX, IDC_TESTMEASUREMENT_LIST, m_TestMeasurementListCtrl);
}
```

4.1.1.1.3.2 CBodystatMFCDlg::OnClickedBtAuthenticateButton Method

Executes the Bluetooth Authenticate button action. Illustrates use of BSAAuthenticateBTDevice (see page 8) from the Bodystat API

C++

```
afx_msg void OnClickedBtAuthenticateButton();
```

C#

```
public afx_msg void OnClickedBtAuthenticateButton();
```

Visual Basic

```
Public Function OnClickedBtAuthenticateButton() As afx_msg void
```


Java

```
public afx_msg void OnClickedBtAuthenticateButton();
```

MATLAB

```
OnClickedBtAuthenticateButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtAuthenticateButton()
{
    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    if(BSAAuthenticateBTDevice(GetSafeHwnd(), m_iTimeout))
    {
        AfxMessageBox(_T("Device authenticated\nClick Device Info button for more information"),
            MB_ICONINFORMATION);
    }
    else
        AfxMessageBox(_T("Error no Bodystat device could be found or authenticated"));
}
```

4.1.1.1.3.3 CBodystatMFCDlg::OnClickedBtAutoSetupButton Method

Executes the Bluetooth Auto Setup button action. Illustrates use of BSAutoSetupBT (see page 8) from the Bodystat API

C++

```
afx_msg void OnClickedBtAutoSetupButton();
```

C#

```
public afx_msg void OnClickedBtAutoSetupButton();
```

Visual Basic

```
Public Function OnClickedBtAutoSetupButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedBtAutoSetupButton();
```

MATLAB

```
OnClickedBtAutoSetupButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtAutoSetupButton()
{
    CString strComPort;

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Perform auto setup. User will be prompted as needed by the API
    // Suitable for supported locales only. Non-supported locales will display messages in
    English
    // For non-supported locales you should use the search and authenticate functions directly
    (silently, prompting the user yourself as required)
    if(BSAutoSetupBT(strComPort.GetBuffer(256), 256, TRUE, GetSafeHwnd()))
    {
        CString strMsg;
        strMsg.Format(_T("Device setup complete.\nCom Port: %s\nClick Device Info button for more
information"), strComPort);

        UpdateData(TRUE);
        m_strComPort = strComPort; // Assign this port to the edit box on the form
        UpdateData(FALSE);

        AfxMessageBox(strMsg, MB_ICONINFORMATION);
    }
    else
        AfxMessageBox(_T("Error no Bodystat device could be found or authenticated"));

    strComPort.ReleaseBuffer();
}
```

```
}
```

4.1.1.1.3.4 CBodystatMFCDlg::OnClickedBtDeviceInfoButton Method

Executes the Bluetooth Device Info button action. Illustrates use of BSGetBTBodystatDevice (see page 12) from the Bodystat API

C++

```
afx_msg void OnClickedBtDeviceInfoButton();
```

C#

```
public afx_msg void OnClickedBtDeviceInfoButton();
```

Visual Basic

```
Public Function OnClickedBtDeviceInfoButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedBtDeviceInfoButton();
```

MATLAB

```
OnClickedBtDeviceInfoButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtDeviceInfoButton()
{
    CString strDeviceName;
    CString strComPort;
    CString strBDA;

    CWaitCursor wc; // Lengthy operation (if not already searched), show the user we are doing something

    // Retrieve information about the Bodystat device
    BOOL bSuccess = BSGetBTBodystatDevice(strDeviceName.GetBuffer(512), 512,
    strBDA.GetBuffer(18), 18, strComPort.GetBuffer(256), 256, m_iTimeout);
    strDeviceName.ReleaseBuffer();
    strComPort.ReleaseBuffer();
    strBDA.ReleaseBuffer();

    if(bSuccess)
    {
        CString strMsg;
        strMsg.Format(_T("Found Bodystat device.\n\nDevice Name: %s\nBluetooth Address: %s\nCom Port: %s"), strDeviceName, strBDA, strComPort.IsEmpty() ? _T("<requires authentication>") : strComPort);

        UpdateData(TRUE);
        m_strComPort = strComPort; // Assign this port to the edit box on the form
        UpdateData(FALSE);

        AfxMessageBox(strMsg, MB_ICONINFORMATION);
    }
    else
    {
        AfxMessageBox(_T("Error no Bodystat device was found.\nClick Search Devices button to retry."));
    }
}
```

4.1.1.1.3.5 CBodystatMFCDlg::OnClickedBtSearchButton Method

Executes the Bluetooth Search button action. Illustrates use of BSSearchBTDevices (see page 22) from the Bodystat API

C++

```
afx_msg void OnClickedBtSearchButton();
```

C#

```
public afx_msg void OnClickedBtSearchButton();
```

Visual Basic

```
Public Function OnClickedBtSearchButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedBtSearchButton();
```

MATLAB

```
OnClickedBtSearchButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtSearchButton()
{
    int iNewDevices = 0;
    int iAuthenticatedDevices = 0;

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Begin the search...
    if(BSSearchBTDevices(iNewDevices, iAuthenticatedDevices, TRUE, GetSafeHwnd(), m_iTimeout))
    {
        CString strMsg;
        strMsg.Format(_T("Search complete.\nFound %d new and %d previously authenticated Bodystat Bluetooth devices.\nSearch waited approximately %.2fs."), iNewDevices, iAuthenticatedDevices, m_iTimeout * 1.28);
        AfxMessageBox(strMsg, MB_ICONINFORMATION);

        // If no new devices were found, lets try extending the search time for next time around
        // We can wait up to the maximum of value of 45 - about 1 minute
        // Obviously lengthy waits like this need to be in a background thread in a real world example
        if (iNewDevices == 0)
            m_iTimeout = min(m_iTimeout * 2, 45);
    }
    else
        AfxMessageBox(_T("Error searching for Bluetooth devices"));
}
```

4.1.1.1.3.6 CBodystatMFCDlg::OnClickedBtStackInfoButton Method

Executes the Bluetooth Get Stack Info button action. Illustrates use of BSGetBTStackInfo (see page 13) from the Bodystat API

C++

```
afx_msg void OnClickedBtStackInfoButton();
```

C#

```
public afx_msg void OnClickedBtStackInfoButton();
```

Visual Basic

```
Public Function OnClickedBtStackInfoButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedBtStackInfoButton();
```

MATLAB

```
OnClickedBtStackInfoButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtStackInfoButton()
{
    CString strInfo;
    if(BSGetBTStackInfo(strInfo.GetBuffer(1024), 1024))
        AfxMessageBox(strInfo, MB_ICONINFORMATION);
    else
        AfxMessageBox(_T("Error: No compatible Bluetooth Stack is available in this pc.\nEither it is switched off, or it is not recognised by this software.\nOnly Microsoft and Widcomm stacks are supported by Bodystat's SDK.));
```

```
    strInfo.ReleaseBuffer();
}
```

4.1.1.1.3.7 CBodystatMFCDlg::OnClickedBtStatusButton Method

Executes the Bluetooth Get Status button action. Illustrates use of BSIsBtAvailable (see page 16) from the Bodystat API

C++

```
afx_msg void OnClickedBtStatusButton();
```

C#

```
public afx_msg void OnClickedBtStatusButton();
```

Visual Basic

```
Public Function OnClickedBtStatusButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedBtStatusButton();
```

MATLAB

```
OnClickedBtStatusButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtStatusButton()
{
    if(BSIsBtAvailable())
        AfxMessageBox(_T("Supported Bluetooth radio module is available!"), MB_ICONINFORMATION);
    else
        AfxMessageBox(_T("No compatible Bluetooth radio module is available in this pc.\nEither it is switched off, or it is not recognised by this software.\nOnly Microsoft and Widcomm stacks are supported by Bodystat's SDK."));
}
```

4.1.1.1.3.8 CBodystatMFCDlg::OnClickedBtUnauthenticateButton Method

Executes the Bluetooth Unauthenticate button action. Illustrates use of BSUnAuthenticateBTDevices (see page 23) from the Bodystat API

C++

```
afx_msg void OnClickedBtUnauthenticateButton();
```

C#

```
public afx_msg void OnClickedBtUnauthenticateButton();
```

Visual Basic

```
Public Function OnClickedBtUnauthenticateButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedBtUnauthenticateButton();
```

MATLAB

```
OnClickedBtUnauthenticateButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedBtUnauthenticateButton()
{
    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    if(BSUnAuthenticateBTDevices(TRUE, GetSafeHwnd(), m_iTimeout))
    {
        AfxMessageBox(_T("Devices authentication removed"), MB_ICONINFORMATION);
    }
    else
        AfxMessageBox(_T("Error removing authentication"));
}
```

4.1.1.1.3.9 CBodystatMFCDlg::OnClickedComConnectButton Method

Executes the communication Connect button action. Illustrates use of BSOpenComport (see page 16) from the Bodystat API

C++

```
afx_msg void OnClickedComConnectButton();
```

C#

```
public afx_msg void OnClickedComConnectButton();
```

Visual Basic

```
Public Function OnClickedComConnectButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComConnectButton();
```

MATLAB

```
OnClickedComConnectButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedComConnectButton()
{
    UpdateData(TRUE);

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
            AfxMessageBox(_T("Connected to device successfully."), MB_ICONINFORMATION);
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
    else
        AfxMessageBox(_T("Error opening com port."));
}
```

4.1.1.1.3.10 CBodystatMFCDlg::OnClickedComReadCalibrationTimeButton Method

Executes the communication Read Calibration Time button action. Illustrates use of BSReadCalibrationTime (see page 17) from the Bodystat API

C++

```
afx_msg void OnClickedComReadCalibrationTimeButton();
```

C#

```
public afx_msg void OnClickedComReadCalibrationTimeButton();
```

Visual Basic

```
Public Function OnClickedComReadCalibrationTimeButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadCalibrationTimeButton();
```

MATLAB

```
OnClickedComReadCalibrationTimeButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedComReadCalibrationTimeButton()
```

```

{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            __time64_t time;

            // Read the calibration date from the device
            bsErr = BSReadCalibrationTime(&time);
            if(bsErr == Bodystat::NoError)
            {
                COleDateTime dtTime(time);
                CString strMsg;
                strMsg.Format(_T("Queried device successfully\nCalibration Date: %s"),
                    dtTime.Format(VAR_DATEVALUEONLY)); // Only interested in the date part of calibration date
                AfxMessageBox(strMsg, MB_ICONINFORMATION);
            }
            else
                AfxMessageBox(_T("Error querying calibration date/time from device.\nNote that legacy
                    devices do store their calibration date electronically."));
        }
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
    else
        AfxMessageBox(_T("Error opening com port."));
}

```

4.1.1.1.3.11 CBodystatMFCDlg::OnClickedComReadCurrentTimeButton Method

Executes the communication Read Current Time button action. Illustrates use of BSReadCurrentTime (see page 17) from the Bodystat API

C++

```
afx_msg void OnClickedComReadCurrentTimeButton();
```

C#

```
public afx_msg void OnClickedComReadCurrentTimeButton();
```

Visual Basic

```
Public Function OnClickedComReadCurrentTimeButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadCurrentTimeButton();
```

MATLAB

```
OnClickedComReadCurrentTimeButton
```

Body Source

```

void CBodystatMFCDlg::OnClickedComReadCurrentTimeButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {

```

```

// Connect to the device
bsErr = BSConnect();
if(bsErr == Bodystat::NoError)
{
    __time64_t time;

    // Read the current time from the internal clock of the device
    bsErr = BSReadCurrentTime(&time);
    if(bsErr == Bodystat::NoError)
    {
        COleDateTime dtTime(time); // Convert it to an MFC COleDateTime for easy display
        CString strMsg;
        strMsg.Format(_T("Queried device successfully\nCurrent Date/Time: %s"),
dtTime.Format());
        AfxMessageBox(strMsg, MB_ICONINFORMATION);
    }
    else
        AfxMessageBox(_T("Error querying current date/time from device.\nNot all models have an
inbuilt real time clock.));
    else
        AfxMessageBox(_T("Error connecting to device.));

    BSCloseComport();
}
else
    AfxMessageBox(_T("Error opening com port.));
}

```

4.1.1.1.3.12 CBodystatMFCDlg::OnClickedComReadModelVersionButton Method

Executes the communication Read Model Version button action. Illustrates use of BSReadModelVersion (see page 18) from the Bodystat API

C++

```
afx_msg void OnClickedComReadModelVersionButton();
```

C#

```
public afx_msg void OnClickedComReadModelVersionButton();
```

Visual Basic

```
Public Function OnClickedComReadModelVersionButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadModelVersionButton();
```

MATLAB

```
OnClickedComReadModelVersionButton
```

Body Source

```

void CBodystatMFCDlg::OnClickedComReadModelVersionButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            BSDeviceModel iModel;
            BYTE byMajor;
            BYTE byMinor;
            BYTE byPsoc2Version;
            BYTE byEepromVersion;

```

```

BYTE byBluetoothInfo;

// Read Model & Version information from the device
bsErr = BSReadModelVersion(&iModel, &byMajor, &byMinor, &byPsoc2Version,
&byEepromVersion, &byBluetoothInfo);
if(bsErr == Bodystat::NoError)
{
    CString strMsg;
    strMsg.Format(_T("Queried device successfully\n\nDevice Type: %d\nFirmware Version:
%d.%d.%d.%d"), iModel, byMajor, byMinor, byPsoc2Version, byEepromVersion);
    AfxMessageBox(strMsg, MB_ICONINFORMATION);
}
else
    AfxMessageBox(_T("Error querying model information from device."));
}
else
    AfxMessageBox(_T("Error connecting to device."));

BSCloseComport();
}
else
    AfxMessageBox(_T("Error opening com port."));
}

```

4.1.1.1.3.13 CBodystatMFCDlg::OnClickedComReadPrinterAddrButton Method

Executes the communication Read Printer Address button action. Illustrates use of BSReadPrinterAddress (see page 19) from the Bodystat API

C++

```
afx_msg void OnClickedComReadPrinterAddrButton();
```

C#

```
public afx_msg void OnClickedComReadPrinterAddrButton();
```

Visual Basic

```
Public Function OnClickedComReadPrinterAddrButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadPrinterAddrButton();
```

MATLAB

```
OnClickedComReadPrinterAddrButton
```

Body Source

```

void CBodystatMFCDlg::OnClickedComReadPrinterAddrButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            CString strPrinterAddr;

            // Read the address of the printer the device is paired to for direct printing
            bsErr = BSReadPrinterAddress(strPrinterAddr.GetBuffer(16), 16);
            strPrinterAddr.ReleaseBuffer();
            if(bsErr == Bodystat::NoError)
            {
                CString strMsg;
                strMsg.Format(_T("Queried device successfully\nPrinter Address: %s"), strPrinterAddr);
                AfxMessageBox(strMsg, MB_ICONINFORMATION);
            }
        }
    }
}

```



```

    }
    else
        AfxMessageBox(_T("Error querying printer address from device.\nNot all models support
direct printing."));
    }
    else
        AfxMessageBox(_T("Error connecting to device."));

    BSCloseComport();
}
else
    AfxMessageBox(_T("Error opening com port."));
}

```

4.1.1.1.3.14 CBodystatMFCDlg::OnClickedComReadProtocolInfoButton Method

Executes the communication ReadProtocolInfo button action. Illustrates use of BSReadProtocolInfo (see page 19) from the Bodystat API

C++

```
afx_msg void OnClickedComReadProtocolInfoButton();
```

C#

```
public afx_msg void OnClickedComReadProtocolInfoButton();
```

Visual Basic

```
Public Function OnClickedComReadProtocolInfoButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadProtocolInfoButton();
```

MATLAB

```
OnClickedComReadProtocolInfoButton
```

Body Source

```

void CBodystatMFCDlg::OnClickedComReadProtocolInfoButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            BYTE byProtocolVersion;
            BYTE byDataVersion;
            BYTE byAuxInfo;

            // Read Protocol Information from the device
            bsErr = BSReadProtocolInfo(&byProtocolVersion, &byDataVersion, &byAuxInfo);
            if(bsErr == Bodystat::NoError)
            {
                CString strMsg;
                strMsg.Format(_T("Queried device successfully\nProtocol Version: %d\nData Version:
%d\nAux Info: %d"), byProtocolVersion, byDataVersion, byAuxInfo);
                AfxMessageBox(strMsg, MB_ICONINFORMATION);
            }
            else
                AfxMessageBox(_T("Error querying protocol information from device."));
        }
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
}

```

```

    }
    else
        AfxMessageBox(_T("Error opening com port."));
}

```

4.1.1.1.3.15 CBodystatMFCDlg::OnClickedComReadSerialNumberButton Method

Executes the communication Read Serial Number button action. Illustrates use of BSReadSerialNumber (see page 20) from the Bodystat API

C++

```
afx_msg void OnClickedComReadSerialNumberButton();
```

C#

```
public afx_msg void OnClickedComReadSerialNumberButton();
```

Visual Basic

```
Public Function OnClickedComReadSerialNumberButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadSerialNumberButton();
```

MATLAB

```
OnClickedComReadSerialNumberButton
```

Body Source

```

void CBodystatMFCDlg::OnClickedComReadSerialNumberButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            ULONG ulSerialNumber;

            // Read the unique serial number from the device
            bsErr = BSReadSerialNumber(&ulSerialNumber);
            if(bsErr == Bodystat::NoError)
            {
                CString strMsg;
                strMsg.Format(_T("Queried device successfully\nSerial Number: %d"), ulSerialNumber);
                AfxMessageBox(strMsg, MB_ICONINFORMATION);
            }
            else
                AfxMessageBox(_T("Error querying serial number from device."));
        }
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
    else
        AfxMessageBox(_T("Error opening com port."));
}

```

4.1.1.1.3.16 CBodystatMFCDlg::OnClickedComReadStoredTestDataButton Method

Executes the communication Read Stored Test Data button action. Illustrates use of BSReadStoredTestData (see page 21) from the Bodystat API to download data

C++

```
afx_msg void OnClickedComReadStoredTestDataButton();
```

C#

```
public afx_msg void OnClickedComReadStoredTestDataButton();
```

Visual Basic

```
Public Function OnClickedComReadStoredTestDataButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadStoredTestDataButton();
```

MATLAB

```
OnClickedComReadStoredTestDataButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedComReadStoredTestDataButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            // Set the array size in the raw data strucutre (SDK expects a fixed value here
            currently)
            m_bsRawData.iRecordArraySize = Bodystat::BS_RAWDATA_ARRAYSIZE;

            // Download test measurement data from the device
            bsErr = BSReadStoredTestData(&m_bsRawData);
            if(bsErr == Bodystat::NoError)
            {
                CString strMsg;
                strMsg.Format(_T("Download complete\nRetrieved %d records\nFirst Test Number: %d"),
                m_bsRawData.iTotalNumRecs, m_bsRawData.ulFirstTestNum);

                // Parse the data and place it in the measurement list control
                PopulateMeasurementList(&m_bsRawData);

                AfxMessageBox(strMsg, MB_ICONINFORMATION);
            }
            else
                AfxMessageBox(_T("Error downloading stored test data from the device."));
        }
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
    else
        AfxMessageBox(_T("Error opening com port."));
}
```

4.1.1.1.3.17 CBodystatMFCDlg::OnClickedComReadTestMeasurementButton Method

Executes the communication Read Test Measurement button action. Performs an immediate test measurement on the device using the BSReadTestBodystat (see page 21) function from the Bodystat API

C++

```
afx_msg void OnClickedComReadTestMeasurementButton();
```

C#

```
public afx_msg void OnClickedComReadTestMeasurementButton();
```

Visual Basic

```
Public Function OnClickedComReadTestMeasurementButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComReadTestMeasurementButton();
```

MATLAB

```
OnClickedComReadTestMeasurementButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedComReadTestMeasurementButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            UINT iz5, iz50, iz100, iz200, ir50;
            float fx50, fPA50;

            // Perform the test measurement on the device and get the raw electrical bio-impedance
            results
            bsErr = BSReadTestBodystat(&iz5, &iz50, &iz100, &iz200, &ir50, &fx50, &fPA50);
            if(bsErr == Bodystat::NoError)
            {
                CString strMsg;
                strMsg.Format(_T("Measurement complete\nImpedance 5kHz: %d\nImpedance 50kHz: %d\nImpedance 100kHz: %d\nImpedance 200kHz: %d\nResistance 50kHz: %d\nReactance 50kHz: %.1f\nPhase Angle 50kHz: %.1f"),
                    iz5, iz50, iz100, iz200, ir50, fx50, fPA50);
                AfxMessageBox(strMsg, MB_ICONINFORMATION);
            }
            else
                AfxMessageBox(_T("Error performing test measurement on the device."));
        }
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
    else
        AfxMessageBox(_T("Error opening com port."));
}
```

4.1.1.1.3.18 CBodystatMFCDlg::OnClickedComWriteCurrentTimeButton Method

Executes the communication Write Current Time button action. Illustrates use of BSWriteCurrentTime (see page 24) from the Bodystat API

C++

```
afx_msg void OnClickedComWriteCurrentTimeButton();
```

C#

```
public afx_msg void OnClickedComWriteCurrentTimeButton();
```

Visual Basic

```
Public Function OnClickedComWriteCurrentTimeButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedComWriteCurrentTimeButton();
```

MATLAB

```
OnClickedComWriteCurrentTimeButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedComWriteCurrentTimeButton()
{
    UpdateData(TRUE); // Grab the com port value from the edit box of the form

    CWaitCursor wc; // Lengthy operation, show the user we are doing something

    // Open a com port to the device
    BSError bsErr = BSOpenComport(m_strComPort.GetBuffer(10), 10);
    if(bsErr == Bodystat::NoError)
    {
        // Connect to the device
        bsErr = BSConnect();
        if(bsErr == Bodystat::NoError)
        {
            // Write the current local PC time to the internal clock of the device
            bsErr = BSWriteCurrentTime();
            if(bsErr == Bodystat::NoError)
            {
                AfxMessageBox(_T("Device updated successfully\nCurrent date/time updated."),
                MB_ICONINFORMATION);
            }
            else
                AfxMessageBox(_T("Error writing current date/time to device.\nNot all models have an
                inbuilt real time clock."));
        }
        else
            AfxMessageBox(_T("Error connecting to device."));

        BSCloseComport();
    }
    else
        AfxMessageBox(_T("Error opening com port."));
}
```

4.1.1.1.3.19 CBodystatMFCDlg::OnClickedResultsButton Method

Executes the Results button action. Displays the results dialog with the selected test measurement or a sample result if none is selected.

C++

```
afx_msg void OnClickedResultsButton();
```

C#

```
public afx_msg void OnClickedResultsButton();
```

Visual Basic

```
Public Function OnClickedResultsButton() As afx_msg void
```

Java

```
public afx_msg void OnClickedResultsButton();
```

MATLAB

```
OnClickedResultsButton
```

Body Source

```
void CBodystatMFCDlg::OnClickedResultsButton()
{
    ShowResultsDialog(); // Display the selected result
}
```

4.1.1.1.3.20 CBodystatMFCDlg::OnDbclkTestMeasurementList Method

Executes when an item in the measurement list is double clicked. Displays the results dialog with the selected test measurement (or a sample result if none is selected.)

C++

```
afx_msg void OnDbclkTestMeasurementList(NMHDR * pNMHDR, LRESULT * pResult);
```

C#

```
public afx_msg void OnDbclkTestMeasurementList(ref NMHDR pNMHDR, ref LRESULT pResult);
```

Visual Basic

```
Public Function OnDbclkTestMeasurementList(ByRef pNMHDR As NMHDR, ByRef pResult As LRESULT) As afx_msg void
```

Java

```
public afx_msg void OnDbclkTestMeasurementList(NMHDR pNMHDR, LRESULT pResult);
```

MATLAB

```
OnDbclkTestMeasurementList
```

Body Source

```
void CBodystatMFCDlg::OnDbclkTestMeasurementList(NMHDR *pNMHDR, LRESULT *pResult)
{
    ShowResultsDialog();    // Display the selected result
    *pResult = 0;
}
```

4.1.1.1.3.21 CBodystatMFCDlg::OnInitDialog Method

C++

```
virtual BOOL OnInitDialog();
```

C#

```
protected BOOL OnInitDialog();
```

Visual Basic

```
Protected Function OnInitDialog() As BOOL
```

Java

```
protected BOOL OnInitDialog();
```

MATLAB

```
OnInitDialog
```

Description

Generated message map functions

CBodystatMFCDlg (see page 156) message handlers

Body Source

```
BOOL CBodystatMFCDlg::OnInitDialog()
{
    CDialogEx::OnInitDialog();

    // Add "About..." menu item to system menu.

    // IDM_ABOUTBOX must be in the system command range.
    ASSERT((IDM_ABOUTBOX & 0xFFF0) == IDM_ABOUTBOX);
    ASSERT(IDM_ABOUTBOX < 0xF000);

    CMenu* pSysMenu = GetSystemMenu(FALSE);
    if (pSysMenu != NULL)
    {
        BOOL bNameValid;
```

```

CString strMenuTitle;

// Add About menu
bNameValid = strMenuTitle.LoadString(IDS_ABOUTBOX);
ASSERT(bNameValid);
if (!strMenuTitle.IsEmpty())
{
    pSysMenu->AppendMenu(MF_SEPARATOR);
    pSysMenu->AppendMenu(MF_STRING, IDM_ABOUTBOX, strMenuTitle);    // Add about menu item to
the main file menu
}

// Set the icon for this dialog. The framework does this automatically
// when the application's main window is not a dialog
SetIcon(m_hIcon, TRUE);    // Set big icon
SetIcon(m_hIcon, FALSE);    // Set small icon

// Add Vista Shield icon to some buttons which may require elevated user privileges
Button_SetElevationRequiredState(m_btnAutoSetup, TRUE);
Button_SetElevationRequiredState(m_btnAuthenticate, TRUE);
Button_SetElevationRequiredState(m_btnUnauthenticate, TRUE);

// Add columns to the test measurement list ready for download of data
int j=0;
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Test No"), 0, 50);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Test Date/Time"), 0, 130);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Model"), 0, 140);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Serial No"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Gender"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Age"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Height"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Weight"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Activity"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Waist"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Hip"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Imp 5 kHz"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Imp 50 kHz"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Imp 100 kHz"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Imp 200 kHz"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Res 50 kHz"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Reac 50 kHz"), 0, 75);
m_TestMeasurementListCtrl.InsertColumn(j++, _T("Phase 50 kHz"), 0, 80);
m_TestMeasurementListCtrl.SetExtendedStyle( m_TestMeasurementListCtrl.GetExtendedStyle() |
LVS_EX_FULLROWSELECT );

return TRUE;    // return TRUE unless you set the focus to a control
}

```

4.1.1.1.3.22 CBodystatMFCDlg::OnPaint Method

C++

```
afx_msg void OnPaint();
```

C#

```
protected afx_msg void OnPaint();
```

Visual Basic

```
Protected Function OnPaint() As afx_msg void
```

Java

```
protected afx_msg void OnPaint();
```

MATLAB

```
OnPaint
```

Description

If you add a minimize button to your dialog, you will need the code below to draw the icon. For MFC applications using the document/view model, this is automatically done for you by the framework.

Body Source

```

void CBodystatMFCDlg::OnPaint()
{
    if (IsIconic())
    {
        CPaintDC dc(this); // device context for painting

        SendMessage(WM_ICONERASEBKGND, reinterpret_cast<WPARAM>(dc.GetSafeHdc()), 0);

        // Center icon in client rectangle
        int cxIcon = GetSystemMetrics(SM_CXICON);
        int cyIcon = GetSystemMetrics(SM_CYICON);
        CRect rect;
        GetClientRect(&rect);
        int x = (rect.Width() - cxIcon + 1) / 2;
        int y = (rect.Height() - cyIcon + 1) / 2;

        // Draw the icon
        dc.DrawIcon(x, y, m_hIcon);
    }
    else
    {
        CDialogEx::OnPaint();
    }
}

```

4.1.1.1.3.23 CBodystatMFCDlg::OnQueryDragIcon Method**C++**

```
afx_msg HCURSOR OnQueryDragIcon();
```

C#

```
protected afx_msg HCURSOR OnQueryDragIcon();
```

Visual Basic

```
Protected Function OnQueryDragIcon() As afx_msg HCURSOR
```

Java

```
protected afx_msg HCURSOR OnQueryDragIcon();
```

MATLAB

```
OnQueryDragIcon
```

Description

The system calls this function to obtain the cursor to display while the user drags the minimized window.

Body Source

```

HCURSOR CBodystatMFCDlg::OnQueryDragIcon()
{
    return static_cast<HCURSOR>(m_hIcon);
}

```

4.1.1.1.3.24 CBodystatMFCDlg::OnSysCommand Method**C++**

```
afx_msg void OnSysCommand(UINT nID, LPARAM lParam);
```

C#

```
protected afx_msg void OnSysCommand(UINT nID, LPARAM lParam);
```

Visual Basic

```
Protected Function OnSysCommand(nID As UINT, lParam As LPARAM) As afx_msg void
```

Java

```
protected afx_msg void OnSysCommand(UINT nID, LPARAM lParam);
```


MATLAB

OnSysCommand

Description

This is OnSysCommand, a member of class CBodystatMFCDlg.

Body Source

```
void CBodystatMFCDlg::OnSysCommand(UINT nID, LPARAM lParam)
{
    if ((nID & 0xFFF0) == IDM_ABOUTBOX)
    {
        CAboutDlg dlgAbout;
        dlgAbout.DoModal();
    }
    else
    {
        CDialogEx::OnSysCommand(nID, lParam);
    }
}
```

4.1.1.1.3.25 CBodystatMFCDlg::PopulateMeasurementList Method

Populate the measurement list control with test measurement data from the API. Illustrates iteration through a BSRawData (see page 40) structure, working with data in a BSMeasurement (see page 36) structure and helper functions from the Bodystat API like BSGetDeviceModelName (see page 14)

C++

```
void PopulateMeasurementList(BSRawData * pRawData);
```

C#

```
public PopulateMeasurementList(ref BSRawData pRawData);
```

Visual Basic

```
Public Sub PopulateMeasurementList(ByRef pRawData As BSRawData)
```

Java

```
public PopulateMeasurementList(BSRawData pRawData);
```

MATLAB

PopulateMeasurementList

Parameters

Parameters	Description
BSRawData * pRawData	[in] RawData to be parsed.

Body Source

```
void CBodystatMFCDlg::PopulateMeasurementList(BSRawData *pRawData)
{
    // Clear any previous items in the list
    m_TestMeasurementListCtrl.DeleteAllItems();

    // Check we have some raw data to work with
    if(pRawData)
    {
        CString strItem;
        int nItem;

        // Process each record in the structure
        for(int i=0; i<pRawData->iTotalNumRecs; i++)
        {
            int j=1; // Column number

            // Test no
            strItem.Format(_T("%d"), i + pRawData->ulFirstTestNum);
            nItem = m_TestMeasurementListCtrl.InsertItem(i, strItem);
        }
    }
}
```

```

m_TestMeasurementListCtrl.SetItemData(nItem, (DWORD_PTR)&(pRawData->record[i]));

// Test date/time
COleDateTime dtTestDate(pRawData->record[i].tTestDate);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, dtTestDate.Format());

// Device model
BSGetDeviceModelName(pRawData->record[i].iDeviceModel, strItem.GetBuffer(32), 32);
strItem.ReleaseBuffer();
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Device serial number
strItem.Format(_T("%d"), pRawData->record[i].ulDeviceSerialNumber);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Gender
m_TestMeasurementListCtrl.SetItemText(nItem, j++, pRawData->record[i].iGender ==
Bodystat::BSMale ? _T("Male") : _T("Female"));

// Age
strItem.Format(_T("%d"), pRawData->record[i].iAge);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Height
strItem.Format(_T("%d cm"), pRawData->record[i].iHeight);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Weight
strItem.Format(_T("%.1f kg"), pRawData->record[i].fWeight);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Activity
switch(pRawData->record[i].iActivity)
{
case 1: strItem = _T("Very Low"); break;
case 2: strItem = _T("Low/Medium"); break;
case 3: strItem = _T("Medium"); break;
case 4: strItem = _T("Medium/High"); break;
case 5: strItem = _T("Very High"); break;
default: strItem = _T("Unknown"); break;
}
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Waist
strItem.Format(_T("%d cm"), pRawData->record[i].iWaist);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Hip
strItem.Format(_T("%d cm"), pRawData->record[i].iHip);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Z5 (Impedance at 5 kHz)
strItem.Format(_T("%d ohm"), pRawData->record[i].iZ_5kHz);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Z50 (Impedance at 50 kHz)
strItem.Format(_T("%d ohm"), pRawData->record[i].iZ_50kHz);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Z100 (Impedance at 100 kHz)
strItem.Format(_T("%d ohm"), pRawData->record[i].iZ_100kHz);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// Z200 (Impedance at 200 kHz)
strItem.Format(_T("%d ohm"), pRawData->record[i].iZ_200kHz);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// R50 (Resistance at 50 kHz)
strItem.Format(_T("%d ohm"), pRawData->record[i].iR_50kHz);
m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

// X50 (Reactance at 50 kHz)

```

```

    strItem.Format(_T("%.1f ohm"), pRawData->record[i].fX_50kHz);
    m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);

    // PA50 (Phase Angle at 50 kHz)
    strItem.Format(_T("%.1f deg"), pRawData->record[i].fPA_50kHz);
    m_TestMeasurementListCtrl.SetItemText(nItem, j++, strItem);
}
}
}

```

4.1.1.1.3.26 CBodystatMFCDlg::ShowResultsDialog Method

Displays the results dialog with the selected test measurement (or a sample result if none is selected.)

C++

```
void ShowResultsDialog();
```

C#

```
public ShowResultsDialog();
```

Visual Basic

```
Public Sub ShowResultsDialog()
```

Java

```
public ShowResultsDialog();
```

MATLAB

```
ShowResultsDialog
```

Body Source

```

void CBodystatMFCDlg::ShowResultsDialog()
{
    CBodystatResultsDlg resultsDlg; // Initialise the results dialog

    // Populate the results dialog with the details of the selected measurement
    // If no item is selected, we simply show the results dialog with its default values
    int nItem = m_TestMeasurementListCtrl.GetNextItem( -1, LVNI_SELECTED );
    if(nItem >= 0)
    {
        BSMeasurement *pBSM = (BSMeasurement *)m_TestMeasurementListCtrl.GetItemData(nItem);
        if(pBSM)
        {
            resultsDlg.m_iModel = pBSM->iDeviceModel;
            resultsDlg.m_dtTestDate = pBSM->tTestDate;
            resultsDlg.m_iGender = pBSM->iGender == BSMale ? 0 : 1;
            resultsDlg.m_iAge = pBSM->iAge;
            resultsDlg.m_iHeight = pBSM->iHeight;
            resultsDlg.m_fWeight = pBSM->fWeight;
            resultsDlg.m_iActivityLevel = pBSM->iActivity;
            resultsDlg.m_iWaist = pBSM->iWaist;
            resultsDlg.m_iHip = pBSM->iHip;

            resultsDlg.m_bModeBodyComposition = TRUE;
            resultsDlg.m_bModeHydration = TRUE;

            resultsDlg.m_iImpedance5K = pBSM->iZ_5kHz;
            resultsDlg.m_iImpedance50K = pBSM->iZ_50kHz;
            resultsDlg.m_iImpedance100K = pBSM->iZ_100kHz;
            resultsDlg.m_iImpedance200K = pBSM->iZ_200kHz;
            resultsDlg.m_iResistance50K = pBSM->iR_50kHz;
            resultsDlg.m_fReactance50K = pBSM->fX_50kHz;
            resultsDlg.m_fPhaseAngle50K = pBSM->fPA_50kHz;
        }
    }
    else
        AfxMessageBox(_T("No measurement selected, an example measurement will be shown instead."), MB_ICONINFORMATION);
}


```

```
// Show the results dialog
resultsDlg.DoModal();
}
```

4.1.1.2 CBodystatResultsDlg Class

CBodystatResultsDlg. Dialog for viewing the Bodystat test measurement results and normals.

Class Hierarchy

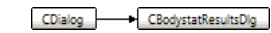

C++

```
class CBodystatResultsDlg : public CDialog;
```

C#

```
public class CBodystatResultsDlg : CDialog;
```

Visual Basic

```
Public Class CBodystatResultsDlg
Inherits CDialog
```

Java

```
public class CBodystatResultsDlg CDialog;
```

MATLAB

```
CBodystatResultsDlg
```

File









```
BodystatResultsDlg.h
```

Methods








	Name	Description
	~CBodystatResultsDlg (see page 186)	This is ~CBodystatResultsDlg, a member of class CBodystatResultsDlg.
	CBodystatResultsDlg (see page 186)	CBodystatResultsDlg Constructor.

CBodystatResultsDlg Data Members

	Name	Description
	m_bModeBodyComposition (see page 187)	This is m_bModeBodyComposition, a member of class CBodystatResultsDlg.
	m_bModeHydration (see page 188)	This is m_bModeHydration, a member of class CBodystatResultsDlg.
	m_dtTestDate (see page 188)	This is m_dtTestDate, a member of class CBodystatResultsDlg.
	m_fPhaseAngle50K (see page 188)	This is m_fPhaseAngle50K, a member of class CBodystatResultsDlg.
	m_fReactance50K (see page 189)	This is m_fReactance50K, a member of class CBodystatResultsDlg.
	m_fWeight (see page 189)	This is m_fWeight, a member of class CBodystatResultsDlg.
	m_iActivityLevel (see page 189)	This is m_iActivityLevel, a member of class CBodystatResultsDlg.
	m_iAge (see page 190)	This is m_iAge, a member of class CBodystatResultsDlg.
	m_iGender (see page 190)	This is m_iGender, a member of class CBodystatResultsDlg.
	m_iHeight (see page 190)	This is m_iHeight, a member of class CBodystatResultsDlg.
	m_iHip (see page 191)	This is m_iHip, a member of class CBodystatResultsDlg.
	m_ilmpedance100K (see page 191)	This is m_ilmpedance100K, a member of class CBodystatResultsDlg.

	m_iImpedance200K (see page 191)	This is m_iImpedance200K, a member of class CBodystatResultsDlg.
	m_iImpedance50K (see page 192)	This is m_iImpedance50K, a member of class CBodystatResultsDlg.
	m_iImpedance5K (see page 192)	This is m_iImpedance5K, a member of class CBodystatResultsDlg.
	m_iModel (see page 192)	This is m_iModel, a member of class CBodystatResultsDlg.
	m_iResistance50K (see page 193)	This is m_iResistance50K, a member of class CBodystatResultsDlg.
	m_iWaist (see page 193)	This is m_iWaist, a member of class CBodystatResultsDlg.
	m_ModelComboCtrl (see page 193)	This is m_ModelComboCtrl, a member of class CBodystatResultsDlg.
	m_ResultsListCtrl (see page 194)	This is m_ResultsListCtrl, a member of class CBodystatResultsDlg.

CBodystatResultsDlg Methods

	Name	Description
	CopyToClipboard (see page 194)	Copies measurement parameters, results and normals to the clipboard. Useful if you need to make note of a particular test set for diagnostics.
	DoDataExchange (see page 196)	DDX/DDV support Bind the form controls to member variables and set basic input range limits
	GetSystemErrorMessage (see page 197)	Gets a system error message.
	OnClickedCalculateResultsButton (see page 198)	Executes the Calculate Results button action. Uses the Bodystat API to calculate the results and normals for this test measurement.
	OnClickedCopyToClipboardButton (see page 200)	Executes the Copy To Clipboard button action.
	OnInitDialog (see page 200)	CBodystatResultsDlg message handlers
	PopulateListCtrl (see page 201)	Populate the list control with results and normals.

4.1.1.2.1 CBodystatResultsDlg::~CBodystatResultsDlg Destructor

C++

```
virtual ~CBodystatResultsDlg();
```

C#

```
public ~CBodystatResultsDlg();
```

Visual Basic

```
Public Sub Finalize()
```

Java

```
public finalize();
```

MATLAB

```
~CBodystatResultsDlg
```

Description

This is ~CBodystatResultsDlg, a member of class CBodystatResultsDlg.

Body Source

```
CBodystatResultsDlg::~CBodystatResultsDlg()
{
}
```

4.1.1.2.2 CBodystatResultsDlg::CBodystatResultsDlg Constructor

CBodystatResultsDlg Constructor.

C++

```
CBodystatResultsDlg(CWnd* pParent = NULL);
```

C#

```
public CBodystatResultsDlg(ref CWnd pParent);
```

Visual Basic

```
Public Sub New(ByRef pParent As CWnd = NULL)
```

Java

```
public CBodystatResultsDlg(CWnd pParent);
```

MATLAB

```
CBodystatResultsDlg
```

Parameters

Parameters	Description
CWnd* pParent = NULL	=NULL

Description

standard constructor

Body Source

```
CBodystatResultsDlg::CBodystatResultsDlg(CWnd* pParent /*=NULL*/)
: CDialog(CBodystatResultsDlg::IDD, pParent)
{
    // Setup the form with some sensible default values
    m_iModel = (int)BSQuadScan_BT;
    m_iGender = 0;
    m_iAge = 30;
    m_iHeight = 165;
    m_fWeight = 85.0f;
    m_iActivityLevel = 3;
    m_iWaist = 85;
    m_iHip = 90;

    m_bModeBodyComposition = TRUE;
    m_bModeHydration = TRUE;

    m_iImpedance5K = 500;
    m_iImpedance50K = 500;
    m_iImpedance100K = 500;
    m_iImpedance200K = 500;
    m_iResistance50K = 500;
    m_fReactance50K = 52.0;
    m_fPhaseAngle50K = (float)6.1;
}
```

4.1.1.2.3 CBodystatResultsDlg Data Members

4.1.1.2.3.1 CBodystatResultsDlg::m_bModeBodyComposition Data Member

C++

```
BOOL m_bModeBodyComposition;
```

C#

```
public BOOL m_bModeBodyComposition;
```

Visual Basic

```
Public m_bModeBodyComposition As BOOL
```

Java

```
public BOOL m_bModeBodyComposition;
```

MATLAB

m_bModeBodyComposition

Description

This is m_bModeBodyComposition, a member of class CBodystatResultsDlg.

4.1.1.2.3.2 CBodystatResultsDlg::m_bModeHydration Data Member**C++**

```
BOOL m_bModeHydration;
```

C#

```
public BOOL m_bModeHydration;
```

Visual Basic

```
Public m_bModeHydration As BOOL
```

Java

```
public BOOL m_bModeHydration;
```

MATLAB

m_bModeHydration

Description

This is m_bModeHydration, a member of class CBodystatResultsDlg.

4.1.1.2.3.3 CBodystatResultsDlg::m_dtTestDate Data Member**C++**

```
ColeDateTime m_dtTestDate;
```

C#

```
public ColeDateTime m_dtTestDate;
```

Visual Basic

```
Public m_dtTestDate As ColeDateTime
```

Java

```
public ColeDateTime m_dtTestDate;
```

MATLAB

m_dtTestDate

Description

This is m_dtTestDate, a member of class CBodystatResultsDlg.

4.1.1.2.3.4 CBodystatResultsDlg::m_fPhaseAngle50K Data Member**C++**

```
float m_fPhaseAngle50K;
```

C#

```
public float m_fPhaseAngle50K;
```

Visual Basic

```
Public m_fPhaseAngle50K As float
```

Java

```
public float m_fPhaseAngle50K;
```

MATLAB

```
m_fPhaseAngle50K
```

Description

This is m_fPhaseAngle50K, a member of class CBodystatResultsDlg.

4.1.1.2.3.5 CBodystatResultsDlg::m_fReactance50K Data Member**C++**

```
float m_fReactance50K;
```

C#

```
public float m_fReactance50K;
```

Visual Basic

```
Public m_fReactance50K As float
```

Java

```
public float m_fReactance50K;
```

MATLAB

```
m_fReactance50K
```

Description

This is m_fReactance50K, a member of class CBodystatResultsDlg.

4.1.1.2.3.6 CBodystatResultsDlg::m_fWeight Data Member**C++**

```
float m_fWeight;
```

C#

```
public float m_fWeight;
```

Visual Basic

```
Public m_fWeight As float
```

Java

```
public float m_fWeight;
```

MATLAB

```
m_fWeight
```

Description

This is m_fWeight, a member of class CBodystatResultsDlg.

4.1.1.2.3.7 CBodystatResultsDlg::m_iActivityLevel Data Member**C++**

```
int m_iActivityLevel;
```

C#

```
public int m_iActivityLevel;
```

Visual Basic

```
Public m_iActivityLevel As Integer
```

Java

```
public int m_iActivityLevel;
```


MATLAB

```
m_iActivityLevel
```

Description

This is m_iActivityLevel, a member of class CBodystatResultsDlg.

4.1.1.2.3.8 CBodystatResultsDlg::m_iAge Data Member**C++**

```
int m_iAge;
```

C#

```
public int m_iAge;
```

Visual Basic

```
Public m_iAge As Integer
```

Java

```
public int m_iAge;
```

MATLAB

```
m_iAge
```

Description

This is m_iAge, a member of class CBodystatResultsDlg.

4.1.1.2.3.9 CBodystatResultsDlg::m_iGender Data Member**C++**

```
int m_iGender;
```

C#

```
public int m_iGender;
```

Visual Basic

```
Public m_iGender As Integer
```

Java

```
public int m_iGender;
```

MATLAB

```
m_iGender
```

Description

This is m_iGender, a member of class CBodystatResultsDlg.

4.1.1.2.3.10 CBodystatResultsDlg::m_iHeight Data Member**C++**

```
int m_iHeight;
```

C#

```
public int m_iHeight;
```

Visual Basic

```
Public m_iHeight As Integer
```

Java

```
public int m_iHeight;
```

MATLAB

m_iHeight

Description

This is m_iHeight, a member of class CBodystatResultsDlg.

4.1.1.2.3.11 CBodystatResultsDlg::m_iHip Data Member**C++**

```
int m_iHip;
```

C#

```
public int m_iHip;
```

Visual Basic

```
Public m_iHip As Integer
```

Java

```
public int m_iHip;
```

MATLAB

m_iHip

Description

This is m_iHip, a member of class CBodystatResultsDlg.

4.1.1.2.3.12 CBodystatResultsDlg::m_iImpedance100K Data Member**C++**

```
int m_iImpedance100K;
```

C#

```
public int m_iImpedance100K;
```

Visual Basic

```
Public m_iImpedance100K As Integer
```

Java

```
public int m_iImpedance100K;
```

MATLAB

m_iImpedance100K

Description

This is m_iImpedance100K, a member of class CBodystatResultsDlg.

4.1.1.2.3.13 CBodystatResultsDlg::m_iImpedance200K Data Member**C++**

```
int m_iImpedance200K;
```

C#

```
public int m_iImpedance200K;
```

Visual Basic

```
Public m_iImpedance200K As Integer
```

Java

```
public int m_iImpedance200K;
```

MATLAB

```
m_iImpedance200K
```

Description

This is m_iImpedance200K, a member of class CBodystatResultsDlg.

4.1.1.2.3.14 CBodystatResultsDlg::m_iImpedance50K Data Member**C++**

```
int m_iImpedance50K;
```

C#

```
public int m_iImpedance50K;
```

Visual Basic

```
Public m_iImpedance50K As Integer
```

Java

```
public int m_iImpedance50K;
```

MATLAB

```
m_iImpedance50K
```

Description

This is m_iImpedance50K, a member of class CBodystatResultsDlg.

4.1.1.2.3.15 CBodystatResultsDlg::m_iImpedance5K Data Member**C++**

```
int m_iImpedance5K;
```

C#

```
public int m_iImpedance5K;
```

Visual Basic

```
Public m_iImpedance5K As Integer
```

Java

```
public int m_iImpedance5K;
```

MATLAB

```
m_iImpedance5K
```

Description

This is m_iImpedance5K, a member of class CBodystatResultsDlg.

4.1.1.2.3.16 CBodystatResultsDlg::m_iModel Data Member**C++**

```
int m_iModel;
```

C#

```
public int m_iModel;
```

Visual Basic

```
Public m_iModel As Integer
```

Java

```
public int m_iModel;
```

MATLAB

```
m_iModel
```

Description

This is m_iModel, a member of class CBodystatResultsDlg.

4.1.1.2.3.17 CBodystatResultsDlg::m_iResistance50K Data Member**C++**

```
int m_iResistance50K;
```

C#

```
public int m_iResistance50K;
```

Visual Basic

```
Public m_iResistance50K As Integer
```

Java

```
public int m_iResistance50K;
```

MATLAB

```
m_iResistance50K
```

Description

This is m_iResistance50K, a member of class CBodystatResultsDlg.

4.1.1.2.3.18 CBodystatResultsDlg::m_iWaist Data Member**C++**

```
int m_iWaist;
```

C#

```
public int m_iWaist;
```

Visual Basic

```
Public m_iWaist As Integer
```

Java

```
public int m_iWaist;
```

MATLAB

```
m_iWaist
```

Description

This is m_iWaist, a member of class CBodystatResultsDlg.

4.1.1.2.3.19 CBodystatResultsDlg::m_ModelComboCtrl Data Member**C++**

```
CComboBox m_ModelComboCtrl;
```

C#

```
protected CComboBox m_ModelComboCtrl;
```

Visual Basic

```
Protected m_ModelComboCtrl As CComboBox
```

Java

```
protected CComboBox m_ModelComboCtrl;
```

MATLAB

```
m_ModelComboCtrl
```

Description

This is m_ModelComboCtrl, a member of class CBodystatResultsDlg.

4.1.1.2.3.20 CBodystatResultsDlg::m_ResultsListCtrl Data Member**C++**

```
CListCtrl m_ResultsListCtrl;
```

C#

```
protected CListCtrl m_ResultsListCtrl;
```

Visual Basic

```
Protected m_ResultsListCtrl As CListCtrl
```

Java

```
protected CListCtrl m_ResultsListCtrl;
```

MATLAB

```
m_ResultsListCtrl
```

Description

This is m_ResultsListCtrl, a member of class CBodystatResultsDlg.

4.1.1.2.4 CBodystatResultsDlg Methods**4.1.1.2.4.1 CBodystatResultsDlg::CopyToClipboard Method**

Copies measurement parameters, results and normals to the clipboard. Useful if you need to make note of a particular test set for diagnostics.

C++

```
void CopyToClipboard();
```

C#

```
protected CopyToClipboard();
```

Visual Basic

```
Protected Sub CopyToClipboard()
```

Java

```
protected CopyToClipboard();
```

MATLAB

```
CopyToClipboard
```

Body Source

```
void CBodystatResultsDlg::CopyToClipboard()
{
    CString strErrorMsg, strBuffer, strTemp;

    // Copy measurement input parameters to the clipboard
    strBuffer += _T("Bodystat Data\r\n");
    strBuffer += _T("-----\r\n");

    CString strModel;
    BSGetDeviceModelName((BSDeviceModel)m_iModel, strModel.GetBuffer(32), 32);
    strModel.ReleaseBuffer();
    strTemp.Format(_T("Model:\t%s\r\n"), strModel);    strBuffer += strTemp;
```

```

    strTemp.Format(_T("Test Date/Time:\t%s\r\n"), m_dtTestDate.m_status == COleDateTime::null
? m_dtTestDate.Format() : _T("Unknown")); strBuffer += strTemp;
    strTemp.Format(_T("Gender:\t%s\r\n"), m_iGender == 0 ? _T("Male") : _T("Female"));
strBuffer += strTemp;
    strTemp.Format(_T("Age:\t%d\tyears\r\n"), m_iAge); strBuffer += strTemp;
    strTemp.Format(_T("Height:\t%d\tcm\r\n"), m_iHeight); strBuffer += strTemp;
    strTemp.Format(_T("Weight:\t%.1f\tkg\r\n"), m_fWeight); strBuffer += strTemp;
    strTemp.Format(_T("Activity Level:\t%d\r\n"), m_iActivityLevel); strBuffer += strTemp;
    strTemp.Format(_T("Waist:\t%d\tcm\r\n"), m_iWaist); strBuffer += strTemp;
    strTemp.Format(_T("Hip:\t%d\tcm\r\n"), m_iHip); strBuffer += strTemp;

    strTemp.Format(_T("Body Composition Mode:\t%s\r\n"), m_bModeBodyComposition ? _T("Enabled")
: _T("Disabled")); strBuffer += strTemp;
    strTemp.Format(_T("Hydration Mode:\t%s\r\n"), m_bModeHydration ? _T("Enabled") :
_T("Disabled")); strBuffer += strTemp;

    strTemp.Format(_T("Impedance 5 kHz:\t%d\ttohmr\n"), m_iImpedance5K); strBuffer +=
strTemp;
    strTemp.Format(_T("Impedance 50 kHz:\t%d\ttohmr\n"), m_iImpedance50K); strBuffer +=
strTemp;
    strTemp.Format(_T("Impedance 100 kHz:\t%d\ttohmr\n"), m_iImpedance100K); strBuffer +=
strTemp;
    strTemp.Format(_T("Impedance 200 kHz:\t%d\ttohmr\n"), m_iImpedance200K); strBuffer +=
strTemp;
    strTemp.Format(_T("Resistance 50 kHz:\t%d\ttohmr\n"), m_iResistance50K); strBuffer +=
strTemp;
    strTemp.Format(_T("Reactance 50 kHz:\t%.1f\ttohmr\n"), m_fReactance50K); strBuffer +=
strTemp;
    strTemp.Format(_T("Phase Angle 50 kHz:\t%.1f\t°\r\n"), m_fPhaseAngle50K); strBuffer +=
strTemp;

    strBuffer += _T("\r\n");

    // Copy the results and normals to the clipboard from the results list control
    int iRowCount = m_ResultsListCtrl.GetItemRowCount();
    int iColCount = 0;
    CHeaderCtrl* pHeaderCtrl = m_ResultsListCtrl.GetHeaderCtrl();
    if (pHeaderCtrl != NULL)
    {
        HDITEM hdi;
        const int HDR_TEXT_SIZE = 256;
        TCHAR lpBuffer[HDR_TEXT_SIZE];
        hdi.mask = HDI_TEXT;
        hdi.pszText = lpBuffer;
        hdi.cchTextMax = HDR_TEXT_SIZE;

        // Copy column titles
        for(int iDashed=0; iDashed < 2; iDashed++) // Iterate through column titles
            twice, once for the title, once to underline it with dashes
        {
            iColCount = pHeaderCtrl->GetItemCount();
            for(int i=0; i < iColCount; i++)
            {
                if(pHeaderCtrl->GetItem(i, &hdi))
                {
                    if(hdi.mask & HDI_TEXT && hdi.pszText != LPSTR_TEXTCALLBACK)
                    {
                        strTemp = CString(lpBuffer);
                        if(iDashed == 0)
                            strBuffer += strTemp; // Use the title
                        else
                            strBuffer += CString(_T('-'), strTemp.GetLength()); // Underline the title with
dashes
                    }
                }
            }
            strBuffer += '\t';
        }
        strBuffer += "\r\n";
    }
}

```

```

// Build clipboard text from the list control contents (results and normals)
for (int i=0; i < iRowCount; i++)
{
    for (int j=0; j < iColCount; j++)
        strBuffer += m_ResultsListCtrl.GetItemText(i, j) + '\t';

    strBuffer += "\r\n";
}

// Prepare the clipboard
if (!OpenClipboard())
{
    AfxMessageBox(_T("The clipboard is temporarily unavailable"));
    return;
}
if (!EmptyClipboard())
{
    CloseClipboard();
    AfxMessageBox(_T("The clipboard cannot be emptied"));
    return;
}

// Allocate clipboard buffer
HGLOBAL hGlobal;
hGlobal = GlobalAlloc(GHND|GMEM_SHARE, (strBuffer.GetLength()+1)* sizeof (TCHAR));
if (!hGlobal)
{
    CloseClipboard();
    AfxMessageBox(CString("Error allocating memory for clipboard data: ") +
        GetSystemErrorMessage(GetLastError()));
    return;
}

// Prevent buffer from moving around and copy the text to it
TCHAR * szClipboard = (TCHAR *)GlobalLock(hGlobal);
if(szClipboard)
{
    StringCchCopy(szClipboard, (strBuffer.GetLength()+1)* sizeof (TCHAR), strBuffer);
    GlobalUnlock(hGlobal);
}

// Determine clipboard format (Unicode if available)
#ifdef UNICODE
    UINT uFmt = CF_UNICODETEXT;
#else
    UINT uFmt = CF_TEXT;
#endif

// Assign clipboard buffer to the clipboard
if (!SetClipboardData(uFmt, hGlobal))
{
    AfxMessageBox(CString("Error copying to the clipboard: ") +
        GetSystemErrorMessage(GetLastError()));
}
CloseClipboard();
}

```

4.1.1.2.4.2 CBodystatResultsDlg::DoDataExchange Method

C++

```
virtual void DoDataExchange(CDataExchange* pDX);
```

C#

```
protected DoDataExchange(ref CDataExchange pDX);
```

Visual Basic

```
Protected Sub DoDataExchange(ByRef pDX As CDataExchange)
```

Java

```
protected DoDataExchange(CDataExchange pDX);
```

MATLAB

DoDataExchange

Description

DDX/DDV support

Bind the form controls to member variables and set basic input range limits

Body Source

```

void CBodystatResultsDlg::DoDataExchange(CDataExchange* pDX)
{
    CDialog::DoDataExchange(pDX);

    DDX_Control(pDX, IDC_MODEL_COMBO, m_ModelComboCtrl);
    DDX_CBIndex(pDX, IDC_MODEL_COMBO, m_iModel);

    DDX_Radio(pDX, IDC_GENDER_MALE_RADIO, m_iGender);
    DDX_Text(pDX, IDC_AGE_EDIT, m_iAge);
    DDV_MinMaxInt(pDX, m_iAge, 1, 99);
    DDX_Text(pDX, IDC_HEIGHT_EDIT, m_iHeight);
    DDV_MinMaxInt(pDX, m_iHeight, 30, 229);
    DDX_Text(pDX, IDC_WEIGHT_EDIT, m_fWeight);
    DDV_MinMaxFloat(pDX, m_fWeight, 0.5, 300.0);
    DDX_CBIndex(pDX, IDC_ACTIVITY_COMBO, m_iActivityLevel);
    DDX_Text(pDX, IDC_WAIST_EDIT, m_iWaist);
    DDV_MinMaxInt(pDX, m_iWaist, 5, 648);
    DDX_Text(pDX, IDC_HIP_EDIT, m_iHip);
    DDV_MinMaxInt(pDX, m_iHip, 5, 648);
    DDX_Check(pDX, IDC_MODE_BODYCCOMP_CHECK, m_bModeBodyComposition);
    DDX_Check(pDX, IDC_MODE_HYDRATION_CHECK, m_bModeHydration);

    DDX_Text(pDX, IDC_IMP_5K_EDIT, m_iImpedance5K);
    DDV_MinMaxInt(pDX, m_iImpedance5K, 0, 2000);
    DDX_Text(pDX, IDC_IMP_50K_EDIT, m_iImpedance50K);
    DDV_MinMaxInt(pDX, m_iImpedance50K, 0, 2000);
    DDX_Text(pDX, IDC_IMP_100K_EDIT, m_iImpedance100K);
    DDV_MinMaxInt(pDX, m_iImpedance100K, 0, 2000);
    DDX_Text(pDX, IDC_IMP_200K_EDIT, m_iImpedance200K);
    DDV_MinMaxInt(pDX, m_iImpedance200K, 0, 2000);
    DDX_Text(pDX, IDC_RES_50K_EDIT, m_iResistance50K);
    DDV_MinMaxInt(pDX, m_iResistance50K, 0, 2000);
    DDX_Text(pDX, IDC_REAC_50K_EDIT, m_fReactance50K);
    DDV_MinMaxFloat(pDX, m_fReactance50K, 0, 200);
    DDX_Text(pDX, IDC_PHASEANGLE_50K_EDIT, m_fPhaseAngle50K);
    DDV_MinMaxFloat(pDX, m_fPhaseAngle50K, 0, 15);

    DDX_Control(pDX, IDC_RESULTS_LIST, m_ResultsListCtrl);
}

```

4.1.1.2.4.3 CBodystatResultsDlg::GetSystemErrorMessage Method

Gets a system error message.

C++

```
CString GetSystemErrorMessage(DWORD dwError);
```

C#

```
protected CString GetSystemErrorMessage(DWORD dwError);
```

Visual Basic

```
Protected Function GetSystemErrorMessage(dwError As DWORD) As CString
```

Java

```
protected CString GetSystemErrorMessage(DWORD dwError);
```

MATLAB

GetSystemErrorMessage

Parameters

Parameters	Description
DWORD dwError	The system error code to look up - typically from GetLastError()

Returns

The system error message as a string.

Body Source

```
CString CBodystatResultsDlg::GetSystemErrorMessage( DWORD dwError )
{
    LPTSTR lpMsgBuf;

    FormatMessage(
        FORMAT_MESSAGE_ALLOCATE_BUFFER |
        FORMAT_MESSAGE_FROM_SYSTEM |
        FORMAT_MESSAGE_IGNORE_INSERTS,
        NULL,
        dwError,
        MAKELANGID(LANG_NEUTRAL, SUBLANG_DEFAULT),
        (LPTSTR) &lpMsgBuf,
        0, NULL );

    CString strErrorMsg;
    strErrorMsg = lpMsgBuf;

    return strErrorMsg;
}
```

4.1.1.2.4.4 CBodystatResultsDlg::OnClickedCalculateResultsButton Method

Executes the Calculate Results button action. Uses the Bodystat API to calculate the results and normals for this test measurement.

C++

```
afx_msg void OnClickedCalculateResultsButton();
```

C#

```
protected afx_msg void OnClickedCalculateResultsButton();
```

Visual Basic

```
Protected Function OnClickedCalculateResultsButton() As afx_msg void
```

Java

```
protected afx_msg void OnClickedCalculateResultsButton();
```

MATLAB

```
OnClickedCalculateResultsButton
```

Body Source

```
void CBodystatResultsDlg::OnClickedCalculateResultsButton()
{
    UpdateData(TRUE);          // Update member variables from form data

    m_ResultsListCtrl.DeleteAllItems(); // Clear the results list

    if ((BSDeviceModel)m_iModel == BSDeviceModel::BSQuadScanTouch || (BSDeviceModel)m_iModel
    == BSDeviceModel::BSMultiScanTouch)
    {
        AfxMessageBox(_T("The selected model is not supported by this SDK sample."),
        MB_ICONINFORMATION);
        return;
    }

    //
```

```

// Prepare input structure for the Bodystat SDK (using member values from form data)
// More likely this would be obtained directly from a BSMeasurement downloaded from the
unit
//
BSMeasurement bm;
bm.iDeviceModel = (BSDeviceModel)m_iModel; // Use our combo selection for Bodystat model
type (its a direct mapping)
bm.tTestDate = m_dtTestDate; // Not relevant for performing calculations, but might be
copied to clipboard later
bm.iGender = m_iGender == 0 ? Bodystat::BSMale : Bodystat::BSFemale; // Convert our
gender radio buttons into Bodystat gender type
bm.iAge = m_iAge; // Age in years
bm.iHeight = m_iHeight; // Height in centimeters
bm.fWeight = m_fWeight; // Weight in kilograms
bm.iActivity = m_iActivityLevel; // Physical activity level (range 1 to 5) 1=very low
2=low/medium 3=medium 4=medium/high 5=very high
bm.iWaist = m_iWaist; // Waist in centimeters
bm.iHip = m_iHip; // Hip in centimeters
bm.iZ_5kHz = m_iImpedance5K; // Impedance 5 KHz (in ohms) (MDD/Quad only)
bm.iZ_50kHz = m_iImpedance50K; // Impedance 50 KHz (in ohms)
bm.iZ_100kHz = m_iImpedance100K; // Impedance 100 KHz (in ohms) (Quad only)
bm.iZ_200kHz = m_iImpedance200K; // Impedance 200 KHz (in ohms) (Quad only)
bm.iR_50kHz = m_iResistance50K; // Resistance 50 kHz (in ohms) (MDD/Quad only)
bm.fX_50kHz = m_fReactance50K; // Reactance 50 kHz (in ohms) (MDD/Quad only)
bm.fPA_50kHz = m_fPhaseAngle50K; // Phase angle 50 kHz (in degrees) (MDD/Quad only)
bm.iFrequencies = 0; // Extended QST/MST frequencies not supported in this sample
bm.pMultifreqData = NULL; // Extended QST/MST frequencies not supported in this sample

// Choose analysis mode based upon selected analysis checkboxes
// These options might normally be stored as program setting rather than asked of the user
each time (like the installation mode of the device)
// The analysis mode impacts the equations used behind some results in certain scenarios
BSAnalysisMode iAnalysisMode;
if(m_bModeBodyComposition && m_bModeHydration)
    iAnalysisMode = BSAnalysisModeBoth;
else if(m_bModeBodyComposition)
    iAnalysisMode = BSAnalysisModeBC;
else if(m_bModeHydration)
    iAnalysisMode = BSAnalysisModeHN;
else
    iAnalysisMode = BSAnalysisModeNone;

//
// Get the results from the Bodystat SDK
//
BSError bseStatus;
BSResults br;
bseStatus = BSCalculateResults(&bm, &br, iAnalysisMode); // Results stored in br structure

// Check for any errors from the calculation function
if(bseStatus != Bodystat::NoError)
{
    CString strMsg;
    strMsg.Format(_T("Error calculating results.\nError Code: %d"), (int)bseStatus);
    AfxMessageBox(strMsg, MB_ICONEXCLAMATION);

    return; // Nothing to display
}

//
// Get the normal ranges from the Bodystat SDK
//
BSNormals bn;
bseStatus = BSCalculateNormals(&bm, &br, &bn); // Normals stored in bn structure

// Check for any errors from the normals function
if(bseStatus != Bodystat::NoError)
{
    CString strMsg;
    strMsg.Format(_T("Error calculating normals.\nError Code: %d"), (int)bseStatus);
    AfxMessageBox(strMsg, MB_ICONEXCLAMATION);
}

```

```

    return; // Nothing to display
}

// Populate the list control with the calculated results and calculated normals
PopulateListCtrl(bm, br, bn);

}

```

4.1.1.2.4.5 CBodystatResultsDlg::OnClickedCopyToClipboardButton Method

Executes the Copy To Clipboard button action.

C++

```
afx_msg void OnClickedCopyToClipboardButton();
```

C#

```
protected afx_msg void OnClickedCopyToClipboardButton();
```

Visual Basic

```
Protected Function OnClickedCopyToClipboardButton() As afx_msg void
```

Java

```
protected afx_msg void OnClickedCopyToClipboardButton();
```

MATLAB

```
OnClickedCopyToClipboardButton
```

Body Source

```

void CBodystatResultsDlg::OnClickedCopyToClipboardButton()
{
    CopyToClipboard();
}

```

4.1.1.2.4.6 CBodystatResultsDlg::OnInitDialog Method

C++

```
virtual BOOL OnInitDialog();
```

C#

```
public BOOL OnInitDialog();
```

Visual Basic

```
Public Function OnInitDialog() As BOOL
```

Java

```
public BOOL OnInitDialog();
```

MATLAB

```
OnInitDialog
```

Description

CBodystatResultsDlg (see page 185) message handlers

Body Source

```

BOOL CBodystatResultsDlg::OnInitDialog()
{
    CDialog::OnInitDialog();

    // Setup the results list control with three columns ready for the results and normal
    values
    m_ResultsListCtrl.InsertColumn(0, _T("Result"), 0, 250);
    m_ResultsListCtrl.InsertColumn(1, _T("Value"), 0, 110);
    m_ResultsListCtrl.InsertColumn(2, _T("Normal"), 0, 140);
    m_ResultsListCtrl.SetExtendedStyle( m_ResultsListCtrl.GetExtendedStyle() |
    LVS_EX_FULLROWSELECT );
}

```

```
//
// Populate model combo with all possible models
//
int nItem = m_ModelComboCtrl.InsertString(0, _T("Unknown"));
CString strModel, strComboItem;
BSError bsErr = NoError;
int iModel = 1;
do
{
    bsErr = BSGetDeviceModelName((BSDeviceModel)iModel, strModel.GetBuffer(32), 32);
    strModel.ReleaseBuffer();
    if(bsErr == NoError)
    {
        strComboItem.Format(_T("%s (ID:%d)", strModel, iModel); // Append model id number to
the string (for illustration and to make it unique for the combo)
        nItem = m_ModelComboCtrl.InsertString(iModel, strComboItem);
    }
    iModel++;
}while(bsErr == NoError); // Stop when we run out of valid model ids

// Select the desired model
m_ModelComboCtrl.SetCurSel(m_iModel);

return TRUE; // return TRUE unless you set the focus to a control
// EXCEPTION: OCX Property Pages should return FALSE
}
```

4.1.1.2.4.7 CBodystatResultsDlg::PopulateListCtrl Method

Populate the list control with results and normals.

C++

```
void PopulateListCtrl(const BSMeasurement &bm, const BSResults &br, const BSNormals &bn);
```

C#

```
protected PopulateListCtrl(BSMeasurement &bm, BSResults &br, BSNormals &bn);
```

Visual Basic

```
Protected Sub PopulateListCtrl(bm As BSMeasurement &, br As BSResults &, bn As BSNormals &)
```

Java

```
protected PopulateListCtrl(final BSMeasurement &bm, final BSResults &br, final BSNormals
& bn);
```

MATLAB

```
PopulateListCtrl
```

Parameters

Parameters	Description
const BSMeasurement &bm	The measurement parameters for this test.
const BSResults &br	The calculated results for this test from the SDK.
const BSNormals &bn	The calculated normals for this test from the SDK.

Body Source

```
void CBodystatResultsDlg::PopulateListCtrl( const BSMeasurement &bm, const BSResults &br,
const BSNormals &bn )
{
    CString strResult, strNormal;
    int nItem;

    // Fat %
    nItem = m_ResultsListCtrl.InsertItem(0, _T("Fat percentage"));
    strResult.Format(_T("%.1f"), br.fFatPerc); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
    strNormal.Format(_T("%d - %d"), bn.iFatPerc_L, bn.iFatPerc_H);
    m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);
}
```

```

// Fat kg
nItem = m_ResultsListCtrl.InsertItem(1, _T("Fat (in kg)"));
strResult.Format(_T("%.1f"), br.fFatKg); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%d - %d"), bn.iFatKg_L, bn.iFatKg_H);
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Lean %
nItem = m_ResultsListCtrl.InsertItem(2, _T("Lean percentage"));
strResult.Format(_T("%.1f"), br.fLeanPerc); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%d - %d"), bn.iLeanPerc_L, bn.iLeanPerc_H);
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Lean kg
nItem = m_ResultsListCtrl.InsertItem(3, _T("Lean (in kg)"));
strResult.Format(_T("%.1f"), br.fLeanKg); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%d - %d"), bn.iLeanKg_L, bn.iLeanKg_H);
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Total Weight
nItem = m_ResultsListCtrl.InsertItem(4, _T("Total Weight (in kg)"));
strResult.Format(_T("%.1f"), br.fTotalWeight); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
// Total weight normal range is calculated using composition or BMI (whichever is the
greater). We indicate which method was used in brackets after the value.
strNormal.Format(_T("%d - %d (%s)"), bn.iTotalWeight_L, bn.iTotalWeight_H,
bn.iTotalWeightMethod == 0 ? _T("Composition") : _T("BMI"));
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Dry lean weight
nItem = m_ResultsListCtrl.InsertItem(5, _T("Dry Lean Weight (in kg)"));
strResult.Format(_T("%.1f"), br.fDryLW); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// TBW %
nItem = m_ResultsListCtrl.InsertItem(6, _T("Total Body Water percentage"));
strResult.Format(_T("%.1f"), br.fTBWPerc); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%d - %d"), bn.iTBWPerc_L, bn.iTBWPerc_H);
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// TBW lt
nItem = m_ResultsListCtrl.InsertItem(7, _T("Total Body Water (in litres)"));
strResult.Format(_T("%.1f"), br.fTBW); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%d - %d"), bn.iTBW_L, bn.iTBW_H);
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// ECW %
nItem = m_ResultsListCtrl.InsertItem(8, _T("Extra Cellular Water percentage"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.1f"),
br.fECWPerc) : strResult = _T("(MDD/Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strNormal.Format(_T("%d"),
bn.iECWPerc_Norm) : strNormal = _T("(MDD/Quad Only)");
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// ECW lt
nItem = m_ResultsListCtrl.InsertItem(9, _T("Extra Cellular Water (in litres)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.1f"), br.fECW)
: strResult = _T("(MDD/Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1, strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// ICW %
nItem = m_ResultsListCtrl.InsertItem(10, _T("Intra Cellular Water percentage"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||

```

```

BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.1f"),
br.fICWPerc) : strResult = _T("(MDD/Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strNormal.Format(_T("%d"),
bn.iICWPerc_Norm) : strNormal = _T("(MDD/Quad Only)");
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// ICW lt
nItem = m_ResultsListCtrl.InsertItem(11, _T("Intra Cellular Water (in litres)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.1f"), br.fICW)
: strResult = _T("(MDD/Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1, strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Body Cell Mass
nItem = m_ResultsListCtrl.InsertItem(12, _T("Body Cell Mass"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ? strResult.Format(_T("%.1f"),
br.fBCM) : strResult = _T("(Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// 3rd space water
nItem = m_ResultsListCtrl.InsertItem(13, _T("3rd space water (in litres)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ? strResult.Format(_T("%.1f"),
br.fThirdSpace) : strResult = _T("(Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Nutrition index
nItem = m_ResultsListCtrl.InsertItem(14, _T("Nutrition index"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ? strResult.Format(_T("%.2f"),
br.fNutrition) : strResult = _T("(Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ? strNormal.Format(_T("%.2f"),
bn.fNutrition_Norm) : strNormal = _T("(Quad Only)"); m_ResultsListCtrl.SetItemText(nItem,
2, strNormal);

// Prediction Marker (TM) - formerly known as Illness Marker (TM)
nItem = m_ResultsListCtrl.InsertItem(15, _T("Prediction Marker (TM)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ? strResult.Format(_T("%.3f"),
br.fIllness) : strResult = _T("(Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ? strNormal.Format(_T("%.3f -
%.3f"), bn.fIllness_L, bn.fIllness_H) : strNormal = _T("(Quad Only)");
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// BMR kcal
nItem = m_ResultsListCtrl.InsertItem(16, _T("Basal Metabolic Rate (in kcal)"));
strResult.Format(_T("%.0f"), br.fBMR); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// BMR kcal/kg
nItem = m_ResultsListCtrl.InsertItem(17, _T("Basal Metabolic Rate per kilogram (in
kcal/kg)"));
strResult.Format(_T("%.1f"), br.fBMRkg); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// EAR
nItem = m_ResultsListCtrl.InsertItem(18, _T("Estimated Average Requirement (in kcal)"));
strResult.Format(_T("%.0f"), br.fEstAvg); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("N/A")); m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// BMI
nItem = m_ResultsListCtrl.InsertItem(19, _T("Body Mass Index (BMI)"));
strResult.Format(_T("%.1f"), br.fBMI); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%d - %d"), bn.iBMI_L, bn.iBMI_H);
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

```

```

// BFMI
nItem = m_ResultsListCtrl.InsertItem(20, _T("Body Fat Mass Index (BFMI)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.1f"), br.fBFMI)
: strResult = _T("(MDD/Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1, strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strNormal.Format(_T("%d - %d"),
bn.iBFMI_L, bn.iBFMI_H) : strNormal = _T("(MDD/Quad Only)");
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// FFMI
nItem = m_ResultsListCtrl.InsertItem(21, _T("Fat Free Mass Index (FFMI)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.1f"), br.fFFMI)
: strResult = _T("(MDD/Quad Only)"); m_ResultsListCtrl.SetItemText(nItem, 1, strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSQuadScan_Family ||
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strNormal.Format(_T("%d - %d"),
bn.iFFMI_L, bn.iFFMI_H) : strNormal = _T("(MDD/Quad Only)");
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);

// Waist/Hip
nItem = m_ResultsListCtrl.InsertItem(22, _T("Waist/Hip ratio"));
strResult.Format(_T("%.2f"), br.fWaistHip); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
strNormal.Format(_T("%.2f"), bn.fWaistHip_Norm); m_ResultsListCtrl.SetItemText(nItem, 2,
strNormal);

// Wellness marker
nItem = m_ResultsListCtrl.InsertItem(23, _T("Wellness Marker (TM)"));
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strResult.Format(_T("%.3f"),
br.fWellness) : strResult = _T("(MDD Only)"); m_ResultsListCtrl.SetItemText(nItem, 1,
strResult);
BSGetDeviceFamily(bm.iDeviceModel) == BSMDD_Family ? strNormal.Format(_T("%.3f - %.3f"),
bn.fWellness_L, bn.fWellness_H) : strNormal = _T("(MDD Only)");
m_ResultsListCtrl.SetItemText(nItem, 2, strNormal);
}

```

4.2 Bodystat C# WinForms Sample

A sample application written in C# .NET based upon WinForms within the Microsoft .NET Framework. Solutions are provided for Visual Studio 2008 and 2010. You must have Visual Studio 2008 or 2010 installed to work with this sample.

Namespaces

Name	Description
BodystatCSWinForms (see page 205)	This is the main namespace for our sample project.

Files

Name	Description
MainForm.cs (see page 216)	<p>Bodystat VB.Net Sample (C) 2010-2012 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p>

ResultsForm.cs (see page 216)	<p>Bodystat C# WinForms Sample (C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p>
-------------------------------	---



4.2.1 BodystatCSWinForms Namespace

This is the main namespace for our sample project.

Module

Bodystat C# WinForms Sample (see page 204)



Classes

	Name	Description
	MainForm (see page 205)	Main sample dialog. Provides buttons utilising the various functions of the Bodystat API. A list control also shows data (test measurements) downloaded from the device.
	ResultsForm (see page 209)	Results sample dialog. Displays the various parameters relating to a measurement (allowing direct input by the user), then utilises the Bodystat API to calculate detailed results and normals.

4.2.1.1 Classes

The following table lists classes in this documentation.


Classes

	Name	Description
	MainForm (see page 205)	Main sample dialog. Provides buttons utilising the various functions of the Bodystat API. A list control also shows data (test measurements) downloaded from the device.
	ResultsForm (see page 209)	Results sample dialog. Displays the various parameters relating to a measurement (allowing direct input by the user), then utilises the Bodystat API to calculate detailed results and normals.

4.2.1.1.1 MainForm Class

Main sample dialog. Provides buttons utilising the various functions of the Bodystat API. A list control also shows data (test measurements) downloaded from the device.

Class Hierarchy

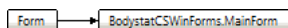

C++

```
class MainForm : public Form;
```

C#

```
public class MainForm : Form;
```

Visual Basic

```
Public Class MainForm
```


Inherits Form

Java

`public class MainForm Form;`

MATLAB

MainForm

File

MainForm.cs (🔗 see page 216)

Methods

	Name	Description
🔗	MainForm (🔗 see page 206)	Prepare the controls on the sample dialog

MainForm Methods

	Name	Description
🔗	PopulateMeasurementList (🔗 see page 207)	Populate the measurement list control with test measurement data from the API. Illustrates iteration through a BSRawData structure, working with data in a BSMeasurement structure and helper functions from the Bodystat API like BSGetDeviceModelName

4.2.1.1.1.1 **MainForm.MainForm Constructor**

Prepare the controls on the sample dialog

C++

`MainForm();`

C#

`public MainForm();`

Visual Basic

`Public Sub New()`

Java

`public MainForm();`

MATLAB

MainForm

Body Source

```
public MainForm()
{
    InitializeComponent();

    // Add columns to the test measurement list ready for download of data
    this.lstviewTestMeasurementList.Columns.Add("Test No", 50);
    this.lstviewTestMeasurementList.Columns.Add("Test Date/Time", 130);
    this.lstviewTestMeasurementList.Columns.Add("Model", 140);
    this.lstviewTestMeasurementList.Columns.Add("Serial No", 75);
    this.lstviewTestMeasurementList.Columns.Add("Gender", 75);
    this.lstviewTestMeasurementList.Columns.Add("Age", 75);
    this.lstviewTestMeasurementList.Columns.Add("Height", 75);
    this.lstviewTestMeasurementList.Columns.Add("Weight", 75);
    this.lstviewTestMeasurementList.Columns.Add("Activity", 75);
    this.lstviewTestMeasurementList.Columns.Add("Waist", 75);
    this.lstviewTestMeasurementList.Columns.Add("Hip", 75);
    this.lstviewTestMeasurementList.Columns.Add("Imp 5 kHz", 75);
    this.lstviewTestMeasurementList.Columns.Add("Imp 50 kHz", 75);
    this.lstviewTestMeasurementList.Columns.Add("Imp 100 kHz", 75);
    this.lstviewTestMeasurementList.Columns.Add("Imp 200 kHz", 75);
    this.lstviewTestMeasurementList.Columns.Add("Res 50 kHz", 75);
    this.lstviewTestMeasurementList.Columns.Add("Reac 50 kHz", 75);
    this.lstviewTestMeasurementList.Columns.Add("Phase 50 kHz", 80);
```

```
}
```

4.2.1.1.1.2 MainForm Methods

4.2.1.1.1.2.1 MainForm.PopulateMeasurementList Method

Populate the measurement list control with test measurement data from the API. Illustrates iteration through a BSRawData structure, working with data in a BSMeasurement structure and helper functions from the Bodystat API like BSGetDeviceModelName

C++

```
void PopulateMeasurementList(BodystatAPI.BodystatAPI.BSRawData ^ pRawData);
```

C#

```
public void PopulateMeasurementList(ref BodystatAPI.BodystatAPI.BSRawData pRawData);
```

Visual Basic

```
Public Function PopulateMeasurementList(ByRef pRawData As
BodystatAPI.BodystatAPI.BSRawData) As void
```

Java

```
public void PopulateMeasurementList(BodystatAPI.BodystatAPI.BSRawData pRawData);
```

MATLAB

```
PopulateMeasurementList
```

Parameters

Parameters	Description
ref BodystatAPI.BodystatAPI.BSRawData pRawData	[in] RawData to be parsed.

Body Source

```
public void PopulateMeasurementList(ref BodystatAPI.BodystatAPI.BSRawData pRawData)
{
    // Clear any previous items in the list
    this.lstviewTestMeasurementList.Items.Clear();

    string strItem = null;

    System.Windows.Forms.ListViewItem lvi = null;
    int i = 0;

    // Process each record in the structure
    for (i = 0; i <= pRawData.iTotalNumRecs - 1; i++) {
        // Test no
        strItem = string.Format("{0}", i + pRawData.ulFirstTestNum);
        lvi = lstviewTestMeasurementList.Items.Insert(i, strItem);
        lvi.Tag = pRawData.record[i];
        // Assign this BSMeasurement to the listview item (so we can reference it later)

        // Test date/time
        lvi.SubItems.Add(pRawData.record[i].tTestDate.ToString());

        // Device model
        BodystatAPI.BodystatAPI.BSGetDeviceModelName(pRawData.record[i].iDeviceModel, out
strItem);
        lvi.SubItems.Add(strItem);

        // Device serial number
        strItem = string.Format("{0}", pRawData.record[i].ulDeviceSerialNumber);
        lvi.SubItems.Add(strItem);

        // Gender
        lvi.SubItems.Add((pRawData.record[i].iGender == BSGender.BSMale ? "Male" : "Female"));

        // Age
```

```
strItem = string.Format("{0}", pRawData.record[i].iAge);
lvi.SubItems.Add(strItem);

// Height
strItem = string.Format("{0} cm", pRawData.record[i].iHeight);
lvi.SubItems.Add(strItem);

// Weight
strItem = string.Format("{0:F1} kg", pRawData.record[i].fWeight);
lvi.SubItems.Add(strItem);

// Activity
switch (pRawData.record[i].iActivity) {
    case 1:
        strItem = "Very Low";
        break;
    case 2:
        strItem = "Low/Medium";
        break;
    case 3:
        strItem = "Medium";
        break;
    case 4:
        strItem = "Medium/High";
        break;
    case 5:
        strItem = "Very High";
        break;
    default:
        strItem = "Unknown";
        break;
}
lvi.SubItems.Add(strItem);

// Waist
strItem = string.Format("{0} cm", pRawData.record[i].iWaist);
lvi.SubItems.Add(strItem);

// Hip
strItem = string.Format("{0} cm", pRawData.record[i].iHip);
lvi.SubItems.Add(strItem);

// Z5 (Impedance at 5 kHz)
strItem = string.Format("{0} ohm", pRawData.record[i].iZ_5kHz);
lvi.SubItems.Add(strItem);

// Z50 (Impedance at 50 kHz)
strItem = string.Format("{0} ohm", pRawData.record[i].iZ_50kHz);
lvi.SubItems.Add(strItem);

// Z100 (Impedance at 100 kHz)
strItem = string.Format("{0} ohm", pRawData.record[i].iZ_100kHz);
lvi.SubItems.Add(strItem);

// Z200 (Impedance at 200 kHz)
strItem = string.Format("{0} ohm", pRawData.record[i].iZ_200kHz);
lvi.SubItems.Add(strItem);

// R50 (Resistance at 50 kHz)
strItem = string.Format("{0} ohm", pRawData.record[i].iR_50kHz);
lvi.SubItems.Add(strItem);


// X50 (Reactance at 50 kHz)
strItem = string.Format("{0:F1} ohm", pRawData.record[i].fX_50kHz);
lvi.SubItems.Add(strItem);

// PA50 (Phase Angle at 50 kHz)
strItem = string.Format("{0:F1} deg", pRawData.record[i].fPA_50kHz);
lvi.SubItems.Add(strItem);
}
}
```

4.2.1.1.2 ResultsForm Class

Results sample dialog. Displays the various parameters relating to a measurement (allowing direct input by the user), then utilises the Bodystat API to calculate detailed results and normals.

Class Hierarchy

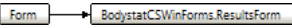

C++

```
class ResultsForm : public Form;
```

C#

```
public class ResultsForm : Form;
```

Visual Basic

```
Public Class ResultsForm
Inherits Form
```

Java

```
public class ResultsForm Form;
```

MATLAB

```
ResultsForm
```

File

ResultsForm.cs (see page 216)

Methods

	Name	Description
	ResultsForm (see page 209)	Initialise the dialog controls with their content

ResultsForm Methods

	Name	Description
	InitialiseFromMeasurement (see page 210)	Initialise the dialog controls with a values from an actual measurement result
	InitialiseSampleMeasurement (see page 211)	Initialise the dialog controls with a sample measurement result
	PopulateListCtrl (see page 212)	Populate the list control with results and normals.
	UseBSMeasurement (see page 215)	Allows the caller to specify a downloaded measurement structure to be use when initialising the dialog. Must be called before the form is loaded (displayed)

4.2.1.1.2.1 ResultsForm.ResultsForm Constructor

Initialise the dialog controls with their content

C++

```
ResultsForm() ;
```

C#

```
public ResultsForm() ;
```

Visual Basic

```
Public Sub New()
```

Java

```
public ResultsForm() ;
```

MATLAB

ResultsForm

Body Source

```

public ResultsForm()
{
    InitializeComponent();

    //
    // Populate model combo with all possible models
    //
    this.cboModelCombo.Items.Add("Unknown");
    string strModel;
    string strComboItem = null;
    BodystatAPI.BSError bsErr = BodystatAPI.BSError.NoError;
    int iModel = 1;
    do
    {
        bsErr = BodystatAPI.BodystatAPI.BSGetDeviceModelName((BSDeviceModel)iModel, out
strModel);

        if ((bsErr == BodystatAPI.BSError.NoError))
        {
            strComboItem = string.Format("{0} (ID:{1})", strModel, iModel);
            // Append model id number to the string (for illustration and to make it unique
for the combo)
            this.cboModelCombo.Items.Add(strComboItem);

        }
        iModel += 1;
    } while ((bsErr == BodystatAPI.BSError.NoError));
    // Stop when we run out of valid model ids

    // Setup the results list control with three columns ready for the results and normal
values
    lstviewResultsList.Columns.Add("Result", 250);
    lstviewResultsList.Columns.Add("Value", 110);
    lstviewResultsList.Columns.Add("Normal", 140);
}

```

4.2.1.1.2.2 ResultsForm Methods**4.2.1.1.2.2.1 ResultsForm.InitialiseFromMeasurement Method**

Initialise the dialog controls with a values from an actual measurement result

C++

```
void InitialiseFromMeasurement();
```

C#

```
public void InitialiseFromMeasurement();
```

Visual Basic

```
Public Function InitialiseFromMeasurement() As void
```

Java

```
public void InitialiseFromMeasurement();
```

MATLAB

InitialiseFromMeasurement

Body Source

```

public void InitialiseFromMeasurement()
{
    this.cboModelCombo.SelectedIndex = (int)bsmMeasurement.iDeviceModel;
    //Me.dtTestDate = bsmMeasurement.tTestDate      ' Not shown on dialog presently
    this.rdoGenderMaleRadio.Checked = bsmMeasurement.iGender == BodystatAPI.BSGender.BSMale;
}

```

```

        this.rdoGenderFemaleRadio.Checked = bsmMeasurement.iGender ==
BodystatAPI.BSGender.BSFemale;
        this.txtAgeEdit.Text = bsmMeasurement.iAge.ToString();
        this.txtHeightEdit.Text = bsmMeasurement.iHeight.ToString();
        this.txtWeightEdit.Text = bsmMeasurement.fWeight.ToString();
        this.cboActivityCombo.SelectedIndex = bsmMeasurement.iActivity;
        this.txtWaistEdit.Text = bsmMeasurement.iWaist.ToString();
        this.txtHipEdit.Text = bsmMeasurement.iHip.ToString();

        this.chkModeBodyCompCheck.Checked = true;
        this.chkModeHydrationCheck.Checked = true;

        this.txtImp5kEdit.Text = bsmMeasurement.iZ_5kHz.ToString();
        this.txtImp50kEdit.Text = bsmMeasurement.iZ_50kHz.ToString();
        this.txtImp100kEdit.Text = bsmMeasurement.iZ_100kHz.ToString();
        this.txtImp200kEdit.Text = bsmMeasurement.iZ_200kHz.ToString();
        this.txtRes50kEdit.Text = bsmMeasurement.iR_50kHz.ToString();
        this.txtReac50kEdit.Text = bsmMeasurement.fX_50kHz.ToString();
        this.txtPhaseAngle50kEdit.Text = bsmMeasurement.fPA_50kHz.ToString();
    }

```

4.2.1.1.2.2.2 ResultsForm.InitialiseSampleMeasurement Method

Initialise the dialog controls with a sample measurement result

C++

```
void InitialiseSampleMeasurement();
```

C#

```
public void InitialiseSampleMeasurement();
```

Visual Basic

```
Public Function InitialiseSampleMeasurement() As void
```

Java

```
public void InitialiseSampleMeasurement();
```

MATLAB

```
InitialiseSampleMeasurement
```

Body Source

```

public void InitialiseSampleMeasurement()
{
    // Setup the form with some sensible default values
    this.cboModelCombo.SelectedIndex = (int)BodystatAPI.BSDeviceModel.BSQuadScan_BT;
    //Me.dtTestDate = DateTime.Now() ' Not shown on dialog presently
    this.rdoGenderMaleRadio.Checked = true;
    this.rdoGenderFemaleRadio.Checked = false;
    this.txtAgeEdit.Text = "30";
    this.txtHeightEdit.Text = "165";
    this.txtWeightEdit.Text = "85.0";
    this.cboActivityCombo.SelectedIndex = 3;
    this.txtWaistEdit.Text = "85";
    this.txtHipEdit.Text = "90";

    this.chkModeBodyCompCheck.Checked = true;
    this.chkModeHydrationCheck.Checked = true;

    this.txtImp5kEdit.Text = "500";
    this.txtImp50kEdit.Text = "500";
    this.txtImp100kEdit.Text = "500";
    this.txtImp200kEdit.Text = "500";
    this.txtRes50kEdit.Text = "500";
    this.txtReac50kEdit.Text = "52.0";
    this.txtPhaseAngle50kEdit.Text = "6.1";
}

```

4.2.1.1.2.2.3 ResultsForm.PopulateListCtrl Method

Populate the list control with results and normals.

C++

```
void PopulateListCtrl(BodystatAPI.BodystatAPI.BSMeasurement ^ bm,
BodystatAPI.BodystatAPI.BSResults ^ br, BodystatAPI.BodystatAPI.BSNormals ^ bn);
```

C#

```
public void PopulateListCtrl(ref BodystatAPI.BodystatAPI.BSMeasurement bm, ref
BodystatAPI.BodystatAPI.BSResults br, ref BodystatAPI.BodystatAPI.BSNormals bn);
```

Visual Basic

```
Public Function PopulateListCtrl(ByRef bm As BodystatAPI.BodystatAPI.BSMeasurement, ByRef
br As BodystatAPI.BodystatAPI.BSResults, ByRef bn As BodystatAPI.BodystatAPI.BSNormals) As
void
```

Java

```
public void PopulateListCtrl(BodystatAPI.BodystatAPI.BSMeasurement bm,
BodystatAPI.BodystatAPI.BSResults br, BodystatAPI.BodystatAPI.BSNormals bn);
```

MATLAB

PopulateListCtrl

Parameters

Parameters	Description
ref BodystatAPI.BodystatAPI.BSMeasurement bm	The measurement parameters for this test.
ref BodystatAPI.BodystatAPI.BSResults br	The calculated results for this test from the SDK.
ref BodystatAPI.BodystatAPI.BSNormals bn	The calculated normals for this test from the SDK.

Body Source

```
public void PopulateListCtrl(ref BodystatAPI.BodystatAPI.BSMeasurement bm, ref
BodystatAPI.BodystatAPI.BSResults br, ref BodystatAPI.BodystatAPI.BSNormals bn)
{
    ListViewItem lvi = default(ListViewItem);

    // Fat %
    lvi = lstviewResultsList.Items.Add("Fat percentage");
    lvi.SubItems.Add(string.Format("{0:F1}", br.FatPerc));
    lvi.SubItems.Add(string.Format("{0} - {1}", bn.FatPerc_L, bn.FatPerc_H));

    // Fat kg
    lvi = lstviewResultsList.Items.Add("Fat (in kg)");
    lvi.SubItems.Add(string.Format("{0:F1}", br.FatKg));
    lvi.SubItems.Add(string.Format("{0} - {1}", bn.FatKg_L, bn.FatKg_H));

    // Lean %
    lvi = lstviewResultsList.Items.Add("Lean percentage");
    lvi.SubItems.Add(string.Format("{0:F1}", br.LeanPerc));
    lvi.SubItems.Add(string.Format("{0} - {1}", bn.LeanPerc_L, bn.LeanPerc_H));

    // Lean kg
    lvi = lstviewResultsList.Items.Add("Lean (in kg)");
    lvi.SubItems.Add(string.Format("{0:F1}", br.LeanKg));
    lvi.SubItems.Add(string.Format("{0} - {1}", bn.LeanKg_L, bn.LeanKg_H));

    // Total Weight
    lvi = lstviewResultsList.Items.Add("Total Weight (in kg)");
    lvi.SubItems.Add(string.Format("{0:F1}", br.TotalWeight));
    // Total weight normal range is calculated using composition or BMI (whichever is the
    // greater). We indicate which method was used in brackets after the value.
    lvi.SubItems.Add(string.Format("{0} - {1} ({2})", bn.TotalWeight_L, bn.TotalWeight_H,
    (bn.TotalWeightMethod == 0 ? "Composition" : "BMI")));

    // Dry lean weight
    lvi = lstviewResultsList.Items.Add("Dry Lean Weight (in kg)");
```

```

lvi.SubItems.Add(string.Format("{0:F1}", br.DryLW));
lvi.SubItems.Add("N/A");

// TBW %
lvi = lstviewResultsList.Items.Add("Total Body Water percentage");
lvi.SubItems.Add(string.Format("{0:F1}", br.TBWPerC));
lvi.SubItems.Add(string.Format("{0} - {1}", bn.TBWPerC_L, bn.TBWPerC_H));

// TBW lt
lvi = lstviewResultsList.Items.Add("Total Body Water (in litres)");
lvi.SubItems.Add(string.Format("{0:F1}", br.TBW));
lvi.SubItems.Add(string.Format("{0} - {1}", bn.TBW_L, bn.TBW_H));

// ECW %
lvi = lstviewResultsList.Items.Add("Extra Cellular Water percentage");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0:F1}", br.ECWPerC) : "(MDD/Quad
Only)");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0}", bn.ECWPerC_Norm) :
"(MDD/Quad Only)");

// ECW lt
lvi = lstviewResultsList.Items.Add("Extra Cellular Water (in litres)");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0:F1}", br.ECW) : "(Quad Only)");
lvi.SubItems.Add("N/A");

// ICW %
lvi = lstviewResultsList.Items.Add("Intra Cellular Water percentage");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0:F1}", br.ICWPerC) : "(MDD/Quad
Only)");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0}", bn.ICWPerC_Norm) :
"(MDD/Quad Only)");

// ICW lt
lvi = lstviewResultsList.Items.Add("Intra Cellular Water (in litres)");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0:F1}", br.ICW) : "(MDD/Quad
Only)");
lvi.SubItems.Add("N/A");

// Body Cell Mass
lvi = lstviewResultsList.Items.Add("Body Cell Mass");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ? string.Format("{0:F1}", br.BCM) : "(Quad
Only)");
lvi.SubItems.Add("N/A");

// 3rd space water
lvi = lstviewResultsList.Items.Add("3rd space water (in litres)");
lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ? string.Format("{0:F1}", br.ThirdSpace) :
"(Quad Only)");
lvi.SubItems.Add("N/A");

// Nutrition index
lvi = lstviewResultsList.Items.Add("Nutrition index");

```



```

        lvi.SubItems.Add(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ? string.Format("{0:F2}", br.Nutrition) :
"(Quad Only)");
        lvi.SubItems.Add(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ? string.Format("{0:F2}", bn.Nutrition_Norm) :
"(Quad Only)");

        // Prediction Marker (TM) - formerly known as Illness Marker (TM)
        lvi = lstviewResultsList.Items.Add("Prediction Marker (TM)");
        lvi.SubItems.Add(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ? string.Format("{0:F3}", br.Illness) : "(Quad
Only)");
        lvi.SubItems.Add(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ? string.Format("{0:F3} - {1:F3}",
bn.Illness_L, bn.Illness_H) : "(Quad Only)");

        // BMR kcal
        lvi = lstviewResultsList.Items.Add("Basal Metabolic Rate (in kcal)");
        lvi.SubItems.Add(string.Format("{0:F0}", br.BMR));
        lvi.SubItems.Add("N/A");

        // BMR kcal/kg
        lvi = lstviewResultsList.Items.Add("Basal Metabolic Rate per kilogram (in kcal/kg)");
        lvi.SubItems.Add(string.Format("{0:F1}", br.BMRkg));
        lvi.SubItems.Add("N/A");

        // EAR
        lvi = lstviewResultsList.Items.Add("Estimated Average Requirement (in kcal)");
        lvi.SubItems.Add(string.Format("{0:F0}", br.EstAvg));
        lvi.SubItems.Add("N/A");

        // BMI
        lvi = lstviewResultsList.Items.Add("Body Mass Index (BMI)");
        lvi.SubItems.Add(string.Format("{0:F1}", br.BMI));
        lvi.SubItems.Add(string.Format("{0} - {1}", bn.BMI_L, bn.BMI_H));

        // BFMI
        lvi = lstviewResultsList.Items.Add("Body Fat Mass Index (BFMI)");
        lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0:F1}", br.BFMI) : "(MDD/Quad
Only)");
        lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0} - {1}", bn.BFMI_L, bn.BFMI_H)
: "(MDD/Quad Only)");

        // FFMI
        lvi = lstviewResultsList.Items.Add("Fat Free Mass Index (FFMI)");
        lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0:F1}", br.FFMI) : "(MDD/Quad
Only)");
        lvi.SubItems.Add((BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSQuadScan_Family ||
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family) ? string.Format("{0} - {1}", bn.FFMI_L, bn.FFMI_H)
: "(MDD/Quad Only)");

        // Waist/Hip
        lvi = lstviewResultsList.Items.Add("Waist/Hip ratio");
        lvi.SubItems.Add(string.Format("{0:F2}", br.WaistHip));
        lvi.SubItems.Add(string.Format("{0:F2}", bn.WaistHip_Norm));

        // Wellness marker
        lvi = lstviewResultsList.Items.Add("Wellness Marker (TM)");
        lvi.SubItems.Add(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family ? string.Format("{0:F3}", br.Wellness) : "(MDD
Only)");

```

```
lvi.SubItems.Add(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) ==
BodystatAPI.BSDeviceFamily.BSMDD_Family ? string.Format("{0:F3} - {0:F3}", bn.Wellness_L,
bn.Wellness_H) : "(MDD Only)");
}
```

4.2.1.1.2.2.4 ResultsForm.UseBSMeasurement Method

Allows the caller to specify a downloaded measurement structure to be use when initialising the dialog. Must be called before the form is loaded (displayed)

C++

```
void UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm);
```

C#

```
public void UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm);
```

Visual Basic

```
Public Function UseBSMeasurement(bsm As BodystatAPI.BodystatAPI.BSMeasurement) As void
```

Java

```
public void UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm);
```

MATLAB

```
UseBSMeasurement
```

Body Source

```
public void UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm)
{
    bsmMeasurement = bsm;
    bUseActualMeasurment = true;
}
```

4.2.2 Files

The following table lists files in this documentation.

Files

Name	Description
MainForm.cs (see page 216)	<p>Bodystat VB.Net Sample (C) 2010-2012 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p>
ResultsForm.cs (see page 216)	<p>Bodystat C# WinForms Sample (C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved. Use subject to license.</p> <p>Web: http://www.bodystat.com Email: info@bodystat.com Tel: +44 (0)1624 629 571 Fax: +44 (0)1624 611 544</p>

Module

Bodystat C# WinForms Sample (see page 204)

4.2.2.1 MainForm.cs

Bodystat VB.Net Sample

(C) 2010-2012 Bodystat (Isle of Man) Ltd. All rights reserved.

Use subject to license.

Web: <http://www.bodystat.com>

Email: info@bodystat.com

Tel: +44 (0)1624 629 571

Fax: +44 (0)1624 611 544

Namespaces

Name	Description
BodystatCSWinForms (see page 205)	This is the main namespace for our sample project.

4.2.2.2 ResultsForm.cs

Bodystat C# WinForms Sample

(C) 2010-2014 Bodystat (Isle of Man) Ltd. All rights reserved.

Use subject to license.

Web: <http://www.bodystat.com>

Email: info@bodystat.com

Tel: +44 (0)1624 629 571

Fax: +44 (0)1624 611 544

Namespaces

Name	Description
BodystatCSWinForms (see page 205)	This is the main namespace for our sample project.

4.3 Bodystat Visual Basic .NET Sample

A sample application written in Visual Basic .NET based upon the Microsoft .NET Framework. Solutions are provided for Visual Basic .NET 2008 and 2010. You must have Visual Basic .NET 2008 or 2010 installed to work with this sample.

Namespaces

Name	Description
BodystatVBNetWinFormSample (see page 217)	This is the main namespace for our sample project.

4.3.1 BodystatVBNetWinFormSample Namespace

This is the main namespace for our sample project.

Module

Bodystat Visual Basic .NET Sample (🔗 see page 216)

Classes

	Name	Description
🔗	FrmBodystatMain (🔗 see page 217)	Main sample dialog. Provides buttons utilising the various functions of the Bodystat API. A list control also shows data (test measurements) downloaded from the device.
🔗	FrmBodystatResults (🔗 see page 237)	Results sample dialog. Displays the various parameters relating to a measurement (allowing direct input by the user), then utilises the Bodystat API to calculate detailed results and normals.

4.3.1.1 Classes

The following table lists classes in this documentation.

Classes

	Name	Description
🔗	FrmBodystatMain (🔗 see page 217)	Main sample dialog. Provides buttons utilising the various functions of the Bodystat API. A list control also shows data (test measurements) downloaded from the device.
🔗	FrmBodystatResults (🔗 see page 237)	Results sample dialog. Displays the various parameters relating to a measurement (allowing direct input by the user), then utilises the Bodystat API to calculate detailed results and normals.

4.3.1.1.1 FrmBodystatMain Class

Main sample dialog. Provides buttons utilising the various functions of the Bodystat API. A list control also shows data (test measurements) downloaded from the device.

Class Hierarchy

BodystatVBNetWinFormSample.FrmBodystatMain

C++

```
class FrmBodystatMain;
```

C#

```
public class FrmBodystatMain;
```

Visual Basic

```
Public Class FrmBodystatMain
```

Java

```
public class FrmBodystatMain;
```


MATLAB

```
FrmBodystatMain
```








File

```
FrmBodystatMain.vb
```

FrmBodystatMain Fields

	Name	Description
	iTimeout (see page 219)	Bluetooth timeout multiplier. Initial value of 7 (equates to about 10 seconds). This is incremented upon each Bluetooth search. A maximum value of 45 is permitted (about 1 minute).

FrmBodystatMain Methods

	Name	Description
	btnBtAuthenticateButton_Click (see page 219)	Executes the Bluetooth Authenticate button action. Illustrates use of BSAAuthenticateBTDevice from the Bodystat API
	btnBtAutoSetupButton_Click (see page 220)	Executes the Bluetooth Auto Setup button action. Illustrates use of BSAutoSetupBT from the Bodystat API
	btnBtDeviceInfoButton_Click (see page 220)	Executes the Bluetooth Device Info button action. Illustrates use of BSGetBTBodystatDevice from the Bodystat API
	btnBtSearchButton_Click (see page 221)	Executes the Bluetooth Search button action. Illustrates use of BSSearchBTDevices from the Bodystat API
	btnBtStackInfoButton_Click (see page 222)	Executes the Bluetooth Get Stack Info button action. Illustrates use of BSGetBTStackInfo from the Bodystat API
	btnBtStatusButton_Click (see page 222)	Executes the Bluetooth Get Status button action. Illustrates use of BSIsBTAvailable from the Bodystat API
	btnBtUnauthenticateButton_Click (see page 223)	Executes the Bluetooth Unauthenticate button action. Illustrates use of BSUnAuthenticateBTDevices from the Bodystat API
	btnComConnectButton_Click (see page 223)	Executes the communication Connect button action. Illustrates use of BSOpenComport from the Bodystat API
	btnComReadCalibrationTimeButton_Click (see page 224)	Executes the communication Read Calibration Time button action. Illustrates use of BSReadCalibrationTime from the Bodystat API
	btnComReadCurrentTimeButton_Click (see page 225)	Executes the communication Read Current Time button action. Illustrates use of BSReadCurrentTime from the Bodystat API
	btnComReadModelVersionButton_Click (see page 226)	Executes the communication Read Model Version button action. Illustrates use of BSReadModelVersion from the Bodystat API
	btnComReadPrinterAddrButton_Click (see page 227)	Executes the communication Read Printer Address button action. Illustrates use of BSReadPrinterAddress from the Bodystat API
	btnComReadProtocolInfoButton_Click (see page 228)	Executes the communication Read ProtocolInfo button action. Illustrates use of BSReadProtocolInfo from the Bodystat API
	btnComReadSerialNumberButton_Click (see page 229)	Executes the communication Read Serial Number button action. Illustrates use of BSReadSerialNumber from the Bodystat API
	btnComReadStoredTestDataButton_Click (see page 230)	Executes the communication Read Stored Test Data button action. Illustrates use of BSReadStoredTestData from the Bodystat API to download data
	btnComReadTestMeasurementButton_Click (see page 231)	Executes the communication Read Test Measurement button action. Performs an immediate test measurement on the device using the BSReadTestBodystat function from the Bodystat API
	btnComWriteCurrentTimeButton_Click (see page 232)	Executes the communication Write Current Time button action. Illustrates use of BSWriteCurrentTime from the Bodystat API
	btnResultsButton_Click (see page 233)	Executes the View Results button action. Opens the results form ready for more detailed analysis and normals
	FrmBodystatMain_Load (see page 234)	Prepare the controls on the sample dialog
	IstviewTestMeasurementList_DoubleClick (see page 235)	Executes the Measurements lists double click action. Opens the results form ready for more detailed analysis and normals
	PopulateMeasurementList (see page 235)	Populate the measurement list control with test measurement data from the API. Illustrates iteration through a BSRawData structure, working with data in a BSMeasurement structure and helper functions from the Bodystat API like BSGetDeviceModelName

4.3.1.1.1.1 FrmBodystatMain Fields

4.3.1.1.1.1.1 FrmBodystatMain.iTimeout Field

C++

```
UInt16 iTimeout;
```

C#

```
public UInt16 iTimeout;
```

Visual Basic

```
Public iTimeout As UInt16
```

Java

```
public UInt16 iTimeout;
```

MATLAB

```
iTimeout
```

Description

Bluetooth timeout multiplier. Initial value of 7 (equates to about 10 seconds). This is incremented upon each Bluetooth search. A maximum value of 45 is permitted (about 1 minute).

4.3.1.1.1.2 FrmBodystatMain Methods

4.3.1.1.1.2.1 FrmBodystatMain.btnBtAuthenticateButton_Click Method

Executes the Bluetooth Authenticate button action. Illustrates use of BSAAuthenticateBTDevice from the Bodystat API

C++

```
btnBtAuthenticateButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtAuthenticateButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtAuthenticateButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnBtAuthenticateButton.Click
```

Java

```
private btnBtAuthenticateButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtAuthenticateButton_Click
```

Body Source

```
Private Sub btnBtAuthenticateButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnBtAuthenticateButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing something

    If (BodystatAPI.BodystatAPI.BSAAuthenticateBTDevice(Me.Handle, iTimeout)) Then
        MsgBox("Device authenticated\nClick Device Info button for more information",
        MsgBoxStyle.Information)
    Else
        MsgBox("Error no Bodystat device could be found or authenticated",
        MsgBoxStyle.Information)
    End If

    Me.Cursor = Cursors.Default
End Sub
```

4.3.1.1.1.2.2 FrmBodystatMain.btnBtAutoSetupButton_Click Method

Executes the Bluetooth Auto Setup button action. Illustrates use of BSAutoSetupBT from the Bodystat API

C++

```
btnBtAutoSetupButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtAutoSetupButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtAutoSetupButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnBtAutoSetupButton.Click
```

Java

```
private btnBtAutoSetupButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtAutoSetupButton_Click
```

Body Source

```
Private Sub btnBtAutoSetupButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnBtAutoSetupButton.Click
    Dim strComPort As New StringBuilder(256)

    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing something

    ' Perform auto setup. User will be prompted as needed by the API
    ' Suitable for supported locales only. Non-supported locales will display messages in English
    ' For non-supported locales you should use the search and authenticate functions directly (silently, prompting the user yourself as required)
    If (BodystatAPI.BodystatAPI.BSAutoSetupBT(strComPort, strComPort.Capacity, True, Me.Handle)) Then
        Dim strMsg As String
        strMsg = String.Format("Device setup complete." & vbCrLf & "Com Port: {0}" & vbCrLf & "Click Device Info button for more information", strComPort.ToString())

        Me.txtComPortEdit.Text = strComPort.ToString() ' Assign this port to the edit box on the form

        MsgBox(strMsg, MsgBoxStyle.Information)
    Else
        MsgBox("Error no Bodystat device could be found or authenticated")
    End If

    Me.Cursor = Cursors.Default
End Sub
```

4.3.1.1.1.2.3 FrmBodystatMain.btnBtDeviceInfoButton_Click Method

Executes the Bluetooth Device Info button action. Illustrates use of BSGetBTBodystatDevice from the Bodystat API

C++

```
btnBtDeviceInfoButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtDeviceInfoButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtDeviceInfoButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnBtDeviceInfoButton.Click
```

Java

```
private btnBtDeviceInfoButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtDeviceInfoButton_Click
```

Body Source

```
Private Sub btnBtDeviceInfoButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtDeviceInfoButton.Click
    Dim strDeviceName As New StringBuilder(512)
    Dim strComPort As New StringBuilder(18)
    Dim strBDA As New StringBuilder(256)

    Me.Cursor = Cursors.WaitCursor ' Lengthy operation (if not already searched), show the
user we are doing something

    ' Retrieve information about the Bodystat device
    Dim bSuccess As Boolean
    bSuccess = BodystatAPI.BodystatAPI.BSGetBTBodystatDevice(strDeviceName,
strDeviceName.Capacity, strBDA, strBDA.Capacity, strComPort, strComPort.Capacity, iTimeout)

    If (bSuccess) Then
        Dim strMsg As String
        strMsg = String.Format("Found Bodystat device." & vbCrLf & vbCrLf & "Device Name:
{0}" & vbCrLf & "Bluetooth Address: {1}" & vbCrLf & "Com Port: {2}", strDeviceName, strBDA,
IIf(strComPort.Length = 0, "<requires authentication>", strComPort))

        Me.txtComPortEdit.Text = strComPort.ToString() ' Assign this port to the
edit box on the form

        MsgBox(strMsg, MsgBoxStyle.Information)
    Else
        MsgBox("Error no Bodystat device was found." & vbCrLf & "Click Search Devices
button to retry.")
    End If

    Me.Cursor = Cursors.Default
End Sub
```

4.3.1.1.1.2.4 FrmBodystatMain.btnBtSearchButton_Click Method

Executes the Bluetooth Search button action. Illustrates use of BSSearchBTDevices from the Bodystat API

C++

```
btnBtSearchButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtSearchButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtSearchButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtSearchButton.Click
```

Java

```
private btnBtSearchButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtSearchButton_Click
```

Body Source

```
Private Sub btnBtSearchButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtSearchButton.Click
    Dim iNewDevices As Integer = 0
    Dim iAuthenticatedDevices As Integer = 0

    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    ' Begin the search...
    If (BodystatAPI.BodystatAPI.BSSearchBTDevices(iNewDevices, iAuthenticatedDevices, True,
Me.Handle, iTimeout)) Then
```



```

        MsgBox(String.Format("Search complete." & vbCrLf & "Found {0} new and {1}
previously authenticated Bodystat Bluetooth devices." & vbCrLf & "Search waited
approximately {2}s.", iNewDevices, iAuthenticatedDevices, iTimeout * 1.28),
MsgBoxStyle.Information)
    Else
        MsgBox("Error searching for Bluetooth devices")
    End If

    ' If no new devices were found, lets try extending the search time for next time around
    ' We can wait up to the maximum of value of 45 - about 1 minute
    ' Obviously lengthy waits like this need to be in a background thread in a real world
example
    If (iNewDevices = 0) Then
        iTimeout = Math.Min(iTimeout * 2, 45)
    End If

    Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.5 FrmBodystatMain.btnBtStackInfoButton_Click Method

Executes the Bluetooth Get Stack Info button action. Illustrates use of BSGetBTStackInfo from the Bodystat API

C++

```
btnBtStackInfoButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtStackInfoButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtStackInfoButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtStackInfoButton.Click
```

Java

```
private btnBtStackInfoButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtStackInfoButton_Click
```

Body Source

```

Private Sub btnBtStackInfoButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtStackInfoButton.Click
    Dim strInfo As New StringBuilder(1024)
    If (BodystatAPI.BodystatAPI.BSGetBTStackInfo(strInfo, strInfo.Capacity)) Then
        MsgBox(strInfo.ToString(), MsgBoxStyle.Information)
    Else
        MsgBox("Error: No compatible Bluetooth Stack is available in this pc." & vbCrLf &
"Either it is switched off, or it is not recognised by this software." & vbCrLf & "Only
Microsoft and Widcomm stacks are supported by Bodystat's SDK.")
    End If
End Sub

```

4.3.1.1.2.6 FrmBodystatMain.btnBtStatusButton_Click Method

Executes the Bluetooth Get Status button action. Illustrates use of BSIsBTAvailable from the Bodystat API

C++

```
btnBtStatusButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtStatusButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtStatusButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtStatusButton.Click
```

Java

```
private btnBtStatusButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtStatusButton_Click
```

Body Source

```
Private Sub btnBtStatusButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtStatusButton.Click
    If (BodystatAPI.BodystatAPI.BSIsBTAvailable()) Then
        MsgBox("Supported Bluetooth radio module is available!", MsgBoxStyle.Information)
    Else
        MsgBox("No compatible Bluetooth radio module is available in this pc." & vbCrLf &
"Either it is switched off, or it is not recognised by this software." & vbCrLf & "Only
Microsoft and Widcomm stacks are supported by Bodystat's SDK.")
    End If
End Sub
```

4.3.1.1.1.2.7 FrmBodystatMain.btnBtUnauthenticateButton_Click Method

Executes the Bluetooth Unauthenticate button action. Illustrates use of BSUnAuthenticateBTDevices from the Bodystat API

C++

```
btnBtUnauthenticateButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnBtUnauthenticateButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnBtUnauthenticateButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtUnauthenticateButton.Click
```

Java

```
private btnBtUnauthenticateButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnBtUnauthenticateButton_Click
```

Body Source

```
Private Sub btnBtUnauthenticateButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnBtUnauthenticateButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    If (BodystatAPI.BodystatAPI.BSUnAuthenticateBTDevices(True, Me.Handle, iTimeout)) Then
        MsgBox("Devices authentication removed", MsgBoxStyle.Information)
    Else
        MsgBox("Error removing authentication")
    End If

    Me.Cursor = Cursors.Default
End Sub
```

4.3.1.1.1.2.8 FrmBodystatMain.btnComConnectButton_Click Method

Executes the communication Connect button action. Illustrates use of BSOpenComport from the Bodystat API

C++

```
btnComConnectButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComConnectButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComConnectButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComConnectButton.Click
```

Java

```
private btnComConnectButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComConnectButton_Click
```

Body Source

```
Private Sub btnComConnectButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComConnectButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
    If (bsErr = BSError.NoError) Then

        bsErr = BodystatAPI.BodystatAPI.BSConnect()
        If (bsErr = BSError.NoError) Then
            MsgBox("Connected to device successfully.", MsgBoxStyle.Information)
        Else
            MsgBox("Error connecting to device.")
        End If

        ' Close the com port
        BodystatAPI.BodystatAPI.BSCloseComport()
    Else
        MsgBox("Error opening com port.")
    End If

    Me.Cursor = Cursors.Default
End Sub
```

4.3.1.1.2.9 FrmBodystatMain.btnComReadCalibrationTimeButton_Click Method

Executes the communication Read Calibration Time button action. Illustrates use of BSReadCalibrationTime from the Bodystat API

C++

```
btnComReadCalibrationTimeButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadCalibrationTimeButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadCalibrationTimeButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadCalibrationTimeButton.Click
```

Java

```
private btnComReadCalibrationTimeButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadCalibrationTimeButton_Click
```

Body Source

```
Private Sub btnComReadCalibrationTimeButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadCalibrationTimeButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
    If (bsErr = BSError.NoError) Then
```

```

' Connect to the device
bsErr = BodystatAPI.BodystatAPI.BSConnect()
If (bsErr = BSError.NoError) Then

    Dim dateTime As Date

    ' Read the calibration date from the device
    bsErr = BodystatAPI.BodystatAPI.BSReadCalibrationTime(dateTime)
    If (bsErr = BSError.NoError) Then

        Dim strMsg As String
        strMsg = String.Format("Queried device successfully" & vbCrLf &
"Calibration Date/Time: {0}", dateTime.ToString())
        MsgBox(strMsg, MsgBoxStyle.Information)

    Else
        MsgBox("Error querying calibration date/time from device." & vbCrLf & "Note
that legacy devices do store their calibration date electronically.")
    End If

Else
    MsgBox("Error connecting to device.")
End If

' Close the com port
BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.10 FrmBodystatMain.btnComReadCurrentTimeButton_Click Method

Executes the communication Read Current Time button action. Illustrates use of BSReadCurrentTime from the Bodystat API

C++

```
btnComReadCurrentTimeButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadCurrentTimeButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadCurrentTimeButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadCurrentTimeButton.Click
```

Java

```
private btnComReadCurrentTimeButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadCurrentTimeButton_Click
```

Body Source

```

Private Sub btnComReadCurrentTimeButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadCurrentTimeButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
    If (bsErr = BSError.NoError) Then

        ' Connect to the device

```

```

bsErr = BodystatAPI.BodystatAPI.BSConnect()
If (bsErr = BSError.NoError) Then

    Dim dateTime As Date

    ' Read the current time from the internal clock of the device
    bsErr = BodystatAPI.BodystatAPI.BSReadCurrentTime(dateTime, -1)
    If (bsErr = BSError.NoError) Then

        Dim strMsg As String
        strMsg = String.Format("Queried device successfully" & vbCrLf & "Current
Date/Time: {0}", dateTime.ToString())
        MsgBox(strMsg, MsgBoxStyle.Information)

    Else
        MsgBox("Error querying current date/time from device." & vbCrLf & "Not all
models have an inbuilt real time clock.")
    End If

Else
    MsgBox("Error connecting to device.")
End If

' Close the com port
BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.11 FrmBodystatMain.btnComReadModelVersionButton_Click Method

Executes the communication Read Model Version button action. Illustrates use of BSReadModelVersion from the Bodystat API

C++

```
btnComReadModelVersionButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadModelVersionButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadModelVersionButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadModelVersionButton.Click
```

Java

```
private btnComReadModelVersionButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadModelVersionButton_Click
```

Body Source

```

Private Sub btnComReadModelVersionButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadModelVersionButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
    If (bsErr = BSError.NoError) Then

        ' Connect to the device
        bsErr = BodystatAPI.BodystatAPI.BSConnect()

```

```

    If (bsErr = BSError.NoError) Then

        Dim iModel As BSDeviceModel
        Dim byMajor As Byte
        Dim byMinor As Byte
        Dim byPsoc2Version As Byte
        Dim byEepromVersion As Byte
        Dim byBluetoothInfo As Byte

        ' Read Model & Version information from the device
        bsErr = BodystatAPI.BodystatAPI.BSReadModelVersion(iModel, byMajor, byMinor,
        byPsoc2Version, byEepromVersion, byBluetoothInfo)
        If (bsErr = BSError.NoError) Then

            Dim strMsg As String
            strMsg = String.Format("Queried device successfully" & vbCrLf & vbCrLf &
            "Device Type: {0}" & vbCrLf & "Firmware Version: {1}.{2}.{3}.{4}", iModel, byMajor,
            byMinor, byPsoc2Version, byEepromVersion)
            MsgBox(strMsg, MsgBoxStyle.Information)

        Else
            MsgBox("Error querying model information from device.")
        End If

    Else
        MsgBox("Error connecting to device.")
    End If

    ' Close the com port
    BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.12 FrmBodystatMain.btnComReadPrinterAddrButton_Click Method

Executes the communication Read Printer Address button action. Illustrates use of BSReadPrinterAddress from the Bodystat API

C++

```
btnComReadPrinterAddrButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadPrinterAddrButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadPrinterAddrButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadPrinterAddrButton.Click
```

Java

```
private btnComReadPrinterAddrButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadPrinterAddrButton_Click
```

Body Source

```

Private Sub btnComReadPrinterAddrButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadPrinterAddrButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
    something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError

```

```

bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
If (bsErr = BSError.NoError) Then

    ' Connect to the device
    bsErr = BodystatAPI.BodystatAPI.BSConnect()
    If (bsErr = BSError.NoError) Then

        Dim strPrinterAddr As New StringBuilder(16)

        ' Read the address of the printer the device is paired to for direct printing
        bsErr = BodystatAPI.BodystatAPI.BSReadPrinterAddress(strPrinterAddr,
strPrinterAddr.Capacity)
        If (bsErr = BSError.NoError) Then

            Dim strMsg As String
            strMsg = String.Format("Queried device successfully" & vbCrLf & "Printer
Address: {0}", strPrinterAddr.ToString())
            MsgBox(strMsg, MsgBoxStyle.Information)

        Else
            MsgBox("Error querying printer address from device." & vbCrLf & "Not all
models support direct printing.")
        End If

    Else
        MsgBox("Error connecting to device.")
    End If

    ' Close the com port
    BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.13 FrmBodystatMain.btnComReadProtocolInfoButton_Click Method

Executes the communication ReadProtocolInfo button action. Illustrates use of BSReadProtocolInfo from the Bodystat API

C++

```
btnComReadProtocolInfoButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadProtocolInfoButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadProtocolInfoButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadProtocolInfoButton.Click
```

Java

```
private btnComReadProtocolInfoButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadProtocolInfoButton_Click
```

Body Source

```

Private Sub btnComReadProtocolInfoButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadProtocolInfoButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)

```

```

If (bsErr = BSError.NoError) Then
    ' Connect to the device
    bsErr = BodystatAPI.BodystatAPI.BSConnect()
    If (bsErr = BSError.NoError) Then

        Dim byProtocolVersion As Byte
        Dim byDataVersion As Byte
        Dim byAuxInfo As Byte

        ' Read Protocol Information from the device
        bsErr = BodystatAPI.BodystatAPI.BSReadProtocolInfo(byProtocolVersion,
byDataVersion, byAuxInfo)
        If (bsErr = BSError.NoError) Then

            Dim strMsg As String
            strMsg = String.Format("Queried device successfully" & vbCrLf & "Version:
{0}" & vbCrLf & "Data Version: {1}" & vbCrLf & "Aux Info: {2}", byProtocolVersion,
byDataVersion, byAuxInfo)
            MsgBox(strMsg, MsgBoxStyle.Information)

        Else
            MsgBox("Error querying protocol information from device.")
        End If

    Else
        MsgBox("Error connecting to device.")
    End If

    ' Close the com port
    BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.14 FrmBodystatMain.btnComReadSerialNumberButton_Click Method

Executes the communication Read Serial Number button action. Illustrates use of BSReadSerialNumber from the Bodystat API

C++

```
btnComReadSerialNumberButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadSerialNumberButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadSerialNumberButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadSerialNumberButton.Click
```

Java

```
private btnComReadSerialNumberButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadSerialNumberButton_Click
```

Body Source

```

Private Sub btnComReadSerialNumberButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadSerialNumberButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device

```



```

Dim bsErr As BSError
bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
If (bsErr = BSError.NoError) Then

    ' Connect to the device
    bsErr = BodystatAPI.BodystatAPI.BSConnect()
    If (bsErr = BSError.NoError) Then

        Dim ulSerialNumber As UInteger

        ' Read the unique serial number from the device
        bsErr = BodystatAPI.BodystatAPI.BSReadSerialNumber(ulSerialNumber)
        If (bsErr = BSError.NoError) Then

            Dim strMsg As String
            strMsg = String.Format("Queried device successfully" & vbCrLf & "Serial
Number: {0}", ulSerialNumber)
            MsgBox(strMsg, MsgBoxStyle.Information)

        Else
            MsgBox("Error querying serial number from device.")
        End If

    Else
        MsgBox("Error connecting to device.")
    End If

    ' Close the com port
    BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.15 FrmBodystatMain.btnComReadStoredTestDataButton_Click Method

Executes the communication Read Stored Test Data button action. Illustrates use of BSReadStoredTestData from the Bodystat API to download data

C++

```
btnComReadStoredTestDataButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadStoredTestDataButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadStoredTestDataButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnComReadStoredTestDataButton.Click
```

Java

```
private btnComReadStoredTestDataButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadStoredTestDataButton_Click
```

Body Source

```

Private Sub btnComReadStoredTestDataButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnComReadStoredTestDataButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)

```

```

If (bsErr = BSError.NoError) Then
    ' Connect to the device
    bsErr = BodystatAPI.BodystatAPI.BSConnect()
    If (bsErr = BSError.NoError) Then

        ' Need to allocate our raw data structure with precisely 100 records
        Dim m_bsRawData As New BodystatAPI.BodystatAPI.BSRawData(99)

        ' Download test measurement data from the device
        bsErr = BodystatAPI.BodystatAPI.BSReadStoredTestData(m_bsRawData)
        If (bsErr = BSError.NoError) Then

            Dim strMsg As String
            strMsg = String.Format("Download complete" & vbCrLf & "Retrieved {0}
records" & vbCrLf & "First Test Number: {1}", m_bsRawData.iTotalNumRecs,
m_bsRawData.ulFirstTestNum)

            ' Parse the data and place it in the measurement list control
            PopulateMeasurementList(m_bsRawData)

            MsgBox(strMsg, MsgBoxStyle.Information)

        Else
            MsgBox("Error downloading stored test data from the device.")
        End If

    Else
        MsgBox("Error connecting to device.")
    End If

    ' Close the com port
    BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.16 FrmBodystatMain.btnComReadTestMeasurementButton_Click Method

Executes the communication Read Test Measurement button action. Performs an immediate test measurement on the device using the BSReadTestBodystat function from the Bodystat API

C++

```
btnComReadTestMeasurementButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComReadTestMeasurementButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComReadTestMeasurementButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadTestMeasurementButton.Click
```

Java

```
private btnComReadTestMeasurementButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComReadTestMeasurementButton_Click
```

Body Source

```

Private Sub btnComReadTestMeasurementButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComReadTestMeasurementButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

```

```

' Open a com port to the device
Dim bsErr As BSError
bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
If (bsErr = BSError.NoError) Then

    ' Connect to the device
    bsErr = BodystatAPI.BodystatAPI.BSConnect()
    If (bsErr = BSError.NoError) Then

        Dim iz5, iz50, iz100, iz200, iR50 As Integer
        Dim fX50, fPA50 As Single

        ' Perform the test measurement on the device and get the raw electrical
        bio-impedance results
        bsErr = BodystatAPI.BodystatAPI.BSReadTestBodystat(iz5, iz50, iz100, iz200,
        iR50, fX50, fPA50)
        If (bsErr = BSError.NoError) Then

            Dim strMsg As String
            strMsg = String.Format("Measurement complete" & vbCrLf & _
            "Impedance 5kHz: {0}" & vbCrLf & _
            "Impedance 50kHz: {1}" & vbCrLf & _
            "Impedance 100kHz: {2}" & vbCrLf & _
            "Impedance 200kHz: {3}" & vbCrLf & _
            "Resistance 50kHz: {4}" & vbCrLf & _
            "Reactance 50kHz: {5:F1}" & vbCrLf & _
            "Phase Angle 50kHz: {6:F1}", _
            iz5, iz50, iz100, iz200, iR50, fX50, fPA50)

            MsgBox(strMsg, MsgBoxStyle.Information)

        Else
            MsgBox("Error performing test measurement on the device.")
        End If

    Else
        MsgBox("Error connecting to device.")
    End If

    ' Close the com port
    BodystatAPI.BodystatAPI.BSCloseComport()
Else
    MsgBox("Error opening com port.")
End If

Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.2.17 FrmBodystatMain.btnComWriteCurrentTimeButton_Click Method

Executes the communication Write Current Time button action. Illustrates use of BSWriteCurrentTime from the Bodystat API

C++

```
btnComWriteCurrentTimeButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnComWriteCurrentTimeButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnComWriteCurrentTimeButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnComWriteCurrentTimeButton.Click
```

Java

```
private btnComWriteCurrentTimeButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnComWriteCurrentTimeButton_Click
```

Body Source

```

Private Sub btnComWriteCurrentTimeButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnComWriteCurrentTimeButton.Click
    Me.Cursor = Cursors.WaitCursor ' Lengthy operation, show the user we are doing
something

    Dim strComPort As New StringBuilder(10)
    strComPort.Append(txtComPortEdit.Text)

    ' Open a com port to the device
    Dim bsErr As BSError
    bsErr = BodystatAPI.BodystatAPI.BSOpenComport(strComPort, strComPort.Capacity)
    If (bsErr = BSError.NoError) Then

        ' Connect to the device
        bsErr = BodystatAPI.BodystatAPI.BSConnect()
        If (bsErr = BSError.NoError) Then

            ' Write a specific time to the device
            'Dim dtNow As DateTime = DateTime.Now().ToUniversalTime()
            'bsErr = BodystatAPI.BodystatAPI.BSWriteCurrentTimeAs(dtNow)

            ' Write the current local PC time to the internal clock of the device
            bsErr = BodystatAPI.BodystatAPI.BSWriteCurrentTime()
            If (bsErr = BSError.NoError) Then
                MsgBox("Device updated successfully" & vbCrLf & "Current date/time
updated.", MsgBoxStyle.Information)
            Else
                MsgBox("Error writing current date/time to device." & vbCrLf & "Not all
models have an inbuilt real time clock.")
            End If

        Else
            MsgBox("Error connecting to device.")
        End If

        ' Close the com port
        BodystatAPI.BodystatAPI.BSCloseComport()
    Else
        MsgBox("Error opening com port.")
    End If

    Me.Cursor = Cursors.Default
End Sub

```

4.3.1.1.1.2.18 FrmBodystatMain.btnResultsButton_Click Method

Executes the View Results button action. Opens the results form ready for more detailed analysis and normals

C++

```
btnResultsButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnResultsButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnResultsButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnResultsButton.Click
```

Java

```
private btnResultsButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnResultsButton_Click
```

Body Source

```

Private Sub btnResultsButton_Click(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles btnResultsButton.Click

```

```

Dim frmResults As New FrmBodystatResults()      ' Initialise the results dialog

' Populate the results dialog with the details of the selected measurement
' If no item is selected, we simply show the results dialog with its default values
If (lstviewTestMeasurementList.SelectedItems.Count() = 1) Then
    Dim pBSM As BodystatAPI.BodystatAPI.BSMMeasurement =
lstviewTestMeasurementList.SelectedItems.Item(0).Tag
    frmResults.UseBSMeasurement(pBSM)
Else
    MsgBox("No measurement selected, an example measurement will be shown instead.",
MsgBoxStyle.Information)
End If

' Show the results dialog
frmResults.ShowDialog()

End Sub

```

4.3.1.1.2.19 FrmBodystatMain.FrmBodystatMain_Load Method

Prepare the controls on the sample dialog

C++

```
FrmBodystatMain_Load(System.Object sender, System.EventArgs e);
```

C#

```
private FrmBodystatMain_Load(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub FrmBodystatMain_Load(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles MyBase.Load
```

Java

```
private FrmBodystatMain_Load(System.Object sender, System.EventArgs e);
```

MATLAB

```
FrmBodystatMain_Load
```

Body Source

```
Private Sub FrmBodystatMain_Load(ByVal sender As System.Object, ByVal e As
System.EventArgs) Handles MyBase.Load
```

```

' Add columns to the test measurement list ready for download of data
Me.lstviewTestMeasurementList.Columns.Add("Test No", 50)
Me.lstviewTestMeasurementList.Columns.Add("Test Date/Time", 130)
Me.lstviewTestMeasurementList.Columns.Add("Model", 140)
Me.lstviewTestMeasurementList.Columns.Add("Serial No", 75)
Me.lstviewTestMeasurementList.Columns.Add("Gender", 75)
Me.lstviewTestMeasurementList.Columns.Add("Age", 75)
Me.lstviewTestMeasurementList.Columns.Add("Height", 75)
Me.lstviewTestMeasurementList.Columns.Add("Weight", 75)
Me.lstviewTestMeasurementList.Columns.Add("Activity", 75)
Me.lstviewTestMeasurementList.Columns.Add("Waist", 75)
Me.lstviewTestMeasurementList.Columns.Add("Hip", 75)
Me.lstviewTestMeasurementList.Columns.Add("Imp 5 kHz", 75)
Me.lstviewTestMeasurementList.Columns.Add("Imp 50 kHz", 75)
Me.lstviewTestMeasurementList.Columns.Add("Imp 100 kHz", 75)
Me.lstviewTestMeasurementList.Columns.Add("Imp 200 kHz", 75)
Me.lstviewTestMeasurementList.Columns.Add("Res 50 kHz", 75)
Me.lstviewTestMeasurementList.Columns.Add("Reac 50 kHz", 75)
Me.lstviewTestMeasurementList.Columns.Add("Phase 50 kHz", 80)

```

```

iTimeout = 7      ' Initial Bluetooth timeout multiplier of 7 (equates to about 10
seconds). This is incremented upon each Bluetooth search. A maximum value of 45 is
permitted (about 1 minute).
End Sub

```

4.3.1.1.1.2.20 FrmBodystatMain.lstviewTestMeasurementList_DoubleClick Method

Executes the Measurements lists double click action. Opens the results form ready for more detailed analysis and normals

C++

```
lstviewTestMeasurementList_DoubleClick(System.Object sender, System.EventArgs e);
```

C#

```
private lstviewTestMeasurementList_DoubleClick(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub lstviewTestMeasurementList_DoubleClick(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles lstviewTestMeasurementList.DoubleClick
```

Java

```
private lstviewTestMeasurementList_DoubleClick(System.Object sender, System.EventArgs e);
```

MATLAB

```
lstviewTestMeasurementList_DoubleClick
```

Body Source

```
Private Sub lstviewTestMeasurementList_DoubleClick(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles lstviewTestMeasurementList.DoubleClick
    btnResultsButton_Click(sender, e) ' Display the selected result
End Sub
```

4.3.1.1.1.2.21 FrmBodystatMain.PopulateMeasurementList Method

Populate the measurement list control with test measurement data from the API. Illustrates iteration through a BSRawData structure, working with data in a BSMeasurement structure and helper functions from the Bodystat API like BSGetDeviceModelName

C++

```
PopulateMeasurementList(BodystatAPI.BodystatAPI.BSRawData ^ pRawData);
```

C#

```
public PopulateMeasurementList(ref BodystatAPI.BodystatAPI.BSRawData pRawData);
```

Visual Basic

```
Public Sub PopulateMeasurementList(ByRef pRawData As BodystatAPI.BodystatAPI.BSRawData)
```

Java

```
public PopulateMeasurementList(BodystatAPI.BodystatAPI.BSRawData pRawData);
```

MATLAB

```
PopulateMeasurementList
```

Parameters

Parameters	Description
ByRef pRawData As BodystatAPI.BodystatAPI.BSRawData	[in] RawData to be parsed.

Body Source

```
Public Sub PopulateMeasurementList(ByRef pRawData As BodystatAPI.BodystatAPI.BSRawData)

    ' Clear any previous items in the list
    Me.lstviewTestMeasurementList.Items.Clear()

    Dim strItem As String
    Dim strBuilder As StringBuilder

    Dim lvi As System.Windows.Forms.ListViewItem
    Dim i As Integer

    ' Process each record in the structure
```

```

For i = 0 To pRawData.iTotalNumRecs - 1

    ' Test no
    strItem = String.Format("{0}", i + pRawData.ulFirstTestNum)
    lvi = lstviewTestMeasurementList.Items.Insert(i, strItem)
    lvi.Tag = pRawData.record(i)      ' Assign this BSMeasurement to the listview item
(so we can reference it later)

    ' Test date/time
    lvi.SubItems.Add(pRawData.record(i).tTestDate.ToString())

    ' Device model
    strBuilder = New StringBuilder(32)
    BodystatAPI.BodystatAPI.BSGetDeviceModelName(pRawData.record(i).iDeviceModel,
strBuilder, strBuilder.Capacity)
    lvi.SubItems.Add(strBuilder.ToString())

    ' Device serial number
    strItem = String.Format("{0}", pRawData.record(i).ulDeviceSerialNumber)
    lvi.SubItems.Add(strItem)

    ' Gender
    lvi.SubItems.Add(If(pRawData.record(i).iGender = BSGender.BSMale, "Male",
"Female"))

    ' Age
    strItem = String.Format("{0}", pRawData.record(i).iAge)
    lvi.SubItems.Add(strItem)

    ' Height
    strItem = String.Format("{0} cm", pRawData.record(i).iHeight)
    lvi.SubItems.Add(strItem)

    ' Weight
    strItem = String.Format("{0:F1} kg", pRawData.record(i).fWeight)
    lvi.SubItems.Add(strItem)

    ' Activity
    Select Case pRawData.record(i).iActivity
        Case 1
            strItem = "Very Low"
        Case 2
            strItem = "Low/Medium"
        Case 3
            strItem = "Medium"
        Case 4
            strItem = "Medium/High"
        Case 5
            strItem = "Very High"
        Case Else
            strItem = "Unknown"
    End Select
    lvi.SubItems.Add(strItem)

    ' Waist
    strItem = String.Format("{0} cm", pRawData.record(i).iWaist)
    lvi.SubItems.Add(strItem)

    ' Hip
    strItem = String.Format("{0} cm", pRawData.record(i).iHip)
    lvi.SubItems.Add(strItem)

    ' Z5 (Impedance at 5 kHz)
    strItem = String.Format("{0} ohm", pRawData.record(i).iZ_5kHz)
    lvi.SubItems.Add(strItem)

    ' Z50 (Impedance at 50 kHz)
    strItem = String.Format("{0} ohm", pRawData.record(i).iZ_50kHz)
    lvi.SubItems.Add(strItem)

    ' Z100 (Impedance at 100 kHz)
    strItem = String.Format("{0} ohm", pRawData.record(i).iZ_100kHz)

```

```
lvi.SubItems.Add(strItem)

' Z200 (Impedance at 200 kHz)
strItem = String.Format("{0} ohm", pRawData.record(i).iZ_200kHz)
lvi.SubItems.Add(strItem)

' R50 (Resistance at 50 kHz)
strItem = String.Format("{0} ohm", pRawData.record(i).iR_50kHz)
lvi.SubItems.Add(strItem)

' X50 (Reactance at 50 kHz)
strItem = String.Format("{0:F1} ohm", pRawData.record(i).fX_50kHz)
lvi.SubItems.Add(strItem)


' PA50 (Phase Angle at 50 kHz)
strItem = String.Format("{0:F1} deg", pRawData.record(i).fPA_50kHz)
lvi.SubItems.Add(strItem)

Next i
End Sub
```

4.3.1.1.2 FrmBodystatResults Class

Results sample dialog. Displays the various parameters relating to a measurement (allowing direct input by the user), then utilises the Bodystat API to calculate detailed results and normals.

Class Hierarchy

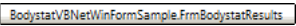

C++

```
class FrmBodystatResults;
```

C#

```
public class FrmBodystatResults;
```

Visual Basic

```
Public Class FrmBodystatResults
```

Java

```
public class FrmBodystatResults;
```

MATLAB

```
FrmBodystatResults
```

File






```
FrmBodystatResults.vb
```

FrmBodystatResults Fields

	Name	Description
	bsmMeasurement (🔗 see page 238)	Base the results dialog on this downloaded measurement
	bUseActualMeasurment (🔗 see page 238)	Determines whether the dialog should use the measurement values or sample values
	dtTestDate (🔗 see page 238)	Test date/time, not presently shown on the form (but may be passed to clipboard data)

FrmBodystatResults Methods

	Name	Description
	btnCalculateResultsButton_Click (🔗 see page 239)	Executes the Calculate Results button action. Uses the Bodystat API to calculate the results and normals for this test measurement.
	btnCopytoclipboardButton_Click (🔗 see page 240)	Implements the Copy To Clipboard button action. Copies the various parameters relating to a measurement directly to the clipboard

	FrmBodystatResults_Load (see page 242)	Initialise the dialog controls with their content
	InitialiseFromMeasurement (see page 242)	Initialise the dialog controls with a values from an actual measurement result
	InitialiseSampleMeasurement (see page 243)	Initialise the dialog controls with a sample measurement result
	PopulateListCtrl (see page 244)	Populate the list control with results and normals.
	UseBSMeasurement (see page 247)	Allows the caller to specify a downloaded measurement structure to be use when initialising the dialog. Must be called before the form is loaded (displayed)

4.3.1.1.2.1 FrmBodystatResults Fields

4.3.1.1.2.1.1 FrmBodystatResults.bsmMeasurement Field

C++

```
BodystatAPI.BodystatAPI.BSMeasurement bsmMeasurement;
```

C#

```
private BodystatAPI.BodystatAPI.BSMeasurement bsmMeasurement;
```

Visual Basic

```
Private bsmMeasurement As BodystatAPI.BodystatAPI.BSMeasurement
```

Java

```
private BodystatAPI.BodystatAPI.BSMeasurement bsmMeasurement;
```

MATLAB

```
bsmMeasurement
```

Description

Base the results dialog on this downloaded measurement

4.3.1.1.2.1.2 FrmBodystatResults.bUseActualMeasurment Field

C++

```
Boolean bUseActualMeasurment = False;
```

C#

```
private Boolean bUseActualMeasurment = False;
```

Visual Basic

```
Private bUseActualMeasurment As Boolean = False
```

Java

```
private Boolean bUseActualMeasurment = False;
```

MATLAB

```
bUseActualMeasurment
```

Description

Determines whether the dialog should use the measurement values or sample values

4.3.1.1.2.1.3 FrmBodystatResults.dtTestDate Field

C++

```
DateTime dtTestDate;
```

C#

```
private DateTime dtTestDate;
```

Visual Basic

```
Private dtTestDate As DateTime
```

Java

```
private DateTime dtTestDate;
```

MATLAB

```
dtTestDate
```

Description

Test date/time, not presently shown on the form (but may be passed to clipboard data)

4.3.1.1.2.2 FrmBodystatResults Methods**4.3.1.1.2.2.1 FrmBodystatResults.btnCalculateResultsButton_Click Method**

Executes the Calculate Results button action. Uses the Bodystat API to calculate the results and normals for this test measurement.

C++

```
btnCalculateResultsButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnCalculateResultsButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnCalculateResultsButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnCalculateResultsButton.Click
```

Java

```
private btnCalculateResultsButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnCalculateResultsButton_Click
```

Body Source

```
Private Sub btnCalculateResultsButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnCalculateResultsButton.Click
```

```
    lstviewResultsList.Items.Clear()           ' Clear the results list

    If Me.cboModelCombo.SelectedIndex = BodystatAPI.BSDeviceModel.BSQuadScanTouch Or
Me.cboModelCombo.SelectedIndex = BodystatAPI.BSDeviceModel.BSMultiScanTouch Then
        MsgBox("The selected model is not supported by this SDK sample.",
MsgBoxStyle.Information)
        Exit Sub
    End If

    '
    ' Prepare input structure for the Bodystat SDK (using member values from form data)
    ' More likely this would be obtained directly from a BSMeasurement downloaded from the
unit
    '
    Dim bm As BodystatAPI.BodystatAPI.BSMeasurement
    bm.iDeviceModel = Me.cboModelCombo.SelectedIndex           ' Use our combo selection for
Bodystat model type (its a direct mapping)
    bm.tTestDate = Me.dtTestDate                               ' Not relevant for performing
calculations, but might be copied to clipboard later
    bm.iGender = IIf(Me.rdoGenderMaleRadio.Checked, BodystatAPI.BSGender.BSMale,
BodystatAPI.BSGender.BSFemale)                               ' Convert our gender radio buttons into Bodystat
gender type
    bm.iAge = Me.txtAgeEdit.Text                               ' Age in years
    bm.iHeight = Me.txtHeightEdit.Text                         ' Height in centimeters
    bm.fWeight = Me.txtWeightEdit.Text                         ' Weight in kilograms
    bm.iActivity = Me.cboActivityCombo.SelectedIndex           ' Physical activity level (range 1
```

```

to 5) 1=very low 2=low/medium 3=medium 4=medium/high 5=very high
    bm.iWaist = Me.txtWaistEdit.Text           ' Waist in centimeters
    bm.iHip = Me.txtHipEdit.Text               ' Hip in centimeters
    bm.iZ_5kHz = Me.txtImp5kEdit.Text          ' Impedance 5 KHz (in ohms)
(MDD/Quad only)
    bm.iZ_50kHz = Me.txtImp50kEdit.Text        ' Impedance 50 KHz (in ohms)
    bm.iZ_100kHz = Me.txtImp100kEdit.Text      ' Impedance 100 KHz (in ohms) (Quad
only)
    bm.iZ_200kHz = Me.txtImp200kEdit.Text      ' Impedance 200 KHz (in ohms) (Quad
only)
    bm.iR_50kHz = Me.txtRes50kEdit.Text        ' Resistance 50 kHz (in ohms)
(MDD/Quad only)
    bm.fX_50kHz = Me.txtReac50kEdit.Text       ' Reactance 50 kHz (in ohms)
(MDD/Quad only)
    bm.fPA_50kHz = Me.txtPhaseAngle50kEdit.Text ' Phase angle 50 kHz (in degrees)
(MDD/Quad only)
    bm.iFrequencies = 0                       ' MST frequencies not supported by
this sample

    ' Choose analysis mode based upon selected analysis checkboxes
    ' These options might normally be stored as program setting rather than asked of the
user each time (like the installation mode of the device)
    ' The analysis mode impacts the equations used behind some results in certain scenarios
    Dim iAnalysisMode As BodystatAPI.BSAnalysisMode
    If (Me.chkModeBodyCompCheck.Checked And Me.chkModeHydrationCheck.Checked) Then
        iAnalysisMode = BodystatAPI.BSAnalysisMode.BSAnalysisModeBoth
    ElseIf (chkModeBodyCompCheck.Checked) Then
        iAnalysisMode = BodystatAPI.BSAnalysisMode.BSAnalysisModeBC
    ElseIf (chkModeHydrationCheck.Checked) Then
        iAnalysisMode = BodystatAPI.BSAnalysisMode.BSAnalysisModeHN
    Else
        iAnalysisMode = BodystatAPI.BSAnalysisMode.BSAnalysisModeNone
    End If

    '
    ' Get the results from the Bodystat SDK
    Dim bseStatus As BodystatAPI.BSError
    Dim br As BodystatAPI.BodystatAPI.BSResults
    bseStatus = BodystatAPI.BodystatAPI.BSCalculateResults(bm, br, iAnalysisMode)
Results stored in br structure

    ' Check for any errors from the calculation function
    If (bseStatus <> BodystatAPI.BSError.NoError) Then
        MsgBox(String.Format("Error calculating results." & vbCrLf & "Error Code: {0}",
bseStatus), MsgBoxStyle.Exclamation)
        Return ' Nothing to display
    End If

    '
    ' Get the normal ranges from the Bodystat SDK
    '
    Dim bn As BodystatAPI.BodystatAPI.BSNormals
    bseStatus = BodystatAPI.BodystatAPI.BSCalculateNormals(bm, br, bn)
in bn structure

    ' Check for any errors from the normals function
    If (bseStatus <> BodystatAPI.BSError.NoError) Then
        MsgBox(String.Format("Error calculating normals." & vbCrLf & "Error Code: {0}",
bseStatus), MsgBoxStyle.Exclamation)
        Return ' Nothing to display
    End If

    ' Populate the list control with the calculated results and calculated normals
    PopulateListCtrl(bm, br, bn)

End Sub

```

4.3.1.1.2.2.2 FrmBodystatResults.btnCopytoclipboardButton_Click Method

Implements the Copy To Clipboard button action. Copies the various parameters relating to a measurement directly to the clipboard

C++

```
btnCopytoclipboardButton_Click(System.Object sender, System.EventArgs e);
```

C#

```
private btnCopytoclipboardButton_Click(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub btnCopytoclipboardButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnCopytoclipboardButton.Click
```

Java

```
private btnCopytoclipboardButton_Click(System.Object sender, System.EventArgs e);
```

MATLAB

```
btnCopytoclipboardButton_Click
```

Body Source

```
Private Sub btnCopytoclipboardButton_Click(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles btnCopytoclipboardButton.Click
```

```
Dim strBuffer As New StringBuilder()
```

```
' Output the measurement input parameters to a buffer ready for the clipboard
strBuffer.Append("Bodystat Data" & vbCrLf)
strBuffer.Append("-----" & vbCrLf)
```

```
Dim strModel As New StringBuilder(32)
BodystatAPI.BodystatAPI.BSGetDeviceModelName(cboModelCombo.SelectedIndex, strModel, strModel.Capacity)
strBuffer.AppendFormat("Model:" & vbTab & "{0}" & vbCrLf, strModel.ToString)

strBuffer.AppendFormat("Test Date/Time:" & vbTab & "{0}" & vbCrLf, dtTestDate.ToString())
strBuffer.AppendFormat("Gender:" & vbTab & "{0}" & vbCrLf, IIf(rdoGenderMaleRadio.Checked, "Male", "Female"))
strBuffer.AppendFormat("Age:" & vbTab & "{0}" & vbTab & "years" & vbCrLf, txtAgeEdit.Text)
strBuffer.AppendFormat("Height:" & vbTab & "{0}" & vbTab & "cm" & vbCrLf, txtHeightEdit.Text)
strBuffer.AppendFormat("Weight:" & vbTab & "{0:F1}" & vbTab & "kg" & vbCrLf, txtWeightEdit.Text)
strBuffer.AppendFormat("Activity Level:" & vbTab & "{0}" & vbCrLf, cboActivityCombo.SelectedItem.ToString())
strBuffer.AppendFormat("Waist:" & vbTab & "{0}" & vbTab & "cm" & vbCrLf, txtWaistEdit.Text)
strBuffer.AppendFormat("Hip:" & vbTab & "{0}" & vbTab & "cm" & vbCrLf, txtHipEdit.Text)

strBuffer.AppendFormat("Body Composition Mode:" & vbTab & "{0}" & vbCrLf, IIf(chkModeBodyCompCheck.Checked, "Enabled", "Disabled"))
strBuffer.AppendFormat("Hydration Mode:" & vbTab & "{0}" & vbCrLf, IIf(chkModeHydrationCheck.Checked, "Enabled", "Disabled"))

strBuffer.AppendFormat("Impedance 5 kHz:" & vbTab & "{0}" & vbTab & "ohm" & vbCrLf, txtImp5kEdit.Text)
strBuffer.AppendFormat("Impedance 50 kHz:" & vbTab & "{0}" & vbTab & "ohm" & vbCrLf, txtImp50kEdit.Text)
strBuffer.AppendFormat("Impedance 100 kHz:" & vbTab & "{0}" & vbTab & "ohm" & vbCrLf, txtImp100kEdit.Text)
strBuffer.AppendFormat("Impedance 200 kHz:" & vbTab & "{0}" & vbTab & "ohm" & vbCrLf, txtImp200kEdit.Text)
strBuffer.AppendFormat("Resistance 50 kHz:" & vbTab & "{0}" & vbTab & "ohm" & vbCrLf, txtRes50kEdit.Text)
strBuffer.AppendFormat("Reactance 50 kHz:" & vbTab & "{0:F1}" & vbTab & "ohm" & vbCrLf, txtReac50kEdit.Text)
strBuffer.AppendFormat("Phase Angle 50 kHz:" & vbTab & "{0:F1}" & vbTab & "°" & vbCrLf, txtPhaseAngle50kEdit.Text)

strBuffer.AppendLine()
```

```

    ' Copy the buffer to the clipboard
    Clipboard.SetText(strBuffer.ToString())
End Sub

```

4.3.1.1.2.2.3 FrmBodystatResults.FrmBodystatResults_Load Method

Initialise the dialog controls with their content

C++

```
FrmBodystatResults_Load(System.Object sender, System.EventArgs e);
```

C#

```
private FrmBodystatResults_Load(System.Object sender, System.EventArgs e);
```

Visual Basic

```
Private Sub FrmBodystatResults_Load(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles MyBase.Load
```

Java

```
private FrmBodystatResults_Load(System.Object sender, System.EventArgs e);
```

MATLAB

```
FrmBodystatResults_Load
```

Body Source

```

Private Sub FrmBodystatResults_Load(ByVal sender As System.Object, ByVal e As System.EventArgs) Handles MyBase.Load
    '
    ' Populate model combo with all possible models
    '
    Me.cboModelCombo.Items.Add("Unknown")
    Dim strModel As New StringBuilder(32)
    Dim strComboItem As String
    Dim bsErr As BodystatAPI.BSError = BodystatAPI.BSError.NoError
    Dim iModel As Integer = 1
    Do
        bsErr = BodystatAPI.BodystatAPI.BSGetDeviceModelName(iModel, strModel, strModel.Capacity)
        If (bsErr = BodystatAPI.BSError.NoError) Then
            strComboItem = String.Format("{0} (ID:{1})", strModel.ToString(), iModel)
            Append model id number to the string (for illustration and to make it unique for the combo)
            Me.cboModelCombo.Items.Add(strComboItem)

            End If
            iModel += 1
        Loop While (bsErr = BodystatAPI.BSError.NoError) ' Stop when we run out of valid model ids

        ' Setup the results list control with three columns ready for the results and normal values
        lstviewResultsList.Columns.Add("Result", 250)
        lstviewResultsList.Columns.Add("Value", 110)
        lstviewResultsList.Columns.Add("Normal", 140)

        ' Populate the measurement control values based on the downloaded measurement or with sample data as appropriate
        If (bUseActualMeasurment) Then
            InitialiseFromMeasurement()
        Else
            InitialiseSampleMeasurement()
        End If
    End Sub

```

4.3.1.1.2.2.4 FrmBodystatResults.InitialiseFromMeasurement Method

Initialise the dialog controls with a values from an actual measurement result

C++

```
InitialiseFromMeasurement();
```

C#

```
public InitialiseFromMeasurement();
```

Visual Basic

```
Public Sub InitialiseFromMeasurement()
```

Java

```
public InitialiseFromMeasurement();
```

MATLAB

```
InitialiseFromMeasurement
```

Body Source

```
Public Sub InitialiseFromMeasurement()

    Me.cboModelCombo.SelectedIndex = bsmMeasurement.iDeviceModel
    'Me.dtTestDate = bsmMeasurement.tTestDate      ' Not shown on dialog presently
    Me.rdoGenderMaleRadio.Checked = bsmMeasurement.iGender = BodystatAPI.BSGender.BSMale
    Me.rdoGenderFemaleRadio.Checked = bsmMeasurement.iGender = BodystatAPI.BSGender.BSFemale
    Me.txtAgeEdit.Text = bsmMeasurement.iAge
    Me.txtHeightEdit.Text = bsmMeasurement.iHeight
    Me.txtWeightEdit.Text = bsmMeasurement.fWeight
    Me.cboActivityCombo.SelectedIndex = bsmMeasurement.iActivity
    Me.txtWaistEdit.Text = bsmMeasurement.iWaist
    Me.txtHipEdit.Text = bsmMeasurement.iHip

    Me.chkModeBodyCompCheck.Checked = True
    Me.chkModeHydrationCheck.Checked = True

    Me.txtImp5kEdit.Text = bsmMeasurement.iZ_5kHz
    Me.txtImp50kEdit.Text = bsmMeasurement.iZ_50kHz
    Me.txtImp100kEdit.Text = bsmMeasurement.iZ_100kHz
    Me.txtImp200kEdit.Text = bsmMeasurement.iZ_200kHz
    Me.txtRes50kEdit.Text = bsmMeasurement.iR_50kHz
    Me.txtReac50kEdit.Text = bsmMeasurement.fX_50kHz
    Me.txtPhaseAngle50kEdit.Text = bsmMeasurement.fPA_50kHz

End Sub
```

4.3.1.1.2.2.5 FrmBodystatResults.InitialiseSampleMeasurement Method

Initialise the dialog controls with a sample measurement result

C++

```
InitialiseSampleMeasurement();
```

C#

```
public InitialiseSampleMeasurement();
```

Visual Basic

```
Public Sub InitialiseSampleMeasurement()
```

Java

```
public InitialiseSampleMeasurement();
```

MATLAB

```
InitialiseSampleMeasurement
```

Body Source

```
Public Sub InitialiseSampleMeasurement()

    ' Setup the form with some sensible default values
    Me.cboModelCombo.SelectedIndex = BodystatAPI.BSDeviceModel.BSQuadScan_BT
```

```

    'Me.dtTestDate = DateTime.Now()           ' Not shown on dialog presently
    Me.rdoGenderMaleRadio.Checked = True
    Me.rdoGenderFemaleRadio.Checked = False
    Me.txtAgeEdit.Text = "30"
    Me.txtHeightEdit.Text = "165"
    Me.txtWeightEdit.Text = "85.0"
    Me.cboActivityCombo.SelectedIndex = 3
    Me.txtWaistEdit.Text = "85"
    Me.txtHipEdit.Text = "90"

    Me.chkModeBodyCompCheck.Checked = True
    Me.chkModeHydrationCheck.Checked = True

    Me.txtImp5kEdit.Text = "500"
    Me.txtImp50kEdit.Text = "500"
    Me.txtImp100kEdit.Text = "500"
    Me.txtImp200kEdit.Text = "500"
    Me.txtRes50kEdit.Text = "500"
    Me.txtReac50kEdit.Text = "52.0"
    Me.txtPhaseAngle50kEdit.Text = "6.1"
End Sub

```

4.3.1.1.2.2.6 FrmBodystatResults.PopulateListCtrl Method

Populate the list control with results and normals.

C++

```

PopulateListCtrl(BodystatAPI.BodystatAPI.BSMeasurement ^ bm,
BodystatAPI.BodystatAPI.BSResults ^ br, BodystatAPI.BodystatAPI.BSNormals ^ bn);

```

C#

```

public PopulateListCtrl(ref BodystatAPI.BodystatAPI.BSMeasurement bm, ref
BodystatAPI.BodystatAPI.BSResults br, ref BodystatAPI.BodystatAPI.BSNormals bn);

```

Visual Basic

```

Public Sub PopulateListCtrl(ByRef bm As BodystatAPI.BodystatAPI.BSMeasurement, ByRef br As
BodystatAPI.BodystatAPI.BSResults, ByRef bn As BodystatAPI.BodystatAPI.BSNormals)

```

Java

```

public PopulateListCtrl(BodystatAPI.BodystatAPI.BSMeasurement bm,
BodystatAPI.BodystatAPI.BSResults br, BodystatAPI.BodystatAPI.BSNormals bn);

```

MATLAB

```

PopulateListCtrl

```

Parameters

Parameters	Description
ByRef bm As BodystatAPI.BodystatAPI.BSMeasurement	The measurement parameters for this test.
ByRef br As BodystatAPI.BodystatAPI.BSResults	The calculated results for this test from the SDK.
ByRef bn As BodystatAPI.BodystatAPI.BSNormals	The calculated normals for this test from the SDK.

Body Source

```

Public Sub PopulateListCtrl(ByRef bm As BodystatAPI.BodystatAPI.BSMeasurement, ByRef br As
BodystatAPI.BodystatAPI.BSResults, ByRef bn As BodystatAPI.BodystatAPI.BSNormals)

    Dim lvi As ListViewItem

    ' Fat %
    lvi = lstviewResultsList.Items.Add("Fat percentage")
    lvi.SubItems.Add(String.Format("{0:F1}", br.fFatPerc))
    lvi.SubItems.Add(String.Format("{0} - {1}", bn.iFatPerc_L, bn.iFatPerc_H))

    ' Fat kg
    lvi = lstviewResultsList.Items.Add("Fat (in kg)")
    lvi.SubItems.Add(String.Format("{0:F1}", br.fFatKg))
    lvi.SubItems.Add(String.Format("{0} - {1}", bn.iFatKg_L, bn.iFatKg_H))

```

```

' Lean %
lvi = lstviewResultsList.Items.Add("Lean percentage")
lvi.SubItems.Add(String.Format("{0:F1}", br.fLeanPerc))
lvi.SubItems.Add(String.Format("{0} - {1}", bn.iLeanPerc_L, bn.iLeanPerc_H))

' Lean kg
lvi = lstviewResultsList.Items.Add("Lean (in kg)")
lvi.SubItems.Add(String.Format("{0:F1}", br.fLeanKg))
lvi.SubItems.Add(String.Format("{0} - {1}", bn.iLeanKg_L, bn.iLeanKg_H))

' Total Weight
lvi = lstviewResultsList.Items.Add("Total Weight (in kg)")
lvi.SubItems.Add(String.Format("{0:F1}", br.fTotalWeight))
' Total weight normal range is calculated using composition or BMI (whichever is the
greater). We indicate which method was used in brackets after the value.
lvi.SubItems.Add(String.Format("{0} - {1} ({2})", bn.iTotalWeight_L, bn.iTotalWeight_H,
IIf(bn.iTotalWeightMethod = 0, "Composition", "BMI")))

' Dry lean weight
lvi = lstviewResultsList.Items.Add("Dry Lean Weight (in kg)")
lvi.SubItems.Add(String.Format("{0:F1}", br.fDryLW))
lvi.SubItems.Add("N/A")

' TBW %
lvi = lstviewResultsList.Items.Add("Total Body Water percentage")
lvi.SubItems.Add(String.Format("{0:F1}", br.fTBWPerc))
lvi.SubItems.Add(String.Format("{0} - {1}", bn.iTBWPerc_L, bn.iTBWPerc_H))

' TBW lt
lvi = lstviewResultsList.Items.Add("Total Body Water (in litres)")
lvi.SubItems.Add(String.Format("{0:F1}", br.fTBW))
lvi.SubItems.Add(String.Format("{0} - {1}", bn.iTBW_L, bn.iTBW_H))

' ECW %
lvi = lstviewResultsList.Items.Add("Extra Cellular Water percentage")
lvi.SubItems.Add(IIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F1}", br.fECWPerc), "(MDD/Quad
Only)"))
lvi.SubItems.Add(IIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0}", bn.iECWPerc_Norm), "(MDD/Quad
Only)"))

' ECW lt
lvi = lstviewResultsList.Items.Add("Extra Cellular Water (in litres)")
lvi.SubItems.Add(IIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F1}", br.fECW), "(MDD/Quad
Only)"))
lvi.SubItems.Add("N/A")

' ICW %
lvi = lstviewResultsList.Items.Add("Intra Cellular Water percentage")
lvi.SubItems.Add(IIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F1}", br.fICWPerc), "(MDD/Quad
Only)"))
lvi.SubItems.Add(IIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0}", bn.iICWPerc_Norm), "(MDD/Quad
Only)"))

' ICW lt
lvi = lstviewResultsList.Items.Add("Intra Cellular Water (in litres)")
lvi.SubItems.Add(IIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or

```



```

BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F1}", br.fICW), "(MDD/Quad
Only)")
    lvi.SubItems.Add("N/A")

    ' Body Cell Mass
    lvi = lstviewResultsList.Items.Add("Body Cell Mass")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family, String.Format("{0:F1}", br.fBCM), "(Quad
Only)")
    lvi.SubItems.Add("N/A")

    ' 3rd space water
    lvi = lstviewResultsList.Items.Add("3rd space water (in litres)")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family, String.Format("{0:F1}", br.fThirdSpace),
"(Quad Only)")
    lvi.SubItems.Add("N/A")

    ' Nutrition index
    lvi = lstviewResultsList.Items.Add("Nutrition index")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family, String.Format("{0:F2}", br.fNutrition),
"(Quad Only)")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family, String.Format("{0:F2}", bn.fNutrition_Norm),
"(Quad Only)")

    ' Prediction marker - formerly known as Illness Marker (TM)
    lvi = lstviewResultsList.Items.Add("Prediction Marker (TM)")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family, String.Format("{0:F3}", br.fIllness), "(Quad
Only)")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family, String.Format("{0:F3} - {1:F3}",
bn.fIllness_L, bn.fIllness_H), "(Quad Only)")

    ' BMR kcal
    lvi = lstviewResultsList.Items.Add("Basal Metabolic Rate (in kcal)")
    lvi.SubItems.Add(String.Format("{0:F0}", br.fBMR))
    lvi.SubItems.Add("N/A")

    ' BMR kcal/kg
    lvi = lstviewResultsList.Items.Add("Basal Metabolic Rate per kilogram (in kcal/kg)")
    lvi.SubItems.Add(String.Format("{0:F1}", br.fBMRkg))
    lvi.SubItems.Add("N/A")

    ' EAR
    lvi = lstviewResultsList.Items.Add("Estimated Average Requirement (in kcal)")
    lvi.SubItems.Add(String.Format("{0:F0}", br.fEstAvg))
    lvi.SubItems.Add("N/A")

    ' BMI
    lvi = lstviewResultsList.Items.Add("Body Mass Index (BMI)")
    lvi.SubItems.Add(String.Format("{0:F1}", br.fBMI))
    lvi.SubItems.Add(String.Format("{0} - {1}", bn.iBMI_L, bn.iBMI_H))

    ' BFMI
    lvi = lstviewResultsList.Items.Add("Body Fat Mass Index (BFMI)")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F1}", br.fBFMI), "(MDD/Quad
Only)")
    lvi.SubItems.Add(IIIf(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0} - {1}", bn.iBFMI_L,
bn.iBFMI_H), "(MDD/Quad Only)")

    ' FFMI
    lvi = lstviewResultsList.Items.Add("Fat Free Mass Index (FFMI)")

```

```

    lvi.SubItems.Add(If(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F1}", br.fFFMI), "(MDD/Quad
Only)"))
    lvi.SubItems.Add(If(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSQuadScan_Family Or
BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0} - {1}", bn.iFFMI_L,
bn.iFFMI_H), "(MDD/Quad Only)"))

    ' Waist/Hip
    lvi = lstviewResultsList.Items.Add("Waist/Hip ratio")
    lvi.SubItems.Add(String.Format("{0:F2}", br.fWaistHip))
    lvi.SubItems.Add(String.Format("{0:F2}", bn.fWaistHip_Norm))

    ' Wellness marker
    lvi = lstviewResultsList.Items.Add("Wellness Marker (TM)")
    lvi.SubItems.Add(If(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F3}", br.fWellness), "(MDD
Only)"))
    lvi.SubItems.Add(If(BodystatAPI.BodystatAPI.BSGetDeviceFamily(bm.iDeviceModel) =
BodystatAPI.BSDeviceFamily.BSMDD_Family, String.Format("{0:F3} - {0:F3}", bn.fWellness_L,
bn.fWellness_H), "(MDD Only)"))
End Sub

```

4.3.1.1.2.2.7 FrmBodystatResults.UseBSMeasurement Method

Allows the caller to specify a downloaded measurement structure to be use when initialising the dialog. Must be called before the form is loaded (displayed)

C++

```
UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm);
```

C#

```
public UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm);
```

Visual Basic

```
Public Sub UseBSMeasurement(ByVal bsm As BodystatAPI.BodystatAPI.BSMeasurement)
```

Java

```
public UseBSMeasurement(BodystatAPI.BodystatAPI.BSMeasurement bsm);
```

MATLAB

```
UseBSMeasurement
```

Body Source

```

Public Sub UseBSMeasurement(ByVal bsm As BodystatAPI.BodystatAPI.BSMeasurement)
    bsmMeasurement = bsm
    bUseActualMeasurment = True
End Sub

```

Index

B

Bodystat 2

Bodystat API 1, 2, 57

Functions 6

Structs, Records, Enums 29

Variables 52

Bodystat API - Main 1

Files 54

Macros 53

Bodystat C# WinForms Sample 204

Files 215

Bodystat MFC Sample 155

Classes 155

Bodystat Visual Basic .NET Sample 216

Bodystat::BS_RAWDATA_ARRAYSIZE 52

Bodystat::BS_RAWDATA_ARRAYSIZE variable 52

Bodystat::BSAnalysisMode 30

Bodystat::BSAnalysisMode enumeration 30

Bodystat::BSAuthenticateBTDevice 8

Bodystat::BSAuthenticateBTDevice function 8

Bodystat::BSAutoSetupBT 8

Bodystat::BSAutoSetupBT function 8

Bodystat::BSCalculateNormals 10

Bodystat::BSCalculateNormals function 10

Bodystat::BSCalculateResults 10

Bodystat::BSCalculateResults function 10

Bodystat::BSCloseComport 11

Bodystat::BSCloseComport function 11

Bodystat::BSConnect 11

Bodystat::BSConnect function 11

Bodystat::BSDeviceFamily 30

Bodystat::BSDeviceFamily enumeration 30

Bodystat::BSDeviceModel 32

Bodystat::BSDeviceModel enumeration 32

Bodystat::BSError 33

Bodystat::BSError enumeration 33

Bodystat::BSGender 35

Bodystat::BSGender enumeration 35

Bodystat::BSGetBTBodystatDevice 12

Bodystat::BSGetBTBodystatDevice function 12

Bodystat::BSGetBTStackInfo 13

Bodystat::BSGetBTStackInfo function 13

Bodystat::BSGetDeviceFamily 13

Bodystat::BSGetDeviceFamily function 13

Bodystat::BSGetDeviceModelName 14

Bodystat::BSGetDeviceModelName function 14

Bodystat::BSGetDeviceModelNameShort 14

Bodystat::BSGetDeviceModelNameShort function 14

Bodystat::BSGetSdkLibraryVersion 15

Bodystat::BSGetSdkLibraryVersion function 15

Bodystat::BSIsBTAvailable 16

Bodystat::BSIsBTAvailable function 16

Bodystat::BSMeasurement 36

Bodystat::BSMeasurement structure 36

Bodystat::BSNormals 38

Bodystat::BSNormals structure 38

Bodystat::BSOpenComport 16

Bodystat::BSOpenComport function 16

Bodystat::BSRawData 40

Bodystat::BSRawData structure 40

Bodystat::BSReadCalibrationTime 17

Bodystat::BSReadCalibrationTime function 17

Bodystat::BSReadCurrentTime 17

Bodystat::BSReadCurrentTime function 17

Bodystat::BSReadModelVersion 18

Bodystat::BSReadModelVersion function 18

Bodystat::BSReadPrinterAddress 19

Bodystat::BSReadPrinterAddress function 19

Bodystat::BSReadProtocolInfo 19

Bodystat::BSReadProtocolInfo function 19

Bodystat::BSReadSerialNumber 20

Bodystat::BSReadSerialNumber function 20

Bodystat::BSReadStoredTestData 21

Bodystat::BSReadStoredTestData function 21

Bodystat::BSReadTestBodystat 21

Bodystat::BSReadTestBodystat function 21

Bodystat::BSResults 41

Bodystat::BSResults structure 41

Bodystat::BSSearchBTDevices 22

Bodystat::BSSearchBTDevices function 22

Bodystat::BSUnAuthenticateBTDevices 23

Bodystat::BSUnAuthenticateBTDevices function 23	BSGetBTStackInfo 121
Bodystat::BSWriteCurrentTime 24	BSGetDeviceFamily 122
Bodystat::BSWriteCurrentTime function 24	BSGetDeviceModelName 123
Bodystat::BSWriteCurrentTimeAs 24	BSGetDeviceModelNameShort 124
Bodystat::BSWriteCurrentTimeAs function 24	BSGetSdkLibraryVersion 124
Bodystat::BSWritePrinterAddress 25	BSIsBTAvailable 125
Bodystat::BSWritePrinterAddress function 25	BSOpenComport 126
Bodystat::ES_RAWDATA_ARRAYSIZE 52	BSReadCalibrationTime 127
Bodystat::ES_RAWDATA_ARRAYSIZE variable 52	BSReadCurrentTime 128, 129
Bodystat::ESCalculateNormals 26	BSReadModelVersion 130
Bodystat::ESCalculateNormals function 26	BSReadPrinterAddress 131
Bodystat::ESCalculateResults 26	BSReadProtocolInfo 132
Bodystat::ESCalculateResults function 26	BSReadSerialNumber 133
Bodystat::ESMeasurement 44	BSReadStoredTestData 133
Bodystat::ESMeasurement structure 44	BSReadTestBodystat 134
Bodystat::ESNormals 45	BSSearchBTDevices 136
Bodystat::ESNormals structure 45	BSUnAuthenticateBTDevices 137
Bodystat::ESRawData 47	BSWriteCurrentTime 138
Bodystat::ESRawData structure 47	BSWriteCurrentTimeAs 138, 139
Bodystat::ESReadStoredTestData 27	BSWritePrinterAddress 140
Bodystat::ESReadStoredTestData function 27	ES_RAWDATA_ARRAYSIZE 115
Bodystat::ESReadTestEquistat 28	ESCalculateNormals 141
Bodystat::ESReadTestEquistat function 28	ESCalculateResults 142
Bodystat::ESResults 48	ESReadStoredTestData 142
Bodystat::ESResults structure 48	ESReadTestEquistat 143
Bodystat::ImpData 49	BodystatAPI.BodystatAPI 81
Bodystat::ImpData::fFrequency 50	BodystatAPI.BodystatAPI.BS_RAWDATA_ARRAYSIZE 114
Bodystat::ImpData::fImpedance 50	BodystatAPI.BodystatAPI.BSAuthenticateBTDevice 115
Bodystat::ImpData::fPhaseAngle 51	BodystatAPI.BodystatAPI.BSAutoSetupBT 116
Bodystat::ImpData::ImZ 51	BodystatAPI.BodystatAPI.BSAutoSetupBT@out string@bool@IntPtr 57
Bodystat::ImpData::ReZ 51	BodystatAPI.BodystatAPI.BSAutoSetupBT@StringBuilder@int @bool@IntPtr 58
BodystatAPI 57	BodystatAPI.BodystatAPI.BSCalculateNormals 117
BodystatAPI class 81	BodystatAPI.BodystatAPI.BSCalculateNormals@BSMeasure ment@BSResults@out BSNormals 59
about BodystatAPI class 81	BodystatAPI.BodystatAPI.BSCalculateResults 118
BS_RAWDATA_ARRAYSIZE 114	BodystatAPI.BodystatAPI.BSCalculateResults@BSMeasure ment@out BSResults@BSAnalysisMode 59
BSAuthenticateBTDevice 115	BodystatAPI.BodystatAPI.BSCloseComport 60, 119
BSAutoSetupBT 116	BodystatAPI.BodystatAPI.BSConnect 60, 120
BSCalculateNormals 117	BodystatAPI.BodystatAPI.BSGetBTBodystatDevice 120
BSCalculateResults 118	
BSCloseComport 119	
BSConnect 120	
BSGetBTBodystatDevice 120	

BodystatAPI.BodystatAPI.BSGetBTBodystatDevice@out string@out string@out string@UInt16 77	BodystatAPI.BodystatAPI.BSMeasurement.pMultifreqData 92
BodystatAPI.BodystatAPI.BSGetBTBodystatDevice@StringBuilder@int@StringBuilder@int@StringBuilder@int@UInt16 78	BodystatAPI.BodystatAPI.BSMeasurement.tTestDate 93
BodystatAPI.BodystatAPI.BSGetBTStackInfo 121	BodystatAPI.BodystatAPI.BSMeasurement.ulDeviceSerialNumber 92
BodystatAPI.BodystatAPI.BSGetBTStackInfo@out string 60	BodystatAPI.BodystatAPI.BSNormals 93
BodystatAPI.BodystatAPI.BSGetBTStackInfo@StringBuilder@int 60	BodystatAPI.BodystatAPI.BSNormals structure 93
BodystatAPI.BodystatAPI.BSGetDeviceFamily 122	BodystatAPI.BodystatAPI.BSNormals.AerobicCapacity_H 62
BodystatAPI.BodystatAPI.BSGetDeviceModelName 123	BodystatAPI.BodystatAPI.BSNormals.AerobicCapacity_L 62
BodystatAPI.BodystatAPI.BSGetDeviceModelName@BSDeviceModel@out string 60	BodystatAPI.BodystatAPI.BSNormals.BFMI_H 63
BodystatAPI.BodystatAPI.BSGetDeviceModelName@BSDeviceModel@StringBuilder@int 61	BodystatAPI.BodystatAPI.BSNormals.BFMI_L 63
BodystatAPI.BodystatAPI.BSGetDeviceModelNameShort 124	BodystatAPI.BodystatAPI.BSNormals.BMI_H 63
BodystatAPI.BodystatAPI.BSGetDeviceModelNameShort@BSDeviceModel@out string 61	BodystatAPI.BodystatAPI.BSNormals.BMI_L 63
BodystatAPI.BodystatAPI.BSGetDeviceModelNameShort@BSDeviceModel@StringBuilder@int 61	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_HDL_H 63
BodystatAPI.BodystatAPI.BSGetSdkLibraryVersion 124	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_HDL_L 63
BodystatAPI.BodystatAPI.BSGetSdkLibraryVersion@out int@out int 62	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_LDL_H 63
BodystatAPI.BodystatAPI.BSIsBTAvailable 62, 125	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_LDL_L 63
BodystatAPI.BodystatAPI.BSMeasurement 84	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Total_H 63
BodystatAPI.BodystatAPI.BSMeasurement.fPA_50kHz 86	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Total_L 63
BodystatAPI.BodystatAPI.BSMeasurement.fWeight 87	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Triglycerides_H 64
BodystatAPI.BodystatAPI.BSMeasurement.fX_50kHz 87	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_Triglycerides_L 64
BodystatAPI.BodystatAPI.BSMeasurement.iActivity 87	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_VLDL_H 64
BodystatAPI.BodystatAPI.BSMeasurement.iAge 88	BodystatAPI.BodystatAPI.BSNormals.Cholesterol_VLDL_L 64
BodystatAPI.BodystatAPI.BSMeasurement.iDeviceModel 88	BodystatAPI.BodystatAPI.BSNormals.ECWPerNorm 64
BodystatAPI.BodystatAPI.BSMeasurement.iFrequencies 88	BodystatAPI.BodystatAPI.BSNormals.FatKg_H 64
BodystatAPI.BodystatAPI.BSMeasurement.iGender 89	BodystatAPI.BodystatAPI.BSNormals.FatKg_L 64
BodystatAPI.BodystatAPI.BSMeasurement.iHeight 89	BodystatAPI.BodystatAPI.BSNormals.FatPerc_H 64
BodystatAPI.BodystatAPI.BSMeasurement.iHip 89	BodystatAPI.BodystatAPI.BSNormals.FatPerc_L 64
BodystatAPI.BodystatAPI.BSMeasurement.iR_50kHz 90	BodystatAPI.BodystatAPI.BSNormals.FFMI_H 65
BodystatAPI.BodystatAPI.BSMeasurement.iWaist 90	BodystatAPI.BodystatAPI.BSNormals.FFMI_L 65
BodystatAPI.BodystatAPI.BSMeasurement.iZ_100kHz 90	BodystatAPI.BodystatAPI.BSNormals.Flexibility_H 65
BodystatAPI.BodystatAPI.BSMeasurement.iZ_200kHz 91	BodystatAPI.BodystatAPI.BSNormals.Flexibility_L 65
BodystatAPI.BodystatAPI.BSMeasurement.iZ_50kHz 91	BodystatAPI.BodystatAPI.BSNormals.GlucoseFasting_H 65
BodystatAPI.BodystatAPI.BSMeasurement.iZ_5kHz 91	BodystatAPI.BodystatAPI.BSNormals.GlucoseFasting_L 65
BodystatAPI.BodystatAPI.BSMeasurement.IITestDate 92	BodystatAPI.BodystatAPI.BSNormals.GlucoseNonFasting_H 65
	BodystatAPI.BodystatAPI.BSNormals.GlucoseNonFasting_L 65
	BodystatAPI.BodystatAPI.BSNormals.GripLeft 65
	BodystatAPI.BodystatAPI.BSNormals.GripRight 65
	BodystatAPI.BodystatAPI.BSNormals.ICWPerNorm 66

BodystatAPI.BodystatAPI.BSNormals.Illness_H 66	BodystatAPI.BodystatAPI.BSReadPrinterAddress@StringBuilder@int 70
BodystatAPI.BodystatAPI.BSNormals.Illness_L 66	
BodystatAPI.BodystatAPI.BSNormals.LeanKg_H 66	BodystatAPI.BodystatAPI.BSReadProtocolInfo 132
BodystatAPI.BodystatAPI.BSNormals.LeanKg_L 66	BodystatAPI.BodystatAPI.BSReadProtocolInfo@out byte@out byte@out byte 71
BodystatAPI.BodystatAPI.BSNormals.LeanPerc_H 66	
BodystatAPI.BodystatAPI.BSNormals.LeanPerc_L 66	BodystatAPI.BodystatAPI.BSReadSerialNumber 133
BodystatAPI.BodystatAPI.BSNormals.Nutrition_Norm 66	BodystatAPI.BodystatAPI.BSReadSerialNumber@out uint 71
BodystatAPI.BodystatAPI.BSNormals.TBW_H 66	BodystatAPI.BodystatAPI.BSReadStoredTestData 133
BodystatAPI.BodystatAPI.BSNormals.TBW_L 66	BodystatAPI.BodystatAPI.BSReadTestBodystat 134
BodystatAPI.BodystatAPI.BSNormals.TBWPerc_H 66	BodystatAPI.BodystatAPI.BSReadTestBodystat@out int@out int@out int@out int@out float@out float 71
BodystatAPI.BodystatAPI.BSNormals.TBWPerc_L 67	
BodystatAPI.BodystatAPI.BSNormals.TotalWeight_H 67	BodystatAPI.BodystatAPI.BSResults 99
BodystatAPI.BodystatAPI.BSNormals.TotalWeight_L 67	BodystatAPI.BodystatAPI.BSResults structure 99
BodystatAPI.BodystatAPI.BSNormals.TotalWeightMethod 67	BodystatAPI.BodystatAPI.BSResults.BCM 72
BodystatAPI.BodystatAPI.BSNormals.WaistHip_Norm 67	BodystatAPI.BodystatAPI.BSResults.BFMI 72
BodystatAPI.BodystatAPI.BSNormals.Wellness_H 67	BodystatAPI.BodystatAPI.BSResults.BMI 72
BodystatAPI.BodystatAPI.BSNormals.Wellness_L 67	BodystatAPI.BodystatAPI.BSResults.BMR 72
BodystatAPI.BodystatAPI.BSOpenComport 126	BodystatAPI.BodystatAPI.BSResults.BMRkg 72
BodystatAPI.BodystatAPI.BSOpenComport@string 67	BodystatAPI.BodystatAPI.BSResults.DryLW 72
BodystatAPI.BodystatAPI.BSOpenComport@StringBuilder@int 68	BodystatAPI.BodystatAPI.BSResults.ECW 73
BodystatAPI.BodystatAPI.BSRawData 96	BodystatAPI.BodystatAPI.BSResults.ECWperc 73
BodystatAPI.BodystatAPI.BSRawData.BSRawData@int 68	BodystatAPI.BodystatAPI.BSResults.EstAvg 73
BodystatAPI.BodystatAPI.BSRawData.iRecordArraySize 97	BodystatAPI.BodystatAPI.BSResults.FatKg 73
BodystatAPI.BodystatAPI.BSRawData.iTotalNumRecs 98	BodystatAPI.BodystatAPI.BSResults.FatPerc 73
BodystatAPI.BodystatAPI.BSRawData.New 97	BodystatAPI.BodystatAPI.BSResults.FFMI 73
BodystatAPI.BodystatAPI.BSRawData.record 98	BodystatAPI.BodystatAPI.BSResults.ICW 73
BodystatAPI.BodystatAPI.BSRawData.ulFirstTestNum 99	BodystatAPI.BodystatAPI.BSResults.ICWperc 73
BodystatAPI.BodystatAPI.BSReadCalibrationTime 127	BodystatAPI.BodystatAPI.BSResults.Illness 73
BodystatAPI.BodystatAPI.BSReadCalibrationTime@out DateTime 68	BodystatAPI.BodystatAPI.BSResults.LeanKg 73
BodystatAPI.BodystatAPI.BSReadCalibrationTime@out Int64 68	BodystatAPI.BodystatAPI.BSResults.LeanPerc 73
BodystatAPI.BodystatAPI.BSReadCurrentTime 128, 129	BodystatAPI.BodystatAPI.BSResults.Nutrition 74
BodystatAPI.BodystatAPI.BSReadCurrentTime@out DateTime@int 69	BodystatAPI.BodystatAPI.BSResults.TBW 74
BodystatAPI.BodystatAPI.BSReadCurrentTime@out Int64@int 69	BodystatAPI.BodystatAPI.BSResults.TBWPerc 74
BodystatAPI.BodystatAPI.BSReadModelVersion 130	BodystatAPI.BodystatAPI.BSResults.ThirdSpace 74
BodystatAPI.BodystatAPI.BSReadModelVersion@out BSDeviceModel@out byte@out byte@out byte@out byte@out byte 70	BodystatAPI.BodystatAPI.BSResults.TotalWeight 74
BodystatAPI.BodystatAPI.BSReadPrinterAddress 131	BodystatAPI.BodystatAPI.BSResults.WaistHip 74
BodystatAPI.BodystatAPI.BSReadPrinterAddress@out string 70	BodystatAPI.BodystatAPI.BSResults.Wellness 74
	BodystatAPI.BodystatAPI.BSSearchBTDevices 136
	BodystatAPI.BodystatAPI.BSSearchBTDevices@out int@out int@bool@IntPtr@UInt16 78
	BodystatAPI.BodystatAPI.BSUnAuthenticateBTDevices 137
	BodystatAPI.BodystatAPI.BSUnAuthenticateBTDevices@bool@IntPtr@UInt16 79

BodystatAPI.BodystatAPI.BSWriteCurrentTime 74, 138	int@out int@out int@out int@out int@out int@out float 76
BodystatAPI.BodystatAPI.BSWriteCurrentTimeAs 138, 139	BodystatAPI.BodystatAPI.ESResults 112
BodystatAPI.BodystatAPI.BSWritePrinterAddress 140	BodystatAPI.BodystatAPI.ESResults structure 112
BodystatAPI.BodystatAPI.BSWritePrinterAddress@string@int 74	BodystatAPI.BodystatAPI.ImpData 114
BodystatAPI.BodystatAPI.BSWritePrinterAddress@StringBuilder@int 75	BodystatAPI.BodystatAPI.ImpData structure 114
BodystatAPI.BodystatAPI.ES_RAWDATA_ARRAYSIZE 115	BodystatAPI.BSAnalysisMode 145
BodystatAPI.BodystatAPI.ESCalculateNormals 141	BodystatAPI.BSAnalysisMode enumeration 145
BodystatAPI.BodystatAPI.ESCalculateNormals@ESMeasurement@ESResults@out ESNormals 75	BodystatAPI.BSDeviceFamily 146
BodystatAPI.BodystatAPI.ESCalculateResults 142	BodystatAPI.BSDeviceFamily enumeration 146
BodystatAPI.BodystatAPI.ESCalculateResults@ESMeasurement@out ESResults 76	BodystatAPI.BSDeviceModel 147
BodystatAPI.BodystatAPI.ESMeasurement 101	BodystatAPI.BSDeviceModel enumeration 147
BodystatAPI.BodystatAPI.ESMeasurement.fPA_50kHz 103	BodystatAPI.BSError 149
BodystatAPI.BodystatAPI.ESMeasurement.iBCScore 103	BodystatAPI.BSError enumeration 149
BodystatAPI.BodystatAPI.ESMeasurement.iDeviceModel 104	BodystatAPI.BSGender 151
BodystatAPI.BodystatAPI.ESMeasurement.iHeight 104	BodystatAPI.BSGender enumeration 151
BodystatAPI.BodystatAPI.ESMeasurement.iZ_140kHz 104	BodystatAPI.BSMeasurement structure 84
BodystatAPI.BodystatAPI.ESMeasurement.iZ_16kHz 105	about BodystatAPI.BSMeasurement structure 84
BodystatAPI.BodystatAPI.ESMeasurement.iZ_200kHz 105	fPA_50kHz 86
BodystatAPI.BodystatAPI.ESMeasurement.iZ_24kHz 105	fWeight 87
BodystatAPI.BodystatAPI.ESMeasurement.iZ_280kHz 106	fX_50kHz 87
BodystatAPI.BodystatAPI.ESMeasurement.iZ_50kHz 106	iActivity 87
BodystatAPI.BodystatAPI.ESMeasurement.iZ_5kHz 106	iAge 88
BodystatAPI.BodystatAPI.ESMeasurement.IITestDate 107	iDeviceModel 88
BodystatAPI.BodystatAPI.ESMeasurement.tTestDate 107	iFrequencies 88
BodystatAPI.BodystatAPI.ESMeasurement.ulDeviceSerialNumber 107	iGender 89
BodystatAPI.BodystatAPI.ESNormals 108	iHeight 89
BodystatAPI.BodystatAPI.ESNormals structure 108	iHip 89
BodystatAPI.BodystatAPI.ESRawData 109	iR_50kHz 90
BodystatAPI.BodystatAPI.ESRawData.ESRawData@int 76	iWaist 90
BodystatAPI.BodystatAPI.ESRawData.iRecordArraySize 111	iZ_100kHz 90
BodystatAPI.BodystatAPI.ESRawData.iTotalNumRecs 111	iZ_200kHz 91
BodystatAPI.BodystatAPI.ESRawData.New 110	iZ_50kHz 91
BodystatAPI.BodystatAPI.ESRawData.record 112	iZ_5kHz 91
BodystatAPI.BodystatAPI.ESRawData.ulFirstTestNum 112	IITestDate 92
BodystatAPI.BodystatAPI.ESReadStoredTestData 142	pMultifreqData 92
BodystatAPI.BodystatAPI.ESReadTestEquistat 143	tTestDate 93
BodystatAPI.BodystatAPI.ESReadTestEquistat@out int@out	ulDeviceSerialNumber 92
	BodystatAPI.BSRawData structure 96
	about BodystatAPI.BSRawData structure 96
	iRecordArraySize 97
	iTotalNumRecs 98
	New 97

record 98	BODYSTATSDK_API 53
ulFirstTestNum 99	BODYSTATSDK_API macro 53
BodystatAPI.ESMeasurement structure 101	BodystatVBNetWinFormSample 217
about BodystatAPI.ESMeasurement structure 101	BodystatVBNetWinFormSample namespace 217
fPA_50kHz 103	Classes 217
iBCScore 103	BodystatVBNetWinFormSample.FrmBodystatMain 217
iDeviceModel 104	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtAuthenticateButton_Click 219
iHeight 104	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtAutoSetupButton_Click 220
iZ_140kHz 104	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtDeviceInfoButton_Click 220
iZ_16kHz 105	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtSearchButton_Click 221
iZ_200kHz 105	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtStackInfoButton_Click 222
iZ_24kHz 105	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtStatusButton_Click 222
iZ_280kHz 106	BodystatVBNetWinFormSample.FrmBodystatMain.btnBtUnauthenticateButton_Click 223
iZ_50kHz 106	BodystatVBNetWinFormSample.FrmBodystatMain.btnComConnectButton_Click 223
iZ_5kHz 106	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadCalibrationTimeButton_Click 224
lITestDate 107	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadCurrentTimeButton_Click 225
tTestDate 107	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadModelVersionButton_Click 226
ulDeviceSerialNumber 107	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadPrinterAddrButton_Click 227
BodystatAPI.ESRawData structure 109	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadProtocolInfoButton_Click 228
about BodystatAPI.ESRawData structure 109	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadSerialNumberButton_Click 229
iRecordArraySize 111	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadStoredTestDataButton_Click 230
iTotalNumRecs 111	BodystatVBNetWinFormSample.FrmBodystatMain.btnComReadTestMeasurementButton_Click 231
New 110	
record 112	
ulFirstTestNum 112	
BodystatCSWinForms 205	
BodystatCSWinForms namespace 205	
Classes 205	
BodystatCSWinForms.MainForm 205	
BodystatCSWinForms.MainForm.MainForm 206	
BodystatCSWinForms.MainForm.PopulateMeasurementList 207	
BodystatCSWinForms.ResultsForm 209	
BodystatCSWinForms.ResultsForm.InitialiseFromMeasurement 210	
BodystatCSWinForms.ResultsForm.InitialiseSampleMeasurement 211	
BodystatCSWinForms.ResultsForm.PopulateListCtrl 212	
BodystatCSWinForms.ResultsForm.ResultsForm 209	
BodystatCSWinForms.ResultsForm.UseBSMeasurement 215	
BodystatSDK.cs 80	
BodystatSDK.h 55	

BodystatVBNetWinFormSample.FrmBodystatMain.btnComWriteCurrentTimeButton_Click 232	BSAnalysisModeBoth enumeration member 30, 145
BodystatVBNetWinFormSample.FrmBodystatMain.btnResultsButton_Click 233	BSAnalysisModeHN enumeration member 30, 145
BodystatVBNetWinFormSample.FrmBodystatMain.FrmBodystatMain_Load 234	BSAnalysisModeNone enumeration member 30, 145
BodystatVBNetWinFormSample.FrmBodystatMain.iTimeout 219	BSDualScan_DIU enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatMain.lstviewTestMeasurementList_DoubleClick 235	BSDualScan_Family enumeration member 30, 146
BodystatVBNetWinFormSample.FrmBodystatMain.PopulateMeasurementList 235	BSDualScan_OPTO enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatResults 237	BSEquistatLite_BT enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatResults.bsmMeasurement 238	BSEquistatLite_Family enumeration member 30, 146
BodystatVBNetWinFormSample.FrmBodystatResults.btnCalculateResultsButton_Click 239	BSEquistatPro_BT enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatResults.btnCopytoclipboardButton_Click 240	BSEquistatPro_Family enumeration member 30, 146
BodystatVBNetWinFormSample.FrmBodystatResults.bUseActualMeasurement 238	BSFemale enumeration member 35, 151
BodystatVBNetWinFormSample.FrmBodystatResults.dtTestDate 238	BSMale enumeration member 35, 151
BodystatVBNetWinFormSample.FrmBodystatResults.FrmBodystatResults_Load 242	BSMDD_BT enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatResults.InitialiseFromMeasurement 242	BSMDD_DIU enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatResults.InitialiseSampleMeasurement 243	BSMDD_Family enumeration member 30, 146
BodystatVBNetWinFormSample.FrmBodystatResults.PopulateListCtrl 244	BSMDD_OPTO enumeration member 32, 147
BodystatVBNetWinFormSample.FrmBodystatResults.UseBSMeasurement 247	BSMultiScan_DIU enumeration member 32, 147
BS1500_BT enumeration member 32, 147	BSMultiScan_Family enumeration member 30, 146
BS1500_DIU enumeration member 32, 147	BSMultiScanTouch enumeration member 32, 147
BS1500_Family enumeration member 30, 146	BSQuadScan_BT enumeration member 32, 147
BS1500_OPTO enumeration member 32, 147	BSQuadScan_DIU enumeration member 32, 147
BSAnalysisModeBC enumeration member 30, 145	BSQuadScan_Family enumeration member 30, 146
	BSQuadScanTouch enumeration member 32, 147
	BSUnknown enumeration member 32, 147
	BSUnknown_Family enumeration member 30, 146
	C
	C# .NET Wrapper 56
	Files 79
	CBodystatMFCDlg 156
	CBodystatMFCDlg class 156
	about CBodystatMFCDlg class 156
	CBodystatMFCDlg 158
	DoDataExchange 165
	m_bsRawData 159
	m_btnAuthenticate 159
	m_btnAutoSetup 159
	m_btnConnect 160
	m_btnDisconnect 160
	m_btnReadCalibrationTime 160
	m_btnReadCurrentTime 161
	m_btnReadModelVersion 161

m_btnReadPrinterAddr 161	CBodystatMFCDlg::m_btnDisconnect 160
m_btnReadProtocolInfo 162	CBodystatMFCDlg::m_btnReadCalibrationTime 160
m_btnReadSerialNumber 162	CBodystatMFCDlg::m_btnReadCurrentTime 161
m_btnReadStoredTestData 162	CBodystatMFCDlg::m_btnReadModelVersion 161
m_btnReadTestMeasurement 162	CBodystatMFCDlg::m_btnReadPrinterAddr 161
m_btnUnauthenticate 163	CBodystatMFCDlg::m_btnReadProtocolInfo 162
m_btnWriteCurrentTime 163	CBodystatMFCDlg::m_btnReadSerialNumber 162
m_hIcon 163	CBodystatMFCDlg::m_btnReadStoredTestData 162
m_iTimeout 164	CBodystatMFCDlg::m_btnReadTestMeasurement 162
m_strComPort 164	CBodystatMFCDlg::m_btnUnauthenticate 163
m_TestMeasurementListCtrl 164	CBodystatMFCDlg::m_btnWriteCurrentTime 163
OnClickedBtAuthenticateButton 165	CBodystatMFCDlg::m_hIcon 163
OnClickedBtAutoSetupButton 166	CBodystatMFCDlg::m_iTimeout 164
OnClickedBtDeviceInfoButton 167	CBodystatMFCDlg::m_strComPort 164
OnClickedBtSearchButton 167	CBodystatMFCDlg::m_TestMeasurementListCtrl 164
OnClickedBtStackInfoButton 168	CBodystatMFCDlg::OnClickedBtAuthenticateButton 165
OnClickedBtStatusButton 169	CBodystatMFCDlg::OnClickedBtAutoSetupButton 166
OnClickedBtUnauthenticateButton 169	CBodystatMFCDlg::OnClickedBtDeviceInfoButton 167
OnClickedComConnectButton 170	CBodystatMFCDlg::OnClickedBtSearchButton 167
OnClickedComReadCalibrationTimeButton 170	CBodystatMFCDlg::OnClickedBtStackInfoButton 168
OnClickedComReadCurrentTimeButton 171	CBodystatMFCDlg::OnClickedBtStatusButton 169
OnClickedComReadModelVersionButton 172	CBodystatMFCDlg::OnClickedBtUnauthenticateButton 169
OnClickedComReadPrinterAddrButton 173	CBodystatMFCDlg::OnClickedComConnectButton 170
OnClickedComReadProtocolInfoButton 174	CBodystatMFCDlg::OnClickedComReadCalibrationTimeButton 170
OnClickedComReadSerialNumberButton 175	CBodystatMFCDlg::OnClickedComReadCurrentTimeButton 171
OnClickedComReadStoredTestDataButton 175	CBodystatMFCDlg::OnClickedComReadModelVersionButton 172
OnClickedComReadTestMeasurementButton 176	CBodystatMFCDlg::OnClickedComReadPrinterAddrButton 173
OnClickedComWriteCurrentTimeButton 177	CBodystatMFCDlg::OnClickedComReadProtocolInfoButton 174
OnClickedResultsButton 178	CBodystatMFCDlg::OnClickedComReadSerialNumberButton 175
OnDbClickTestMeasurementList 179	CBodystatMFCDlg::OnClickedComReadStoredTestDataButton 175
OnInitDialog 179	CBodystatMFCDlg::OnClickedComReadTestMeasurementButton 176
OnPaint 180	CBodystatMFCDlg::OnClickedComWriteCurrentTimeButton 177
OnQueryDragIcon 181	CBodystatMFCDlg::OnClickedResultsButton 178
OnSysCommand 181	CBodystatMFCDlg::OnDbClickTestMeasurementList 179
PopulateMeasurementList 182	CBodystatMFCDlg::OnInitDialog 179
ShowResultsDialog 184	
CBodystatMFCDlg::CBodystatMFCDlg 158	
CBodystatMFCDlg::DoDataExchange 165	
CBodystatMFCDlg::m_bsRawData 159	
CBodystatMFCDlg::m_btnAuthenticate 159	
CBodystatMFCDlg::m_btnAutoSetup 159	
CBodystatMFCDlg::m_btnConnect 160	

- CBodystatMFCDlg::OnPaint 180
 - CBodystatMFCDlg::OnQueryDragIcon 181
 - CBodystatMFCDlg::OnSysCommand 181
 - CBodystatMFCDlg::PopulateMeasurementList 182
 - CBodystatMFCDlg::ShowResultsDialog 184
 - CBodystatResultsDlg 185
 - CBodystatResultsDlg class 185
 - ~CBodystatResultsDlg 186
 - about CBodystatResultsDlg class 185
 - CBodystatResultsDlg 186
 - CopyToClipboard 194
 - DoDataExchange 196
 - GetSystemErrorMessage 197
 - m_bModeBodyComposition 187
 - m_bModeHydration 188
 - m_dtTestDate 188
 - m_fPhaseAngle50K 188
 - m_fReactance50K 189
 - m_fWeight 189
 - m_iActivityLevel 189
 - m_iAge 190
 - m_iGender 190
 - m_iHeight 190
 - m_iHip 191
 - m_iImpedance100K 191
 - m_iImpedance200K 191
 - m_iImpedance50K 192
 - m_iImpedance5K 192
 - m_iModel 192
 - m_iResistance50K 193
 - m_iWaist 193
 - m_ModelComboCtrl 193
 - m_ResultsListCtrl 194
 - OnClickedCalculateResultsButton 198
 - OnClickedCopyToClipboardButton 200
 - OnInitDialog 200
 - PopulateListCtrl 201
 - CBodystatResultsDlg::~CBodystatResultsDlg 186
 - CBodystatResultsDlg::CBodystatResultsDlg 186
 - CBodystatResultsDlg::CopyToClipboard 194
 - CBodystatResultsDlg::DoDataExchange 196
 - CBodystatResultsDlg::GetSystemErrorMessage 197
 - CBodystatResultsDlg::m_bModeBodyComposition 187
 - CBodystatResultsDlg::m_bModeHydration 188
 - CBodystatResultsDlg::m_dtTestDate 188
 - CBodystatResultsDlg::m_fPhaseAngle50K 188
 - CBodystatResultsDlg::m_fReactance50K 189
 - CBodystatResultsDlg::m_fWeight 189
 - CBodystatResultsDlg::m_iActivityLevel 189
 - CBodystatResultsDlg::m_iAge 190
 - CBodystatResultsDlg::m_iGender 190
 - CBodystatResultsDlg::m_iHeight 190
 - CBodystatResultsDlg::m_iHip 191
 - CBodystatResultsDlg::m_iImpedance100K 191
 - CBodystatResultsDlg::m_iImpedance200K 191
 - CBodystatResultsDlg::m_iImpedance50K 192
 - CBodystatResultsDlg::m_iImpedance5K 192
 - CBodystatResultsDlg::m_iModel 192
 - CBodystatResultsDlg::m_iResistance50K 193
 - CBodystatResultsDlg::m_iWaist 193
 - CBodystatResultsDlg::m_ModelComboCtrl 193
 - CBodystatResultsDlg::m_ResultsListCtrl 194
 - CBodystatResultsDlg::OnClickedCalculateResultsButton 198
 - CBodystatResultsDlg::OnClickedCopyToClipboardButton 200
 - CBodystatResultsDlg::OnInitDialog 200
 - CBodystatResultsDlg::PopulateListCtrl 201
 - ComErrorConfiguringPort enumeration member 33, 149
 - ComErrorConnecting enumeration member 33, 149
 - ComErrorOpeningPort enumeration member 33, 149
 - ComErrorProcessingCommand enumeration member 33, 149
 - ComErrorReading enumeration member 33, 149
 - ComErrorUnsupportedCommand enumeration member 33, 149
 - ComErrorUserAborted enumeration member 33, 149
 - ComErrorWriting enumeration member 33, 149
- ## F
- FAQs 152
 - FrmBodystatMain class 217
 - about FrmBodystatMain class 217
 - btnBtAuthenticateButton_Click 219
 - btnBtAutoSetupButton_Click 220
 - btnBtDeviceInfoButton_Click 220
 - btnBtSearchButton_Click 221
 - btnBtStackInfoButton_Click 222

btnBtStatusButton_Click 222
 btnBtUnauthenticateButton_Click 223
 btnComConnectButton_Click 223
 btnComReadCalibrationTimeButton_Click 224
 btnComReadCurrentTimeButton_Click 225
 btnComReadModelVersionButton_Click 226
 btnComReadPrinterAddrButton_Click 227
 btnComReadProtocolInfoButton_Click 228
 btnComReadSerialNumberButton_Click 229
 btnComReadStoredTestDataButton_Click 230
 btnComReadTestMeasurementButton_Click 231
 btnComWriteCurrentTimeButton_Click 232
 btnResultsButton_Click 233
 FrmBodystatMain_Load 234
 iTimeout 219
 IstviewTestMeasurementList_DoubleClick 235
 PopulateMeasurementList 235
 FrmBodystatResults class 237
 about FrmBodystatResults class 237
 bsmMeasurement 238
 btnCalculateResultsButton_Click 239
 btnCopytoclipboardButton_Click 240
 bUseActualMeasurment 238
 dtTestDate 238
 FrmBodystatResults_Load 242
 InitialiseFromMeasurement 242
 InitialiseSampleMeasurement 243
 PopulateListCtrl 244
 UseBSMeasurement 247

I

ImpData structure 49
 about ImpData structure 49
 fFrequency 50
 fImpedance 50
 fPhaseAngle 51
 ImZ 51
 ReZ 51
 InvalidActivity enumeration member 33, 149
 InvalidAge enumeration member 33, 149
 InvalidDeviceFamily enumeration member 33, 149
 InvalidDeviceMode enumeration member 33, 149

InvalidDeviceModel enumeration member 33, 149
 InvalidGender enumeration member 33, 149
 InvalidHeight enumeration member 33, 149
 InvalidHip enumeration member 33, 149
 InvalidImpedanceZ100 enumeration member 33, 149
 InvalidImpedanceZ200 enumeration member 33, 149
 InvalidImpedanceZ5 enumeration member 33, 149
 InvalidImpedanceZ50 enumeration member 33, 149
 InvalidParamSupplied enumeration member 33, 149
 InvalidPhaseAnglePA50 enumeration member 33, 149
 InvalidReactanceX50 enumeration member 33, 149
 InvalidResistanceR50 enumeration member 33, 149
 InvalidWaist enumeration member 33, 149
 InvalidWeight enumeration member 33, 149

M

M_PI 54
 M_PI macro 54
 MainForm class 205
 about MainForm class 205
 MainForm 206
 PopulateMeasurementList 207
 MainForm.cs 216

N

NoError enumeration member 33, 149

R

Redistribution 154
 ResultsForm class 209
 about ResultsForm class 209
 InitialiseFromMeasurement 210
 InitialiseSampleMeasurement 211
 PopulateListCtrl 212
 ResultsForm 209
 UseBSMeasurement 215
 ResultsForm.cs 216

S

Samples 155

V

Visual Basic .NET wrapper 80

Classes 81

Structs, Records, Enums 145